

The image is a cover page for a manual. It is divided into four quadrants. The top-left quadrant is a solid bright blue rectangle. The top-right quadrant features a dark blue background with a close-up, artistic photograph of blue flower petals. The bottom-left quadrant is a solid white rectangle. The bottom-right quadrant features a close-up, artistic photograph of yellow flower petals. The text is positioned in the top-right and bottom-right quadrants.

breco™
Head Massager

breco 1180

Operating Instructions

simply relax

Catalogue

Thank you for using Breo. Please read the manual carefully before using it and keep it properly for future reference.

Catalogue

Product theory and applicable range	2
Product Characteristic	3
Product component scheme	4
Preparation before using	5
Remote control operation	6
Instrument operate flow plan	7
Cautions	8
Maintenance@ Technology Index	9
After-sales service and warranty	10

Product principle /Apply range

Product principle

INOW1180 is a head massager with intelligentized airpressure, through different massaging, heating vibrating , Iknow1180 that adopts computer chip control technology and far infrared can help increase oxygen supply of brain, improve blood circulation of brain, relax the pressure of brain and ease the tiredness of your brain, promote metabolism and keep the health to brain. The device also adopts natural music , relax music and brain wave music to activate the potential power of brain.

Product apply range

Which can help from insomnia, headache , and over used brain. It is used to be the aid of cure, and relax the pressure of your brain, prevent the disease of brain and keep the health to brain.



undergraduate
secondary school students



insomnia



people who overuse brain



accountant



designers Secretary teachers



People who long-time
use brain



drivers




headache

Product Specification

Product Specification


- This product is a multifunctional high-tech device for health care of brain
With the help of computer chip control technique, this product gives an all-way humanization simulated control.
- The product successfully applies the concept of air pressing to the protection and health keeping of the brain. The air pressing can give the brain an overall massage under the control of the minor computer. This process can lessen the stress of the brain, easy the tiredness of the brain.
- With a built-in three-dimensional magnetic field, when you use the product, the constant movement of the field gives an effective massage to the main points of the brain, generating micro-current in blood vessels around the brain and activating the blood circulation.
- This product adopts multi-vibration, the brain can active according to the vibration, it makes the brain obtain effective exercise, ease the tiredness of brain.
- This product can adjust according to the brain's size
- This product including natural music, relax music and brain wave music, it can stimulate the potential power of the brain, improve the potential power.
- This product's timer including 5/10/15 minutes, it 's convenient for customers
- This product has obtained and passed many domestic and international patent

Point function

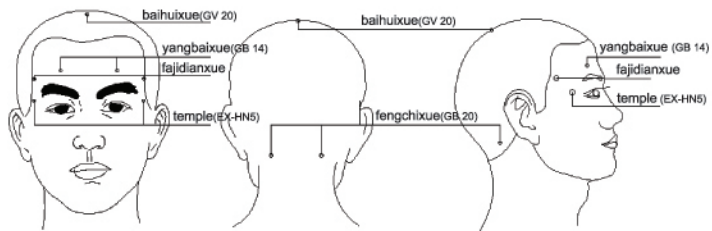


Product introduction/structure

Product component scheme

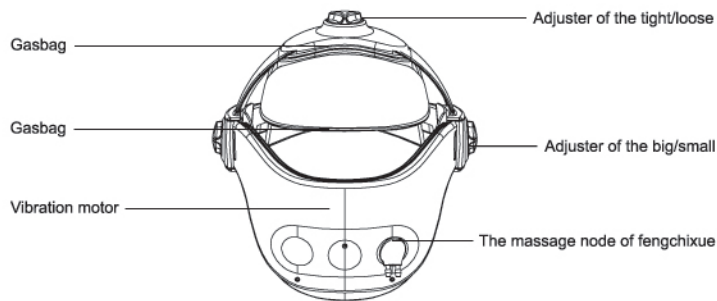

Structure

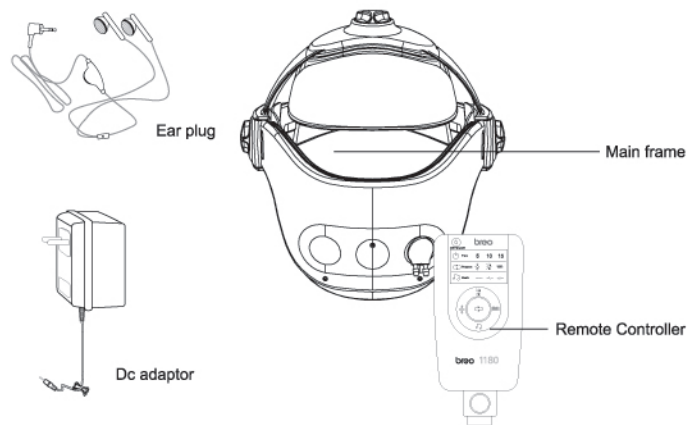


How to use

Preparations before using

Battery installation

1. press the the back of the clip, take out of it form the controller
2. open the battery cover from the back of controller and insert 4pcs AAA batteries properly
3. close the battery cover of the controller
4. check the anode and cathode of battery. Keep the them correctly installed; check the battery cover and keep it tight, the anode and cathode, please refer to the right picture, alkaline battery is suggested.

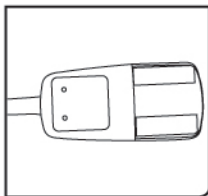


! cautions

check the anode and cathode of battery. Keep the them correctly installed; check the battery cover and keep it tight, the anode and cathode, please refer to the right picture, alkaline battery is suggested.

Adapter installation

1. please use the DC6V adapter that has been prepared by Breeze
2. plug the DC Pin into the DC socket of the controller
3. connect the adaptor to AC power socket properly

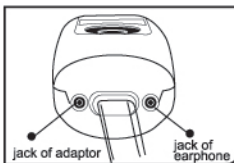

connect the earphone to the socket




the front image effect of wearing the head massager



the side face of wearing the head massager

Remote control operation

Remote control operation


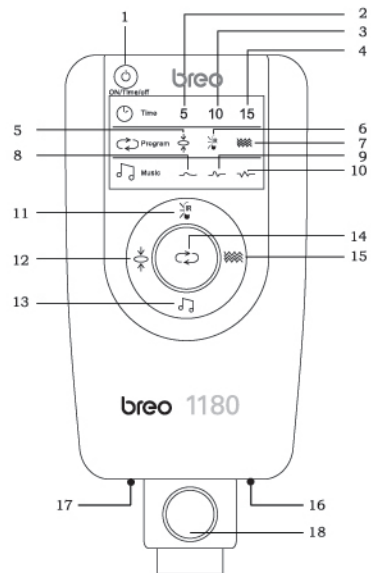


button indication

- | | | |
|---------------------------------|--|---|
| 1: ON/OFF/TIMER | 8: Natural music indicator light | 14: Massaging modes selection |
| 2: 5 Minutes | 9: Relax music indicator light | 15: Vibration modes selection |
| 3: 10 Minutes | 10: a brain wave music indicator light | 16: Earphone plug |
| 4: 15 Minutes | 11: Far infrared selection | 17: DC plug(DC:6V, 1000Ma)
<small>(Adaptor specification: DC6V 1000mA)</small> |
| 5: Air pressing indicator light | 12: Air pressure modes selection | 18: Air valve |
| 6: Far infrared indicator light | 13: Music modes selection
<small>(Press 1 second to select music
Press 3 second to off music)</small> | |
| 7: Vibration indicator light | | |

Instrument operate flow plan

Instrument operate flow plan

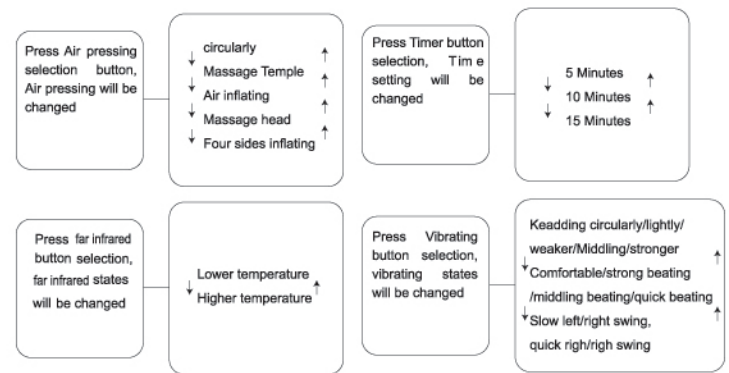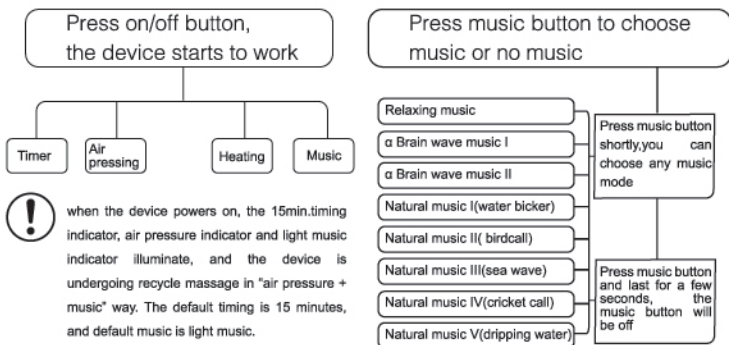

Notice

! Notice

- Please read the manual carefully before operation
- Operation should according to the cautions and instructions
- Please don't use the product when suffering from brain cancer, or have some serious brain diseases.
- Please take off contact lens, barrette and hair accouterment before using
- Please do not wash the head massager by water. Clean them lightly by using alcohol if necessary
- Do not grasp the wire only when using, it's forbid to dismantle the product

! User notice

- Power on after wearing the head massager, avoid to spoil the gasbag
- Please take off contact lens, barrette and hair accouterment before using
- Please do not use the product under high temperature and humidity moist environment
- Please avoid strong sunlight and corrosion gas
- Please do not use this product when outdoor walking or driving
- Please do not drop the main unites to the ground or dip water into it
- Please use AAA batteries or matched DC adaptor, when use outer electrical source, the batteries are not allowed to put into the batteries box. Take the batteries out of the machine when it is power off, avoid the trouble arouse from the liquid of batteries.

How to wear main unit

Maintenance & Technology

If any trouble occur when using the product, please check it according to the below table, if the product still can not get rid of the trouble or can not be operated ideal, please contact the product service center to require maintenance.

⚠ Please do not dismantle or repair this product optionally

malfunction maintenance	possible cause of malfunction	remedy way
can not turn on the power supply	non-damage of wire	check the condition of wire connection
	battery's power is insufficient	if it's true, please change the battery
	does the installation of battery is correct	assemble the anode and cathode of the battery correct, and keep the battery fixed in the battery box
	does the connection of DC adaptor is correct	please assure the correct connection
indicator's light is insufficient	does the battery's power is weak	if it's true, please change the battery
the air pressure and heating intensity is too weak	does the battery's power is weak	if it's true, please change the battery, and keep the battery fixed in the battery box

Technology index

work environment

Temperature: 5°C-40°C
Relative humidity: ≤80%

specification

Air pressure: 860hpa-1060hpa
Intensity of magnetic field: 50mT-250mT
Heating temperature: 30°C-45°C

power supply

DC: 6V
Power: <6W

dimension of figuration

main frame dimension:
230mm×220mm×234mm (length* width*height)
controller dimension:
141mm×68.5mm×44mm (length* width*height)
net weight: 0.625kg

After-sales service and warranty

After-sales service

Power inspection

- Inspect the power; do not start the device for use if the power pressure exceeds the prescribed scope.
- Inspect eyepatch and cable frequently. If any damage or break is found, it shall stop using immediately and send to the professional after-sales service to handle.

Maintenance

- Clean the handle controller when the device is turned off.
- Place the device at safe place after using, and do not put any heavy on it.
- Take measures to prevent humidity, mice and pollution.
- Do not tear down the device without permission when fault occurs but should resort to professional technical personnel.