




BM190 Massage Chair

OPERATION INSTRUCTIONS

SAVE THESE INSTRUCTIONS

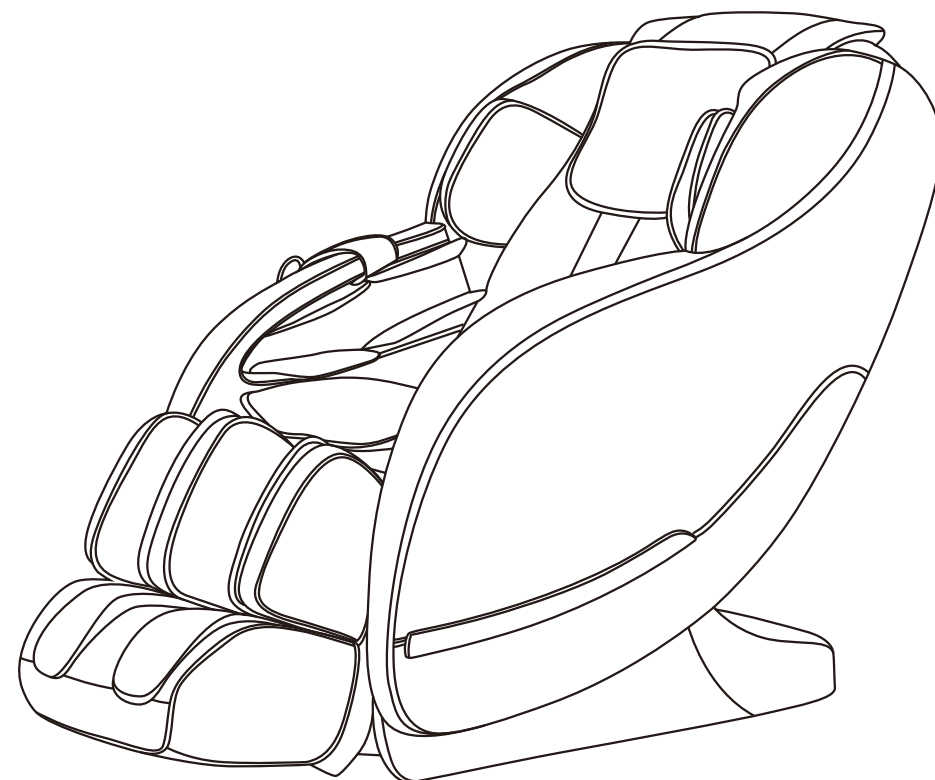

Thanks for purchasing this product. Please read this manual instruction carefully in order to know how to use this product correctly before operation. Please pay more attention to some important safety warning and keep this manual book properly for future use.

Note: Our Company reserves the right to revise the design and description of this product without any further notice. Please refer to the actual color of this product.

CONTENTS

Safety and maintenance	02
Name and function of components	04
Function	05
Method of usage	06
Product specification	13

Safety and maintenance

1. Important Safety Warning.



- Don't allow children to touch moveable part of this product.
- Please use well-grounding power supply suited to this product.
- Please pull out the plug after use or before cleaning to avoid injury or Damage to the product.
- Please operate this product according to this manual instruction.
- Don't use accessories which are not recommended.
- Don't use this product outdoors.
- Please read this manual instruction carefully before operation.
- Any other usage that is not listed in this manual instruction is forbidden.
- 20 minutes of usage each time is recommended.
- Please don't use it if the leather is broken or damaged.
- Please don't use it if the cover or leather is damaged.
- Please don't use it if the radiator hole is covered or blocked.
- Please don't drop anything into this product.
- Please don't fall asleep while using this product.
- Please don't use it if drunk or feel unwell.
- Please don't use it within one hour after having meal.
- Please don't make the massage function too strong for avoiding Injury.

2. Environment for usage.

- Please don't use it under high temperature and moist environment such as bathroom.
- Please don't use it immediately while the environmental temperature changes sharply.
- Please don't use it under heavy-dusty or caustic environment.
- Please don't use it where there is not enough space or not good ventilation.

3. People who are not suitable to use this product.

- People who are suffering from osteoporosis are not suitable to use this product.
- People who are suffering from heart disease or carrying electronic medical equipment such as heart pacemaker are not suitable to use this product.
- People who are suffering from fever are not suitable to use this product.
- Pregnant women or women that are in menses period are not suitable to use this product.
- People who got injured or skin disease are not suitable to use this product.
- Children under 14-year old or people mentally unbalanced without being Supervised are not allowed to use this machine.
- People who are told by doctor to have rest or people who feel unwell are not suitable to use this product.
- People whose bodies are wet are forbidden to use this product.
- The appliance has a heated surface. persons insensitive to Heat must be careful when using the appliance.

4. Safety.

- Check the voltage if it is suitable to the specification of this product.
- Don't use wet hand to pull the plug.
- Don't make water go inside this product to avoid electric shock or cause damage to this product.

Safety and maintenance

- Don't damage the wires or change the circuit of this product.
- Don't use wet cloth to clean the electric parts such as switch and plug.
- Be away from this product under power cut state to avoid injury if the power comeback suddenly.
- Stop using this product while it works abnormal, and please consult local agent immediately.
- Stop using this product if you feel unwell, and please consult health care doctor.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- If the supply cord is damaged, it must be replaced by the manufacturer, its Service agent or similarly qualified persons in order to avoid a hazard.


5. Maintenance.

- This product should be maintained by certain appointed agents, users are prohibited to disassemble or maintain by themselves.
- Please don't forget to cut the power after usage.
- Don't use this product if the socket looses.
- If this product will be left unused for long, please curl up the wires and reserve this product in dry and dustless environment.
- Don't reserve this product under high temperature or near fire and avoid longtime direct sunshine.
- If the detachable cord / supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a Hazard.
- Please clean this product by dry cloth. Don't use thinner, benzene or alcohol.
- The mechanical components of this product are specially designed and made, it is not necessary to be specially maintained.
- Don't use edged res to thrust this product.
- Don't roll or pull this product above uneven ground, it should be lifted before moving.
- Please use it intermittently; don't make this product work continuously for long.


6. Solution of ordinary malfunctions.

- It is nature movement voice if the motor sounds while in usage.
- If the controller can't work normally, please check the plug and socket if they are connected strongly and make sure the switch is turned on.
- If the rated working time is over, the product will automatically activate the power-off button; If this product works continuously for long, the temperature protector will make the product automatically power-off, it should be used again after half an hour's rest.

Nom et fonctions des composantes

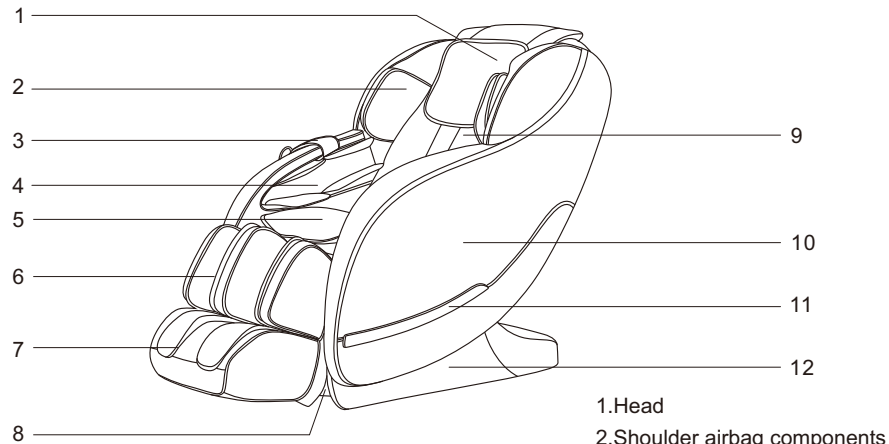

- 1.Head
- 2.Shoulder airbag components
- 3.Arm airbag components
- 4.Seat cushion
- 5.Calf rest
- 6.Foot rest
- 7.Calf rest end sealing plate
- 8.Back cushion
- 9.Armrest
- 10.Armrest decoration
- 12.Side cover

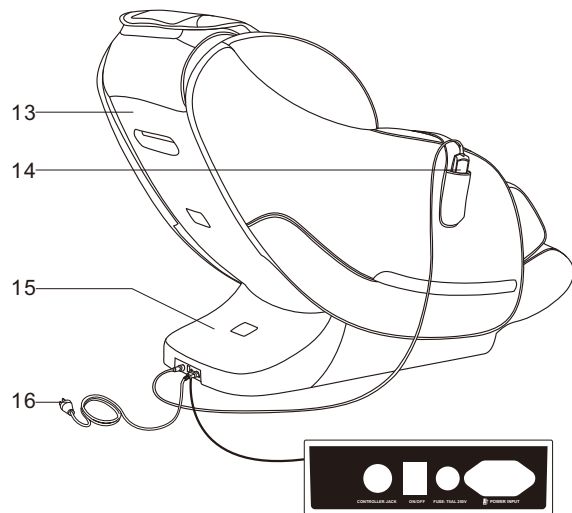

- 13.Back cover
- 14.Remote controller
- 15.PCB box cover
- 16.Power cord
- 17.Controller cord socket
- 18.Switch
- 19.Fuse box
- 20.Power cord socket

Function

- This product is designed with a set of four- wheel driven mechanical massage hands that can move up and down while stretching back and forth.
- Designed with five emulation massage methods of shiatsu, kneading, flapping, knocking, simultaneous kneading and flapping
- Set with six kinds of characterized auto massage function: comfort, ease, relax, neck& shoulder massage, back& waist massage, fast improving experience massage.
- Set with 1 memory function
- Manually select the type of upper body massage with the five massage point options which include fixed position, partial position, full back, neck& shoulder, back& waist. With three level of speeds for adjustable. In state of fixed position, the massage ball can move up and down to adjust the massage point accurately. In the state of tapping, shiatsu and knocking, the width between the two kneading balls is adjustable with three levels: wide, medium and narrow.
- Air pressure massage function with three modes: full body air pressure, arms air pressure, leg air pressure. with 3 intensity levels adjustable.
- Leg designed with roller massage function. At the bottom of foot, two pairs of papillae on toe root rotating massage; on the arch and heel mastoid skin scraping massage.
- The back heating function uses carbon fiber as an infrared heating source.
- The footrest can be extended to suit different leg lengths.
- set the calf and the back of a chair lift regulating function alone. The calf will be auto separate the mechanism when it touch foreign body.
- zero space, reach the maximum space saving, and design with three auto-relining mode.
- L shape back track fitting your body curves , can realize the zero gravity when back down
- The chair will be reset when switch off. The massage hands will be reset as well.

Method of usage

Note: The specific displays refer to the actual object of this product!



06

Method of usage


1. Switch power supply on



2. Start massage

- ① Press power button, start massage.
- ② Automatic shoulder position and massage acupressure point detection, waiting for the completion of testing.
- ③ If there is difference with the massage balls and the actual shoulder position, you can press the up and down button to adjust the ball to a suitable position. Shoulder position has 10 grades adjustable. It will run automatic massage if there is no adjustment on shoulder position within 10 seconds.

Button	Description	Display
"Up" Shoulder Position	Press to move the mechanical hands higher to a suitable position.	
"Down" Shoulder Position	Press to move the mechanical hands lower. to a suitable position	

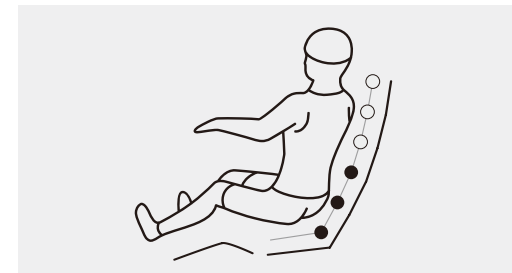

During the detection, the controller indicator light with "PP" blink, press any button of the 6 auto modes (comfort, ease, relax, neck & shoulder massage, back & waist massage, fast improving experience massage), or choose memory program, it will enter to the massage situation which stored on the chair.

④. Choose massage function to massage

- Memory function massage: press memory it will revert to the previous memory massage mode.

Button	Description	Display
Memory program		Previous memory massage mode.

07

Method of usage

· Characterized auto massage

Function	Button	Description	Display
Characterized auto massage		Deeply massage on acupressure points, effectively alleviate pain	A1
		Loosen body physique	A2
		Promote blood circulation	A3
		Focused on neck and shoulder massage	A4
		Focused on back and waist massage	A5
		Fast to experience all the massage function, 5minutes	A0

· Upper body manual massage

Function	Button	Description	Display
Technique		Kneading: speed	
		Tapping: speed, width	
		Kneading & flapping: speed	
		Shiatsu: Speed, width	
		Knocking: speed, width	
Speed		Speed: 3 adjustable levels	Speed
Width		Width: 3 adjustable levels	Width

Method of usage

Function	Button	Description	Display
Back Stretch		The massage hands will stay in one fixed position	H0
		The massage hands will move up and down in a partial position	H1
		The massage hands will move up and down the entire back area	H2
		The massage hands will move up and down on the neck and shoulder area	H3
		The massage hands will move up and down on the back and waist area	H4
Massage position adjust (in fixed position, partial modes you can adjust the massage position)		Massage hand move up	
		Massage hand move down	

· Air pressure function massage






Function	Button	Description	Display
Air pressure mode		Air pressure massage on arm, calf and foot part	
		Air pressure massage on arm	
		Air pressure on calf and foot	
Intensity		Intensity : high, middle, lower.	High Middle Lower

· Other manual function

Function	Button	Description	Display
Foot roller		Sole roller: sole roller switch, with 3 speed levels available	ROLLERS
Heating		Far infrared back heating, has warm feeling when it start to work 3 minutes later	HEATING


Method of usage

- Backrest and calf position adjust

Function	Button	Description
Calf rise		Press and hold this button to slowly lift the calf rest. Release the button to stop. If the calf reached max. Level, it will stop with "di...di" alarming
Calf drop		Press and hold this button to slowly lower the calf rest. If the calf reached max. Level, it will stop with "di...di" alarming
Backrest rise		Press and hold this button to slowly lift the back rest. If the backrest reached max. Level, it will stop with "di...di" alarming
Backrest drop		Press and hold this button to slowly lower the back rest. Release the button to stop. If the backrest reached max. Level, it will stop with "di...di" alarming
Auto lay down mode		With three modes circulation function, when press the button, turn into auto lay down mode.

3. Turn off power, stop massage

- Press power button during massage progress, all massage function will be stopped immediately, meanwhile you can set up memory function, the backrest and foot part will restore. Or switch off the chair when massage time finished, you can press the memory function button to storage the massage situation, backrest and foot part will be restore.
- Memory function: display with "EA" flashing instructions, press the button this time, it start to storage the massage mode doing when stop the chair. The remote screen will display E1. Or re-press the switch button to stop the chair with the massage hand, backrest and foot part restore. Or waiting for 10 minutes with chair auto stopped, the massage hand, backrest, foot part will restore.

Function	Button	Display
Memory Function		E1

- Cut off the power of the whole machine.
Figure (Switch power off)



Method of usage

4. Calf installation

- The install hole of rail shelf (chart 1).
- Place the isolation sleeve on both sides of the calf hook. (Outside) (chart2).

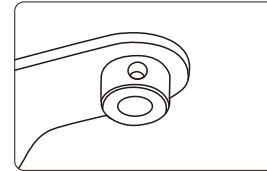


Chart 1

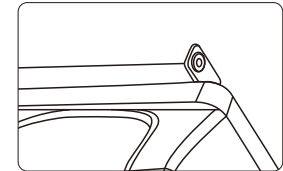


Chart 2

- Place the flat washer into the pin shaft. (chart3).
- Calf hooks install hole and rail shelf install hole should be aligned, insert the pin with a flat washer (both sides are required to be inserted). The pins are aligned with the screw holes for easy installation. (chart4).

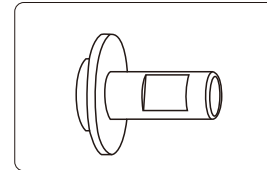


Chart 3

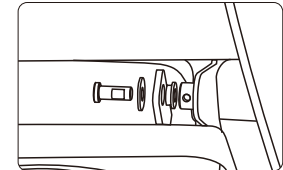


Chart 4

- Use the Allen wrenches lock two hexagon socket head screws. (chart5).

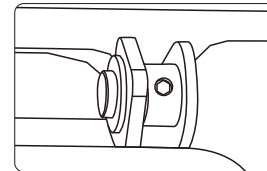


Chart 5

Note

- When the machine under the working state, with lying and calf rest lifting, you could not apply gravity on the calf, to prevent the massage chair does harm to people before dumping (chart1).
- Ensure the mechanism hands not in the seat part, otherwise will cause uncomfortable in the Buttock when you seating in the chair (chart2).

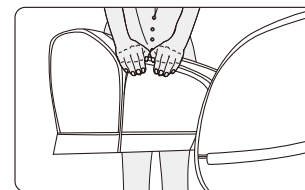


Chart 1


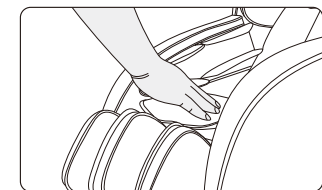


Chart 2

Method of usage


5. Head cushions and backrest pad's dismounting and use manual

- Foldable pillow pad. The users can choice the number of layer to find the best intensity of back massage.
- Separate the head pillow to fall down to the back to reduce the intensity of back massage.
- Backrest cushion and backrest connected by zipper (1), head cushion and backrest pad connected by thread gluing (2)

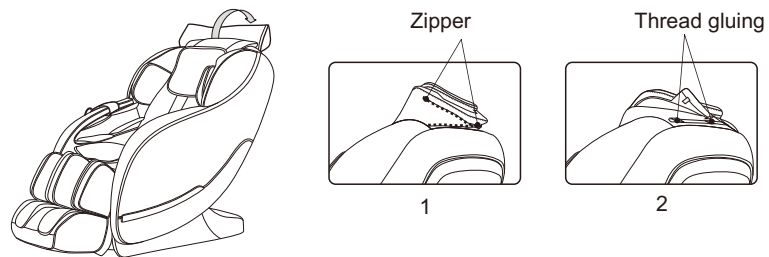

6. Floor Protection

- Put the heavy massage chairs on the wood floor for a longtime may damage the floor, so please place a carpet or other object to avoid this damage.

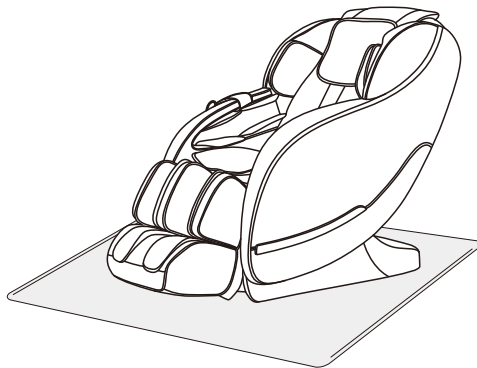


Caution: Do not move the chairs with its wheels on wood floor, rough ground or in a narrow space. Two people are required to lift the chair. During the lifting period, take the calf sealing board and backrest cover grooves as the focal point(Please do not lift armrest Directly.)

Method of usage

7. Method of movement

Press the remote button to make the calf to be lifting to the highest level, then turn off the power switch, unplug the power wire, remote controller wire. Hold a groove in the back plastic cover to tilt the backrest backward to a certain degree (with the gravity center resting on the wheels), push the chair forward or backward with your hands and finally resume the chair to the normal position in a slow and gentle manner.

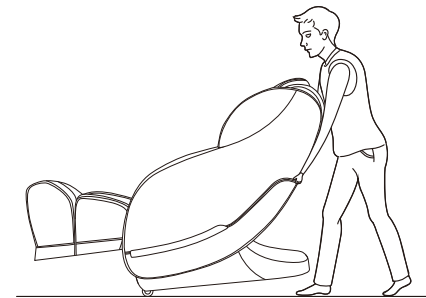

Note: During movement you must turn off the power switch first and then unplug the power line and hand controller connecting line.

Product specification

Model: BM190

Name: Massage Chair

Rated voltage: 110-120V~ 60Hz

220-240V~ 50Hz/60Hz

Rated power input: 150W

Rated time: 20Min

Safety structure: Class I

Material introduction: PVC, PA, steel parts and electric & ele-ctron parts

MATERIAL	PU	LEATHER	CLOTH	WOOD