CruxCare® C1 Home Care Calling System
CruxCare® C1 Hauspflege-Rufsystem
Système d'appel pour soins à domicile CruxCare® C1
CruxCare® C1 Sistema di chiamata per l'assistenza domiciliare
Sistema de llamadas de atención domiciliaria CruxCare® C1
CruxCare® C1 Oproepsysteem voor de thuisverpleging

User Manual
Benutzerhandbuch
Manuel d'utilisation
Manuale utente
Manual del usuario
Gebruikershandleiding

ENGLISH

CruxCare® C1 Home Care Calling System

1. Disclaimer
2. CruxCare C1 system description
3. Instructions for use
4. Accessories
5. Connecting accessories
6. Indication of low batteries
7. How to replace batteries
8. Resetting the C1 pager to factory settings
9. Customer service

Congratulations on your purchase of this CruxCare product. It is created to withstand many years of daily use. It is effortless to use, and we have gone to great lengths to create a CruxCare product portfolio that will allow you to use this product in many situations. Please read this manual carefully to ensure that you will get the most out of the product.

1.Disclaimer

The range of wireless devices may vary and be dependent on your specific environmental conditions.

CruxCare products augment the caregiver's efforts. The use of a CruxCare C1 home care calling system is not a substitute for visual monitoring by a caregiver. Test the device regularly. The CruxCare system relies on wireless radio technology, which is subject to physical and environmental conditions. The pager does not detect when it is out of range. Test it periodically in the setting in which you use it.

2. CruxCare C1 system description

The C1 Home Alert System is a tool that intends to augment caregivers in their daily duties. With a press of a button, a patient can actively ask for help very quickly. Also, with accessories like door sensors and motion sensors, patients can be monitored passively. Connect the pager to up to eight C1-compatible accessories.

3. Instructions for use

Unpacking

Remove the battery strip on the back of the pager. The pager will switch on automatically. The blue light flashes every 5 seconds to indicate that the pager is active.

Switching the pager on and off

Switch the pager on and off by pressing the volume button for 4 seconds. When the pager is switched on, it is immediately active. The blue light flashes every 5 seconds to indicate that the pager is active.

Sound settings

There are two sound types. The red and the yellow zone feature sound 1, the green and the blue zone feature sound 2. This allows the caregiver to distinguish by ear between 2 patients, or 2 locations, or 2 sensors, etc.

Adjusting the volume

The pager features 4 volume levels and a silent mode. Adjust the levels by repeatedly pressing the Volume-button.

Sounds & reset

When a connected device triggers a notification, the pager produces a sound (depending on the zone and volume adjustment), and the respective zone indicator will flash. This sound repeats every 10 seconds, and the zone indicator flashes continuously until the sound is stopped by pressing one of the two buttons.

IMPORTANT: You will never miss a call even if multiple calls are made from more than one device shortly after another because the associated zone indicator will continue to light up until the carer stops the notification. This safety feature prevents an emergency call is missed by the carer in the case of more than one call in a short period.

To stop the sound, press any of the two buttons on the pager.

Status light

The C1 pager features a status light. When the pager is switched on, the blue light will flash every 5 seconds. If desired, this indication can be switched off. Short-pressing the Connect-button will switch the indication on or off.

4. Accessories

Optional compatible accessories (sold separately)

Accessory	Product Code
CruxCare call button (red)	CCCB01R
CruxCare call button (blue)	CCCB02B
Tiiwee PIR motion sensor	TWPIR02
Tiiwee door sensor	TWWS03

CruxCare call buttons red and blue CCCB01R and CCCB02B

Calling the carer

Press the red or blue call button; the button will light up.

Stop the notification sound

To stop (reset) the alarm, press any of the two buttons on the pager.

Installation

This device can be worn around the neck using a lanyard. Alternatively, mount the call button to a fixed position by using double-sided adhesive pads.

Power supply: 1 x CR2450 battery

Dimensions: 5,5 x 5,5 x 2 cm

Housing material: ABS


Radio-frequency: 433 MHz

Tiiwee window & door sensor TWWS03

This sensor can be mounted on doors and windows. When the patient opens a door or window, the sensor sends a wireless signal to the CruxCare C1 Pager. This is helpful when guarding rooms or areas.

Please find below mounting examples for the most common situations in your home:

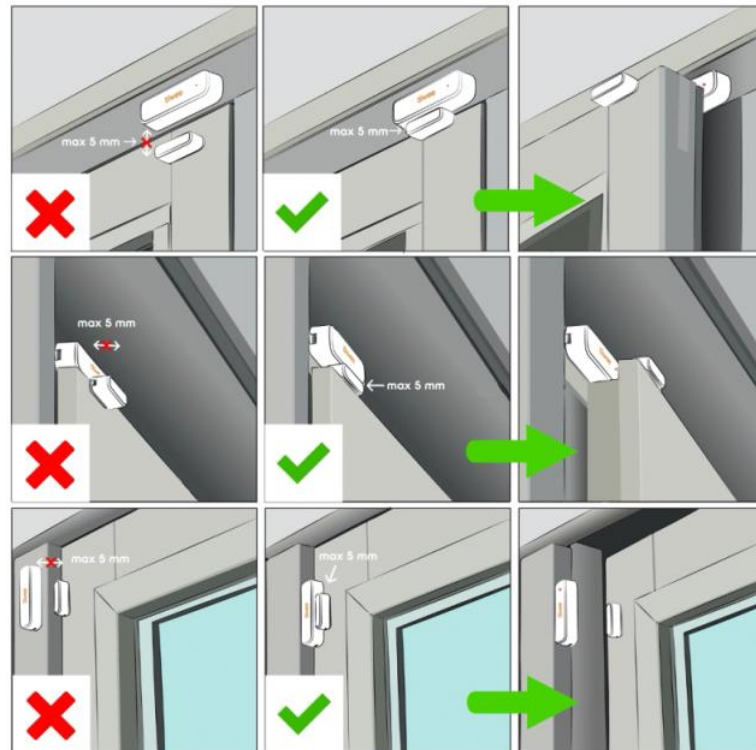Power supply: 1 x AAA alkaline battery

Battery life: 10 - 12 months

Dimensions sensor: 78x24x18mm

Dimensions magnet: 75x10x12mm

Housing material: ABS

Radio-frequency: 433 MHz

Tiiwee PIR motion sensor TWPIR02

This sensor can be used to guard rooms or areas. When the patient moves in the observed area, a signal will be sent to the CruxCare C1 Pager.

Power supply: 3 x AAA alkaline battery or by an external 5V/1A DC power supply connected through the micro-USB port. The batteries can stay inside the sensor while the external power supply is connected and in use. They will work as a power backup in case of 220V power failure. The power adaptor is not included with the product.

Detecting Range: 4-8 m/110 degrees angle

Detector Size: 106x67x45mm

Housing material: ABS

Radio-frequency: 433 MHz

5. Connecting accessories

IMPORTANT: if you have purchased a product containing the pager and an accessory, the accessory is already connected. The pairing procedure is only needed for additional devices or if you want to change the zone for an accessory that is already connected.

Step 1

Enter Programming mode

- Ensure that the pager is switched on.
- Press and hold the Connect-button until you hear a beep.
- The pager immediately goes into pairing mode for the red zone.
- Select the desired zone by short-pressing the Connect-button again.
- The respective zone indicator will be shown.

IMPORTANT: If no button is pressed or no accessory is connected, the pager will automatically exit the programming mode after 20 seconds. If this happens, just enter programming mode again as outlined in Step 1.

Step 2

Pair with compatible devices (use selected CruxCare or Tiiwee accessories only)

- While the pager is in programming mode (respective zone indicator light is on, step 1), trigger the device you would like to connect.
- You hear two short beeps, and the zone indicator light will flash twice to acknowledge successful pairing.
- You can now continue to add more devices to the same or other zones.
- Exit the programming mode by short pressing the volume button.

IMPORTANT: If you hear only one beep during pairing, it means that the accessory was already paired to the pager previously. If you want to change the settings for this accessory, reset the pager and connect the accessory to the zone you wish to.

Step 3

Test the pager

- Trigger the connected accessory
- The pager will sound
- The respective zone indicator light will flash
- To stop the sound, press any of the two buttons on the pager

IMPORTANT: it is always advised to test the operating range within your building or work environment.

6. Indication of low batteries

Product	Low battery indication
CruxCare C1 Pager	Two consecutive beeps every 20 seconds accompanied by a flash of all LEDs
CruxCare call buttons red & blue (CCC01R CCC02B)	No low battery indication. Replace batteries every six months
Tiiwee window & door sensor (TWWS03)	LED light of the sensor is blinking or is continuous lit
Tiiwee PIR motion sensor (TWPIR02)	The indicator light in the front lights up continuously

7. How to replace batteries

CruxCare C1 pager

Take the pager out of the silicone sleeve. If you used the screw to secure the rear panel, remove the screw at the bottom. Slide off the back panel and replace the batteries with 3x AA-size 1.5V alkaline batteries.

IMPORTANT: do not use rechargeable batteries. Only use alkaline batteries.

CruxCare call buttons

Open the button at the site where the little cutout is visible at the rim of the call button. This can be done with the tip of a knife or a screwdriver. You will now see the battery. Replace the battery with a new CR2450 (3V) battery. Please note the correct polarity. Close the call button by placing the cutouts opposite each other and pressing the two halves together.

Tiiwee window & door sensor

Remove the screw at the bottom of the sensor. Now slide the sensor upwards and remove it from the back panel. Take out the old battery and replace it with a new AAA alkaline battery. Please note the correct polarity of the battery. Slide the sensor on the back panel and secure it with the screw at the sensor's bottom. Test the sensor.

Tiiwee PIR motion sensor

Open the back panel and replace the batteries with 3x AAA batteries.

8. Resetting the C1 pager to factory settings

Use this procedure when you want to assign an accessory to a different zone.

- Remove one of the 3 x AA batteries
- Press and hold the Connect-button while inserting the battery again
- You will hear ten beeps, followed by ten short beeps
- Release the button
- After the factory reset, the pager will switch on again automatically

IMPORTANT: After a factory reset, all accessories will need to be connected again. Please refer to chapter 5 for accessory pairing instructions.

IMPORTANT: After the factory reset, the volume setting will be at the highest level.

9. Customer service

If additional support is required, please visit www.tiwee.com. Here you can find supplementary support material as well as our support videos in your language.

You can also contact us at info@tiwee.com. We respond within 24 hours to all questions, but mostly still on the same day. We provide a two years warranty on all our products.

CruxCare® products are manufactured in China.

CruxCare® and tiwee® are registered brands owned by Tiiwee B.V.

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CruxCare® en tiiwee® zijn gedeponeerde merken van Tiiwee B.V.