# Tiiwee Pager

ENGLISH

# Contents

## Introduction

## System Overview

Acoustic vibration receiver

Call button (Transmitter)

## Unpacking and Testing Your Tiiwee Pager

Powering the receiver

Testing the product

Testing the operating range

Adjust the sound volume

## Changing the ring tone for the call button

## Replacing Batteries

## Resetting the Pager to Factory Settings (deleting call button from the receiver)

## Technical Specifications

## Customer Service

## Introduction

Congratulations on your new Tiiwee Pager.

### Working principle

Press the call button, and a wireless signal is sent to the receiver. The receiver will emit a sound, and it will vibrate. The red LED indicator light will also light up.

### Intended use


Use the Tiiwee wireless pager in situations where patients, older adults or disabled people need help. Its purpose is to notify family members, caregivers or neighbours who then can provide support.

### Disclaimer


The Tiiwee pager augments the caregiver's efforts. The use of a Tiiwee pager is not a substitute for visual monitoring by a caregiver. Test the device regularly. The Tiiwee pager relies on wireless radio technology which is subject to physical and environmental conditions. The product does not have an out-of-range function. It should be tested periodically in the setting in which it is used.

## System Overview

### Acoustic vibration receiver
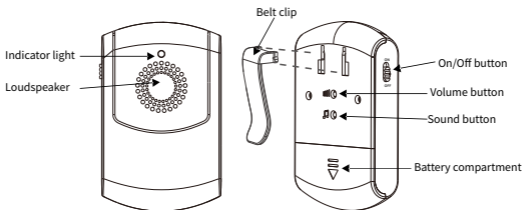
## Call button (Transmitter)
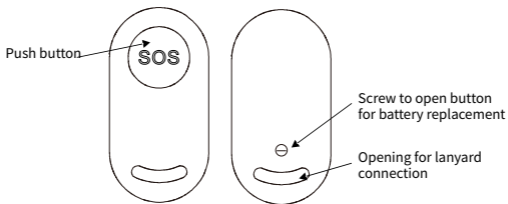
## Unpacking and Testing Your Tiiwee Pager

Unpack all kit components and lay them in front of you.

### Powering the receiver

Open the battery compartment and insert the included AAA batteries. Close the battery compartment.

### Testing the product

- Switch on the receiver by using the On/Off-button.
- Now, press the call button.
- You will hear the "Ding Dong" sound, the receiver will vibrate, and the LED indicator light will flash.

### Testing the product

- The range of the system in an open area is about 75 meters. However, radio signals are affected by physical and environmental conditions like walls, doors and other barriers. Inside homes, the expected real range of the system is about 2-3 floors or 15-20 meters.
- Please always test if the operating range is sufficient for your desired use.

### Adjust the sound volume

The volume of the sound can be adjusted in 4 steps by pressing the Volume button with the tip of a pen.

## Changing the ring tone for the call button

**Note:** The factory setting for the ring tone is "Ding Dong". However, you can choose from 58 melodies.

**IMPORTANT:** In step 2, the call button is eliminated from the receiver and then reconnected in step 3. If no call button is pressed within the 5 seconds as described in step 2, the connection between the call button and the receiver is lost, and you will need to go through the steps 1 to 3 again.

1. Choose your favourite ring tone by pressing the Sound button with the tip of a pen.
2. Press the Volume button for 5 seconds until you hear a low buzzing sound and the LED indicator flashes continuously - you are now in pairing mode.
3. Now, press the call button within 5 seconds. You will hear the selected sound on the receiver. The receiver will vibrate, and the LED indicator light will flash.
4. Your new settings now have been stored.

## Replacing Batteries

The average battery life of the vibration receiver is about three months. The average battery life of the call button is approximately 10-12 months.

Battery type call button: CR2032

Battery type receiver: 2 x 1.5 V AAA

**Note:** Please always use high-quality batteries for the receiver and the call button.

**Note:** The receiver will not lose its setting after battery replacement.

**Note:** The system does not feature a low battery warning. Therefore, we recommend checking the functionality of the pager each time or each day before use. A good practice is to make this part of your standard routine.

## Resetting the Pager to Factory Settings

(deleting call button from the receiver)

- With the tip of a pen press the Sound button for 5 seconds.
- You will hear a low buzzing sound.
- The call button is now disconnected from the receiver.
- The ring tone sound has been reset to the factory setting ("Ding Dong").
- You can now connect the call button again by following the steps outlined in the section "Changing the ring tone for the call button".

## Technical Specifications

Specifications	Vibration (Receiver)	Call button (Transmitter)
Power Supply	3V DC power (2x 1.5V AAA batteries)	CR2032 battery
Wireless Range	75 meters in open area	
RF Frequency	433 MHz	
Sound Volume	adjustable in 4 steps (110dB Max.)	
Number of Melodies	58	
Power Consumption	<0.1W	
Waterproof Rate		IP55
Operation Temperature	-15°C - 55°C	

## Customer Service

Please visit [www.tiivee.com](http://www.tiivee.com) for support.

You can also contact us at [info@tiivee.com](mailto:info@tiivee.com). We respond within 24 hours, but mostly still on the same day. We provide a two years warranty on all our products.

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