



## Shyn Daily Sonic Toothbrush

### IMPORTANT SAFEGUARDS

Your Shyn Daily Sonic Toothbrush provides a safe, highly effective brushing experience for a cleaner mouth, a more confident smile and a healthier you.

#### READ ALL INSTRUCTIONS BEFORE USING

#### DANGER

In order to reduce the risk of electrocution:

1. Don't use while bathing.
2. Don't place or store charging station where it can fall or be pulled into tub or sink.
3. Don't place in or drop the charging station into water or other liquid.
4. Don't reach for a charging station that has fallen into water. Unplug it immediately.

#### WARNING

1. To reduce risk of burns, electrocution, fire, or injury to people:  
This product is not intended for use by anyone (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the product.
3. Use this product only as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if the charging station has been dropped into water.
5. Keep the cord away from heated surfaces.
6. Never use while sleeping or drowsy.
7. Never drop or insert any object into any opening of the product.
8. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
9. Do not plug charger into an outlet with a voltage other than that specified on the bottom of the charger. Do not use this charger with a voltage converter.
10. Do not force plug into an outlet. Discontinue use if plug has been damaged.
11. This product contains no user serviceable parts.
12. Use this product only for its intended use as described in this document or as recommended by your dental professional.
13. This product is a personal care device and is not intended for use on multiple patients in a dental office or institution.

Not intended for children ages 3 and under.


#### MEDICAL WARNINGS

1. Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.
2. Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after 1 week of use.
3. Discontinue use of this product and contact a physician/dentist if discomfort or pain is experienced.
4. This product complies with the safety standards for electromagnetic devices.  
If you have a pacemaker or other implanted device, contact your physician prior to use.
5. Consult your physician if you have other medical concerns.

#### ELECTROMAGNETIC FIELDS

This product complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

### SAVE THESE INSTRUCTIONS

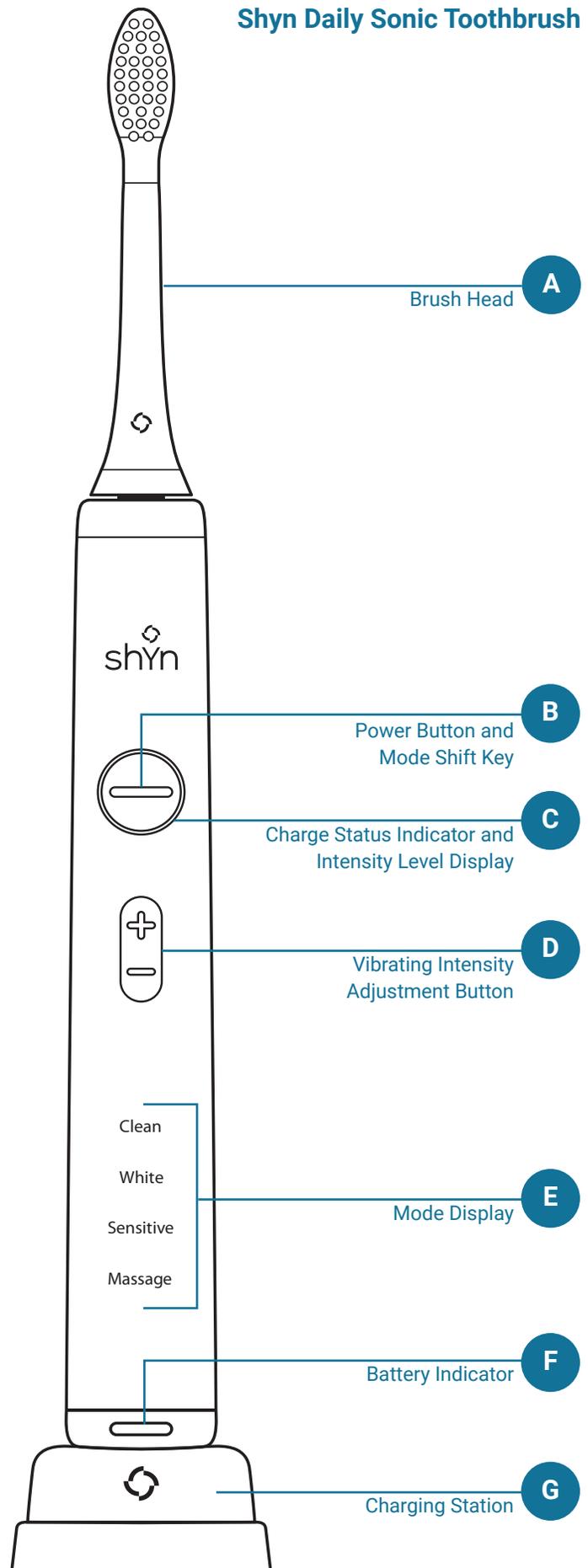

## POWER UP: CHARGING YOUR TOOTHBRUSH

- Before you first use your toothbrush, we recommend charging the brush handle for 16 hours. It is normal for batteries to lose charge while not in use for extended periods of time.
- Plug the charging station into an electrical outlet, then place the toothbrush handle on the charging station. A full charge can take up to 16 hours. The charge status indicator (C) will flash yellow while charging, becoming brighter as the battery level increases. It will turn solid yellow when charging is complete.
- Once fully charged, the battery indicator (F) will turn solid white. The battery should last up to 7 days when brushing for 2 minutes, twice daily.
- When the battery is low, the battery indicator will turn yellow.
- When the battery is critically low, it will flash yellow. Return it to the charging station immediately to avoid damaging the toothbrush.
- If the battery level becomes too low, the toothbrush will perform a safety shutdown. If your toothbrush shuts down, charge it for a full 16 hours before use.
- For maximum battery life, fully charge the battery at least every 6 months when not in use.
- You can't overcharge the battery, but you can save energy. We recommend unplugging the charging station when you're not using it.

CHARGE INDICATOR <b>C</b>	STATUS
Flashing White	Charging
Solid White	Fully Charged

BATTERY INDICATOR <b>F</b>	STATUS
Solid White	Fully Charged
Solid Yellow	Low Battery
Flashing Yellow	Critically Low Battery

## ATTACH & REMOVE THE BRUSH HEAD

- To attach the brush head, align the brush head with the bristles facing the front of the handle, then push it firmly onto the handle until it's completely secured. A small gap should be present between the brush head and handle. **NOTE: the brush handle will only work with a brush head attached.** This prevents the handle from operating accidentally when packed for travel.
- To remove the brush head, simply pull off of the handle.



## HOW TO BRUSH YOUR BEST

- Wet the brush head and apply toothpaste.
- Before turning on your toothbrush, hold the bristles against your teeth and slightly angled toward your gum line. Once your brush is in place, press the power button (B). This prevents toothpaste and water from splattering.
- To brush your teeth, apply light pressure and let the toothbrush do the brushing movement. **Very important: don't scrub or press too hard as this can lead to damage to your teeth and gums.**
- Gently work from tooth to tooth using a small back and forth motion. This helps the longer bristles reach between teeth.
- Dentists recommend spending 30 seconds brushing each quadrant of your mouth, for a total of 2 minutes. To signal the end of each 30-second brushing interval, the brush handle will pulse 2 times. At the end of 2 minutes, the brush handle will stop vibrating and turn off (unless you are using a 3 minute mode, please see Modes section for details.) **NOTE:** If you would like to pause the timer during a brushing cycle, tap the power button (B). When you tap it again, the cycle will start where you left off.
- You can also brush your tongue with the toothbrush turned on or off.

## BRUSH SAFELY WITH A BUILT-IN PRESSURE SENSOR

- As you brush, our built-in pressure-sensing technology will alert you if you are brushing too hard. Brushing too hard can lead to teeth and gum damage. The pressure sensor will alert you by beeping until you brush more softly and safely. This feature is already activated in your brush. **NOTE: the pressure sensor will only work when the brush is set to higher speeds of Clean Mode or the Sensitive Mode.**
- To turn the pressure sensor off, long-press the power button (B) for about 3 seconds. When you hear 2 beeps, the sensor has been de-activated.
- To turn the pressure sensor on, long-press the power button (B) for about 3 seconds. When you hear a beep, the sensor has been activated.

## 4 MODES FOR A CUSTOM CLEAN

Your toothbrush offers 4 modes so you can choose the one that's best for your dental care needs:

- **Clean:** 2 minutes brushing time, 31,000+ vibrations per minute, 5 intensity levels, pressure sensor on intensity levels 3, 4 and 5
- **White:** 2 minutes brushing time, 31,000+ vibrations per minute, 3 intensity levels
- **Sensitive:** 3 minutes brushing time, 31,000+ vibrations per minute, 5 intensity levels, pressure sensor on intensity levels 3, 4 and 5
- **Massage:** 2 minutes brushing time, 31,000+ vibrations per minute, 3 intensity levels
- To change modes, turn on the toothbrush. Within 3 seconds, quick press the power button/mode shift key (B) to set or change the mode. Otherwise, the toothbrush will turn off.
- Your current mode is indicated by the mode display (E).

## INTENSITY SETTINGS FOR TOTAL COMFORT

- Depending on the modes you are using, you can switch between 3 or 5 different intensity settings. This allows your mouth to get used to the sonic vibrations gradually.
- While brushing, the intensity level display (C) is lit lower for low speeds and brighter for high speeds.
- **NOTE:** your toothbrush is equipped with a memory function. It will turn on set to the same mode and speed that you chose when you last used it.

## CUSTOMIZE YOUR CLEAN

Your mouth is unique, and it deserves a brushing experience built around its needs. We offer 3 different brush heads for specific dental care goals—and all of them remove 6x-7x more plaque vs. brushing manually. Learn more at [shyn.com](http://shyn.com).

### Anti-Plaque Brush Head

Impress your dentist. Our high-density bristle design and Advanced Clean tip removes up to 7x more plaque compared to a manual toothbrush.

### Gum Care Brush Head

Go easy on your gums. With extra-soft bristles and a Gentle Clean tip, this brush head is designed to improve gum health and clean effectively along the gum line.

### Whitening Brush Head


Brighten up your look. Highly effective diamond-shaped bristles whiten teeth in just 1 week, while the slim, angled neck offers superior plaque removal.

## DENTAL INFO YOU SHOULD KNOW

- Slight bleeding may occur during the first few days of use. If bleeding persists for more than one week, call your dentist.
- This toothbrush is safe to use on braces, dental restorations, veneers, and implant material.

## TRADE UP EVERY 3 MONTHS

- The "Replace-Me<sup>®</sup>" bristles on your brush head will fade after about 3 months of use to remind you to replace the brush head.
- Dentists recommend replacing your brush head every 3 months for optimal performance.
- Replace sooner if your brush head becomes loose due to damage.
- Use only replacement brush heads identified as compatible with this power toothbrush.
- Every replacement brush head Shyn offers is compatible.

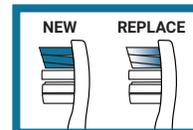

## CLEANER BRUSH, CLEANER MOUTH

- Clean your toothbrush brush after every use.
- Remove the brush head from the handle and rinse both brush head and handle attachment area separately with warm water. Allow both to air dry.
- Handle and charger can be surface-cleaned using a damp cloth. Charger should be unplugged before cleaning.

**DO NOT WASH ANY PARTS IN DISHWASHER.**

## TRAVEL SMART, LOCK YOUR BRUSH

Your toothbrush comes with a handy button lock feature great for travel. Here's how to use it:

**To lock:** while charging the toothbrush, quick press the power button/mode shift key (B) twice in one second. All display lights will flash once to indicate the toothbrush has locked.

**To unlock:** quick press the power button/mode shift key (B) twice in one second. All display lights will flash 3 times to indicate the toothbrush has unlocked. You can unlock while the toothbrush is on or off the charging station.

## STORE IT SO IT LASTS

When not in use, keep the charging station and toothbrush away from heat sources. Store them in a dry place such as a bathroom countertop, cupboard, drawer, or cabinet. Do not wrap the cord around the charging station.

Product is for household use only.

Product has no user serviceable parts.

Product contains Nickel-Metal Hydride (NiMH) batteries.

Please recycle or dispose of properly according to your local regulations.

**Caution!** Opening handle will destroy the product.

## BATTERY REMOVAL PRIOR TO DISPOSAL OF HANDLE

This power toothbrush contains a rechargeable battery that must be disposed of properly according to your local regulations.

Removing the battery will render the device inoperable, as the battery cannot be replaced.

The handle is not intended to be reassembled. To remove battery you will need a Phillips screwdriver, a bladed screwdriver, old towel, and pliers/wire cutters. Make sure the battery power is completely empty before removing it and follow basic safety precautions when following this removal procedure:

1. Using a Phillips head screwdriver, remove the screw from inside the bottom cap. Place a towel (or other item to protect from accidental damage) on a flat, sturdy surface. Then place the handle on the towel. Using a bladed screwdriver, pry the bottom cap from the handle. The cap should come free, but may require prying at multiple areas.
2. While pressing the metal shaft of the handle against a hard surface, use the screwdriver to displace the two locking tabs holding the internal structure in place. Once displaced, the internal structure can easily be removed.
3. Locate the battery on the internal structure and position a screwdriver underneath the battery. Lift and pry the battery out of the plastic housing. The wires connected to the battery will disconnect as the battery is lifted out. (The circuit board breaking during this process is normal.)
4. Place tape over battery ends to avoid any short-circuiting due to energy remaining in the battery.

Recycle or dispose of the battery according to local regulation. Discard the remaining portion of the product as appropriate.

This product is rated IPX7 as tested per IEC standard 60529.

# FOR USA ONLY

## LIMITED 2-YEAR WARRANTY

The manufacturer warrants this product to be free of defects in material and workmanship for 2 years from the date of original purchase at retail.

If the product exhibits such a defect as stipulated above, at its option, the manufacturer will repair or replace the product without cost for parts and labor. The defective product must be shipped to the manufacturer prepaid and insured along with proof of date of original purchase.

The cost to insure and ship the product to the manufacturer is the responsibility of the consumer. The cost to return the product to the consumer is the responsibility of the manufacturer. This warranty does not cover the following: finishes or normal wear and damage caused by tampering, unreasonable use, servicing attempted by or performed by unauthorized service centers, accident, water, dirt, failure to provide reasonable and necessary maintenance including recharging or products that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO 2 YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL THE MANUFACTURER BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you. This warranty gives you specific legal rights and you may have other legal rights which vary from state to state.

To have your defective product repaired or replaced, please call 1-866-532-6041. You will be connected with a customer service representative who will provide you with a return authorization number and address to where you should ship the product. Please insure the product is packaged well prior to shipping. Ship the product prepaid and it is recommended to insure the shipment. Be certain to include a copy of your proof of purchase to verify warranty coverage and a return address and phone number to assist with shipping the product back to you.

## FCC DECLARATION OF CONFORMITY

This device complies with part 18 of the FCC rules.

Operation of this device may cause harmful interference to radios and other radiofrequency devices.

This device should be maintained in accordance with the provided instructions.

If undesirable interference occurs while operating this device, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/television technician for help.



For more information visit [www.shyn.com](http://www.shyn.com)



# FOR CANADA ONLY

## LIMITED 2-YEAR WARRANTY

In the event this product fails to function within the specified 2 year warranty period because of defects in material or workmanship and the consumer returns the product to the manufacturer, at its option, the manufacturer will either repair or replace the product without additional charge to the consumer. This warranty does not cover any product which has been damaged by servicing attempted by or performed by an unauthorized service center, tampering, misuse, dropping or abuse. This warranty excludes products which have been modified or used for commercial purposes.

To the extent allowed by law, the manufacturer shall not be responsible for loss of the product, loss of time, commercial loss, inconvenience, special or consequential damages. The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

To have your defective product repaired or replaced, please call 1-866-532-6041. You will be connected with a customer service representative who will provide you with a return authorization number and address to where you should ship the product. Please insure the product is packaged well prior to shipping. Ship the product prepaid and it is recommended to insure the shipment. Be certain to include a copy of your proof of purchase to verify warranty coverage and a return address and phone number to assist with shipping the product back to you.