



HOME SUITE

Neck & Shoulder Massager

USER MANUAL



NECK & SHOULDER CORDLESS SHIATSU MASSAGER

SAFETY INFORMATION:

To avoid electric shock or malfunction do not use the product in:

- A damp or dusty environment
- In temperatures Higher than 104F
- In an Outdoor Environment

USE WHEN:

Consult a Doctor Before Use If You Are:

- Receiving Medical Treatment
- A Cancer Patient
- Pregnant
- Suffering From Acute Diseases
- A Person with History of Heart Attacks



NOTES BEFORE USE:

- Do not use this product with other therapy equipment or electric blankets
- Do not use the product 30 minutes after eating
- To avoid injury, do not allow children to use product
- Avoid carrying heavy objects on any part of the product, causing mechanical Failure or damage due to fallen objects
- Do not continue to massage a fixed point on the body for more than ten minutes
- This product is designed for one person use at a time.
- Simultaneous use by two or more persons is not recommended
- This product is designed only for those who weigh less than 265lbs

STOP USING WHEN:

In the following circumstances you should immediately stop using the product, Turn off the power and pull the cable from you outlet:

- Accidentally splashing water on the massage heads or handles
- Related accessories are damaged or internal parts are exposed
- Feel unusual or severe pain during massage
- A sudden power failure
- During Lightning



NECK & SHOULDER CORDLESS SHIATSU MASSAGER

ASSEMBLY AND REPAIR:

- *Do not dismantle the protective covering or touch the internal parts of the machine to*
- *Avoid personal injury due to electric shock or mechanical failure*
- *The product can only be repaired by professional service personnel – Do not*


ACCESSIBILITY:

- *Do not pull off adapter with wet hands to avoid electric shock hazard*
- *When you pull off adapter, grasp the whole adapter instead of just the wire*
- *The timing function can automatically stop the massage cushion, please use the Power switch after each use*
- *Adapter shall not be packed by the product or other heavy objects*
- *If the adapter is damaged, immediately stop using the product, turn off the power Switch, pull the adapter out from the outlet and contact a qualified service Personnel*
- *When it is not in use for a long time, turn off the power switch*




PRODUCT FEATURES AND EFFICIENCY:

Neck & Shoulder Shiatsu Massager can provide an in-depth massage for the neck, back, shoulders, legs, and waist. The heating treatment can relax your body, eliminate fatigue, relieve pressure, and promote blood circulation. Compact, stylish design not only meets household demand, but also considers the needs of the driver, so that you can enjoy the fun of massage in the car to achieve the perfect massage

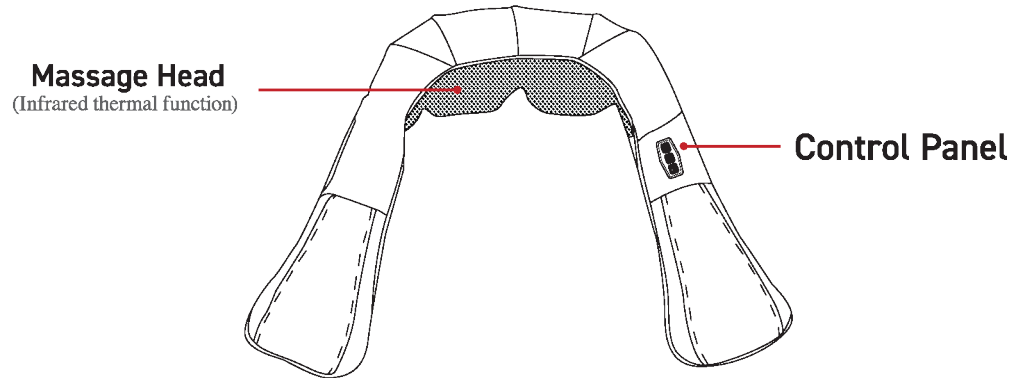

NECK & SHOULDER CORDLESS SHIATSU MASSAGER

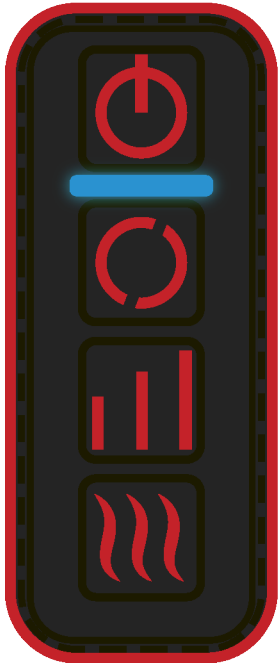

Cordless Battery Operated Indicator Shows Charge

Multi Directional, Rotating Shiatsu Massage

3 Intensities - Light, Medium & Heavy Massage

Heated Massage for Therapeutic Blood Circulation



Products: Neck & Shoulder Shiatsu Massager

Power: 24W

Adapter Input: 110V-240V 50/60HZ

Product Input: 12.6V,1000MA

Safety Class: II

G.W: 1.6KGS

N.W: 1.5KGS

Addapter Output: 12.6V,1000MA



www.zarifausa.com



461 W CENTURY DR
#3
MURRAY UT
84123

