# 4>ZARIFA USA



## **Z-SMART PLUS**

Massage Chair

**USER MANUAL** 

## **THANK YOU!**

## Our Loyal Customer

Our goal is to provide a graceful living experience through drug-free pain relief. Our products are made with the customer in mind. During trying times of an opioid epidemic, we aim to provide the highest quality products using modern methods of traditional practices that achieve beneficial results.

Meeting our customers and hearing their stories is what continually drives us each day. Without customer feedback we would not have the knowledge to continue to innovate our products and perform the research our customers need the most. Thanks to all our customers providing us with their feedback, each of our medical products have been designed with you in mind. We pay the utmost care and attention to detail in order to continue delivering nothing short of the best medical devices possible. We could not be more excited for you to join us as continue our journey towards stopping cycles of chronic pain in your life.



## Massage Chair Contents

#### **Massage Chair**

Thank you for purchasing our product. Please read the safety and operating procedures in this User Manual before using the massager. Keep this booklet in a safe place for future reference.

#### **Contents**

Safety Warnings	03-04
Safety and Maintenance	05-08
Parts and Functions	09-10
Remote Control Introduction	11-12
Operating Instructions	13
1. Plug in and Power on	13
2. Massage Functions and Controls	13-30
3. Quick Access Controls	31-32
4. Voice Commands	33
5. Stop Massage and Power off	34
6. Bluetooth	34
7. Armrest Installation	34-36
8. Footrest Installation	36-38
9. Backrest and Head Pillow Adjustments	38
10. Floor Protection	39
11. Moving the Chair	39
Product Specifications	40



#### **Important Safety Warnings**

- Do not put hands or fingers near the massage mechanism or rollers while this
  product is running. Touching the mechanism can result in injury from the rollers
  squeezing together.
- Ensure that the unit is properly grounded before operating.
- Do not leave the product unattended. Unplug from outlet when not in use. Do not use attachments not recommended by the manufacturer.
- Use the product only as intended and use as described in this operation manual.
- The recommended usage time is 20 min. Do not use excessively.
- Do not use if chair cover material or leather is broken or damaged.
- Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not massage any area of the body that is swollen, inflamed, or where skin broken.
- Do not use this product while sleeping.
- Do not stand on or inside of this appliance. Use only while seated.



#### Saftey

- Do not clean any parts with a wet cloth. Keep moisture away from electric parts.
- Do not use if the power cord is damaged.
- Do not use the product if it becomes damaged or makes irregular noise.
- Do not let children play with the device.
- If the cord is damaged contact the manufacturer for recommended replacement.

#### **People Not Suitable To Use The Chair**

- People with circulatory problems (such as varicose veins or thrombosis), children, and those advised not to use a massage chair by their doctor should not use this chair.
- If you have any of the following health conditions, please consult your doctor for medical advice before using this chair: serious heart problems, serious diabetes, advanced osteoporosis, serious skin conditions, malignant tumor(s), sensory impairment, pregnancy (or suspected pregnancy) or those who have just given birth, acute gastrointestinal complaints (gastritis, hepatitis, enteritis), joint dysfunction such as rheumatism, hammer toe, gout, thecitis (or suspected thecitis), high fever, spinal abnormalities due to past injury or ailment, the curvature of the spine (scoliosis acute neck sprain (whiplash), or any other medical condition that may contraindicate usage of this product.
- Seek medical advice before using this chair if you have a pacemaker or other
  electronic medical device if you have been restricted to bed rest, if you are planning
  to use this chair to treat a specific ailment, to recover from surgery, or for therapy.
- Elderly people and those with weak bones should seek medical advice before using this chair even if they do not have a specific disease or disorder mentioned above.
- Do not use this chair if you have recently consumed alcohol.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

#### **Environment Usage**

- Do not use if air openings are blocked or obstructed. Inspect regularly to ensure air openings are free from debris such as lint, hair, or foreign objects.
- Never drop or insert anything into the openings of the chair.
- Do not use this product when standing on a damp floor, or when any part of the body is in contact with metal or any similar ground.
- · Do not use this product outdoors.
- Do not use in humid or in high-temperature areas, such as a bathroom.
- Do not use it in a heavily dusty or caustic environment.
- · Use in a well-ventilated environment.

#### **Trouble Shooting**

- If the controller is not working normally check the cable and plug socket. Ensure the plug is properly plugged in and secure. Check to ensure that the cable is undamaged.
- This product is equipped with a thermal protection device. This is an added safety
  feature to protect your product from overheating with usage over the rated operation
  time. To cool turn the product off and let it rest for at least 30 min.

#### Safety and maintenance


#### Maintenance

- This product is designed and built to be maintenance-free. No periodic lubrication is required.
- No user-serviceable parts inside. To be repaired by manufacture certified technician only.
- If the product is not in use for a long time, store it appropriately. It's recommended
  to cover the chair, tuck the power cord away, and store it in an environment that is
  dust and moisture-free.
- Do not store near heat or open flame.
- Do not leave exposed to sunlight. This may cause damage or fading to the chair.
- Do not store in extreme temperatures. Excessive heat or cold may cause damage to components.
- Clean with a dry cloth and/or vacuum with soft brush extension regularly. If the
  product becomes dirty, use a water-based upholstery cleaner only. Treat stains with
  spot cleaner with upholstery shampoo or mild detergent foam.
- Do not roll on uneven or rough surfaces.



#### Warning


• Do not push down on the footrest when the footrest is lifting this may result in personal injury or damage to the chair. Do not sit or stand on footrest.

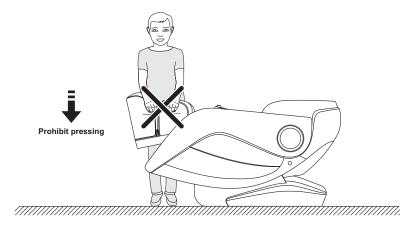



#### Warning

• Ensure massaging hands are lowered before sitting to avoid discomfort.

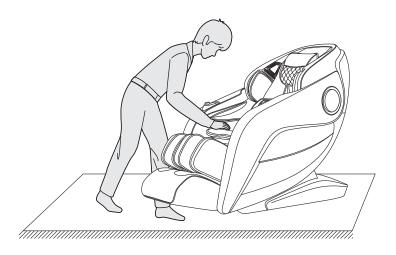


Safety and maintenance



#### Warning


• Do not pull on the remote control cable or power cord. It is recommended that the power cable is unplugged when the chair is not in use.

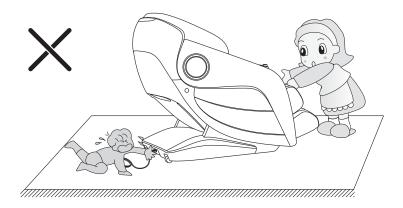




#### Warning

• Inspect the area before enabling zero gravity mode to make sure the area is free of children and pets to avoid injury.

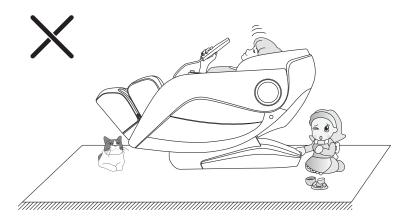

#### Name and function of components

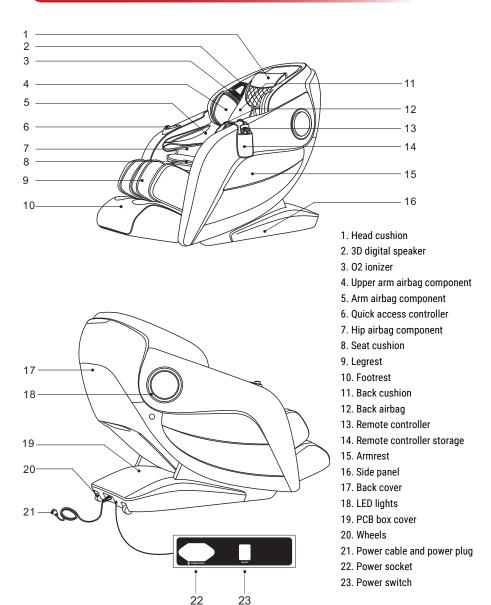


#### Parts and functions

#### **Functions**

- Intelligent, quiet, 3D massage hands that hands combine for a professional massage experience.
- Automatic shoulder detection that can manually be adjusted.
- 8 different intelligent massage modes: fast, comfort, relax, stretch, full body, neck & shoulder, back & waist.
- 3 programmable memory modes.
- Configurable manual massage for the upper body with 5 different modes shiatsu, kneading, flapping, knocking, and simultaneous kneading and flapping. Each technique has five adjustable speed levels and the tapping and shiatsu techniques have five adjustable massage ball positions.
- Air pressure arm massage: upper arm shoulder massage with 4 airbags, Arm sleeve massage with 8 airbags, and 3 levels of air pressure intensity adjustment.
- Lower body air pressure massage functionality: Waist pressure massage with 4 airbags, leg and hip
  pressure massage with 26 airbags, foot air pressure with 18 airbags, and three levels of air pressure
  intensity adjustment.
- Z-Smart carbon fiber back and legrest heat with 3 adjustable levels.
- · Space-saving design that slides the chair forward when reclining.
- Designed with combo "S" and "L" shape track for full body massage from head to toe.
- · Intelligent voice control functionality.
- 3D Best-In-Class Bluetooth Speakers with fast pairing functionality.
- Waist Twisting functionality: stretches and moves the waist using the left and right waist airbags.
- · Quick access controller embedded into right armrest for quick and easy access to adjustments.
- Zero-G experience, the legs extend 135° and the spine stretches along 118° to create the feeling of weightlessness.
- Adjustable leg base with auto-detection of leg length. When adjusted the footrest will automatically detect
  the length of the leg to target the most comfortable massage area.
- Foot massager with two pairs of rollers. One on the heel and arch, and one that rotates on the toe area.
- Built-in LED accent lighting.
- 02 ionizer built into the right armrest.
- Adjustable upper arm: each shoulder and upper arm airbag is independently adjustable. To adjust the
  width, pull inward until you hear a click and release. Continue to pull inward past the 2 settings and release
  to open back up.
- · Automatic power off and shutdown sequence
- Colorful TFT (thin-film transistor) remote controller with adjustable brightness.

#### **Remote Controller Display**

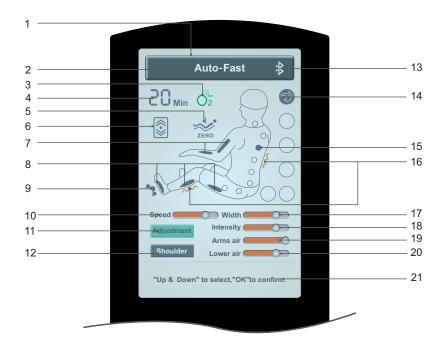


#### Note: Remote pictured may differ from actual unit.

- 1. TFT colorful LCD
- 2. Massage function information display
- 3. 02 ionizer
- 4. Time display
- 5. Zero gravity mode display
- 6. Back adjustment display
- 7. Arm air pressure display
- 8. Leg and foot air pressure display
- 9. Foot roller display
- 10. Speed adjustment
- 11. Adjustment menu

- 12. Shoulder position adjustment
- 13. Bluetooth display
- 14. Massage mode display
- 15. Massage position display
- 16. Back and leg heat display
- 17. Width adjustment
- 18. Massage intensity adjustment
- 19. Upper body air pressure adjustment
- 20. Lower body air pressure adjustment
- 21. Status bar

#### **Button Panel**

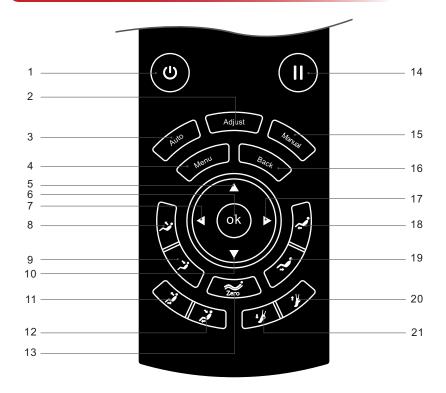

- 1. Power: on/off
- 2. Adjustment: enters adjustment mode
- 3. Auto: enters automatic program selection
- 4. Menu: enter main menu
- 5. Up: guide button goes up
- 6. OK: OK button
- 7. Left: guide button goes left
- 8. Backrest up: backrest goes up
- 9. Backrest down: backrest goes down
- 10. Down: guide button goes down
- 11. Raise up: raises backrest up and lowers footrest down

- 12. Lie down: lowers backrest down and raises footrest up
- 13. Zero gravity: adjust zero gravity function
- 14. Pause: pause at current function
- 15. Manual: enter to manual function selection
- 16. Back: back to function information
- 17. Right: guide button goes right
- 18. Footrest up: adjust footrest up
- 19. Footrest down: adjust footrest rest down
- 20. Footrest shorten: shortens footrest length
- 21. Footrest extend: extends footrest length

#### **Operating Instructions**

#### 1. Plug in and power on

- Plug in the power cable
- Turn on the main power switch



#### 2. Massage function and controls

#### 1. Standby

 When the massage chair is powered on and off it is in standby mode. The set function is the only function that can be performed in standby mode.



#### 2. Start a massage


Press the "power" button, start the massage, the backrest will stay in position and the footrest
will rise slightly. Choose an automatic function mode to start the massage. The massage chair
will automatically shut down if no mode is selected within 5 minutes.






#### Operating instructions

- Choose an automatic massage mode, the backrest and leg rest will adjust after selection to a fixed angle.
- After adjusting the angle the chair will perform a scan for optimal positioning. You can press "OK" to skip this and start the default massage program.

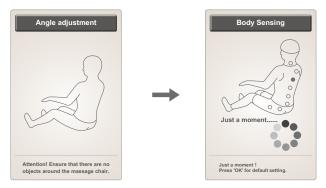


Once the shoulder position is detected the buzzer tone will sound. At this point, if you would like to
adjust the shoulder positioning you can adjust up and down to find the perfect position. The remote
has 11 built-in levels. If no adjustment is made on the shoulder position the chair will automatically
begin massage after 10 seconds.

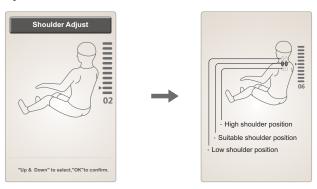


Button	Function		
(OK)	Adjust up to move the massage rollers upward to correct position		
	Adjust up to move the massage rollers downward to correct position		

 After entering the massage modes, the shoulder position can be adjusted using the shoulder adjustment menu.

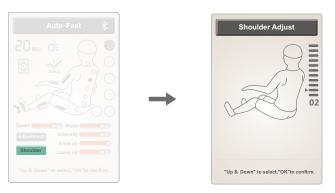

#### 3. Pause button


 When a massage mode is operating, press "pause" to temporarily suspend all massage functions. Press "pause" again to resume. If the function remains paused for over 20 min, then all massage functions will shut down automatically.

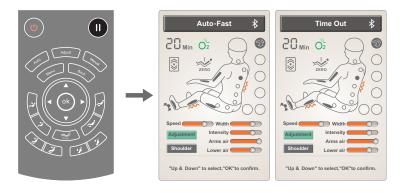

#### Operating instructions


#### 4. Menu

- Press "menu" to navigate to the main menu selections, including automatic mode, manual mode, adjustments, customize, settings, and device information.
- Press "up" or "down" to select a function, then press "OK" to confirm the selection. Press "back" to go back to the previous screen.

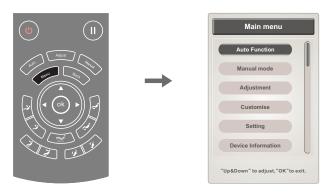

## 5. Automatic function (auto button – pressing the "auto" button will also enable this function.)

- Press "menu" to enter the main menu screen, select the auto function icon, there are 8 automatic modes: fast, comfort, relax, ease ache, fatigue recover, spinal care, pressure relief, waist stretch.
- Press "up" or "down" to select the massage function, press "OK" to confirm and enter the selected massage function. Press "back" to go back to the function setting screen.

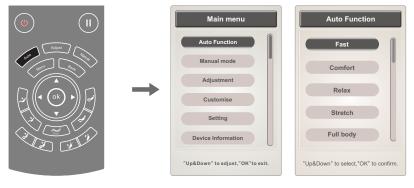


Note: the selected functions will show in the massage function status bar, press "back" to check it directly.


Function	Description		
Fast	Express massage mode for a quick massage.		
Comfort	A full body massage comfortable pressure massage.		
Relax	Kneading and air pressure massage programmed to provide relief from fatigue and stress.		
Stretch	Back and waist targeting massage to loosen stiff muscles and target the lower lumbar and thoracic region.		
Full body	Targeting full body massager to recovery and relief whole body.		
Neck & shoulder	Comfortable neck and shoulder massage, to relieve sore and tired muscles.		
Back & waist	Focused on back and waist to loosen stiff muscles and relieve sore ligaments.		
Full air massage	Focused on upper arm and shoulders, lower arm, knee, leg, and foot to relax the whole body.		

## 6. Manual (manual button – pressing the "manual" button will also enable this function.)

- Press "menu" to enter the main menu screen, select the manual mode icon, including massage mode, massage position, air pressure, etc.
- Press "left" or "right" to switch to the sub menu, then press "up" or "down" to select the massage function, press "back" to return to the status information screen.

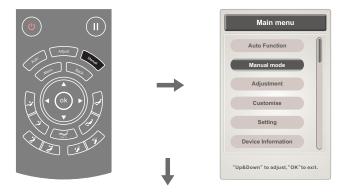


#### Operating instructions

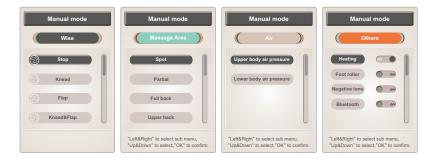

#### 1. Manual mode (manual button - mode)

- Manual mode (mode): including kneading, flap, kneading & flap, shiatsu, tapping, five different massage modes. When in flap, shiatsu, or tapping, the width of massage ball can be adjusted.
- Press "up" or "down" to select the manual function icon, press "back" to confirm starting the selected mode.

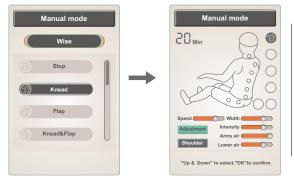


Note: select the stop icon to stop all the manual modes, select the single manual mode icon to start certain massage function alone.

Function	Description	
Kneading	Speed	
Tapping	Speed and width	
Kneading and tapping	Speed	5-level speed adjustable
Shiatsu	Speed and width	5-level width adjustable
Shoulder lifting	Speed	
3D	Speed and width	

#### 2. Manual mode (manual button – massage area)

- Manual mode (massage area): spot, partial, full back, upper back, lower back.
- Press "up" or "down" to select the massage area button, press "OK" to confirm.

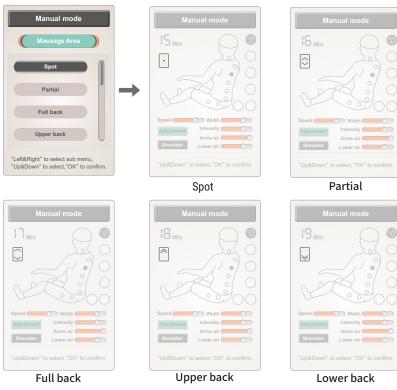


Function	Description		
Spot	Express massage mode for a quick massage.		
Partial	Small area narrow targeting massage	Press "up" or "down" to select the massage area button, press "OK" to confirm.	
Full back	Overall massage of full back		
Upper back	Upper back massage		
Lower back	Lower back massage		

#### Operating instructions

#### 3. Manual mode (manual button – air pressure)


- Manual mode (air pressure): Including upper body, lower body air pressure massage function.
- Press "up" or "down" to choose a different air pressure function, press "OK" to confirm it.

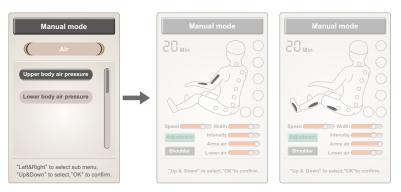

Function	Description		
Upper air	Turn on/off upper air pressure.		
Lower air	Turn on/off lower air pressure.		


#### 4. Manual mode (manual button - others)

- Manual mode (others): heat therapy, foot roller, negative (oxygen) ions, Bluetooth.
- Press "up" or "down" to choose the massage function, press "OK" to toggle on/off the selected function.

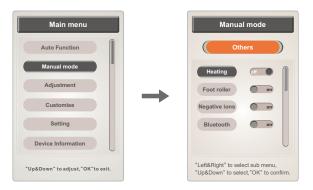


#### Manual mode (other - heat therapy)

- Press "OK" to turn on/off back heating or footrest heating. Once enabled allow 3 min for it to come up to temperature.
- Press "up" or "down" to select other massage functions, press "back" to return to the function selecting screen.

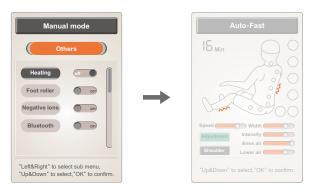

#### Manual mode (other - foot roller)

- Turn the foot roller massage function on/off by pressing the "OK" button.
- Press "up" or "down" to select other massage functions, press "back" to return to the function selecting screen.

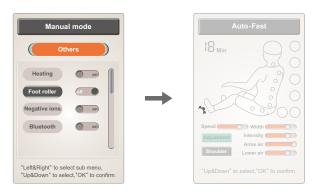

#### Operating instructions

#### Manual mode (other - negative (oxygen) ions)

- Turn the negative oxygen ions function on/off by pressing the "OK" button.
- Press "up" or "down" to select other massage functions, press "back" to return to the function selecting screen.




#### Manual mode (other - Bluetooth)

- Turn the Bluetooth functionality on/off by pressing the "OK" button.
- After the Bluetooth-enabled audio source device (such as mobile phone, tablet, etc.) is paired
  with the massage chair Bluetooth module, the music on the audio source will be wirelessly
  transmitted to the massage chair audio system for music playback via Bluetooth.
- Press "up" or "down" to select other massage functions, press "back" to return to the function selecting screen.




## 7. Adjustment function (adjustment button – you can also directly press the "adjustment button" to enter the adjustment function interface)

- Press "menu" to enter the main menu interface, select the adjustment function button. This
  will adjust massage techniques, massage area, intensity, other (currently timing), and other
  adjustment functions.
- Press the "left" button or "right" button to switch the sub menu function, select the massage function in the menu, press "up" or "Down" to select the function to adjust, press "OK" to confirm, press "back" to return to status screen.

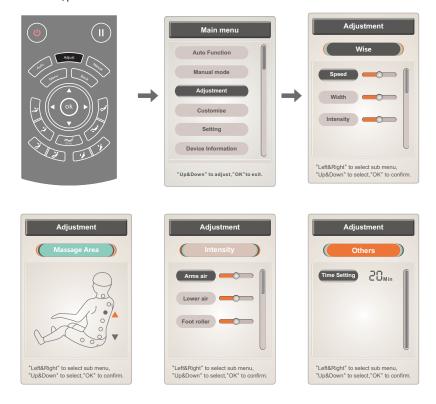

Note: to enter the adjustment menu press the "adjust" button at any point.

#### Operating instructions

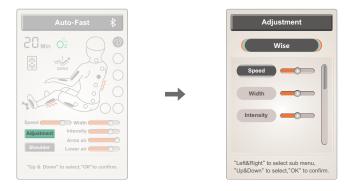

#### 1. Adjust function (adjust button - massage mode)

- Adjustment function (massage mode): adjust massage mode including speed, width, intensity.
- Press the "up" or "down" button to select the speed or width function, and press the "OK" button to adjust.








#### Note:

-Speed: five levels of massage speed adjustment.

-Width: five levels of massage width adjustment when in tapping and shiatsu mode.

#### 2. Adjustment function (adjustment button - massage area)

 Adjustment function (massage area): In the modes of fixed point or partial massage area, press and hold the "up" button to move the massage area up, release to stop. Press and hold the "down" button to move the massage area down, release to stop.

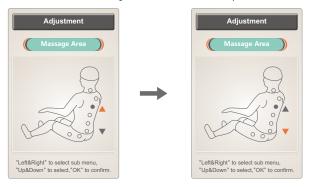

#### 3. Adjustment function (adjustment button - intensity)

- Adjustment function (intensity adjustment): adjusts upper body air pressure, lower body air pressure, foot rollers, with 5 speeds of adjustment.
- Press "up" and "down" to select the air pressure or the wheel roller, and press "OK" to adjust the intensity.








#### 4. Adjustment function (adjustment button – other adjustments)

Adjustment function (other adjustments), Settings for timing currently. Each time you press
the "OK" button it will be adjusted in increments of 5 min. Once the maximum of 40 min is
reached it will return to the minimum amount of 5 minutes. (the cycle will be repeated)




#### 8. Custom

- Press "menu" to select "custom" function under the main menu to enter the custom storage and selection menu.
- Select "save" or "use" with the "left" button or "right" button and select "OK" to confirm selection.
- Select custom Storage I or Custom Storage II with the "Up" or "Down" button.
- Storage: store all the current techniques, pressure, angle, and other massage function settings.
- Call: recall the previously stored techniques, air pressure, angle, and other massage function settings.

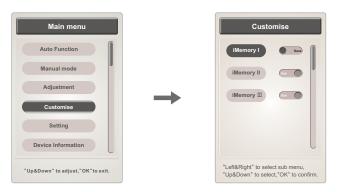


#### 9. Setting

- Press "menu" to enter the main menu screen, select "setting" function button, to enter the settings screen. Adjust language, back light, lock screen, and other settings.
- Press the "left" button or "right" button to switch the sub menu. Press "up" or "down" to select the function, press "back" to return to the status screen.

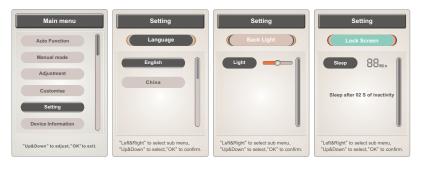

#### 1. Settings (language)

- Settings (language): Sets language to English or Chinese.
- Press "up" or "down" to select the language type, press "OK" to confirm, press "back" to return to the status screen.

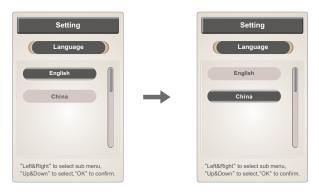


#### Operating instructions

#### 2. Settings (lock screen)

Settings (lock screen): Press "OK" to adjust the lock screen sleep setting time. Each time you
press "OK" the lock screen is adjusted between the following settings: 15 seconds, 30
seconds, 1 minute, 2 minutes, or never lock screen. Press "OK" to enter the loop and press
"back" to return to the status screen.

#### 3. Settings (back light)


Set (back light) brightness, press "enter" to adjust the brightness, a total of 5 level settings.
 Press "back" to return to the status screen.

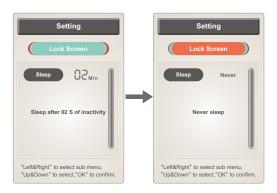



#### 10. Device information

- Press "menu" to enter the main menu interface and select the device information function to display the current model number and other information.
- Press "up" or "down" to select and press "OK" to confirm, press "back" to return to the status screen.

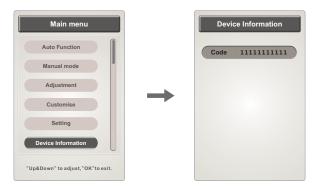

#### 11. Zero gravity button

• Press "zero gravity" to turn on the zero gravity massage function, you can choose "zero gravity I" or "zero gravity II" mode to find the setting that best fits your needs.







Function	Description		
ZEIRO	Zero gravity i: adjust to zero gravity i angle.		
ZERO	Zero gravity ii: adjust to zero gravity ii angle.		

#### 12. Back button

 To return to the main interface screen press this button to display the current massage mode settings and operation status.








#### 13. Position and adjustment buttons

Function	Description
~	Footrest rise button: press and hold this button to slowly adjust the frame up, release to stop rising. After adjusted it will cycle to auto-detect to the length of leg and foot.
<del>~</del>	Lower leg lowering button: press and hold this button to slowly adjust the frame down, release to stop lowering. After adjusted it will cycle to auto-detect to the length of leg and foot.
×	Backrest up button: press and hold this button to slowly raise the backrest, release to stop rising.
2	Backrest down button: press and hold this button to slowly lower the backrest, release to stop lowering.
Ä	Rise button: press and hold this button, the backrest rises slowly, while the footrest frame slowly lowers until the button is released to stop.
7	Lower button: press and hold this button, the backrest lowers slowly, while the footrest fame slowly rises until the button is released to stop.
€ Zero	Zero gravity button: When the button is pressed chair enters zero gravity mode. There are two different zero gravity modes to pick between.
79	Footrest raise button: press and hold this button to manually adjust the footrest to shorten, release to stop.
•	Footrest lower button: press and hold this button to manually adjust the footrest to lengthen, release to stop.

#### 3. Quick access controller

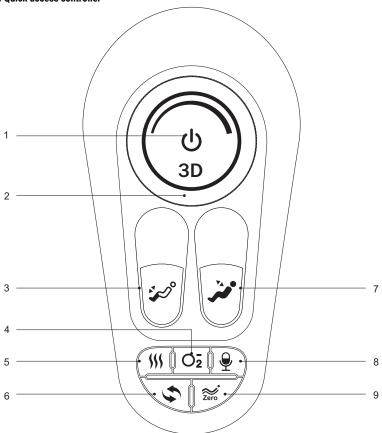

- Power button: press and hold this button to turn all massage functions on or off. It will also pause and ongoing massage program
- 2. 3D intensity button: adjusts the intensity of 3D massage hands.
- Footrest raise/lower button: press and hold this button to manually adjust the footrest to shorten, release to stop.
- Negative oxygen ions button: switch on or off negative O2 ionizer.

- 5. Heat button: turns on or off carbon fiber heating.
- 6. Automatic massage mode cycle: cycles between the automatic massage modes.
- 7. Rise and lower button: adjust the backrest and the footrest frame up or down.
- 8. Voice command: enables voice command.
- Zero gravity button: when the button is pressed chair enters zero gravity mode. There are two different zero gravity modes to switch between.

Function	Description
0	Power button: press 1 second to pause, press 3 seconds to turn on or off. If the functions remains paused for over 20 min, then all massage functions will shut down automatically.
3D	3D intensity button: 3D massage intensity adjustment. Moves the massage rollers inward and outward.
×~~	Footrest raise button: press and hold this button to manually raise footrest, release to stop.
	Footrest lower button: press and hold this button to manually lower footrest, release to stop.
Y4.0	Rise button: press and hold this button, the backrest rises slowly, while the footrest frame slowly lowers until the button is released to stop.
	Lower button: press and hold this button, the backrest lowers slowly, while the footrest fame slowly rises until the button is released to stop.
<b>\$\$\$</b>	Heat button: switches on heating function.
O <sub>2</sub>	Negative oxygen ions button: witch on or off negative oxygen ions function.
<b>9</b>	Voice command: switches on voice command.
\$	Automatic massage mode cycle: cycles between the automatic massage modes.
æro Zero	Zero gravity button: automatically adjusts backrest and footrest into the zero gravity position. There are 2 modes, press again to switch between modes.

#### 4. Voice control instruction


- When the massage chair is turned on or in standby mode, say "hi Alice" to wake up the voice control function.
- After the successful wake-up, the massage chair will reply "I'm here". You now have 6 seconds to respond with one of the voice commands below. If the command is recognized it will answer with the corresponding answer. If the command is not recognized or not heard in 6 seconds, the voice command will automatically turn off and needs to be initialized again.


	Enable voice control			
То а	To awaken say "hi Alice". If successful, the massage chair will reply "I'm here".			
Serial No	Voice control entry Answer entry			
1	Massage on	Answer: ok, massage on.		
2	Massage close	Answer: ok, massage close		
3	Comfort Massage.	Answer: ok, comfort massage		
4	Relax massage.	Answer: ok, relax massage		
5	Full body massage.	Answer: ok,full body massage		
6	Neck and shoulder massage.	Answer: ok,Neck and shoulder massage		
7	Back and waist massage.	Answer: ok,back and waist massage		
8	Stretch massage.	Answer: ok,stretch massage		
9	Fixed massage.	Answer: ok,fixed massage		
10	Change the other mode.	Answer: ok, change the other mode		
11	Go little down.	Answer: ok, go little down		
12	Go little up.	Answer: ok, go little up		
13	Open the air pressure.	Answer: ok, open the air pressure		
14	Close the air pressure.	Answer: ok , close the air pressure		
15	Up the seat position.	Answer: ok, up the seat position		
16	Down the seat position.	Answer: ok, down the seat position		
17	Stop adjustment.	Answer: ok, stop adjustment		

#### Operating instructions

#### 5. Power off, stop massage

- At any time during the massage process press the power key to stop all massage functions.
   The seat and backrest will return to the starting position and all functions of the massage are immediately turned off.
- Turn off the power switch and unplug the machine between uses.

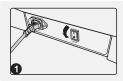

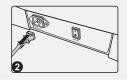




- 1 Power off the main power switch
- 2 Unplug the power cord from the unit
- 3 Unplug from the wall


#### 6. Bluetooth player operation

- Turn on Bluetooth on the massage chair. Press "manual" and navigate by pressing "right" to the others screen. Select "Bluetooth" and press the "OK" button to enable Bluetooth.
- 2. Enter the Bluetooth settings on your smartphone, tablet, or Bluetooth enabled device.
- Pair with the massage chair by selecting the device. Once paired the device will show up under a device id formatted as follows: (IMCM-XXXX).
- 4. Once connected you may now play music from your device.


#### 7. Installing armrest


- Locate the corresponding armrest wiring harness in the center frame section and prepare for connection. (figure 1)
- 2. Match up and connect the wiring to each corresponding plug (follow the markings on each cable). (figure 2).

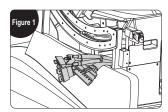

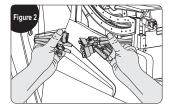



- 3. Connect each plug and the air pipe securely. Once connected tuck excess wire into the center frame. (figure 3)
- 4. Ensure that the air pipe has adequate room and is not being cut off or pinched. (figure 4)

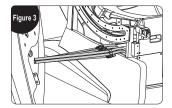

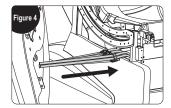



- Hold the handrail with both hands, connect the front end of the armrest into the steel frame, and secure it in place. (figure 5)
- 6. Connect the LED lighting wire at the rear end of the armrest and tuck in the excess wire. Then align the connectors on the rear of the armrest a connect them firmly in place. (note that connectors must be properly spaced for accurate fit) (figure 6)

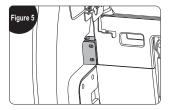

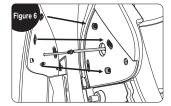



- 7. Secure the rear end of the armrest in place with a M6x20 screw. (figure 7)
- 8. Then secure the upper arm holster and install M5x16 screw. (figure 8)

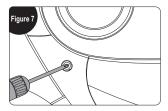

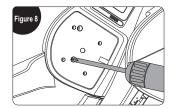


- 9. Align the mounting hardware on the steel frame at the front end of the fixed armrest and secure it with two M5x16 screws. (figure 9)
- 10. Install the screw cover for the rear end screw recess. (figure 10)

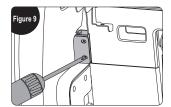

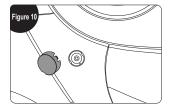


#### 8. Installing footrest


- Remove the tie securing the wiring connection and air pipe connection to the metal for shipping. (figure 1)
- 2. Connect the wiring connector and the air pipe and the leg joint in the lower portion of the frame. Ensure the docking is in place. Then insert the docking cable and the air piper from the lower end of the frame into the lower side of the frame paying attention to the plug. Avoid cutting off the airway when inserting this. (figure 2)

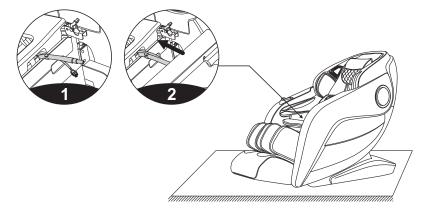


- 3. Insert the rotating inner sleeve and the rotating outer sleeve into the left side of the footrest and the right side of the lower leg (note: the large step must face inward). (figure 3)
- Open the footrest support hinge cover on the left and right mounts of the massage chair. (figure 4)

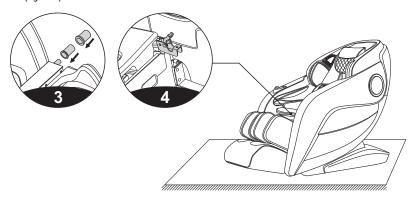

- 5. Hang the lower leg (left and right shaft) into the frame support hinge plat at the same time. (figure 5)
- Cover the rotating shaft cover, press footrest rotating cover, and then tighten footrest rotating cover board with two M5x14-12.9 screws and a flat pad with the supplied hex wrench. (figure 6)

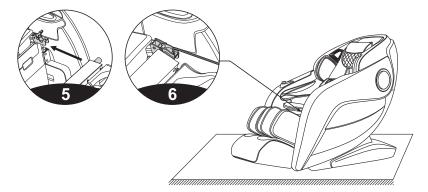



7. After the screw is installed, install the provided footrest cover (figure 7).

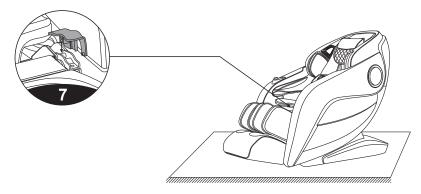

#### 9. Pillow and backrest adjustments

• The head pillow has a three-part (three-layer) fold-down design. Choose the number of layers for the right fit for you. For a gentler massage separate the hook and loop fastener on the headrest to reveal the back pad. The additional head pillow/backrest portion may also be removed by unzipping (1) and/or detaching the hook and loop fastener. (2)

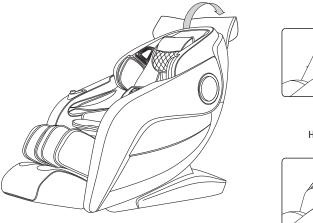


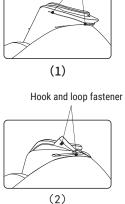

Zipper

#### 10. Floor protection


 Place on carpet or floor protector to avoid damage to surfaces. The massage chair can do damage to sensitive surfaces, such as wood flooring.

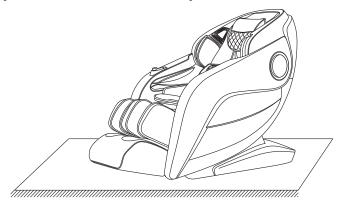

Caution: do not move the chair by the wheels on wood flooring, rough ground, or narrow spaces. Two people are required to lift the massage chair. Remove footrest and armrests.

#### 11. Moving the chair

Before moving detach or secure the power cable. Tilt the chair by the backrest backward until
well balanced at the center of gravity on the wheels. Push the chair forward or backward to
the desired location then set the chair down slowly and carefully.

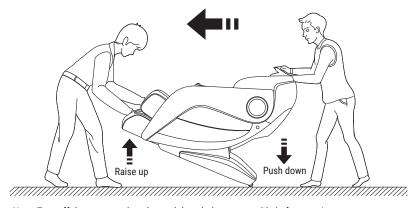

Note: Turn off the power and unplug and detach the power cable before moving.

#### **Product Specifications**

Model: A701 Configure: 2

Name: Massage Chair

Rated voltage: 110-120V~ 60Hz

220-240V~ 50Hz/60Hz ...

Rated power input:220W Rated time: 20Min Safety structure: Class I

Material introduction: PVC, PA, steel parts and electric & electron parts

MATERIAL	PU	LEATHER	CLOTH	WOOD



www.zarifausa.com





