BianStone & Graphene Far Infrared Heating Therapy MagicPad

Product Introduction

This MagicPad Pro adopts a new heat injection technology, natural red needle stone combined with graphene heating cloth to generate heat, human body safe USB-C PD 65W low voltage input, waterproof, fireproof and electric shock-proof. Apply cutting-edge technology PEMF and terahertz, double the energy, give you different warm physical therapy care, and bring you instant relief and relaxationin just 15 minutes. It can be used as a seat cushion or foot cushion, and as a hot compress for shoulders, stomach, and waist; it is a good friend whether you are at home or in the office.


Product material and function


Natural Red BianStone: Bianstone surgery, one of the six medical techniques of Chinese medicine, mentioned in many places in the "Huangdi Neijing", it is called a stone that can cure diseases, using Bianstone as a treatment tool Representative: Chinese medicine Bian Que. After testing by modern instruments, it contains more than 40 kinds of trace elements and minerals that are beneficial to the human body, as well as more than 20 kinds of anti-aging elements such as calcium, magnesium, zinc, chromium, strontium, and Graphene: Nobel Prize in Physics, a revolutionary material for biomedicine in the 21st

tury. Graphene far-infrared rays have the functions of improving blood circulation, ing joint pain, boosting spirit, skin care and beauty. PEMF: World Class Energy Technology. The cells of the human body need energy to achieve optimal functional status. Pulsed magnetic fields are a soothing and energetic method that can help the human body replenish energy, activate the vitality of the Terahertz: It can produce a wave heat effect, increase cell vitality, adjust nerves and

the body, promote body metabolism, improve the human body's magnetic field, open up the meridians in the human body, enhance the human body's mental vitality, and promote the function of sweat glands, thereby accelerating the excretion of waste in the body, beneficial to health.


Product Structure

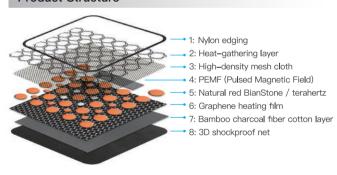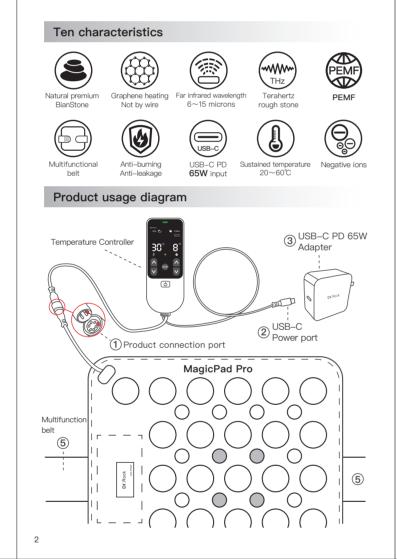
Product connection diagram:

- ① Connect the thermostat six-pin connector (female seat) with the MagicPad Pro six-pin connector (male seat) according to the limit slot and the arrow; (2) Plug the other end of the thermostat USB-C into the power adapter or car charger, and then power on to use;
- 3 Note: The power adapter (included with the product) and the car charger and power bank (please purchase by yourself) must support PD fast charging and the output power of the USB-C single port is ≥65W. If you need to buy, please choose a certified power adapter or car charger, power bank; 4 Please refer to the figure below for vehicle and outdoor use scenarios:

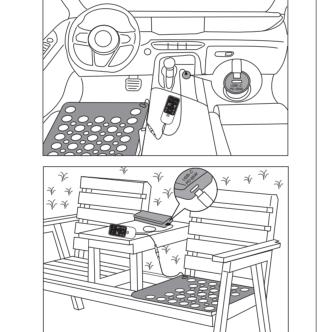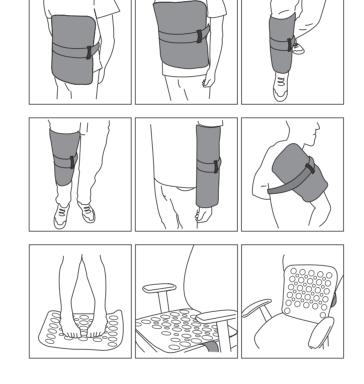4 Multifunctional belt

Controller use and troubleshooting

Timing increase button Timing decrease button - Power switch button default temperature "30°C" is set to "8" hours when it is used for the first time;

Power Indicator

- (PEMF) Night Mode

- Timing unit (hour)

Heating Status Indicator

2. The " icon is for temperature adjustment. Press " or " or " to adjust for temperature rise and fall. Note: After the temperature is set, the displayed temperature will automatically jump to the actual sensing temperature; 3. " $\overset{f c}{f c}$ " is for timing adjustment, press " $ilde{igwedge}$ " or " $ilde{igwedge}$ " key for timing adjustment, the unit is H/hour. Note: After the timer is counted down, it will automatically stop working and shut down when the set time is reached:

4. The icon flashes during adjustment, and the icon is always on for 6 seconds after setting to indicate completion. Note: After the temperature is set, the " 🔆 " icon flashes to indicate heating, and when the set temperature is reached, the " 🔆 " icon disappears to indicate that the heating is stopped: 5. "PEMF" adjustment, the default Night Mode is turned on, press it once for Day Mode, and press it again for off. Note that "PEMF" is turned on for 30 minutes,

and it will be turned off automatically after 30 minutes; 6. The thermostat has a memory function. After the last function adjustment, the default last adjustment parameters will be used again; 7. Temperature mode set up Automatic heating time

safe mode: A.Warm therapy 1 20–45°C Can be set manually for 1–12 hours B.Warm therapy 2 46–50°C Heating for 4 hours, automatically return to 39°C C.Hot therapy 51–60℃ Heating for 4 hours, automatically return to 39℃

5. Thermostat display troubleshooting		
Display	Cause of issue	Troubleshooting
E1	The thermostat is not connected to the product or the connection is in poorcontact; The temperature sensor probe fails.	Check the connection line, whether the connection is normal, and it needs to be re-energized. If the product is still abnormal, please contact after-sales service.
EH	High temperature error or temperature out of protection	Turn off the power, turn on the product again after cooling to confirm, if it shows that the product is abnormal, please contact after–sales

Troubleshooting Suitable for the crowd 1. Check the power supply, whether it is normal EU Abnormal voltage 2. Recommended output power ≥ 65W certified power

1. People with sub-health, low immunity and easy fatigue 2. People with edema, obesity, insomnia, neurasthenia, and poor physical fitness 3. Symptoms such as rheumatism, rheumatoid patients, numbness of hands and feet 4. Cold hands and feet, female physical discomfort, chills, etc.

5. Lumbar and leg pain, shoulder and neck pain, muscle pain and other pain symptoms

Package Contents • MagicPad Pro x1 • Temperature controller x1 • Product storage bag x1

Manual book x1
USB-C PD 65W Adapter x1

Warranty description Dr.Rock provides quality warranty services for hardware products and accessories in the packaging and defects in materials and workmanship. The warranty period of the MagicPad Pro (including temperature controller) is one year

from the date of purchase of the product, and the included storage bag is guaranteed for 6 months. Damages caused by normal wear or aging, acc improper use, alterations, etc. are not covered by this warranty. How to obtain after-sales service: Please call the Dr.Rock authorized brand or go to the authorized retailer's purchase place or online customer service for after-sales service, and please provide the relevant purchase credentials (invoice or purchase record).

Product maintenance methods

You can use a wet towel to wipe the contaminated part when it is not powered on, or use a disinfectant wet towel to disinfect it regularly, and then turn it on for heating after drying.

A Precautions

- 1. Read the manual carefully before use and use it in accordance with the
- 2. Don't let it be used alone by people who cannot take care of themselves, infants or children who are not sensitive to heat.
- 3. It is forbidden to use in water and humid environment 4. Do not dry clean or wash. You can wipe the contaminated part with a wet towel when the power is not turned on, and then turn on the power to heat it after drying.
- 5. Do not insert sharp objects into this product. 6. Do not place bedding and other items in a small area to cover the heating area
- when the product is heated to prevent local overheating and causing malfunctions 7. It is forbidden to place flammable materials such as fireworks and lighters on
- the surface of the product to prevent fire. 8. Please confirm whether the power supply is cut off before going out.

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Using the heating function, the temperature can be adjusted according to the tolerance of the skin. In order to prevent skin burns.if you feel unwell, stop *The temperature controller can set the timing: 1~12 hours *Heating temperature: 20~60°C, adjustable by lifting button *Accessories: 3D shock net, Density mesh cloth, Bamboo charcoal fiber cotton (fireproof)

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adapter or car charger;

Physiotherapy mode

EC Current overcurrent 3. If there is any abnormality after reinserting the power supply and restarting, please contact the after–sales service.

● Daily Warm Compress Therapy: Set the thermostat to 30~40°C, PEMF setting

(use 12Hz day mode during the day, 7.83Hz evening mode at night) for daily use, and recommend heating for 6–10 hours/time.

It is suitable for people with weak and cold constitution, poor microcirculation of the body, long-term fatigue, rheumatic pain, cold hands and feet, and weakened immunity

• Hot Therapy: Set the thermostat The thermostat is set to 50°C-60°C, PEMF is set to

and numbness in hands and feet. It can accelerate blood circulation, sweat and detoxify,

• Natural Therapy: Turn off the power of the thermostat (at room temperature) in hot

*Product name: Dr.Rock BianStone & Graphene Far Infrared Heating Therapy MagicPad Pro

*Physiotherapy mode: Room temperature natural therapy, warm therapy, Hot therapy,

*Heating device: Graphene heating film

weather, or natural care at normal room temperature. Suggested use: Long term. It feels cool and comfortable to use in summer; long-term use can make the body absorb the far-infrared rays that are beneficial to the body, which is conducive to

maintaining and improving the operation of the body's microcirculation system.

*Adapter input voltage: 100-240V~50/60Hz, USB-C PD 20V/3.25A 65W

*Current and voltage: human body safety low voltage maximum 20V/3.25A

PEMF magnetic therapy

*MagicPad Pro input power: USB-C PD 20V/3.25A 65W and above

• Warm Conditioning Therapy: set the thermostat to 41–49°C, PEMF setting (use

12Hz day mode during the day,7.83Hz night mode at night) for daily body

It is suitable for resisting chronic diseases, long-term fatigue, rheumatic pain, dysmenorrhea, cold palace, blood stasis, cold hands and feet, poor body

It is suitable for sports injuries, backache and bone pain, shoulder and neck pair

microcirculation, decreased immunity, edema, obesity, sub-healthy people:

conditioning, and recommend heating for 2-4 hours per time.

12Hz daytime mode, used for hyperthermia conditioning when,

neating use suggestion: 1 hour / each time.

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