

## Charging the device

### Battery Status

Full Charge 

Partial Charge 

Low Charge 

Charge immediately 



The battery icon will indicate low or no power when the device needs to be recharged. When this happens:

1. Insert the charging cord into the power supply connection at the bottom of the device. The charging icon will display on the LCD display screen.
2. Wait until the LCD display screen turns black. This indicates that the device is fully charged.

**NOTE:** DO NOT disconnect the charging cable while the device is charging.

9

## Important Safety Information

**WARNING:** Read this section in its entirety before using the device to avoid injury to yourself and others and to ensure the safety of the device. Follow the following safety precautions whenever operating the device.

The following indicates potentially hazardous situations which, if not avoided, may result in serious injury.

1. Do not use the device on nipples, genitals or any sensitive area such as around the eyes.
2. Do not let the thermodynamic wire touch your skin, always have a tip attached.
3. Only use the device for its intended purposes described in this manual.
4. Never change the tip while the device is turned on. Turn the device off and wait for it to cool completely before doing so.
5. Do not submerge the device in water.
6. Do not use the device near a bath, sink or anything else filled with water.

10

7. Do not continue using the device if it causes pain or adverse skin effects after the first use. Consult a doctor if this persists.
8. The following people should not use the Krasr Blue Light Hair Removal Device:
  - › People with impaired physical, mental or sensory abilities
  - › People with any skin allergies or skin conditions such as dermatitis (eczema) or herpes
  - › People with swelling, sunburns or any other skin injury on the area that will be treated
  - › People recovering from an operation
9. Keep the device away from children at all times. Do not allow children or infants to use the device or come into contact with the device otherwise injury may occur.
10. Do not use the device near flammable substances (alcohol, benzine, arousal sprays, etc.) or explosions and fires may occur.
11. Make sure that the device is turned off at all times when not in use to avoid fires.

11

12. Do not let dirt or any other substances get inside the device.
13. Do not drop the device. If the device sustains damage from a fall, do not use it.
14. Do not attempt to repair, modify or disassemble the device yourself or else injury may occur because of abnormal operation. Contact the store where the device was purchased if you think the device needs to be repaired, or contact Krasr directly
15. Clean the device after each use (only AFTER the device has cooled completely)
16. Properly dispose of the device package after opening it and keep the package and package contents away from children to avoid suffocation.
17. Only use the device with the supplied charger
18. Do not let the device near flammable substances at any time

### SPECIFICATIONS

Weight	Dimensions	Temperature	Humidity	Power Supply
104g	120x40x15mm	5-30 °C	30-80% RH	Input 100-240VAC/50-60Hz 0.2A Output 8.4VDC 0.5A

12

## Warranty Information

The manufacturer (the “warrantor”) warrants to the original purchaser that this equipment will be free from defects under normal and proper operation for a period of one (1) year from the date of purchase. This warranty covers all parts and labor for that one (1) year period. This warranty shall not apply to any equipment that has been subjected to abuse, misuse, and alteration of any type or to any defect or damage caused by repair, replacement, substitution or use with parts other than parts approved by the manufacturer. The manufacturer shall not be liable for the loss of use of this equipment, loss of time or business, inconvenience or any other indirect, consequential, special or incidental damages due to defects covered by this warranty or any implied warranty. To implement this warranty, contact the dealer from whom you purchased the equipment. If necessary, you will be asked to return all, or the appropriate part, of the equipment. If warranted, the repairs will be performed at no charge and the equipment returned to you. Unit owner must prepay all shipping costs when returning items for warranty service. This limited warranty is the only written or express warranty given by the manufacturer. This warranty gives you specific legal rights but you may also have other rights, which may vary from state to state or province to province.

FOR ALL WARRANTY CLAIMS AND OTHER QUESTIONS AND CONCERNS, PLEASE EMAIL [support@krasr.com](mailto:support@krasr.com).

13

## Money Back Guarantee

We’re so confident that our products will help you achieve the results you want that we offer a 30-day money back guarantee on any product. If you are not satisfied with a product for any reason, we will gladly refund your money, no questions asked. Visit

[krasr.com/return-policy](http://krasr.com/return-policy)  
for more details

## Contact Krasr

If you need help choosing a product, are wondering how to use a product, are have problems, or if you have questions about anything else, we’re here and we’re happy to help. If you would like to contact Krasr, please visit [krasr.com](http://krasr.com) or send us an email at

[support@krasr.com](mailto:support@krasr.com)

# KRASR™





## Permanent Hair Removal

# About Krasr Blue Light Hair Removal Device


The Krasr Blue Light Hair Removal Device has been designed to remove hair for long periods of time without the harshness of razors. The thermodynamic wire transmits heat to the root of the hair follicle which not only eliminates existing hair, it slows hair regrowth. The adjustable device strength and interchangeable tips make the device suitable for thin or thick hair and it is perfectly safe to use on many parts of the body, including the face.

**No harsh creams, no chemicals and no lasers. This type of hair removal is kinder to the skin.**


## DEVICE FEATURES

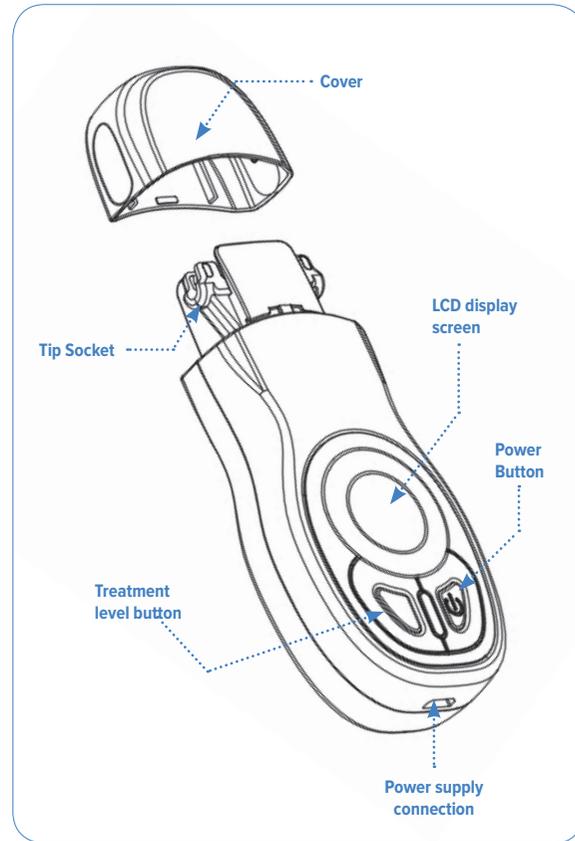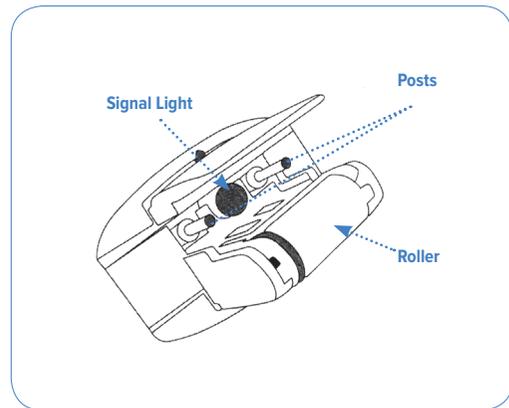


## WHAT'S INCLUDED

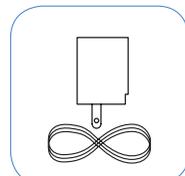


charging cable with outlet adapter

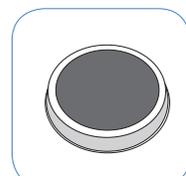

Buffer Pad

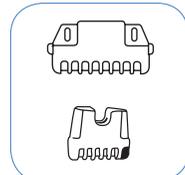

Tips

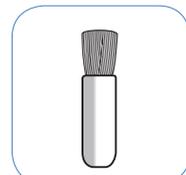

Cleaning Brush

# Using the Krasr Blue Light Hair Removal Device

## INSTRUCTIONS FOR USE

1. Charge the device for 5 hours uninterrupted before using the device for the first time (read section "charging the device" to learn how).
2. Make sure skin is clean and dry
3. Make sure that the device is not plugged in, otherwise it will not work.
4. Remove the device cover (figure 3)
5. Select the narrow tip for small areas like the face, neck and underarms or the wide tip for larger areas like the arms or legs . (figure 4)
6. Insert the tip into the device: Make sure the logos are facing the same direction and line up the posts in the tip socket with the holes in the tip and gently push the adapter into place. (figure 5,6)



Figure 3

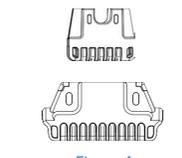


Figure 4



Figure 5


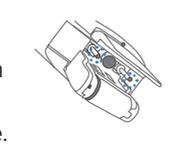


Figure 6



7. Press and hold the on/off button for 2-3 seconds, until the LCD display lights up.

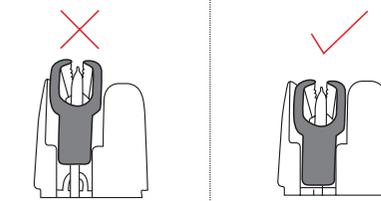



8. Adjust the device intensity: Press once for thin hair, twice for normal hair and three times for thick hair.



Treatment Level Button

- 9.

**NOTE:** The first time you use the device, use it on the lowest setting before increasing it to test the effects on your skin. When the device is turned off and on again it will remain on the last used treatment level.

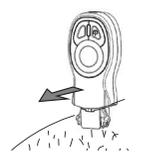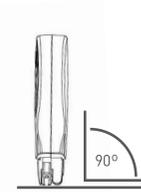

10. Press the adapter gently against the skin on the area that you wish to treat and hold the device at a 90 degree angle with your skin.

11. Gently move the device in one direction over the hair. Repeat 2-4 times for the same area until it has been treated, then move onto another area.

12. Keep the device steady. The blue light will flicker if you are moving the device too slow and you will need to adjust the speed accordingly.

13. **DO NOT** press any buttons on the device while using it

14. A red light will activate if you are not using the device properly. If this happens, restart your movement until the blue light turns on.

**NOTE:** It is normal for skin to feel warm and for you to smell slight burning while the device is working because the device is heating the hair follicles to stop hair regrowth. If the heat is too intense, discontinue use immediately.

15. When you are finished, use the buffer to gently brush the treated area to remove the excess hair that has been removed (move back and forth gently).

16. Apply moisturizer to the treated area to soothe skin.

17. **Wait for the device to cool completely**, then use the cleaning brush to clean the excess hair from the tip.

# Tips for use

1. If treating loose skin, stretch it so that it is tight when being treated. This will ensure best results.
2. When treating the **upper lip**, start inward and move outward and **DO NOT** touch the lip or let the skin bunch up into the tip.
3. To treat the **chin and jawline**, start just below the chin and move upward moving with the curve of the face.
4. To treat sideburns, start from the base of the sideburns and glide upwards.

**NOTE:** It will take 2-3 months, with regular use, for hair regrowth to stop. During the first 4 weeks, use the device 4 times each week. During the second 4 weeks, use the device 3 times per week, and gradually decrease the usage each week from then on. You will notice the hair becoming thinner as you use the device and that less hair will grow, so there will not be a need to use the device as often.