

KRASR™ 3 IN 1 NANO IONIC FACIAL STEAMER



IMPORTANT

Please read all sections of this manual in their entirety before using the Krasr 3-in-1 Nano Ionic Facial Steamer and blackhead extracting tools for the first time and follow all procedures related to safety, device maintenance and skin care before, during and after using the device. This will ensure the safety of yourself and the device and will optimize results.


GETTING STARTED

About the Device

The Krasr Nano Ionic Facial Steamer is a new type of facial steamer that uses nano-scale steam produced from an ultra sonic vaporizer. This type of steam penetrates the skin and opens pores up to 10 times more than steam produced from heating elements on traditional steamers. As a result, blackheads will be easier to remove and you will be able to remove more of them, and skin will be left soft to the touch.

A 5-piece blackhead extractor kit is included with 5 blackhead extractor tools that are made from medical grade ingredients.

Key Features

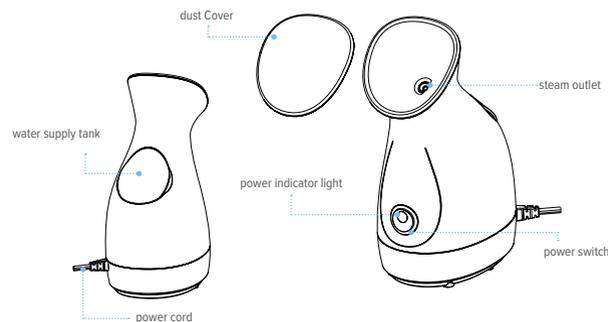

2

What's Included:

BLACKHEAD EXTRACTOR KIT



EYE HOLE



ANGLED



NEEDLE



OVAL



FLAT WIRE

- > (1) Nano-ionic facial steaming device
- > (1) Water pouring cup
- > (1) Power cord and user manual

USING THE DEVICE

Preparing the Steamer

Note: Water may be found in the water supply tank before the first use. This water is leftover from the factory inspection and is normal.

1. Place the steamer on a flat, level and solid surface and keep it there at all times while using it.
2. Remove the water supply tank from the back of the steamer, and twist the cover to remove it.
3. Fill the water supply tank with distilled or purified water.
4. Insert the water supply tank back into the unit. Water is supplied to the main unit through the water supply tank and the water level will decrease when the supply tank is reinserted.

3

5. Insert the electrical plug into a power outlet.
6. Press the power button until the blue light appears. Steam will emit from the steamer after about 40 seconds

Note: A small amount of hot water may be sprayed from the device before steam starts emitting. Be cautious of this. It is not a malfunction.

Using the Steamer

7. Position your face at least 20cm away from the device, with the chin pointed towards the nozzle to ensure that steam covers the entire face.
8. Steam the face for 10 minutes each time.
9. Steam the face before using blackhead extractor tools or before exfoliating.
10. As water in the supply tank decreases, the way that steam emits from the steamer will change. This is normal. If there is no water in the supply tank, no steam will emit and the temperature control will be turned off. Allow the steamer to cool for 10 minutes before refilling the supply tank and using the steamer again.

After using the Steamer

11. Turn off the steamer: Press the button on the steamer so that the blue light goes out.
12. Now is the time to exfoliate or use the blackhead extractor tools.
13. To use the blackhead extractor tools:
 - a. Choose the desired tool based on your specific needs.
 - b. Position the tool so that blackheads are in the centre of the hole.
 - c. Press gently until the blackhead comes out of the pore.
14. When you are completely finished treating your skin, apply moisturizer. This will lock in the moisture from the steamer and give you added benefits, otherwise the moisture will evaporate and leave your skin dry

4

15. Allow the steamer to cool down for at least 10 minutes before putting it away or using it again.

16. Remove the power cord from the outlet.

17. Discard whatever water remains in the supply tank and the main device.
To prevent burning yourself, hold the supply tank and the main device so that your hands are not exposed to the hot water. Do not leave water in the supply tank or the main unit while not in use.

18. Insert the water supply tank back into the unit until it clicks.

19. Store the device only after it has dried sufficiently.

Device Cleaning

1. Fill the water supply tank with cool water, attach the cover to the supply tank and insert it into the back of the steamer.

2. Gently, and without tilting the steamer, shake it horizontally 3-4 times.

3. Discard the water from the main supply tank and the main unit.

4. Repeat steps 1-3 as many times as desired.

5. Using a DRY cloth, wipe the outside of the steamer (only when it has cooled) and the power cord to remove dust.

IMPORTANT SAFETY INFORMATION

WARNING: Read this section in its entirety before using the device to avoid injury to yourself and others and to ensure the safety of the device. Follow the following safety precautions whenever operating the device.

5

⊘ Do not submerge the device in water (Do not use the device near a bath, sink or anything else filled with water), touch it with wet hands or store it in rooms with high levels of humidity (such as a bathroom).

⊘ The device should be used for no more than 10 minutes each time.

⊘ The following people should not use the Krasr Ionic Facial Steamer:

⊘ People with impaired physical, mental or sensory abilities

⊘ People with swelling, sunburns or any other skin injury on the area that will be treated

⊘ Keep the device away from infants at all times. Do not allow children or infants to use the device or come into contact with the device otherwise injury may occur.

⊘ Do not continue to use the device if steam emits in an unusual way.

⊘ Do not use the device near flammable substances (alcohol, benzine, arousal sprays, etc.) or explosions and fires may occur.

⊘ Make sure the device is unplugged when the water tank is being filled or when the device is not being used, otherwise electrocution or other injuries may occur.

⊘ Make sure that the device is turned off at all times when not in use to avoid fires.

⊘ Never use the device if the power cord has become damaged, modified, bent or twisted in any way.

⊘ Do not drop the device. If the device sustains damage from a fall, do not use it.

⊘ Do not insert objects into the steam outlet or block the steam outlet or explosions may occur.

⊘ Do not attempt to repair, modify or disassemble the device yourself or else injury may occur because of abnormal operation. Contact the store where the device was purchased if you think the device needs to be repaired, or contact Krasr directly

6

⊘ Use only distilled or purified water with the device, otherwise equipment malfunction may occur or, in some cases, burning or scalding.

⊘ Always keep face and hands at least 20cm away from the device nozzle, including when turning the device on and off.

⊘ Keep the device on a flat, level and solid surface when using it. The device is not intended for handheld use

⊘ Do not attach, remove or install the water supply tank and do not move the device for at least 10 minutes after each use as steam continues to be emitted and the temperature must cool down.

Clean the device after each use (only AFTER the device has cooled completely) otherwise grease or other residues may accumulate. Do not shake the device while cleaning it.

PRODUCT INFORMATION

Specifications

Model	SK-0806
Power	280W
Frequency	60HZ
Capacity	70ML
Performance Standard	GB4706.1 GB4706.15 GB4343.1 GB17625.1

7

WARRANTY INFORMATION

This limited warranty is the only written or express warranty given by the manufacturer. This warranty gives you specific legal rights but you may also have other rights, which may vary from state to state or province to province.

FOR ALL WARRANTY CLAIMS AND OTHER QUESTIONS AND CONCERNS,
PLEASE EMAIL

support@krasr.com

Contact Krasr

If you need help choosing a product, are wondering how to use a product, are having problems, or if you have questions about anything else, we're here and we're happy to help.

If you would like to contact Krasr, please visit krasr.com or send us an email at

support@krasr.com

KRASR™ NANO IONIC
FACIAL STEAMER