

WARNINGS

1. For sensitive skin, place a piece of voile fabric to avoid excess stimulation.
2. Do not hold on one part of the skin for over 10 seconds.
3. Pregnant women and people with pacemakers are unfit to use this product
4. Make sure that the probe is inserted fully into the device before turning it ON. Turn the intensity down to the lowest setting when you are finished using the device, then turn the device OFF
5. Individuals with freckles might not receive maximum results using this device.

TECHNICAL STANDARD

VOLTAGE	FREQUENCY	POWER
110V, US plug	50-60Hz	10W

IMPORTANT SAFETY NOTICE

1. Don't disassemble the device by yourself. All the repairing should be performed by trained staff.
2. Please don't assemble or operate this device while your hands are wet or when you are near water.
3. If any of the conditions below occur, unplug device immediately:
 - The device touches liquid
 - The device has an abnormal smell, is smoking, or producing
 - loud and strong noises.
 - The cable is damaged
 - The device has physical damage
4. Keep the cable out of the way of others so no one trips on it.

5. Please pull out the power cable after using it for your own safety
6. Don't insert anything into the device's jack, this is a hazard that can potentially cause a fire. If anything gets into the device, contact the manufacturer. Don't try to get it out yourself.
7. Please don't place the device on anything that can't stand steadily.
8. Our company reserves the right to modify the information in this introduction without notification!

CONTACT KRASR

Warranty

Included with the purchase of this device is a 1 year limited warranty from Krasr.

For warranty conditions, please visit krasr.com/warranty

Contact Information

If you would like to contact Krasr, please visit krasr.com or send us an email at support@krasr.com

KRASR RADIANCE



KRASR™


RADIANCE
HIGH FREQUENCY
FACIAL DEVICE

USER MANUAL

IMPORTANT

PLEASE READ ALL INFORMATION IN THIS MANUAL BEFORE USING THIS PRODUCT!

WHAT'S INCLUDED

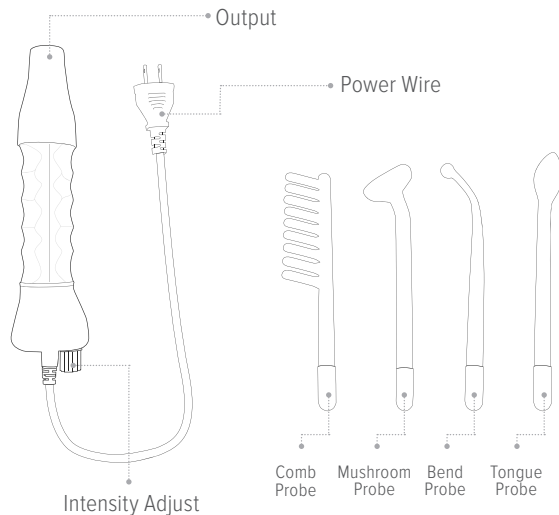

EFFECTS AND BENEFITS

The high-frequency facial device is an important beauty instrument. A high frequency current is produced in the glass probe and when applied to the face, produces an oxygen pocket that has an antibacterial effect. This effect eliminates viruses and bacteria.

1. It helps to calm the nerves and has an analgesia function.
2. It produces the ozone for sterilization, can help to heal wounds and to restrain pores.
3. Accelerates blood circulation and improve metabolism.
4. This device has 2 anti-aging applications, the anti-wrinkle application and the skin-tightening application.

OPERATION PROCEDURES

The Bend Probe: The Bend Probe zaps pimples and reduces their appearance on the skin in a matter of days.

1. Attach the Bend Probe into the device.
2. Make sure the device is on the lowest setting, then turn it ON.
3. Apply the end of the bend probe to pimples and “zap” them individually.
4. Once finished, turn the device down to the lowest setting then turn it OFF.

The Mushroom Probe: Treats broad areas such as cheeks, the jawline, and the forehead. Blood circulation and lymph activity is accelerated to tighten skin and reduce wrinkles and fine lines.

1. Attach the Mushroom Probe into the device.
2. Make sure the device is on the lowest setting, then turn it ON.
3. Apply the end of the Mushroom Probe to the broad areas of the face. This will help those areas of the face get better blood flow. It will also lower the visibility of wrinkles on the face. Keep the device moving along your face. Do not hold on one part of the face for over 10 seconds.
4. Once finished, turn the device down to the lowest setting then turn it OFF.

The Comb Probe: This method is used for stimulating the surface of the skin, reducing folliculitis, and improving hair growth.

1. Attach the Comb Probe into the device.
2. Make sure the device is on the lowest setting, then turn it ON.
3. Move the hair tube in the direction of your hair, brushing gently.
4. Once finished, turn the device down to the lowest setting then turn it OFF.

The Tongue Probe Treats wrinkles and fine lines in smaller areas, more sensitive areas. It is mostly used for treating the wrinkles under the eyes and on the ridges of the nose

1. Attach the Tongue probe into the device.
2. Make sure the device is on the lowest setting, then turn it ON.
3. Apply the end of the probe to under your eyes and the ridges of the nose to make the wrinkles there less visible and to get rid of the fine lines. This probe can also be used for crow's feet.
4. Once finished, turn the device down to the lowest setting then turn it OFF.