



**melodysusie**  
MAKE BEAUTY A NEW SMART



Hot & Cold Dual Facial  
Massager

User Manual

## ● Product Introduction

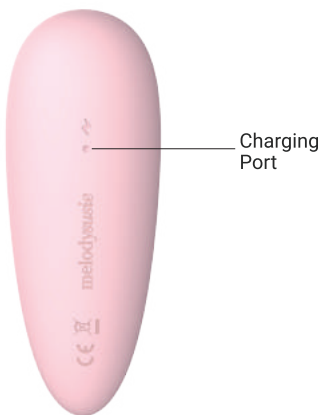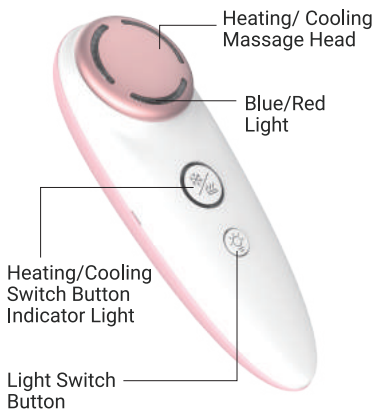

## ● Product Parameters

Model:	MS - 601
Materials:	Silica Gel + ABS + Aluminum alloy
Size:	13.2 x 5 x 4.2 CM
Weight:	165g
Battery:	1000mAH 3.7V
Charge Time:	3H ( DC 5V 2A )
Wavelength:	Red light 620-660nm Blue-ray 460nm
Temperature:	Heating $\leq 42^{\circ}\text{C}$ , Cooling $\geq 15^{\circ}\text{C}$

## ● Charging Method

### 1) USB cable charging

The battery is not fully charged before shipping. Please fully charge the battery before first use.

### 2) You can charge the device at any USB port of Computer or Adapter.

LED lights flicker continuously in the process of charging, and will stop flickering when fully charged.

## ● Instructions

### Cooling Mode







(Figure 1: Standby)



(Figure 2: Cooling)



(Figure 3: Cooling ,  
Blue light off)

- 1) Long press the button  for 1.5 seconds to turn on device, and the indicator light would turn on as shown in Figure 1.
- 2) Press the button  once for Cooling mode, and the blue light will turn on at the same time as shown in Figure 2.
- 3) Press the button  to turn on/off the blue light as shown in Figure 3.
- 4) Long press the button  for 1.5 seconds to turn off device.

## ● Instructions

### Heating Mode







(Figure 4: Standby)



(Figure 5: Heating)



(Figure 6: Heating,  
Red light off)

- 1) Long press the button  for 1.5 seconds to turn on device, and the indicator light would turn on as shown in Figure 4.
- 2) Press the button  twice for Cooling mode, and the red light will turn on at the same time as shown in Figure 5.
- 3) Press the button  turn on/off the red light as show in Figure 6.
- 4) Long press the button  for 1.5 seconds to turn off device.

## ● Instructions

### Red/Blue Light Mode






(Figure 7)



(Figure 8)



(Figure 9)

- 1) Press the button  to turn on the red light as shown in Figure 7.
- 2) Press the button  once more to turn on the blue light as shown in Figure 8.
- 3) Press the button  to turn off the light as shown in Figure 9.

## ● **Cleaning & Maintenance**

- **Please make sure that the device is powered off when not in use.**
- **Please clean the device with cotton pads or soft fabric.**  
DO NOT use diluent, solvent and other corrosive industrial products to clean sonic head in case of damage.
- **Please keep the device in a dry and cool place.**
- **Keep away from children.**