


How to use

1. Slide the Power Switch down to the **ON** position and follow the audio instructions.
2. Press the **First Checks** green button to begin.
3. Once the first checks are complete, select one of the procedure buttons.
4. Pause the audio at any time by pressing the **Pause & Resume** button.
5. Resume the audio by pressing the **Pause & Resume** button again.
NOTE: The audio will resume from the beginning of the sentence in which it was paused.
6. You can restart the audio from the beginning at any time by pressing the procedure button. A different procedure can also be started at any time by pressing the required button.

Power **ON/OFF** switch

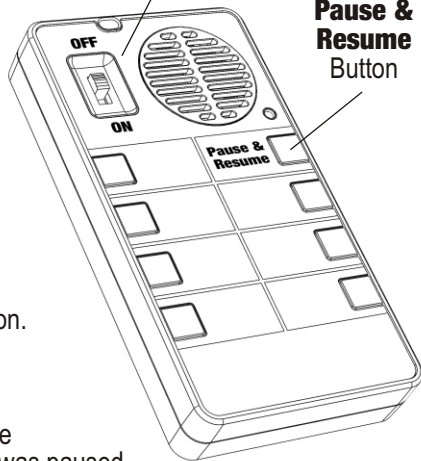


Pause & Resume
Button

FIRST RESPONSE

First Aid Audio Instructions

USER GUIDE

V1.0

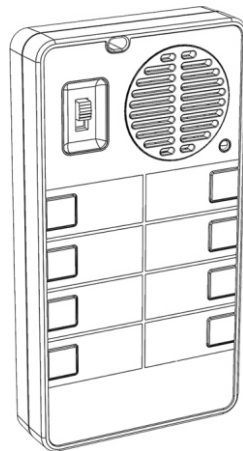

Disclaimer

The procedures and protocols contained on this device are compliant with the Resuscitation Council (UK) guidelines 2015. This product is not intended to replace the need for professional help. The user assumes all risk of liability.

Intelligent**FirstAid**[™]


How to Install Batteries

Slide the power switch to the **OFF** position. The battery compartment is located on the rear of the device. Use a screwdriver to remove the screw. Install new batteries with the correct polarity as below. Refit the cover and the screw.

Requires:
3 x AAA 1.5V
Alkaline Batteries
(not included).

Safety Information

1. Not suitable for children under 3 years. Do not let children play with batteries.
2. Use only new batteries of the stated size and type.
3. DO NOT mix old and new batteries.
4. DO NOT recharge, heat or dispose of batteries in a fire.

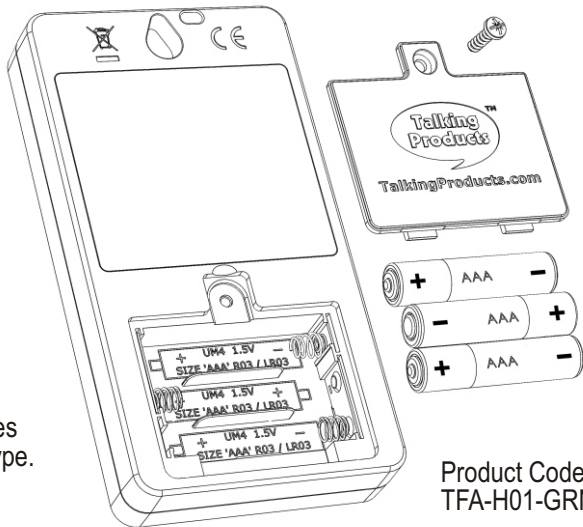

Product Code:
TFA-H01-GRN



Support & Warranty

Please keep this document for future reference. 12 months warranty is provided. Opening the case will invalidate your warranty.

Contact Us

For assistance and technical support.



Talking Products Ltd
Unit C8 The Premier Centre
Abbey Park Industrial Estate
Romsey, Hampshire
SO51 9DG, England

Tel: +44 (0) 1794 278327
Email: sales@TalkingProducts.com

www.TalkingProducts.com