glo

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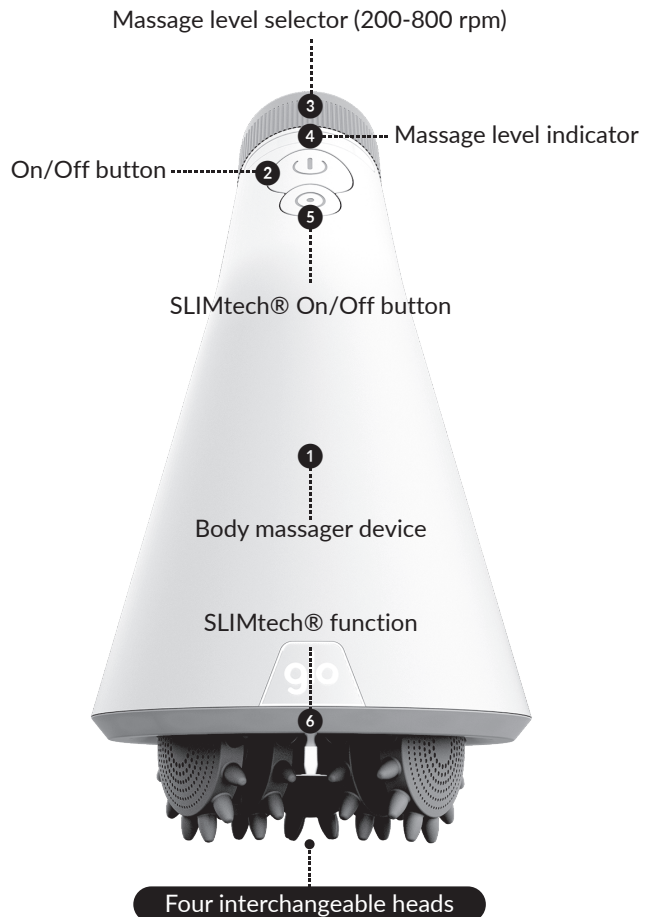
2.USING THE APPLIANCE

2.USING THE APPLIANCE

- Plug the device into the mains.
- Place the selected treatment head.
- Switch on the device by pressing the power button **2**.
- If desired, switch on the SLIMtech® technology by pressing the button **5**.
- Hold the device firmly, preventing it from turning around on itself, as this could end up damaging the cable.
- Select the speed (from 200 rpm to 800 rpm) using the regulator **3**. Start with a gentle speed, and increase according to the sensitivity to massage.
- Pass it slowly in circular movements, avoiding pressing it against the skin with excessive pressure.
- To drain, use the Cellu-Detox head **8**, slide it with upward movements from the ankle to the groin, avoiding touching the bones.
- To change the head, always turn off the device by pressing the power button **2**.

Inserting and removing the head

- To position the head, insert the selected head shaft on the shaft at the base of the device, and secure it by pressing the central button of the head **11**. Pressing or forcing the anchoring system is needless.
- To remove the head, press the central button **11**. Once released, pull it outwards.


Massage Duration

- The massage lasts about 15-20 minutes
- It is recommended to dedicate about 5 minutes per area to be treated. In the case of the legs, spend about 10 minutes on each leg (thigh + buttock).
- If you want to dedicate a longer time to specific areas, we recommend dividing the treatment into two areas--upper and lower--to allow enough time for each and then alternating the days of treating them.

Direction of the massage

The massage should be done in the direction of blood circulation (from the extremities towards the heart), to promote the elimination of toxins through the lymphatic system.

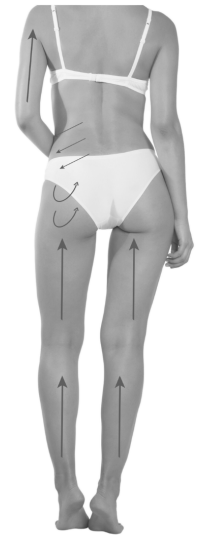Auto-safety shut-off

- The device has a safety system to control the treatment time, which is activated after about 20 minutes of use.
- This sensor is not only activated by time, but also according to the pressure exerted on the skin, the intensity of the treatment, the fat and hardness of the tissue to be massaged and the time of the massage.
- If used too intensely, the safety system may activate before 20 minutes.
- Remember, it is important NOT TO PRESS THE APPLIANCE WITH EXCESSIVE PRESSURE ON THE SKIN.

Recommendations

- It is important to **be consistency**, we recommend use it daily for at least the first month of treatment. Then depending on your personal situation, the treatment can be spaced out from 2 to 3 days a week. Do not stop using it regularly, because the body tends to return to its previous state if it is not taken care of. One of the great advantages of having Glo910+ is that you can enjoy long-lasting results through consistent use.
- It is important to DRINK WATER to help the body drain (at least 1.5 liters/day).
- DO NOT USE THE CELLU-ACTIVE ⑦ AND CELLU-DETOX ⑧ HEADS ON THE BONES.
- AVOID MASSAGING ON THE SPINE.
- In delicate areas, such as the neckline and neck, use the SOFT-MASSAGE head ⑩ with low massage intensity.
- The BODY-BRUSH exfoliating head ⑨ is designed for the body, do not use it on the face, neck or décolleté.

- WHEN MASSAGING THE DECOLLETE, NECK OR TRAPEZIUS, TUCK THE HAIR AWAY TO PREVENT IT FROM ACCIDENTALLY GETTING CAUGHT IN THE DEVICE.
- DO NOT PRESS THE DEVICE WITH EXCESSIVE PRESSURE ON THE SKIN. The motor reaches up to 800 rpm, it has enough power to perform the massage properly.
- At the beginning of the treatment, when using the Cellu-Active head ⑦, bruising may occur, as there may be a large accumulation of encysted cells; this is normal and will progressively improve with the treatment.
- SLIMtech® technology does not emit heat, has no side effects, can be used on the whole body and with all heads.
- It is advisable to change the Cellu-Active treatment head ⑦ when you notice that the tips are rounded due to wear, so as not to lose effectiveness.
- Complete the treatment with Glo cosmetics and nutricosmetics to enhance the result. Discover all Glo products at www.glo910.com

Ending the treatment

- Once the massage is over, turn off the device by pressing the power button ②.
- The device will stop.
- Then, unplug it from the mains.
- Clean it and store it in its bag.
- It is advisable to apply a cream or oil after the massage to moisturize the skin.

3.CLEANING AND UPKEEP

3.CLEANING AND UPKEEP

3.CLEANING AND UPKEEP

3.CLEANING AND UPKEEP

3.CLEANING AND UPKEEP

3.CLEANING AND UPKEEP

3.CLEANING AND UPKEEP

3.CLEANING AND UPKEEP

3.CLEANING AND UPKEEP

3.CLEANING AND UPKEEP

3.CLEANING AND UPKEEP

3.CLEANING AND UPKEEP

- After unplugging the device:
 - Use a dry cloth to clean the exterior.
 - Remove the head to clean any dead skin cells that may remain between the seals. Use the brush provided in the accessories to clean it.
 - Wipe the inside of the base of the device with a dry cloth.
- Do not use solvents, detergents or abrasive products.
- When not in use, store it in the bag.

Cleaning the massage heads

- It is recommended to clean the heads with a small brush or a dry cloth.
- Cleaning Soft-Massage head 10: The sponges can be removed and cleaned with neutral soap and lukewarm water. First gently remove the rubber band holding the sponge, then the sponge. Wait for them to dry completely before replacing them on the head shafts. Place the sponge first and then the rubber.

Upkeep

- When you are not using the device, store it in its bag.
- Do not twist the cable, pull it, or wrap it around the unit to avoid damaging it. Use the Velcro that comes with the device for proper storage. If damaged, replace it immediately by qualified personnel or contact customer service.


Those instructions are also available on our website, www.glo910.com, along with videos and more details on how to get the most out of your Glo910 device. Follow us on Instagram to stay up to date with all the information.

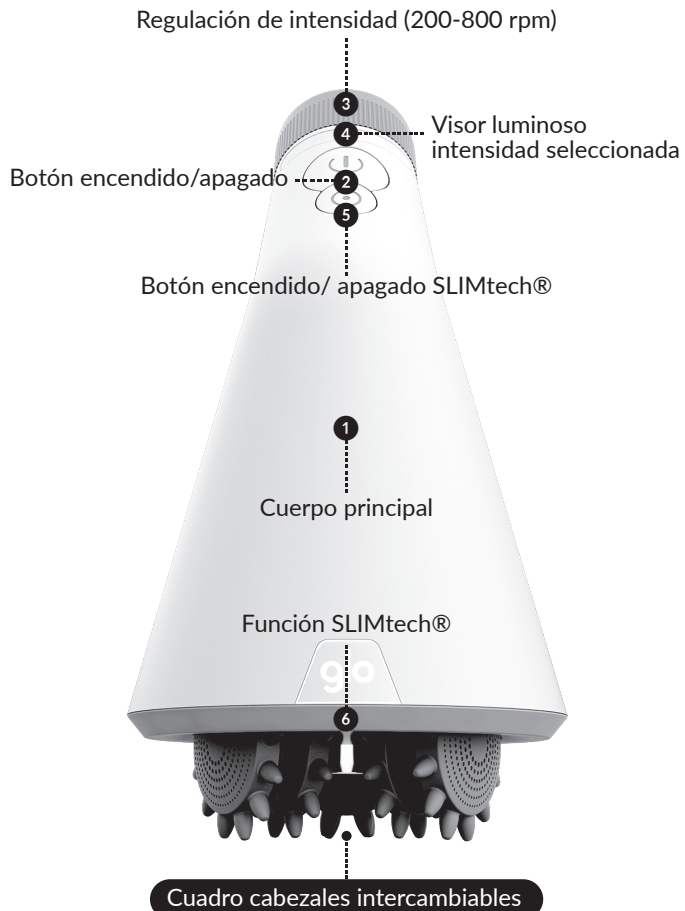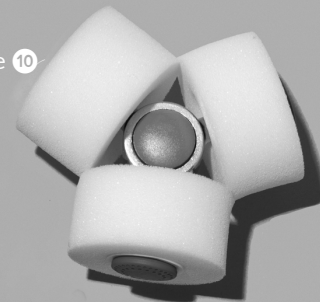
Non-contractual photos: product may differ from photo.

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