



# Nicefeel<sup>®</sup>

# Oral Irrigator

## INSTRUCTION MANUAL


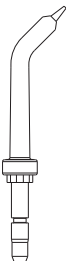
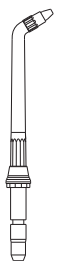


**Please read this manual carefully before use!**

**Keep this manual for reference as it contains important information!**

Model NO: FC168

# Product Structure



				
Ordinary High-Pressure Tip	Periodontal Tip	Orthodontic Tip	Dental Plaque Tip	Tongue Cleaner Tip

- Ordinary High-Pressure Tip: 3 color ring. For hygienic purpose, recommend using by 3 people respectively(family pack).

The other 4 special tips

- Tongue Cleaner Tip x 1
- Periodontal Tip x 1
- Orthodontic Tip x 1
- Dental Plaque Tip x 1


See more instructions on page 10 - 11

## Instructions

Plug the cord into an electrical outlet. If the outlet is controlled by a wall switch, be sure it is turned on.


## Preparing for Reservoir

Remove the lid & Reservoir from the base. Fill the Reservoir with lukewarm water. Replace the reservoir over the base and press down firmly.

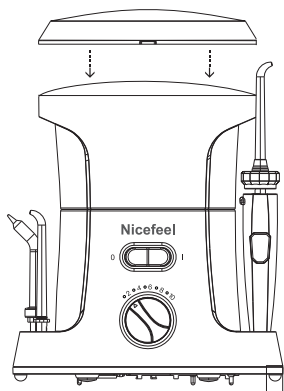


## Inserting and removing tips

Insert the tip into the center of the knob at top of the dental jet handle, the colored ring will be flush with the end of the knob if the tip is correctly locked into place. To remove the tip from the handle, press the jet tip eject button and pull the tip from the handle.

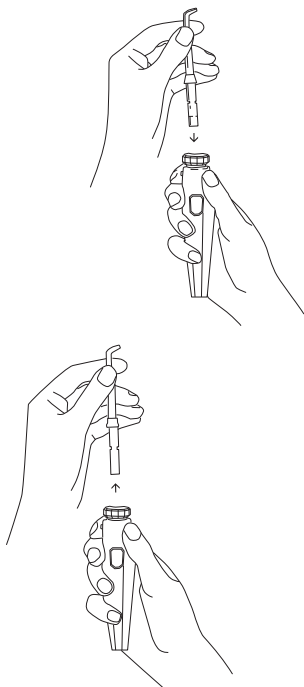

## Adjusting the Pressure Settings

Turn the pressure control dial on the base of the unit to the lowest setting for first-time use. Gradually increase pressure over time to the setting you prefer, or as instructed by your dental professional.

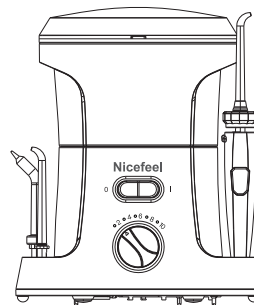


## Using the Oral Irrigator

Lean low over sink and place the jet tip in mouth. Aim the tip toward teeth. With your free hand, turn the on/off button on. Set pressure at desired setting.







Direct the stream by turning the knob at the top of the handle.

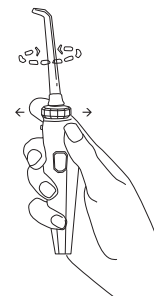

**This appliance should be kept out of reach of children at all times.**

## Recommended technique

Direct the jet stream at 90-degree angle to your gumline, slightly close lips to avoid splashing but allow water to flow freely from mouth into the sink. For best results, start in the molar area (back teeth) working towards the front teeth. Glide jet tip along gumline and pause briefly between teeth. Continue the process until all areas around and between the teeth have been cleaned.

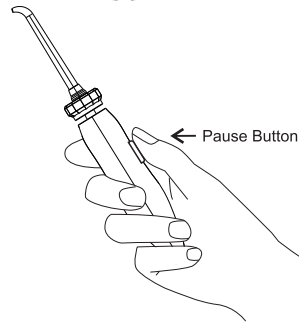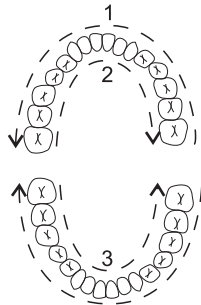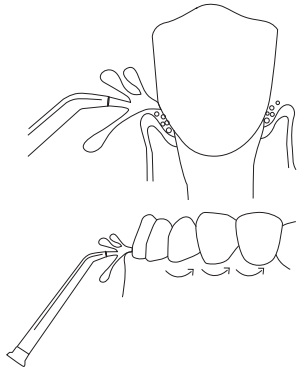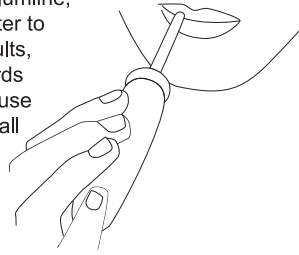

## Pause Control

You can temporarily stop the flow at any time by pressing the pause button on the handle.

## When finished


Turn the unit off. Empty any liquids left in the reservoir. Unplug the unit.

## Special Tips Instruction

### Periodontal Tip


The Periodontal Tip delivers water or anti-bacterial solutions deep into periodontal pockets.

To use the Tooth Pocket Tip, lean low over your sink and place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline. Turn the unit on to the lowest pressure setting and clean gently following gumline, inserting tip in between teeth. Continue tracing along the gum-line until you finish cleaning.

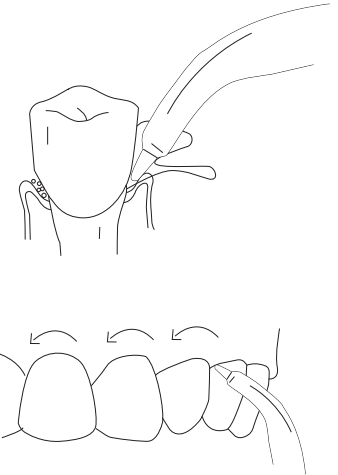

### Tongue Cleaner

Lean over the sink. Place the tongue cleaner on your tongue and start the unit - use the lowest pressure setting, increase as comfortable. Gently place the tongue cleaner in the center/midline of your tongue about half way back and with light pressure pull forward.

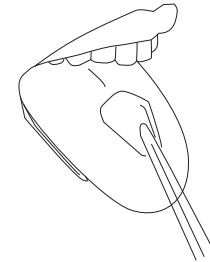




## Orthodontic Tip

The special Orthodontic Tip is uniquely designed to simultaneously brush and rinse hard-to-reach areas around braces and other dental work.


To use the Orthodontic Tip, lean low over the sink and use the brush bristles to clean around dental work. Water will flow through the brush to rinse away any debris that is removed.

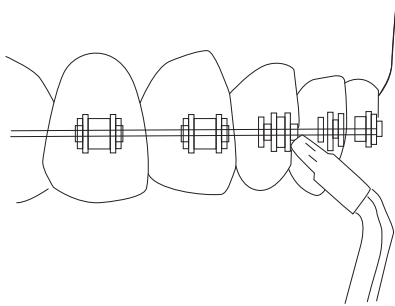

## Dental Plaque Tip

Direct the Jet Stream at 90-degree angle to your gumline. Close lips slightly to avoid splashing but allow water to flow freely from mouth into the sink. Start in the molar (back teeth) area working toward the front teeth. Place the tip close to the teeth so that the bristles are gently touching the teeth.

Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth. Continue until all areas around and between the teeth both inside and outside, top and bottom that have been cleaned.

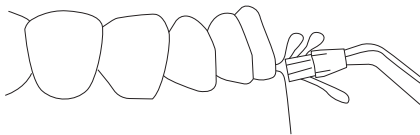

## Cleaning

After using any special solution in the reservoir, take a moment after every use to rinse the unit.

Partially fill the reservoir with warm tap water, aim the jet tip into your sink, turn on the unit and let it run until the reservoir is empty. This will help keep the tip from clogging and help your unit to last longer.

Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser.

Before cleaning, please unplug from the electrical outlet.

Before exposing the dental jet to freezing temperatures, remove the reservoir and run the unit until completely empty.

## Service Maintenance

1. After using the Nicefeel Oral Irrigator, be sure to unplug the power cord, and the tank upside down on a flat surface, when the water is completely dry, put it into a cool, dry place together with the cleaned and dry host. It should be placed in the box or similar storage when not in use for a long time.
2. When a machine malfunctions, you can contact the seller to inquiry related maintenance methods.
3. Do not disassemble or attempt to resolve this product, resulting in any problem, not within the Company's warranty.
4. When the power cord is damaged, go to our designated service point or related maintenance shop for maintenance. Do not replace the power cord or plug yourself.
5. If the product is out of the warranty period, we will provide repairing service as well, we will charge the cost of materials appropriately and a small amount of labor cost.

## Common problems and solutions

Problem	Cause	Solution
Jet Tips and Tongue Cleaner Leakage between the top and handle knob	Tip not fully engaged	Remove tip and reinstall
Periodontal Tip Tear in soft tip	Tip is worn	Replace with new tip
Oral Irrigator Base Hose to the handle leaks	Damage to the hose	Can be replaced by the seller for a self- service repair kit
Reservoir leaks	Missing water tank valve or water tank valve is upside down	Replace Reservoir valve (contact the seller); Turn the valve right side up
Inadequate pressure	Reservoir not positioned properly	Re-position water tan, pr- ess reservoir down firmly
Unit doesn't start to work	Motor not running	Be sure outlet is working, test with another applia- nce (check electrical reset button and/or wall switch)

## Warranty

We warrant to the original purchase / owner of this new product that it is free from defects in materials and workmanship for one year from date of purchase. We will repair or replace any part of the product, which in our opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase and was used according to these instructions. This limited warranty excludes the standard Jet Tips, Periodontal Tip, Tongue Cleaner Tip or Dental Plaque Tip attachments.

**Thank you again for choosing the Nicefeel Oral Irrigator!**


# Thank you for choosing Nicefeel Oral Irrigator!

Welcome! By selecting Nicefeel Oral Irrigator you are on your way to better gum health and a healthier smile!


With the combination of water pressure and pulsation, Nicefeel Oral Irrigator removes harmful bacteria deep between teeth and below gumline where brushing won't reach.

## Traditional Clean

Traditional brushing can't easily reach hidden food debris and harmful bacteria

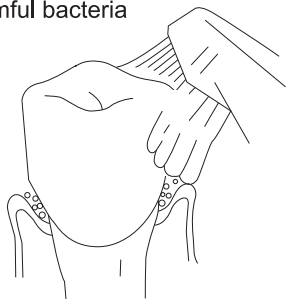


Flossing has limited access below the gumline

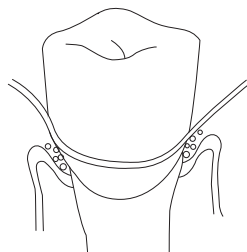

## Nicefeel Oral Irrigator Clean

Removes trapped food debris and harmful bacteria with high frequency pulse water floss

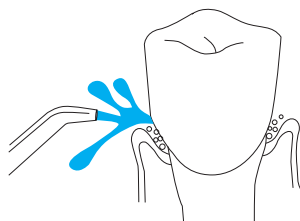

## Content

Important Safeguards	3
Product Structure	5
Getting started	7
Professional Tips Functions	8
Cleaning and Maintenance	10
Troubleshooting Guide	13
Warranty	14

### Accessories:

Main Unit (including Reservoir)	x	1
Instruction Manual	x	1

### Tips

Ordinary High-Pressure Tip	x	3
Periodontal Tip	x	1
Orthodontic Tip	x	1
Dental Plaque Tip	x	1
Tongue Cleaner	x	1

# ◆ IMPORTANT SAFEGUARDS

When using electrical products, especially

**When children are present, the basic safety precautions below should always be followed.**

## **Danger - To reduce the risk of electrocution:**

- ◆ A product should never be left unattended when plugged in. Always unplug product after use. Do not handle plug with wet hands.
- ◆ Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- ◆ Do not use while bathing. Never use while sleeping or drowsy.
- ◆ Do not immerse in water or other liquid.
- ◆ Do not reach for a product that has fallen into any liquid, unplug immediately.
- ◆ Check the power supply cord for damage before the first use and during the life of the product.
- ◆ Power cord is unchangeable, this appliance is to be disused if the cord is damaged.
- ◆ Do not use attachments not recommended by us. Only use attachments (including Jet Tips, Tooth Pocket Tips, Tongue Cleaner Tips and Orthodontic Tips that are recommended by us.
- ◆ Fill reservoir with water or other recommended solutions only.
- ◆ Do not drop or insert any foreign object into any opening or hose.
- ◆ Keep electrical cord away from heated surfaces.
- ◆ Do not operate where oxygen or aerosol sprays are being used.
- ◆ For indoor use only.

## SAVE THESE INSTRUCTIONS

## **WARNING - To reduce the risk of burns, electric shock, fire, serious injury:**

- ◆ Do not leave the appliance when plugged in. Unplug from outlet when not in use.
- ◆ Do not use this product if it has a damaged cord or plug and if it is not working properly. Such as if it has been dropped or damaged, or dropped into any liquid. Stop using the product immediately!
- ◆ Closely instruct and supervise children and individuals with special needs in the proper use of this product. Do not allow to be used as a toy.
- ◆ Do not direct water under the tongue, into the ear, nose, or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas. Pressure generated by this product may cause serious damage to these parts. See "Operating Instructions" for correct usage.
- ◆ Use as described in this manual only. Use our recommended attachments only.
- ◆ Do not use while wearing any oral jewelry, remove prior to use.
- ◆ Do not use if you have an open wound on your tongue or in your mouth.
- ◆ If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- ◆ Do not pull or carry by cord, do not use cord as handle, do not close a door on cord, or pull cord around sharp edges or corners. Do not place / run appliance over cord. Keep cord away from heated surfaces.
- ◆ Do not unplug by pulling on cord. To unplug, grasp the plug, not the cord.
- ◆ Do not handle plug or appliance with wet hands.
- ◆ Turn off the product before unplugging.
- ◆ It's normal phenomenon that caused gum bleeding when you use the water flosser in the first time, which depended on different gum bearing ability.