

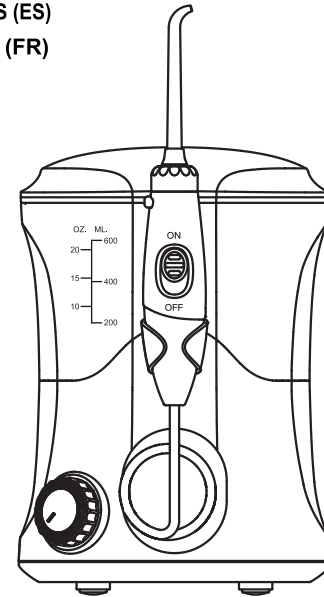
Oral Irrigator

- INSTRUCTION (EN)
- ISTRUZIONI (IT)
- ANLEITUNG (DE)
- INSTRUCCIONES (ES)
- INSTRUCTION (FR)

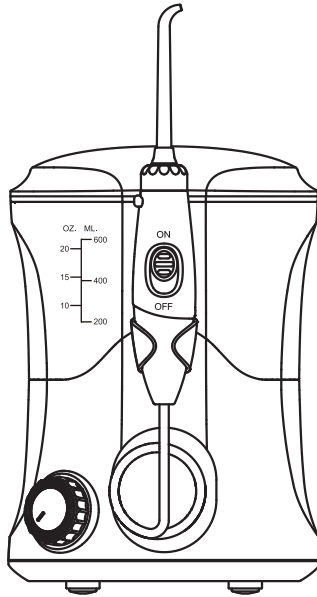Oral Irrigator Usage Manual

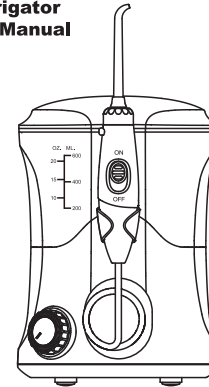Model no.: FC169

After-sales services

Thank you for purchasing the Oral Irrigator. Please note the following:


1. At the time of purchase, please retain the effective sales document issued by the retailer to serve as a maintenance document.
2. A two year product malfunction warranty for non-human malfunctions is provided for this product in accordance with the relevant state provisions. If the product was purchased via Amazon, documentation is not required, and only the order form number needs to be provided. In order for us to provide corresponding after-sales. Free maintenance (based on the invoice date for purchase of the device);
3. Where the following circumstances arise during the maintenance period, costs and services fees shall be charged based on circumstances;
 - The customer causes damage due to inappropriate usage or maintenance;
 - Repairers who are not covered by third-party repair cause damage during repairs;
 - Lack of a relevant sales document or more than two years have elapsed since the date of manufacture;
 - Damage caused during usage for commercial purposes;
Damage caused by force majeure events (such as storms, earthquake or fire).

4. Please carefully fill in the following form, and deliver the device to the company's designated after-sale repairer for repair work.


Maintenance and repair record

Full name of customer		Contact number	
Product model		Product colour	
Purchase model		Purchase platform	
Contact address			
Maintenance date			

Product structure

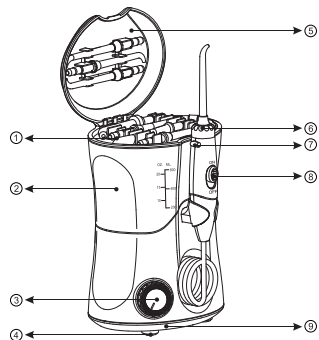1. Reservoir
2. Water tank
3. On/off switch/ pressure adjustment button
4. Base suction cup
5. Water tank cover
6. 360-degree jet tip dial
7. Jet tip disassembly button
8. Handle switch
9. Base

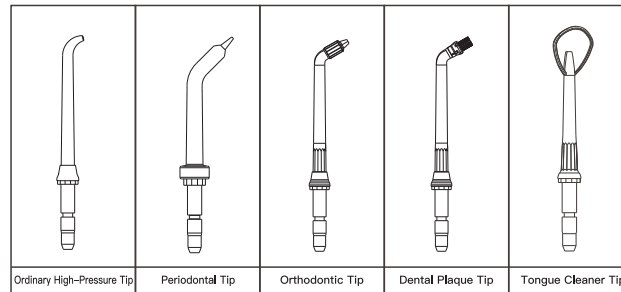- Ordinary High-Pressure Tip
- Periodontal Tip
- Orthodontic Tip
- Dental Plaque Tip
- Tongue Cleaner Tip

- The ordinary high-pressure tip comes with three color rings, for three different people to use (household pack).


There are four other special-use jet tips (Tongue Cleaner Tip, Periodontal Tip, Dental Plaque Tip, Orthodontic Tip)

- Periodontal Tip x 1
- Dental Plaque Tip x 1
- Orthodontic Tip x 1
- Tongue Cleaner Tip x 1


For usage instructions, refer to pages 7 to 8 of this instruction manual.

Preparation for usage


Insert the power plug of the water floss device into the power supply.

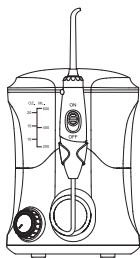Once the water tank is filled with water, place the water tank back into the device, as well as slightly press on top. Ensure that the rubber seal of the water tank base is aligned with the rubber seal of the device base.

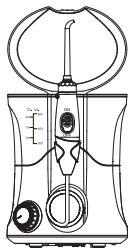
Insertion and removal of the spray tube

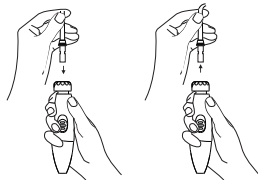When inserting the jet tip into the handle, you should hear a crisp clicking sound. When removing the jet tip from the handle, you only need to lightly press the disassembly button of the jet tip, and remove the jet tip from the handle.


Water pressure adjustment

When you hear a click while turning the on-off switch, this indicates that the device has been switched on. Turning in a clock-wise direction will increase the water pressure. If the device is switched on for more than 3 minutes, the device will automatically turn off. For the next time switching on, it is necessary to rotate the on/off switchback until a clicking sound is heard, then turn it in a clockwise direction again.

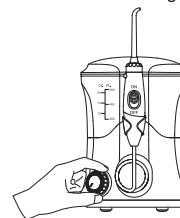Using the Oral Irrigator


Place the jet tip in the oral cavity and direct it towards the teeth. The other hand switches on the device, and sets the pressure at the required value (it is recommended to use the minimum pressure for your first use). The direction of the water flow can be controlled by rotating the 360 degrees button on the top of the handle.

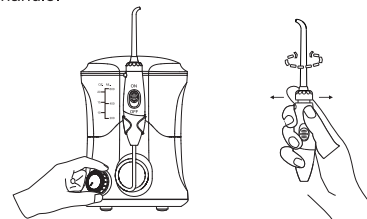Usage Technique

When the water flow and the gums are at a 90 degree angle, slightly close your mouth in order to avoid the water spraying out, yet let the water flow from the corner of the mouth into the sink. In order to achieve better cleaning results, please start flossing from the molars to the incisors. Let the jet tip lightly shift along the gums, and pause briefly between teeth. When water flossing teeth, let the water flow pass through the gaps between the teeth. Repeat the same methx

ods with the each side of the teeth, until the spaces between the teeth, the fronts and backs of the teeth, and the tops and bottoms of the teeth are washed clean.

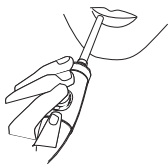
Pause control

You can push the switch on the handle at any time to pause the Oral Irrigator.

Following completion of closing

Switch off the power source, remove the plug, and wipe away any remaining water in the water tank.

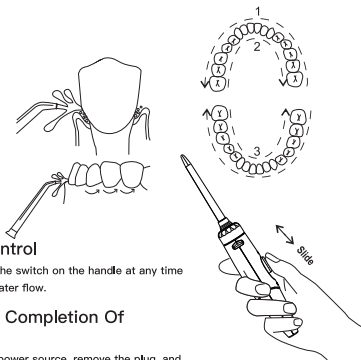Pause Control

You can push the switch on the handle at any time to pause the water flow.

Following Completion Of Closing

Switch off the power source, remove the plug, and wipe away any remaining water in the water tank.

Introduction to the jet tips


Periodontal Tip

The periodontal tip can direct water or liquid oral medicines in specific locations beneath the gum lining.

Ordinary High-Pressure Tip

Flosses deep locations within the oral cavity with a 360 degree range, supplementally covering places which are hard to clean, food scraps and harmful bacteria located deep within dental cracks or beneath the gum. This provides you with a pleasant dental experience.

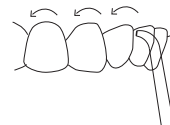
Tongue Cleaner Tip

Place the tongue scraper on the tongue before switching on the Oral Irrigator. Lightly place the tongue scraper on the centre-rear part of the tongue, lightly pulling it forwards while keeping it in contact with the surface of the tongue.

Dental Plaque Tip


Dental plaque is a group of bacteria that stick to the surface of the teeth yet has not mineralized. The jet tip is very effective at preventing dental plaque. The special-use dental plaque tip is for cleaning areas which is difficult to reach in braces and other dental devices.

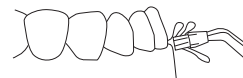Orthodontic Tip

The special Orthodontic Tip is uniquely designed to simultaneously brush and rinse hard-to-reach areas around braces and other dental work.

To use the Orthodontic Tip, lean low over the sink and use the brush bristles to clean around dental work. Water will flow through the brush to rinse away any debris that is removed.

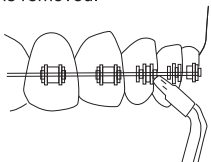Washing:

If any special liquid medicines are used in water tank, please take a while to wash out the device.

Pour some warm water into the water tank, and direct the jet tip at the water sink. Switch on the device and let it work until the water tank is empty. This can help avoid blockage of the jet tip, and extend the usage life of the device.

Use a cloth and neutral non-abrasive cleaning agent to clean the product.

Prior to the cleaning, please ensure that the device is not connected to a power source.

Maintenance and Upkeep

- Each time after using the Oral Irrigator please remove the power plug, and place the water tank on a flat surface. Once the water is fully drained or wiped dry, place it together with the cleaned device in a dark, dry place. When out of use for long periods of time, place it in the package or in a storage box.
- When the device malfunction, proceed to a retailer for advice or call the after-sale services for relevant repair method.
- Please do not disassemble or attempt to dismantle the product. The company does not provide a warranty for any malfunctions caused by that.
- When the power line is damaged, please proceed to the company's designated repairer or qualified repair store for repairs. It is strictly forbidden to change the power line or plug independently.
- When the Oral Irrigator's warranty expires, the company will continue to provide you with services by charging at only materials costs and a small amount of worker services fee.

Common Malfunctions and Handling Methods

Malfunction	Cause	Handling method
Water leakage between the jet tip and handle button	Jet tip hasn't been installed properly	Insert the jet tip into the handle until you hear a clicking sound
Periodontal pocket brush head is split	Jet tip suffers from severe damage	Replace with a new head
Tubing that connects the water flosser base with the handle leaks water	Tube is damaged	Retailer replaces or provides after-sales repair
Water tank leaks	A water tank valve hasn't been installed or the water tank valve has been installed in the wrong direction	Replace the water valve or remove the valve and reinstall it correctly
Pressure is insufficient	Water tank hasn't be installed correctly	Re-install the water tank, and place light downward pressure
Device cannot be switched on	The engine isn't operating	Ensure that the power source is connected, or contact the retailer for repairs

Our company guarantees the original users and owners of this new product that there are no defects with product materials or manufacturing techniques, and the warranty period is one year from the date of purchase.

We thank you again for purchasing the Oral Irrigator!

Contents:

Important safety points: 3
 Product structure: 5
 Product usage instructions: 8
 Special jet tip usage instructions: 11
 Cleaning and maintenance: 13
 Common malfunctions and handling methods: 14

Gentle reminder

Accessories included in this product:

Main device (including water tank) x 1
Instruction manual x 1

Jet Tips

Ordinary High-Pressure Tips x 3
Tongue Cleaner Tip x 1
Periodontal Tip x 1
Dental Plaque Tip x 1
Orthodontic Tip x 1

Thank you for purchasing the Oral Irrigator


The Oral Irrigator uses a pulsing water flow to clean the teeth, and can remove harmful bacteria hidden deep within dental cracks and the gum – places which cannot be cleaned via brushing alone.

Traditional oral hygiene:

Traditional toothbrushes are unable to remove scraps of food or harmful bacteria that are hidden deep.


Dental floss can only reach limited area beneath the gum.

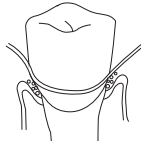The Oral Irrigator

Clean out harmful bacteria hidden in dental cracks and under the gum line.

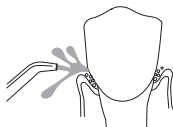♦ IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING

DANGER– To reduce the risk of electrocution:

- ♦ Always unplug this product immediately after using.
- ♦ Do not use while bathing.
- ♦ Do not place or store product where it can fall or be pulled into a tub or sink.
- ♦ Do not place in or drop into water or other liquid.
- ♦ Do not reach for a product that has fallen into water. Unplug immediately.


WARNING– To reduce the risk of burns, electrocution, fire, or injury to persons.

- ♦ A product should never be left unattended when plugged in.
- ♦ Close supervision is necessary when this product is used by, on, or near children or invalids.
- ♦ Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- ♦ Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center for examination and repair.

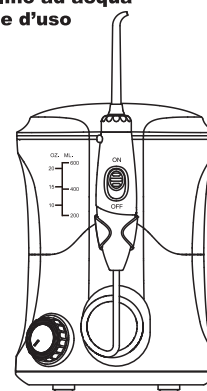- ◆ Keep the cord away from heated surfaces.
- ◆ Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
- ◆ Never use while sleeping or drowsy.
- ◆ Never drop or insert any object into any opening or hose.
- ◆ Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- ◆ Connect this product to a properly grounded outlet only.
- ◆ See Grounding Instructions.
- ◆ Unplug this product before filling. Fill (reservoir) with water only unless otherwise specified by manufacturer. Do not overfill (or specify filling instructions).
- ◆ It's normal phenomenon that caused gum bleeding when you use the water flosser in the first time, which depended on different gingival sensitivity. It is recommended to start at a lower level to use the water flosser.

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