
Nicefeel

# Professional Oral Irrigator Instruction Manual

ENGLISH .....	02
DEUTSCH .....	16
FRANÇAIS .....	30
ITALIANO .....	44
ESPAÑOL .....	58
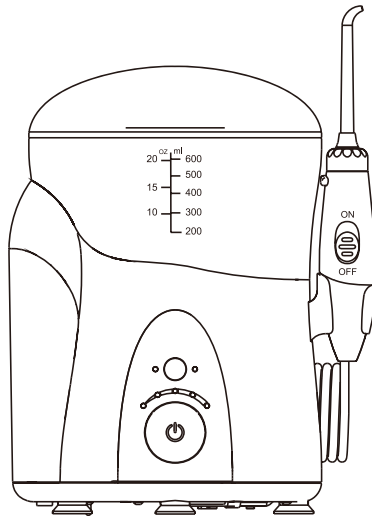
## Thank you for choosing NICEFEEL Oral Irrigator!



You are on your way to better gum health and a healthier smile!

With the combination of water pressure and pulsation, NICEFEEL Oral Irrigator removes harmful bacteria deep between teeth and below gum line where brushing won't reach.

### Traditional Clean


Traditional brushing can't easily reach hidden food debris and harmful bacteria

Flossing has limited access below the gum line
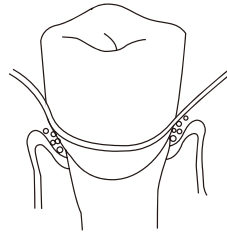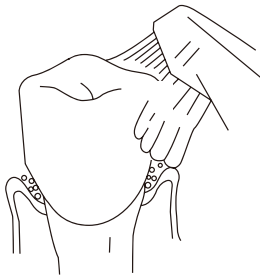
### NICEFEEL Oral Irrigator

Removes trapped food debris and harmful bacteria with high frequency pulse water floss
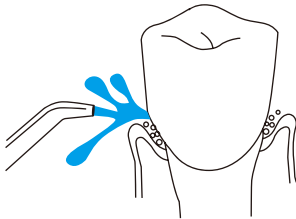
## Content

Important Safeguards.....	04
Product Structure.....	06
How to Use.....	08
Special Tips Instruction.....	12
Cleaning.....	14
Service Maintenance.....	14
Common Problems and Solutions.....	15

## Accessories:

Main Unit

Instruction Manual

## Tips

### **Standard Tips \*3**

Periodontal Tip \*1

Orthodontic Tip \*1

Dental Plaque Tip \*1

Tongue Cleaner \*1

## IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

### Danger – To reduce the risk of electrocution:


- A product should never be left unattended when plugged in. Always unplug products after use. Do not handle plug with wet hands.
- Do not place or store product where it can fall or be pulled into water, shower stall or sink.
- Do not use while bathing. Never use while sleeping or drowsy.
- Do not immerse in water or other liquid.
- Do not reach for a product that has fallen into any liquid, unplug immediately.
- Check the power supply cord for damage before the first use and during the life of the product.
- Power cord is unchangeable, this appliance is to be disused if the cord is damaged.
- Do not use attachments not recommended by us. Only use attachments (including Jet Tip, Tooth Pocket Tip, Tongue Cleaner Tip and Orthodontic Tip that are recommended by us).
- Fill reservoir with water or other recommended solutions only.
- Do not drop or insert any foreign object into any opening or hose.
- Keep electrical cord away from heated surfaces.
- Do not operate where oxygen or aerosol sprays are being used.
- For indoor use only.

## SAVE THESE INSTRUCTIONS

## WARNING – To reduce the risk of burns, electric shock, fire, serious injury:

- Do not leave the appliance when plugged in. Unplug from outlet when not in use.
- Do not use this product if it has a damaged cord or plug and if it is not working properly. Such as if it has been dropped or damaged, or dropped into any liquid. Stop using the product immediately!
- Closely instruct and supervise children and individuals with special needs in the proper use of this product. Do not allow to be used as a toy.
- Do not direct water under the tongue, into the ear, nose, or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas. Pressure generated by this product may cause serious damage to these parts. See “Operating Instructions” for correct usage.
- Use as described in this manual only. Use our recommended attachments only.
- Do not use while wearing any oral jewelry, remove prior to use.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- Do not pull or carry by cord, do not use cord as handle, do not close a door on cord, or pull cord around sharp edges or corners. Do not place / run appliance over cord. Keep cord away from heated surface.
- Do not unplug by pulling on cord. To unplug, grasp the plug, not the cord.
- Do not handle plug or appliance with wet hands.
- Turn off the product before unplugging.
- It's normal phenomenon that caused gum bleeding when you use the water flosser in the first time, which depended on different gingival sensitivity. It is recommended to start at a lower level to use the water flosser.

## PRODUCT STRUCTURE

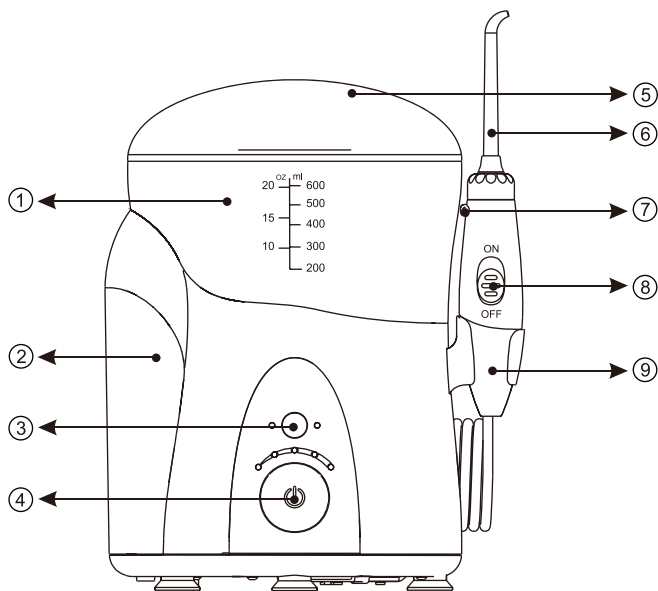1. Water Tank

2. Storage Box

3. Mode Button

4. Power Button/Pressure Control Button


5. Cover

6. Tip

7. Release Button

8. Pause Button

9. Handle

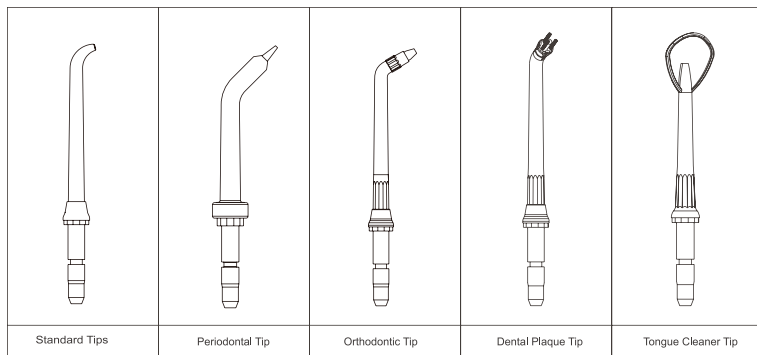**Standard Tips:** 3 color ring.

For hygienic purpose, recommend using by 3 people respectively (family pack).

## The other 4 special tips:


- Periodontal Tip \*1
- Orthodontic Tip \*1
- Dental Plaque Tip \*1
- Tongue Cleaner \*1

## HOW TO USE

Plug into socket. If the socket is controlled by a wall switch, make sure the switch is turned on.

### Preparing for water tank


Remove the water tank from the base, fill the water tank with water, put the water tank back to the base and press down firmly and put on the lid.

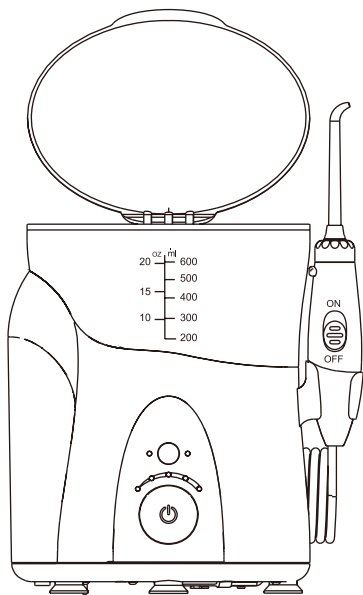
## Inserting and removing tips

Insert the tip into the center of the knob at top of the dental jet handle, the colored ring will be flush with the end of the knob if the tip is correctly locked into place. To remove the tip from the handle, press the release button and pull the tip from the handle.


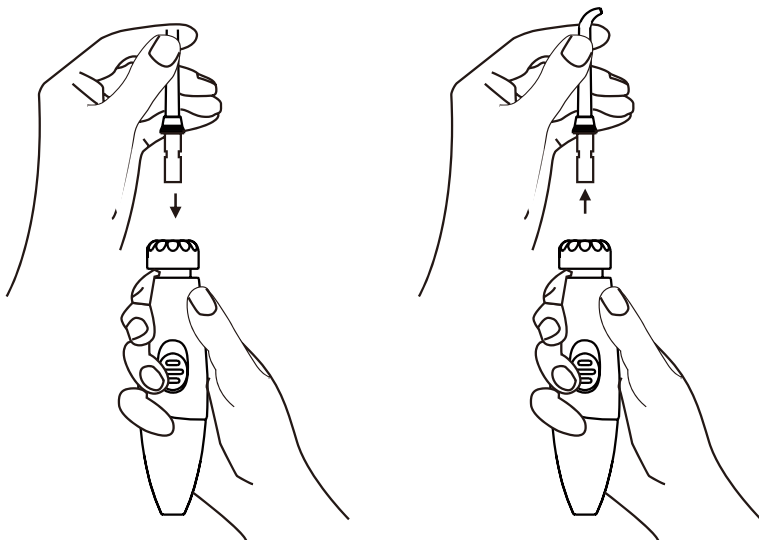
## Adjusting the pressure settings

Turn the pressure control dial on the base of the unit to the lowest setting for first-time use. Gradually increase pressure over time to the setting you prefer, or as instructed by your dental professional.

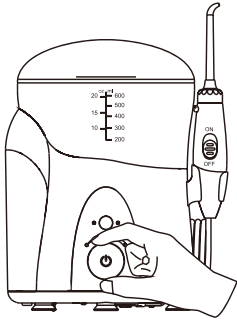
## Using the oral irrigator

Lean low over sink and place the jet tip in mouth. Aim the tip toward teeth. With your free hand, turn the on/off button on. Set pressure at desired setting. Direct the stream by turning the knob at the top of the handle.
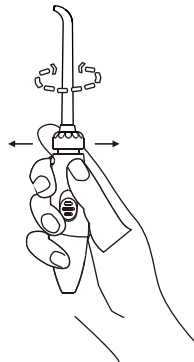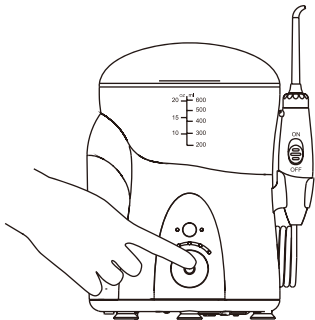
## Recommended technique

Direct the jet stream at 90-degree angle to your gum line, slightly close lips to avoid splashing but allow water to flow freely from mouth into the sink. For best results, start in the molar area (back teeth) working towards the front teeth. Glide jet tip along gum line and pause briefly between teeth. Continue the process until all areas around and between the teeth have been cleaned.

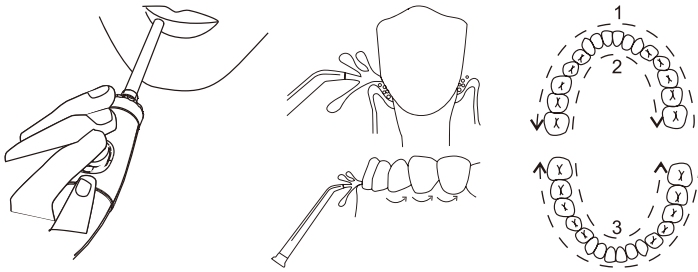
## Pause control

You can temporarily stop the flow at any time by slide down button on the handle.

## When finished

Turn the unit off, pour out the liquid in the water tank and unplug the power plug.

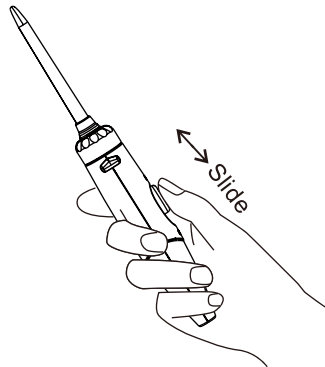
## SPECIAL TIPS INSTRUCTION

### Periodontal Tip


The periodontal tip delivers water or anti-bacterial solutions deep into periodontal pockets.

To use the tooth pocket tip, lean low over your sink and place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gum line. Turn the unit on to the lowest pressure setting and clean gently following gum line, inserting tip in between teeth. Continue tracing along the gum line until you finish cleaning.

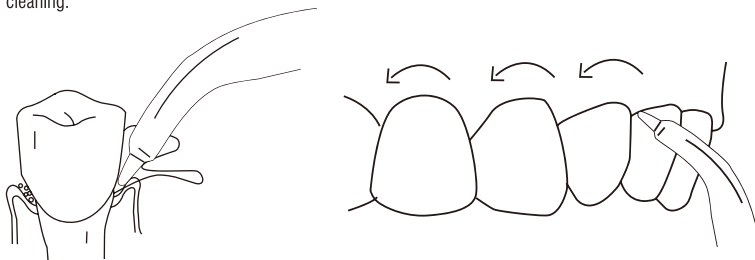
### Tongue Cleaner

Lean over the sink, Place the tongue cleaner on your tongue and start the unit use the lowest pressure setting, increase as comfortable. Gently place the tongue cleaner in the center/midline of your tongue about half way back and with light pressure pull forward.
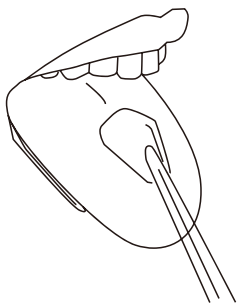
## Orthodontic Tip


The special Orthodontic Tip is uniquely designed to simultaneously brush and rinse hard-to-reach areas around braces and other dental work. To use the Orthodontic Tip, lean low over the sink and use the brush bristles to clean around dental work, water will flow through the brush to rinse away any debris that is removed.
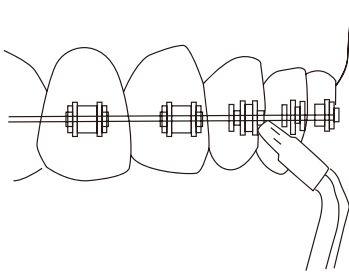
## Dental Plaque Tip

Direct the Jet Stream at 90-degree angle to your gum line, Close lips slightly to avoid splashing but allow water to flow freely from mouth into the sink. Start in the molar (back teeth) area working toward the front teeth, place the tip close to the teeth so that the bristles are gently touching the teeth.

Gently glide tip along gum line, parsing briefly between teeth to gently brush and allow the water to flow between the teeth. Continue until all areas around and between the teeth both inside and outside, top and bottom that have been cleaned.
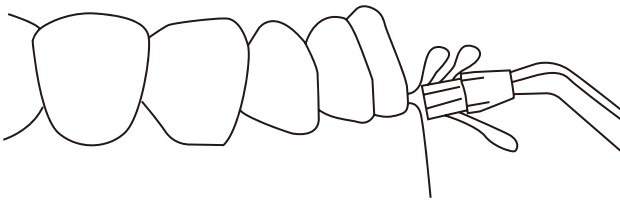
## CLEANING

After using any special solution in the water tank, take a moment after every use to rinse the unit.

Partially fill the water tank with water, aim the jet tip into your sink, turn on the unit and let it run until the water tank is empty, This will help keep the tip from clogging and help your unit to last longer.

Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser.

Before cleaning, please unplug from the electrical outlet.

Before exposing the dental jet to freezing temperatures, remove the water tank and run the unit until completely empty.

## SERVICE MAINTENANCE

1. After using the oral irrigator, be sure to unplug the power cord, and the tank upside down on a flat surface, when the water is completely dry, put it into a cool, dry place together with the cleaned and dry host, it should be placed in the box or similar storage when not in use for a long time.
2. When a machine malfunctions, you can contact the seller to inquiry related maintenance methods.
3. Do not disassemble or attempt to resolve this product. If resulting in any problem, the product will no longer meet the warranty policy.
4. When the power cord is damaged, go to our designated service point or related maintenance shop for maintenance. Do not replace the power cord or plug by yourself.
5. If the product is out of the warranty period, we will provide repairing service as well, we will charge the cost of materials appropriately and small amount of labor cost.

## COMMON PROBLEMS AND SOLUTIONS

Problem	Cause	Solution
Jet tips and tongue cleaner leakage between the top and handle	Tip not fully engaged	Remove tip and reinstall
Periodontal tip tear in soft tip	Tip is damaged	Replace with new tip
Oral irrigator base hose to the handle leaks	The hose is damaged	Contact the seller for service
Water tank leaks	Missing water tank valve or water tank valve is upside down	Reinstall the tank valve
Inadequate pressure	Water tank not positioned properly	Make sure that the water tank is in close contact with the main unit
Without water flow out	Motor not running	Contact the seller for service

We warrant to the original purchaser/owner of this new product that it is free for defects in materials and workmanship for 1 year from the date of purchase.

Thank you again for choosing NICEFEEL oral irrigator!

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Fly Cat Electrical CO.,Ltd