CASADA
HEALTHCARE

time to relax

MINIWELL TWIST

CASADA[®]
HEALTHCARE

since
2000

< timetorelax

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since
2000

< timetorelax

Since being established in 2000, the name CASADA stands for high quality products on the health-spa and fitness market. Today, CASADA-equipment is available in 37 countries world-wide.

Constant further development, an unmistakable design and a focused orientation towards a better quality of life are what characterize CASADA-products, as much as their unbeatable and reliable quality.

Congratulations!

With the purchase of this massage device you have shown awareness of your health.

In order for you to enjoy the advantages of your massage device in the long term, we ask you to read and adhere to the instruction manual carefully.

We wish you a lot of fun with your personal Miniwell Twist.

Table of contents

Safety instructions	10
Inside the package	10
Start up	10
Technical data Miniwell Twist.....	11
Equipment / functions.....	11
Optional	12
Control panel	12
Troubleshooting	12
Storage / care.....	12
Conformity declaration	13
Warranty conditions.....	13

Safety instructions

Please carefully read the following instructions before using this massage device to ensure its proper function and optimal effect. Please retain these instructions for use!

- This massage device complies with the recognised principles of technology and current safety regulations.
- This equipment is maintenance free. Any necessary repairs must be performed by an authorised professional. Improper use and unauthorised repairs are prohibited for safety reasons and will void the warranty.
- Never touch the power plug with wet hands.
- Please avoid the equipment coming into contact with water, high temperatures and direct sunlight.
- Be sure to unplug the equipment from the power point during extended periods of non-use to eliminate the risk of a short circuit.
- Do not use damaged plugs, cords or loose sockets. Damaged plugs or cords must be replaced by the manufacturer, service representative or qualified personnel. In the event of malfunction unplug the equipment immediately. We accept no liability for damages resulting from wrongful or improper use.
- In order to prevent excessive stress on muscles and nerves, the recommended massage time of 30 minutes should not be exceeded.
- Never use pointy or sharp objects.
- Any massage – also massage by hand – has to be avoided during pregnancy and/or in case of one or more complaints in the area of the massage: recent injuries, thrombotic illnesses, inflammations and swelling of any kind as well as cancer. Consultation of your physician is recommended prior to any treatment of afflictions and ailments.
- Do not use the device while driving.
- If you depend on electrical aids such as a pace maker, please consult your physician prior to usage in any case.

Inside the package

- Miniwell Twist
- Adapter
- Car adapter

Start up

Please connect the adapter with the device and subsequently plug the power plug into the power socket. Now turn on the device by the control panel.

With the car adapter the Miniwell Twist can also be used when travelling (not to be used when driving!). Simply plug the car adapter in the cigarette lighter.

Technical data Miniwell Twist


Measurements device:	35 x 23 x 15 cm
Weight:	2 kg
Power source:	
USA	Input: 100-240V~50/60HZ, 600mA / Output: DC 12 V 2000mA
CANADA	Input: 100-240V~50/60HZ, 600mA / Output: DC 12 V 2000mA
AUSTRALIA	Input: 100-240V~50/60HZ, 600mA / Output: DC 12 V 2000mA
EU	Input: 100-240V~50/60HZ, 600mA / Output: DC 12 V 2000mA
Rated power:	24 Watt
Certificates:	    

Equipment / Functions

The Miniwell Twist's integrated deep warmth intensifies the massage. It permeates the massage heads and quickly becomes noticeable. During the warm up period, the circular movements of the massage heads relax tense muscles. The rotating movement of the massage heads loosen up the tissues.

Thanks to its ergonomic form the Miniwell Twist can give an effective massage to nearly all body parts-be it neck, shoulder, back, legs, feet or stomach. Thanks to the optimal placing of the four massage heads with their relaxing circular motion, a soothing comforting massage is guaranteed every time.

Regular massage with the Miniwell Twist stimulates the functions of the nerves and activates the metabolism. Blood circulation is improved, even the most stubborn muscular tension is gently massaged away. In certain cases improved blood circulation and higher metabolism can even prevent blood clots and malfunctions of the internal organs.

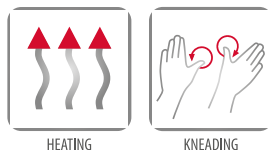Anyone who has had the opportunity to enjoy a Shiatsu-Massage knows how beneficial it is. Regardless of the purpose of the massage, be it as a preventive measure against possible ailments or simply to relax and improve the quality of your life, a good, effective massage is an unforgettable experience.

The central philosophy of the Shiatsu Kneading is that the human being is a combination of body and spirit. This massage technique combines ancient traditional knowledge with modern treatments. The Shiatsu-Massage relaxes and regenerates mind and body giving you a new lease of life; relieving stress and at the same time giving you a feeling of invigorating well-being.

Optional

The additional active bands not only ensure the optimal placement of the Miniwell Twist. They also bring further advantages: You can regulate the pressure on the affected area to your own requirements and you can enjoy your Miniwell twist massage without having to look for something to lean up against.

Control panel

[1]

On / Off

Use this button in order to start the Miniwell Twist. The function Kneading/Shiatsu and heating starts.

In order to turn off Miniwell Twist, press the button again.

[1]

Troubleshooting

Error

Adjustment

The massage device cannot operate even though it has been connected to the power supply.

- Check the connection of the adapter.
- The device could be overheated after prolonged use. It will automatically cool down to resume its normal condition. Please turn off the power and allow the device to cool off for 30 to 50 minutes.

Storage / Care

Regular care will extend the life of your Miniwell Twist whilst also ensuring safe operation at all times.

1. Recommended storage: dry, free of dust.
2. Keep liquids and open fire away from the Miniwell Twist.
3. Clean the equipment using a dry, soft cloth.
4. Always unplug the Miniwell Twist when not in use.

Conformity declaration

Corresponds to the European guidelines EMC 2014/30/EU, LVD 2014/35/EU, RoHs.

Warranty conditions

The dealer grants the statutory warranty period on the product described overleaf. The warranty applies from the date of purchase. The date of purchase is demonstrated by the sales receipt.

The dealer is to repair or replace products registered within Germany which are found to be defective free of charge. Parts subject to wear and tear such as cloth covers are not covered. To assert your claim, simply return the defective product with the sales receipt before the end of the warranty period.

The warranty claim does not apply if it is established that a defect has occurred e.g. through external influence or as a result of repair or alteration not carried out by the manufacturer or an authorised dealer.

The warranty granted by the seller is limited to repair or replacement of the product. Under this warranty, the manufacturer or seller has no further liability and is not responsible for damage caused by disregard of the operating instructions and/or improper use of the product.

If inspection of the product by the seller reveals that the warranty claim relates to a defect not covered by the warranty or that the warranty period has expired, the cost of the inspection and repair are to be borne by the customer.

Casada International GmbH
Obermeiers Feld 3
33104 Paderborn, Germany
info@casada.com
www.casada.com

Copyright Images and texts are subject to the copyright of Casada International GmbH and may not be used for other purposes without express confirmation.

Copyright © 2022 Casada International GmbH. All rights reserved.

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Copyright © 2022 CASADA International GmbH.
Copyright © 2022 Casada International GmbH. Todos los derechos reservados.
Copyright © 2022 Casada International GmbH. Med förbehåll för alla rättigheter.