METHODE
BRIGITTE
KETTNER

Brigitte Kettner

EYE REJUVENATOR

English/Français/Español

Eye Rejuvenator

Thank you for choosing Methode Brigitte Kettner.

Driven by her passion for a homeopathic approach to healthy skin, Ms. Brigitte Kettner, a well-known German biologist and skincare specialist, developed Methode Brigitte Kettner. MBK has been successfully sold throughout Europe for over three decades and continues to gain worldwide success by **offering a full range of targeted and effective skincare products along with supporting hand-held devices** that are designed and developed in Germany.

Please carefully read the following safety precautions and keep them for future reference.


EYE REJUVENATOR

By using microcurrent technology, MBK's Eye Rejuvenator helps minimize the appearance of fine lines and wrinkles around the eye area. It also helps minimize the appearance of under-eye puffiness and sagging.

Contents:


- Eye Rejuvenator Device
- User Manual
- AC Adapter

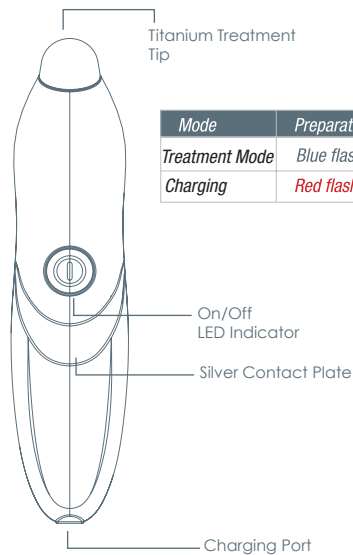Mode	Preparation	Ready
Treatment Mode	Blue flashes	Blue light
Charging	Red flashes	Red light

AC Adapter

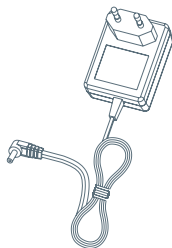IMPORTANT SAFETY PRECAUTIONS:

READ ALL INSTRUCTIONS BEFORE USING

When using electrical appliances, basic safety precautions should be followed, including the following:

- Do not charge this device with any other charger, Only use the charger that was supplied.
- Do not place or store the device near water, such as a sink or tub.
- This device should be kept out of the reach of children.
- Keep the AC adapter and device dry.
- To clean after use, simply wipe with a dry cloth.

DO NOT GET THIS DEVICE WET.

- Do not operate the device while plugged in.
- Never operate this device if it has a damaged cord or plug.
- Do not wrap the cord around the appliance.
- Remove contact lenses before using.
- If you are experiencing serious skin problems, contact your dermatologist.

DO NOT USE THIS APPLIANCE IF:

- You have recently tanned or sun-bathed.
- You have any type of skin irritation, rash, infection, burn, inflammation, active acne breakouts, rosacea, eczema, etc.
- You are recovering from a cosmetic procedure.
- You are using skin-thinning medications and/or topicals such as, but not limited to: Accutane, Tretinoine, and/or steroids

USAGE GUIDE:

1. Charge your Eye Rejuvenator Device for 4-5 hours before using it for the first time.

*-Charge by inserting the AC adapter into the charging socket (refer to diagram on page 4) and plug into a wall socket. The red light will flash to indicate the battery is charging. When fully charged, the red light will stay on. Once charged, unplug the cable; do not keep the charging cable connected to a power source when not charging the device. Once the battery begins to run low again, the red light will flash repeatedly to alert you it is time for a re-charge. *NEVER OPERATE THE DEVICE WHILE IT IS PLUGGED IN.*


2. Turn on by holding down the on/off button for 0.5 seconds. Once the blue indicator stops flashing, the device is ready to use. **Note: The treatment tip will be very warm.**

3. Make sure skin is freshly cleansed and toned, then apply an eye serum or preferred product around the eye area—eyelid, corner of the eye, & under the eye.

4. Hold the device, making sure that your hand is making direct **contact with the silver contact plate (this maintains the microcurrent connection)** and slowly and gently glide the device from the outer corner to inner corner of your eye, pausing for 3-5 seconds on acupoints (refer to diagram). When the device beeps, switch to the other eye. The device will turn off by itself after treatment.

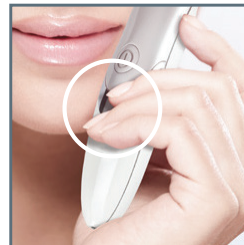5. Follow up by proceeding with your regular facial regimen—such as additional facial serums, moisturizer, and daily SPF.

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