

---

METHODE  
BRIGITTE  
KETTNER

---

## **Introducing GLO!**

Thank you for choosing this MBK facial tool. Please read this instruction manual carefully before using this unit and keep it for future reference.

### **I. Safety Instructions**

1. This device should not be used in the following: children, pregnant women, patients with epilepsy, tumors and cancer patients, patients with cardiac pacemakers, implanted with cardiac defibrillators, and other devices implanted inside the body. If you have any medical questions, such as serious illness, epilepsy, etc., please consult your doctor before use.
2. People who are allergic to light or taking drugs that cause light sensitivity are prohibited to use this device.
3. Do not look directly at the light source, which may cause eye fatigue.
4. Do not use red light for a long time on the thyroid (neck position). Direct stimulation of red light may cause thyroid hyperplasia.
5. Do not use micro-current on the chest, this behavior may cause arrhythmia.
6. To avoid electric leakage, burns, fire, or damage, do not immerse the device in water.
7. Please don't use this device if you have eczema.
8. Please wipe the device with anti-bacterial before using for the first time.
9. If you have sensitive skin, please test an area on your arm first.
10. When the power consumption is only 20 percent left, the red-light indicator will blink twice every 10 seconds to indicate that the device needs to be charged.
11. This product is charged by the USB cable. The red-light blinks during charging. It will go solid when the device is fully charged.
12. The product will automatically shut down when the battery has been drained.

### **Usage Instructions:**

1. Hold on/off button down for 3 seconds to turn on the device. Please do not use the unit while plugged in for charging.
2. Please apply cleanser for deep cleansing mode.
3. Please adjust the vibration to the lowest level when using it for the first time, and then adjust the settings higher if wanted.
4. The device uses microcurrent technology. You may or may not feel the actual microcurrent stimulation. This is normal.

### **II. Device Introduction**

This multifunctional device is used for deep cleansing, product penetration and skin tightening. It integrates positive ion, negative ion, heat, vibration, red light, yellow light and microcurrent (EMS).

Please use cleanser when using cleansing side and a conductive product when using the microcurrent/LED mode.

**Functions:**

**Positive ion:** The positive ion mode helps remove residue and aids in deep cleansing.

**Negative ion:** The negative ion mode helps with product penetration.

**Warm:** The heat helps to open the pores for better product penetration and removal of residue on the skin.


**Vibration:** This mode will help release residue in the pores. It stimulates (awakens) the skin and aids in relaxing and toning.

**Red light:** 630nm-635nm, red light penetrates to the dermis, promoting growth of collagen to maintain the elasticity and water content of skin.

**Yellow light:** 590nm-595nm, yellow light can accelerate the blood circulation, activate cells, promote cell growth, and improve pigmentation.

**EMS:** Electric Muscle Stimulation (microcurrent), the device mimics brainwave signals to stimulate the muscles for a toning effect.

**III. Product Diagram**

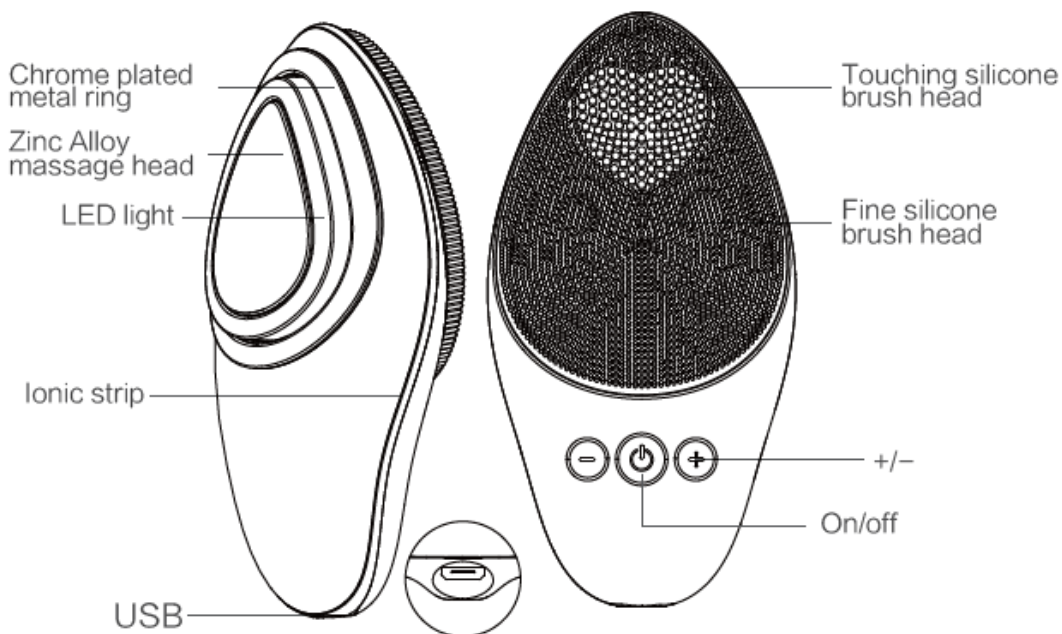

**IV. Operation and Usage**

1. The cleansing vibration has six intensities. Press '+' and '-' key to adjust the vibration. Every press, the red light will flick moving from lowest level and to highest level.

2. Each time you press on/off button it switches the mode. The red-light flashes once.
3. Fingers must be touching the silver ionic strip around the edge to complete the electric circle.

**Cleansing Mode with Silicone Brush Head:**

Please remove base makeup before using the device in cleansing. Apply cleanser to face making sure there is plenty of moisture (water) on the face. The device glides over the face better with moisture. When you turn the device on it will default to level 1 vibration. Chose the suitable vibration using the '+'/'-'key. Cleanse the face for about 2 minutes using circular movements.




**Cleansing Mode with Smooth Metal Side: positive ion, vibration, heat, and red light.**

Take the magnetic ring off the unit, cover with a cotton pad saturated with toner or makeup remover then place ring back on unit. A long press on the on/off key turns the unit on. Press + or – for intensity level. Make sure fingers are touching the ionic strip to close current. Move around the face in circular motions for deep cleansing.




**Infusing Mode with Smooth Metal Side: negative ion, vibration, heat, and yellow light. This mode is for product penetration.**

Apply skin care products on the face. A long press the on/off key turns on the device. Then two short presses on the on/off will put the unit in infusion mode. Press + or – for intensity level. Make sure fingers are touching the ionic strip to close current. Use upward and outward movements.

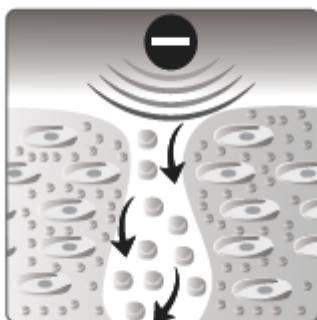

## **Firming Mode with Smooth Metal Side: EMS & red light. This mode is to firm and stimulate/massage the skin.**

Apply skin care products on the face. A long press the on/off key turns on the device. Then three short presses on the on/off will put the unit in firming mode. Press + or – for intensity level. Make sure fingers are touching the ionic strip to close current. Use upward and outward movements. There is no vibration in this mode so you will not hear anything.

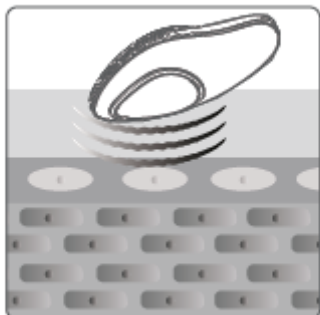

## **V. Troubleshooting**

<b>Breaks</b>	<b>Possible Reason</b>	<b>Solution</b>
Unable to turn on	Not long press on/off key for 3 seconds	Be sure to press the on/off key for 3 seconds
	Under charging	Remove USB, re-start the device
	Low battery	Connect to USB and charge
Unable to charge	USB does not connect to the adapter	Make sure USB cable is fully connected to the adapter
	USB does not insert to the device	Fully insert USB cable to the device
	USB cable broken	Replace another USB cable
Strong tingling feel	Too high level of EMS	Adjust to lower EMS intensity
	Lack of makeup remover or toner	Replenish makeup remover or toner
No vibration	Not touching the ionic stripe	Make sure to touch the ionic stripe while using
	Dry skin	Apply skincare products to face

## **VI. Cleaning and maintenance**

Please clean the brush after using. Please don't apply alcohol, nail remover and the other corrosive solution to wipe the machine. It may cause product failure, cracking parts or discoloration, etc. Please charge this device frequently to maintain the life of the product.

## **VII. Specifications**

Product Name	GLO
--------------	-----

Produce Model	MBKCB22
Vibration Frequency	5500 rpm
Charging Voltage	DC 5V/1A
Charging Time	3 hours
Working Time	1 hour
Lithium Battery	3.7V 500mAh
Vibration Intensity	6 levels
Produce Dimension	120*62*41mm
Net Weight	127g
Package	Brush, USB cable, manual

**This device has a 2-year warranty!**