

Electrical Foot Warmer Heating Pad Instructions

Thank you for choosing the Foot Warmer. Our goal is to provide the best products to enhance your life.

This is an electrical product. Please ensure you read and understand the instructions below before using your Foot Warmer. Failure to do so may result in electrical shock or fire causing serious bodily injury.

WARNING!

RISK OF ELECTRICAL SHOCK

The product is designed specifically for humans.

WARNING: This is an electrical product. Please read all directions before use. Failure to do so may result in electrical shock or fire causing serious bodily injury.

1. Read all instructions before using this product.
2. FOR INDOOR USE ONLY
3. DO NOT use if any visible damage is apparent (cracks, tears or exposed wires and heat element).
4. Must instantly turn off the Foot Warmer and Unplug when overheat or charred flavor in use.
5. Serious heart patients, those that depend on a pace maker and pregnant women should check with their doctor before use.
6. DO NOT use for unconscious and sensitive to heat people.
7. DO NOT make any modifications to switch button ,wires and heat element.
8. DO NOT use in sleeping;Do not put the Foot Warmer under in mattress or quilt when use.
9. DO NOT use cord as handle.
- 10.DO NOT puncture or use any sharp object to secure the Foot Warmer.
- 11.DO NOT use another electronic heating device in conjunction with your Foot Warmer.
- 12.DO NOT tuck the heating pad when use.
- 13.DO NOT Sit and lie on the Foot Warmer heating pad.
- 14.DO NOT Heap and Cover anything on the Foot Warmer when use.
- 15.DO NOT USE for Children.
- 16.Unplug when not in use.
- 17.DO NOT USE WET.
- 18.DO NOT dry clean this Foot Warmer. Dry cleaning solvents may damage the wires.
- 19.DO NOT use bleach or cleaning fluids.
- 20.DO NOT use a wringer.
- 21.DO NOT iron your Foot Warmer.
- 22.DO NOT use mothballs or sprays of any kind.
- 23.DO NOT reconnect the heated product or use it until it is COMPLETELY dry.
- 24.DO NOT use analgesic cream in combination with your Foot Warmer.

Function:

The fast heating technology provides you unmatched warmth and pain relief to cold feet, saving feet from poor circulation, arthritis and cold floors; Besides, it's perfect for treating sore muscles in your back, shoulders, abdomen, legs, and arms; promote blood flow, or keep you warm.

6 Temperature Settings: LED controller adjusts heat easily with 6 temperature settings to accommodate your comfort level (40°C/104°F to 60°C/140°F).

TIMER:4 level settings ranging from 30~120min,conserves energy and prevents burns caused by excessive heating.

6 HEAT SETTINGS:

Warm-temperature range (~40°C/104°F)

Low-temperature range (~45°C/112°F)

Medium Low -temperature range (~48°C/119°F)

Medium-temperature range (~52°C/126°F)

Medium High-temperature range (~56°C/133°F)

High-temperature range (~60°C/140°F)

Features:

ON/OFF Push Switch;IC controller box with LED indication;

6 temperature levels to adjust(40°C/104°F to 60°C/140°F);

4 timer levels to adjust(30/60/90/120min.);

Single circuit provides consistent temperature;

Fuse inside of controller;

Hand and Machine Washable.


Foot Warmer Heating Pad Operation:

- 1.Ensure the plug into a safety-tested electrical outlet.
- 2.Turn On:Push Up the Switch key to ON, the Foot Warmer into Default Mode(1/2/3 heat Indicator LED ON,60min Timer Indicator LED ON).
Default Mode:
Heat level: Medium (3 level)—temperature range (~48°C/119°F)
Timer level:60min to auto shut OFF.
- 3.Adjust Temperature (1~6 levels):
In Default Mode,press the heat key to Adjust Temperature for suitable your requirement,that heat LED indication will rages corresponding Temperature level.
- 4.Timer:
In Default Mode,press the Timer key to Set the heat-up time,that timer LED indication will rages corresponding Timer levels(30/60/90/120min.).
The timer LED will also be shifted down one level in every 30min,When the timer is out,all the LED will be off,and the Foot Warmer will be shut off.In Timer Auto OFF Mode press Heat key or Timer key back to Default Mode.
- 5.Turn Off:Push Down the Switch Key to OFF in any Modes
- 6,Clean: Turn off the Foot Warmer, Unplug Foot Warmer from the wall outlet. Remove the controller and Power unit. Foot Warmer can be machine washed and must be dry before reassembling the unit. Do not use cleaners of chemicals.
- 7,Storage: Turn off the Foot Warmer, dial out the Foot Warmer, Store the pad in a dry ,safe, clean and cool place after clean up.

KEEP SHIPPING PACKAGING for use in the event that your unit needs servicing.(If this unit is shipped without the original packing, damage could occur and void the warranty).

PRE-OPERATING INSTRUCTIONS

When removing the unit from the packing case, check for any shipping damage. Please contact your shipper and dealer immediately if the unit bears any sign of damage from mishandling.

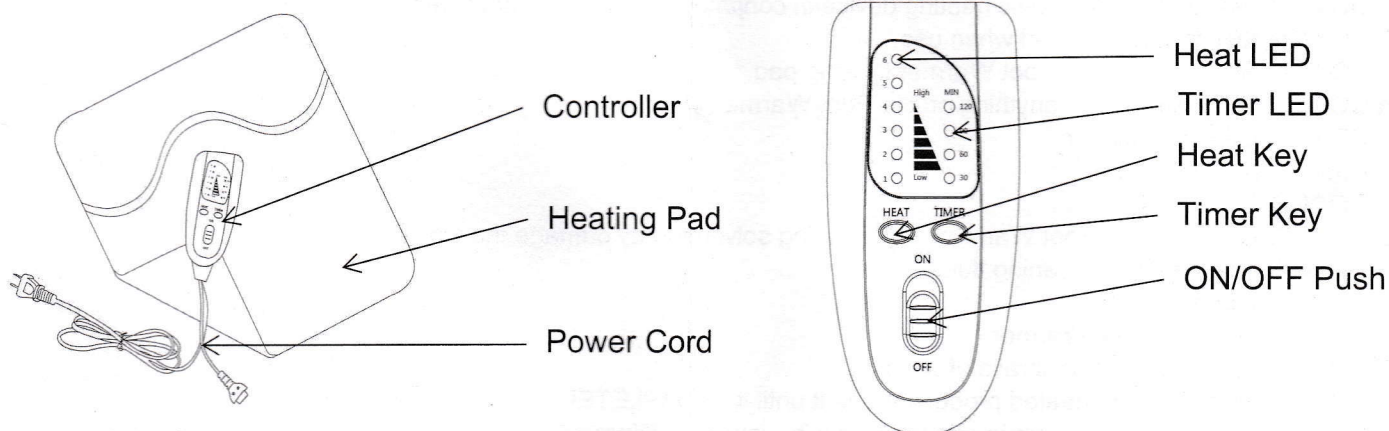

Contents:

One Foot Warmer Heating pad, One controller box with Power Cord,Instructions

Limitations of Liability:

1. Terms of usage is unconditional upon your acceptance, without changes to the terms, conditions or notice listed herein.
2. Proper use including reviewing all safety precautions.
3. No prohibited or unlawful use. The products is intended for humans use only.
4. Our reserves the right to change the terms and conditions under which this product is offered.
5. Buyer assumes all risks and liability from the use of this product. Our is not liable for any damages arising from use or misuse of this product.

Specifications:

Model name:	Foot Warmer Heating pad	
Model No.:	FT-1616	FT-3016
Size:	Small Size 40*40cm(16*16Inch)	Larger Size 76x40cm(30*16Inch)
Power:	50 Watts	120 Watts



MADE IN CHINA