



# USER MANUAL

Inspire Recovery Gun

Model - IRG01

## WARRANTY INFORMATION

### 1-YEAR WARRANTY

If within one year from the date of purchase, this product fails due to a defect in material or workmanship, Pi Fitness and/or the international distributor, dealer or retailer from whom the purchase was made, will repair or replace the product, or necessary components, free of charge.

### THIS WARRANTY EXCLUDES:

- (a ) damage caused by accident, abuse, mishandling, or transport;
- (b) units subjected to unauthorized repair;
- (C) units not used in accordance with INPSPIRE care instructions;
- (d) damage exceeding the cost of the product;
- (e) deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises;
- (f) failure to provide the dated proof of purchase.

Some states and countries do not allow a limitation of damages, so the foregoing limitation may not apply to you. This warranty guarantees specific legal rights, and but other rights may vary from country to country and from state to state. This warranty service is available through the operating policies and procedures at Pi Fitness and all Pi Fitness international distributors. If you are a client, please contact the appropriate Pi Fitness business entity, dealer or retailer from whom you directly purchased the product for warranty and return authorization procedures. This warranty is initiated and executed by Pi Fitness.

## TECHNICAL DETAILS

### AC Charger

Rated input	100-240V~ 50/60Hz, 39W
Rated output	26V/1.5A

### Battery

Rated Voltage	26V/1.5A
Type	Li-ion
Capacity	2000mAh
Working hour	2-3hours
Charging time	90min

### Massager

No-load speed	1200-3200r/min
Power	120W max
Frequency	1 <sup>st</sup> Gear 20Hz, 2 <sup>nd</sup> Gear 30Hz, 3 <sup>rd</sup> Gear 40Hz, 4 <sup>th</sup> Gear 54Hz
Stroke length	12mm
Weight	1.21kg/2.60kg (including charger and box)
Dimension	23*13*5cm

Arm rotator button

Gear indicator

Battery  
indicator

Gear switch

Detachable  
massage head

Ergonomic grip

Main power switch

Charging interface



# CHARGING INTERFACE

## 1. TURN ON/OFF THE MACHINE

- Turn on the main power switch at the bottom of the machine, and the control panel battery indicator.
- Turn off the main power switch, and all functions stop working.

## 2. FUNCTION USE

- When the machine is on, touch the gear switch, and the first gear digit icon lights up, with the machine running at low frequency. Repeat touching the gear switch and the frequency gradually increases.
- Press the gear switch and hold for 3 seconds, Gear lock. Press the gear switch and hold for 3 seconds again, the machine starts.
- When the machine is overloaded, the machine shutdown protection. Please turn off the main power switch and turn on it again.

## 3. ATTACHING THE MESSAGE HEAD

Turn off the machine, apply the attachment tongue side (blue part in Fig.1) into the groove side of the machine (red part in Fig.1). Completion result shows in Fig.2.

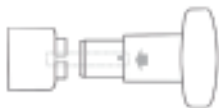


Fig.1 Installation

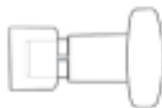


Fig.2 Completion result

## 4. POWER DISPLAY

- The battery displays white light when fully charged, yellow light when less than 30%, red light when less than 20%, and flashes red and motor stops working when less than 10%.
- When the red light flashes, there is a buzzer to indicate charging. A buzzer sounds for 10 seconds a loop, and the battery indicator blinks for 1-2 seconds a loop.

## 5. CHARGING

- When the power connects for charging, the battery indicator shows white light, and the flashing is cycled from strong to weak for 2-3 seconds a loop. When the charging is completed, the battery indicator stops flashing and the white light is fully illuminated.

## 6. ADJUST THE ANGLE OF THE MESSAGE ARM

- Turn off the machine, press the message arm adjustment button, adjust the angle of the message arm to the appropriate position, and release the adjustment button.
- There are three angles for choices: 0°, 45°, 90°. Adjust to the specific position, and the message arm is automatically locked.



Fig.1 Message arm angle

## MAINTENANCE

- The batteries in packs are not designated for exchange or removal by the users, If problems with the power supply should occur, please contact the service center.
- If the device malfunctions, please do not try to repair or to fix it by taking the device apart yourself. This leads to the automatic expiration of warranty.

## WARNINGS AND CAUTION

- Before connecting the device to the power supply, please ensure that the supply voltage stated on the rating plate is compatible with the unit.
- Please charge the new batteries before first use.
- The device should only be used with the provided, certified, INSPIRE charger.
- Do not adjust the speed/frequency or pull out the battery while the INSPIRE is ON.
- Do not put the device in water or in other liquids.
- Do not use the device for more than 30 minutes at a time. Wait 5 minutes prior to using the device again following a 30-min session with the device.
- Do not apply the INSPIRE on your head or near your genitals.
- Do not use the device to supplement or replace any medical treatment. Please ask your doctor before using the device if you have any unexplained pain, if you are receiving medical treatment, or using medical equipment.

## DISCLAIMER

PLEASE CONTACT YOUR DOCTOR BEFORE USING THE INSPIRE IF YOU HAVE ANY PHYSICAL HEALTH CONCERNS OR IF YOU HAVE ANY OF THE FOLLOWING CONDITIONS:

Pregnancy, complications with diabetes, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD's, and metal pins or plates.