

EYE MASSAGER - SRMY01



Please read this manual carefully before use.

Reminder

Thank you for purchasing the eye massager. Please read all the contents of this manual carefully before use and keep it in order to better use this product and warranty service.

The product pictures and related product features shown in this manual and the interpretation of this term are for reference. All matters are subject to the actual product.

Warm + high frequency vibration

constant temperature warm (104 °F ± 5 °F)+ about 160 times /s high-frequency vibration massage for 3 minutes to smooth the eye skin, dilute dark circles, it is a rescue tool for all the night crowd.


EMS micro current

104 ° F ± 5 ° F constant temperature warm+ about 160 times / second high frequency vibration massage + micro current lifting and tightening. Use with skin care products/creams to make the skin more delicate and firm.

Packing list

Product + dust protection cover + USB charging line + manual book

Product description

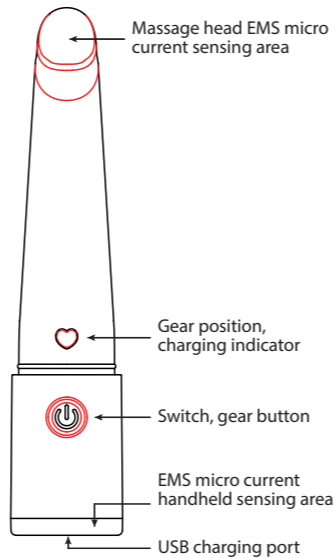

Product use function introduction

- 1.Press and hold the button for about 2 seconds, the product is turned on, the LED starts red lighting, the massage head starts to have heating + vibration.
- 2.Press the button again, the LED light changed to blue color, the massage head start to have heating + vibration + EMS micro current operation.
- 3.The device can be used together with skin care products to relieve eye fatigue and tighten the eye area skin. It is recommended to massage each part for 3 minutes.
- 4.Wipe the massage head with a damp cloth or a wet tissue after use, and then put on the protective case.

Attention

- 1.This product is suitable for the eyes and face areas, CAN NOT be used while charging.
- 2.Use clean hands to apply skin care products to your face before using device. Do not dip the product into skin care products directly. Use a damp cloth or a wet tissue to clean the device.
- 3.After using, cover the protective case, put it in a dry and ventilated place, and can not close to fire or high temperature environment.
- 4.This product can not be used for kids. It is not a toy, keep out of reach of children and babies.

- 5.Don't disassemble and modify this product without permission. The product damage caused by the purchaser's own behavior won't be covered by the warranty service.
- 6.Don't wipe the product with paint thinner, benzine or nail polish, as this may damage the product or cause cracks or fading of the parts.

Product charging

1.Battery reminder:

When the battery is low, the indicator light flashes red.

2.Charging method:

Plug one end of the USB charging cable into the bottom of the product, and the other end into the USB charger or the USB port of the computer. When charging, the indicator light flashes red; when fully charged, the indicator light is on 5 minutes and then out.

3.Charging time:

This product's charging time about 2 hours - 2.5 hours.



Product parameters

Name	Eye massager
Model	SRMY01
Battery capacity	350mAh
Charging time	About 2.5h
Working time	About 60min
Product size	ϕ26.5*L129.5mm
Net weight	150g
Rated input	5VDC/0.5A
Power	≤2W
Vibration frequency	About 180 times/S



<https://facevital.com>



info@facevital.com