
---

# Ultrasonic Skin Scrubber

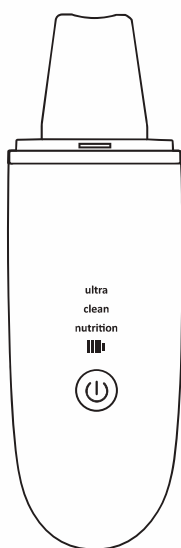User Manual – English

EN

---

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DE

## **ATTENTION**

- 1.** Do not exceed the scope of use described in this manual.
- 2.** Pregnant women should consult their doctor before use.
- 3.** Do not use near eyeballs, inner ear, joints or bone areas. Do not use near metal prosthesis. Do not use this device with or near medical machines such as cardiac pacemakers, oxygenators, electrocardiographs etc.
- 4.** Do not use in areas of vein or arterial pathology. Do not use if you are a patient suffering from Thrombophlebitis or Blood Coagulation conditions (precedents) or related issues. Do not use on skin ulcers, inflammations or broken skin. Do not use on areas with skin allergies or on sunburned skin.
- 5.** Avoid physical impact. Be careful not to drop the device.
- 6.** Only use this device according to the manual. Do not disassemble, transform, retrofit or do maintenance on this device.
- 7.** Please store the device in a clean, dry and stable place. Do not store in excessive heat or cold. Keep away from fire and stoves, ovens and microwave ovens.
- 8.** The charging base is NOT waterproof. Please keep the charging base dry and clean.
- 9.** Unplug the device after charging is finished.
- 10.** Please keep the device away from children.

# ***INTRODUCTION***

Thank you very much for purchasing our products and your confidence in us; we hope our product can bring you a perfect experience. Please read the User Manual and all related information carefully and keep it for future reference.

## ***Product Specification***

Net Weight: 7.2 oz / 205 grams

Measurement: 6.9 x 0.8 x 2 in / 175MM \* 20MM \* 52MM

Working Environment: 32-95 Fahrenheit / 0-35°C

Input: AC 100-240V, 50/60Hz; Output: DC-5V, 1A  
(Rechargeable)

Power: 2.2W

Battery capacity: 600mAh

Ultrasonic frequency: 24kHz

Blue Light Wavelength: 465±3nm

Waterproof Level: IPX5 (Device only, charging base is not waterproof)

## ***Accessories***

1. Main Device
2. Spin USB charging cable
3. Charging Base
4. User Manual
5. Gift Box

## Features

**1.** Standby Mode: Press the button once to turn on the device and enter STANDBY MODE. You can see the lights for “ultra”, “clean” and “nutrition” flash.

**2.** Continuous Wave Deep Cleansing Mode: Press the button twice to enter this mode. You can hear continuous sonic vibrations and see the “ultra” indication light is always on.

**3.** Lift Up Mode (Pulsed Wave Mode): Press the button a third time to enter this mode. You can hear pulsed sonic vibrations and see the “ultra” indication light flash.


**4.** High Level Deep Cleansing Mode: Press the button a fourth time to enter this mode. You can hear continuous sonic vibrations and see the “clean” indication light is always on.

**5.** Sonophoresis Mode: Press the button a fifth time to enter this mode. You can hear continuous sonic vibrations and feel ions. The “nutrition” light is always on.

**6.** Standby Mode again: Press the button a sixth time to enter the Standby Mode again. No functions now. Press again and functions circle.

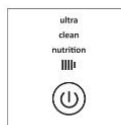**7.** Turn off the device. Long press the button for 1s to power off the device.

**8.** The device is set to work 10 mins and then turn off automatically.

## How to Charge

**1.** When battery is low, the RED indication light lights up and device turns off automatically.

**2.** Put the device on the charging base and make sure it is well connected. Then charge it with the USB cable. During charging the RED light flashes. When device is fully charged, the light turns BLUE and is always on.

During charging,  
Red light flashes.

When charging is  
finished, Blue light  
always on.

## Function Induction

**1.** Ultrasonic deep cleansing technology – The Ultrasonic Skin Scrubber produces 24,000Hz per second of ultrasonic waves that effectively loosen clogged pores and exfoliate the skin to remove dead skin cells painlessly.

**2.** Firms and lifts skin – Micro-current pulse technology and intermittent ultrasonic help boost cell collagen regeneration and activates cell vitality, ensuring your skin remains moisturized and healthy.

**3.** Advanced ion technology – The positive ion helps removes dirt and deeply cleanses the skin. It accelerates the speed of skin cell metabolism, softens old keratin, restores skin elasticity and helps keep skin feeling tight and smooth. The negative ion nutrition leading-in function helps break down and refine larger molecules such as collagen to massage and support your skin to easily absorb nutrients.

**4.** Bio Blue Light – wavelength  $465\pm 3\text{nm}$ : diminishes inflammation, kills bacterium and propionibacterium acnes, removes acne and scars etc.

## ***Deep Cleansing - Overview***



- 1.** Apply skin cleansing serum/oil to your face, keep your face always moist.
- 2.** Turn on the device. Choose the first “Continuous Wave Deep Cleansing Mode” or third “High Level Deep Cleansing Mode”. Hold the unit with the control panel away from the face. Cleanse the face according to the diagram on the bottom of this page. Treatment time is typically 3-5 mins per session – no more than 5 times per week. We recommend starting with the first “Continuous Wave Deep Cleansing Mode”.
- 3.** Wash your face with water after treatment.

## ***Deep Cleansing - Device Usage***

We recommend using the device according to this diagram, as it encourages the body's natural lymphatic drainage.

- 1.** Glide the device from the neck upwards on the face in horizontal strokes. Repeat on both sides of the face.
- 2.** Methodically work over all areas of the face as shown in the illustration.
- 3.** Change the direction of the device and glide it in vertical strokes starting from the neck area. Work from the jaw to the hairline.
- 4.** Cleanse around the curve of the nostrils using the curved inner edge of the blade.
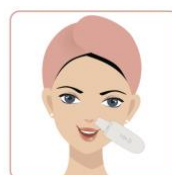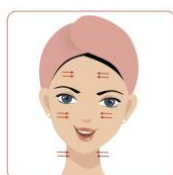
## ***Tips for Effective Cleansing***

Use the scooped out middle section of the blade to cleanse the chin area under the lip. Use vertical strokes.


Select the first “Ultrasonic Deep Cleansing Mode” to work around the eyebrow area. Avoid the delicate under eye area when using the third “High Level Deep Cleansing Mode”.

## ***Lift Up & Sonophoresis***

- 1.** Apply your favorite serum or moisturizer to your face. Hold the device so that the blade is directly placed (flat) on against the skin.
- 2.** Turn on the device and select the second “Lift Up Mode” or the fourth “Sonophoresis Mode”. We recommend starting with the “Lift Up Mode” first.
- 3.** Using gentle pressure, glide the device over the surface of the skin using the flat side of the blade. For best results, work slowly moving the blade backwards at all times.
- 4.** 3-5 minutes is enough for each treatment. Use 3-5 times per week.
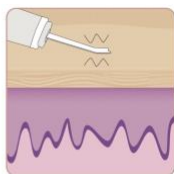
## ***Tips for Lift Up & Sonophoresis***

Begin on the neck area. With gentle pressure, glide the blade over the surface of the skin, using backward movement, working from the center of the face outwards.

Continue upwards from the jawline to the hairline. Glide over the face until all product is absorbed.

## ***Maintenance***

- 1.** After each use, cleanse the device with a wet cloth and then store it in a cool and dry place.
- 2.** Before cleansing the device, make sure it is turned off, and do not connect it to the power supply.
- 3.** Device is for individual, personal use only. Do not share the device with others.
- 4.** Keep the User Manual for future reference.


## ***Failure***

If the device does not work properly, please check the following:

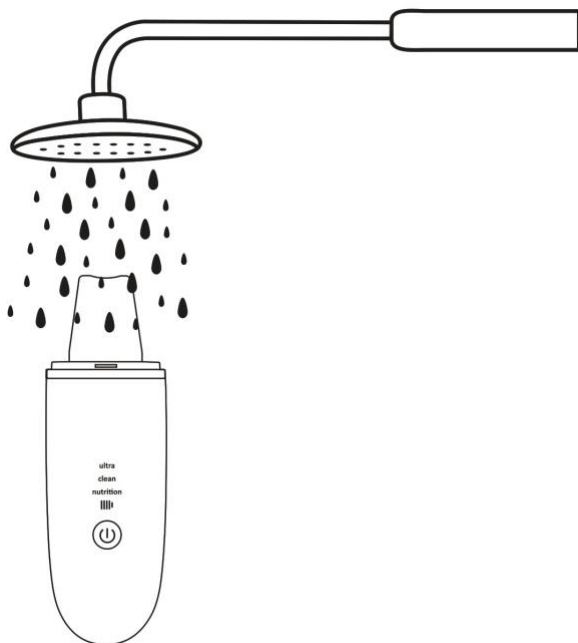- 1.** Make sure the device is fully charged.
- 2.** Make sure to press the “Power” button and that the working indication lights are on – lit up.
- 3.** If there is no reaction after device is turned on, the battery may be low. Please re-charge the device and try again.
- 4.** If there is no action for 2 mins, the device will turn off automatically to protect circuits.

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