

full body

# massage mat with heat

READ AND SAVE THESE INSTRUCTIONS



## SPECIFICATIONS

Part #: 702378

Name: Full Length Body Massage Mat with Heat

Input Voltage: 110 – 240 V 50Hz/ 60Hz

Output Voltage: 12V 1.2A

Rated Power: 14.4W





## IMPORTANT

As with all electric devices, basic care and attention need to be followed:

1. Always unplug the power cord from the power source immediately after use and before attempting to clean the unit.
2. DO NOT operate this appliance if it has a damaged power cord or electric plug; if it is not working properly; or if it has been dropped or exposed to water or steam.
3. Never attempt to repair or disassemble this device.
4. Keep power cord and appliance away from heated surfaces. Excessive heating can occur and cause fire, electric shock or injury to person.
5. Never leave an appliance unattended when connected to power source. Always unplug the power cord from power outlet when it is not in use .
6. NEVER drop or insert objects into any opening.
7. This device is for home use only. Do not use outdoors.
8. DO NOT immerse this device in water or any other liquid.
9. To disconnect – remove cord for the power source.
10. Only use for its intended use.

## FEATURES

- 9 Vibrating Motors
- 7 Vibration Modes
- ON/Off Heat function
- Targeted massage for:
  - Neck & Shoulders
  - Lumbar (upper/lower back)
  - Leg
  - Full Body Massage

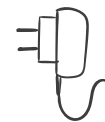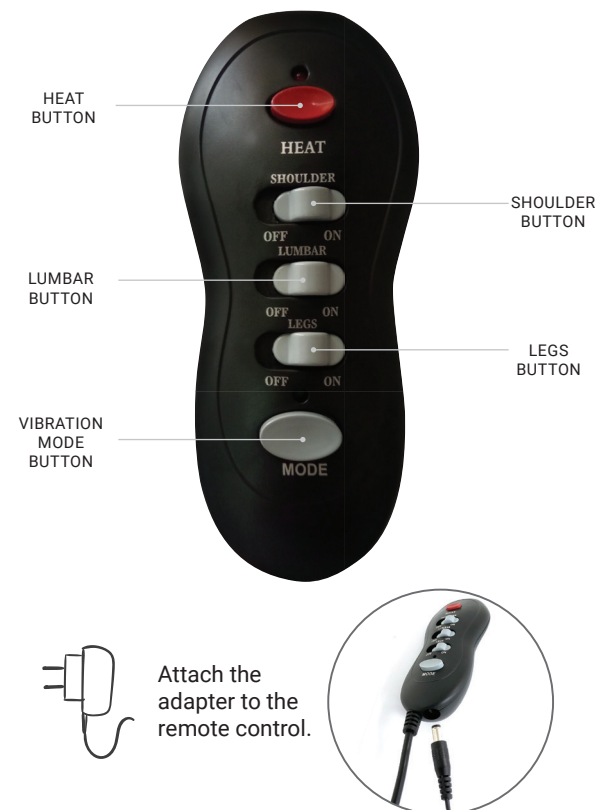

Attach the adapter to the remote control.

## OPERATION

Place massage mat on chair, bed or floor and plug unit into power source. Don't forget to attach the adapter to the remote control.

### Targeted Massage

- To target the neck and shoulders, slide the SHOULDER button to the ON position. To stop shoulder/neck massage, slide button to the OFF position.
- To target the lumbar region, slide the LUMBAR button to the ON position. To stop, slide button to the OFF position.
- To target the legs, slide the LEGS button to the ON position. To stop leg massage, slide button to the OFF position.

### Full Body Massage

- For a full body massage slide the SHOULDER, LUMBAR AND LEG buttons to the ON position. To stop, slide all buttons to the OFF position.

### Heat

- To turn on heat function, press the HEAT button once. To turn off, press again.

## VIBRATION MODES

There are 7 different vibration modes that cycle through automatically when the unit is operating.

Press the MODE button to skip to the next vibration style.

## CLEANING AND MAINTENANCE

1. Unplug the unit.
2. Use a clean, slightly damp cloth and wipe down the surface. (Do not submerge in water.)
3. Wipe clean and allow to dry. Do not fold the unit if damp.

## STORAGE

- Place massager in a safe, dry and cool place.
- Avoid contact with sharp edges or pointed objects which might cut or puncture the surface.
- To avoid breakage, do not wrap the power cord around the unit.

## CAUTION

- Keep device away from children and/or pets.
- Do not use on open sores or broken bones.
- Those with sensory deficiencies should be supervised during use.
- Consult a doctor prior to using if wearing a pacemaker or portable EKG monitor.

## WARRANTY

Relaxus Products Ltd. warrants this product to be free from manufacturer's defect for one year from the date of purchase. Original purchase receipt is required for any warranty claim.

Warranty excludes:

1. Damage caused by improper or incorrect installation or use.
2. Breakage or damage from mishandling or accidents.
3. Use of product for purposes for which it was not intended, or commercial use.

Relaxus Products Ltd. retains discretion to repair, replace parts or supply a new product of the same or similar model. Warranty is not transferable. Any associated shipping costs are not included in warranty coverage.



Relaxus Products Ltd.  
1590 Powell St., Vancouver, BC  
V5L 1H3, Canada  
[relaxus.com](http://relaxus.com)