

breo

breo® iDream 3

EYE & HEAD MASSAGER



User Manual

Thank you for purchasing iDream3.
Please read this manual carefully before using the product,
and keep it properly for future reference.

1. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
2. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
3. The battery must be removed from the appliance before it is scrapped;
4. The appliance must be disconnected from the supply mains when removing the battery;
5. The battery is to be disposed of safely.
6. WARNING: For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.
7. The appliance is only to be used with the power supply unit provided with the appliance.
8. This appliance contains batteries that are only replaceable by skilled persons.

changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

Contents

Introduction	2
Inside the Package	3
Function Deployment	4
Preparation	5
Operation	6
Precautions	8
Troubleshooting & Technical Parameters	9

Breo WEEE-Reg-Nr.DE51664358

Introduction, scope and groups of application

Introduction

IDream3 in virtue of humanized design and adjustable knobs, the product is applicable for heads in various circumferences. The massager performs kneading, vibrating and hot compress massage against special acupuncture points on head and eyes to assist in supplying more oxygen to brain and eyes, improve blood circulation of brain and eyes, and relieve pressure and fatigue of brain and eyes for promoting sleep and maintaining health. While performing massage, the product enables you to enjoy music for relaxing both your brain and soul.

Scope of application

Kneading and vibration massaging the specific acupuncture points on head and eye, so as to alleviate strained feeling and stress, release fatigue, channel massage improves sleeping, which can make you feel fresh, calm and relax completely.

Inside the Package

Your massager is composed of the main unit, remote, power adaptor, user manual, earphone, battery.

Massager



Main Unit



Remote

Basic components



Adaptor



Earphone




User Manual



Battery (AAAx3)

Function Deployment

Main unit and Remote

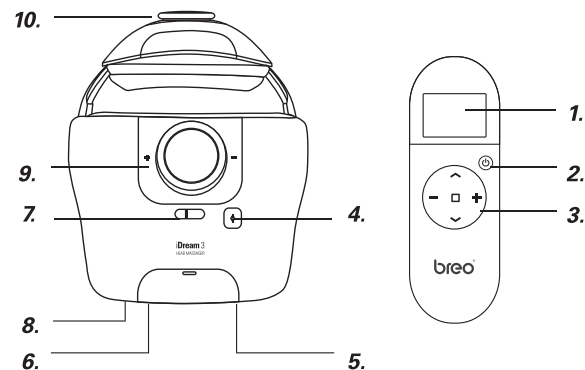

- | | |
|-------------------------|---------------------------------|
| 1. LCD Display | 6. Earphone Jack |
| 2. ON/Off Key | 7. Power Switch Key |
| 3. Operating Buttons | 8. Encoding Joint Operation Key |
| 4. Air Discharge Button | 9. Size Adjusting Knob |
| 5. DC Power Interface | 10. Tightness Adjusting Knob |

LCD Display:

From top to bottom: time/air pressure/hot compress/music

Operating Buttons:

Clockwise: Up key/Mode+/Down key/Mode-, OK key in the center



Power Switch Key:

It's a must to push the Power Switch Key to be OFF after using.

Function Deployment

Attention :

If there was problem of receiving signal between the main unit and remote, you can perform encoding joint operation:

1. Turn off the main unit and remote, press and hold **【8】** encoding operation key of the main unit, then turn on the main unit until the LED blue indicator that on front of the main unit flashing.
2. Release **【8】** encoding operation key of the main unit, synchronously press and hold  key of the remote and key  for 2-3 seconds, the LCD display of the remote flashes, in the meantime, the blue LED indicator of the main unit turns off, which means it successes for the code verification. (The operation needs to be finished within 6 seconds)
3. Turn off the main unit. When turn on the product again, the remote works.

Please make sure that one main unit is corresponding to one remote. The operation of encoding joint with multiply remote controller to one main unit is prohibited.

Preparation Before Use

1. Please check the machine and accessories whether complete.
2. Please check and make sure there is enough power.
3. Please check and make sure there is no damage of the adaptor.


Operation

Charging The Unit

1. Open **【5】** DC Power Interface, and insert the adaptor.
2. During charging, the power indicator on the front of the main unit indicate red and keep flashing, when it fully charged, it indicates blue.

Attention:

1. Keep the massager in off state while charging the main unit. Charging time is roughly 3-4 hours.
2. The main unit doesn't work while charging.
3. If voltage of the battery is lower than 3.2 ± 0.1 V, the product will automatically power off. If LED indicator of the main unit flashes, and all functions shut down, please recharge the main unit.

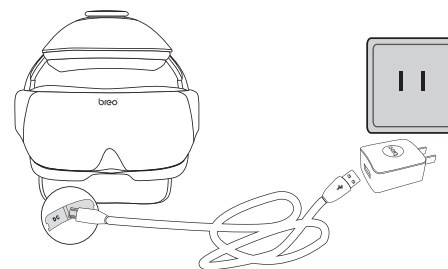

Power on

1. Remove the cover of battery pack on back of the remote controller, place 3 pieces of AAA batteries with correct pole, and then replace the battery pack cover; (Picture1)
2. The main unit of the product is powered by high performance environmental friendly lithium battery.

Operation

How to wear

1. Adjust the top and two sides to a proper state, and then to wear on the main unit.
2. Meanwhile adjust the top to a proper distance to head.

Operation

Power On/Off


Power-on — Push the power protection key on the main unit to Level "ON"

- Press and hold  Button
- Preset air pressure + hot compress for 15 minutes' massage.

(Default air pressure: Circulation; Default hot compress: High temperature; Music: Bird music)

Power-off — Short press  Button — Push **【7】** Power Switch Key on the main unit to Level "OFF" (picture 2)

If switchover of time, air pressure, music and hot compress mode is needed, please have it done by pressing the related functional keys manually.

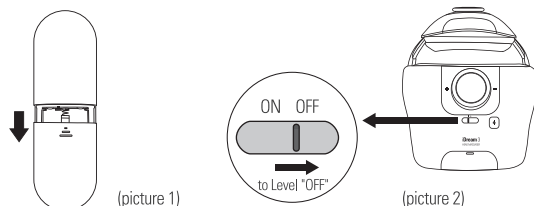




Manual Mode

Short press  Button, choose function mode





- Short press  Button
- Short press  Button to choose

Operation





Time

- Short press  Button to select "Time", short press  Button
- Short press  Button for time switchover among 15minutes/10 minutes/5 minutes





Air pressure

- Short press  Button to select "Air pressure", short press  Button
- Short press  Button to select five modes of air pressure (The default is circulation)
- Off/On, short press  Button

Heat

- Short press  Button to select "Hot press", short press  Button
- Short press  Button: Low temperature/High temperature
- Off/On, short press  Button

Music

- Short press  Button to select "Music", short press  Button
- Short press  Button to select Next
- Short press  Button to select previous

Precautions

1. Do not use the product during driving.
2. Prior to use, please remove any jewelry or object that may affect massage from the massage position.
3. Do not use the product in bathroom or other damp environments.
4. Keep the product away from water and fire sources and corrosive environments.
5. Please observe safety regulations when using the product in specific environments.
6. Do not use the product if the massage position has undergone operation.
7. Do not use the product if the massage position has skin injury, redness and inflammation, congestion or other discomfort symptoms.
8. Persons with heart disease or history of heart disease, or with acute disease, skin disease, malignant tumor or physical abnormality should not use the product.
9. Persons with bleeding tendency or other blood diseases, or with osteoporosis, soft-tissue disorders or apoplexia should not use the product.
10. Persons with slow action, psychosomatic disorders or disturbance of intelligence should not use the product; please follow the doctor's advice if any of these persons intends to use the product.
11. Stop using the product in case of any discomfort during use.
12. Women who are pregnant or in menstrual period should not use the massager.
13. Children under six years of age are forbidden to use the product; children between 6-12 years of age must use the product under the guidance and care of their parents or other persons with supervision capability.
14. Use alcohol cotton to clean the external appearance; avoid using water for cleaning.
15. Do not use any sharp objects to scratch or stab the product.
16. The product is an electronic wellness product for home use rather than a medical device. It is not intended for diagnosis, therapy or other commercial purposes.
17. Use the charger supplied by our company to charge the product, and pay attention to electrical safety.
18. Cut off power supply immediately in case of any abnormality during use. Do not attempt to disassemble or repair the product.
19. The product has a heating surface; those who are insensitive to heat must take care.

Operation steps for battery removal

1. Cut off the external power supply, and remove the product housing;
2. Disconnect the power supply to the battery;
3. Remove the battery;
4. Place the waste battery at a fixed recycling station;
5. Please observe local environmental protection laws and waste disposal guidelines and protect the environment.

Notes

Remove the battery out of the appliance before scrapping the appliance;

The appliance must be disconnected from power supply before removing the battery;

The battery should be disposed of safely.



Except professionals, no one shall attempt to remove the battery!

Troubleshooting

If there is any problem when you use the device, please follow the instructions below to check. Do not disassemble or repair the device yourself.

Fault	Possible Reasons	Solutions
Unable to boot.	If power switch key is turned on.	Push this key to Level ON.
	If the lithium-battery of the main unit is low.	If yes, please charge the main unit at once.
	If the battery in remote is low.	Place new battery into the remote.
The air pressure or vibration is over low.	If the lithium-battery of the main unit is low.	If yes, please charge the main unit at once.

Technical parameters

Product name	Head Massager
Product model	iDream 3
Rated voltage	5V \equiv
Product weight	1125g
Main unit size	L282 × W206 × H254(mm)
Remote size	L132 × W44 × H24(mm)

Shenzhen Breo Technology Co., LTD

19/F, Hisense Southern Building, No.1777 Chuangye Road, Yuehai Street,
Nanshan District, Shenzhen 518054, China
support@breocare.com

CE RoHS FCC

Storage

1. If the storage period exceeds 6 months, please check and recharge it.
2. Do not store close to floor, wall or roof.
3. Keep in good ventilation condition and away from direct sunshine and corrosive gases.