




## CONTENTS

Precuations	3
Caution	4
Necessary knowledge	5
Product Configuration	6
Features	7
Product Operation	8
How to use	9
Opening Guide	10
Instructions	11
Charging	12
Charging instructions	13
Troubleshooting	14
Symbol description	16
FCC Statement	17
Product parameters	18

## PRECAUTIONS

1. Please read the instructions carefully, do not use this product on the same part continuously for a long time, to prevent skin damage and muscle pain caused by excessive massage.
2. Do not use this product during precision operation or driving.
3. Please avoid areas that has undergone surgery, or has skin damage, redness, inflammation, congestion, etc..
4. People with heart disease or medical history, acute disease patients, skin disease patients, malignant tumor patients , people with blood diseases such as bleeding tendency, people with osteoporosis, soft tissue diseases, stroke and other body abnormalities should not use.
5. Do not use it if you have slow movement, physical and mental discomfort or intellectual disability.Please follow the doctor's advice.
6. It is not recommended to use the massager during pregnancy and menstrual period. Please follow the doctor's advice if you need.
7. When storing the product, please put it in a place where children are not easy to reach.
8. Children under the age of 6 should not use it, and children aged 6-12 must be under the supervision of their parents or use under the guidance and supervision of competent persons.

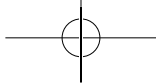

## CAUTION

1. This product is a household electronic health care product and cannot replace medical equipment. It is not suitable for diagnostic, therapeutic and other commercial uses.
2. Before use, please remove the jewelry on the massaged area or other items that may affect the massage.
3. Please stop using if there is pain or discomfort in the massage area.
4. If the product is abnormal, please shut down immediately, cut off the power supply and contact customer service, please do not disassemble and repair by yourself.
5. When using in other specific environments, please comply with safety regulations.
6. Please keep this product away from fire sources and corrosive environments.
7. It must only be supplied at safety extra low voltage corresponding to the marking on the appliance.

## NECESSARY KNOWLEDGE

### About frequency of use

1. Use as needed, within 10 minutes at a time, it is recommended to use it at an appropriate interval reuse.
2. If overused, it may cause physical discomfort and skin injury.
3. Therefore, when you use it for the first time, you need to observe your physical condition and skin changes, and try to shorten usage time.


### Caution before use

1. Use as needed, within 10 minutes at a time, it is recommended to use it at an appropriate interval reuse.
2. If overused, it may cause physical discomfort and skin injury.
3. Therefore, when you use it for the first time, you need to observe your physical condition and skin changes, and try to shorten usage time.

### Wash about the state after use

1. The massaged area may become red and hot for a short time, which is a normal phenomenon. After a while will recover naturally.
2. When the hair is accidentally involved in the product, please do not pull it hard, just turn off the product and take it off.

**PRODUCT CONFIGURATION**

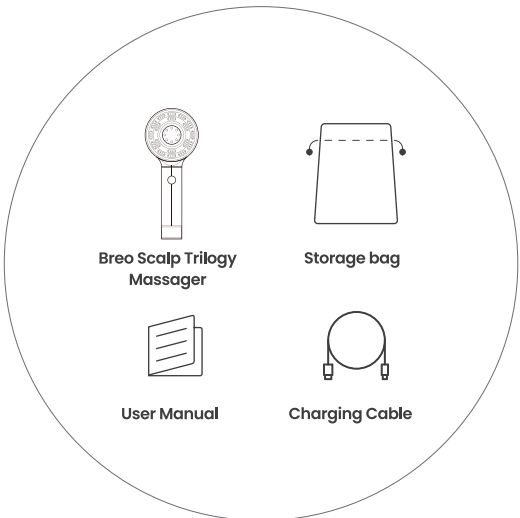

**Breo Scalp Trilogy  
Massager**


**Storage bag**

**User Manual**

**Charging Cable**

6

**FEATURES**

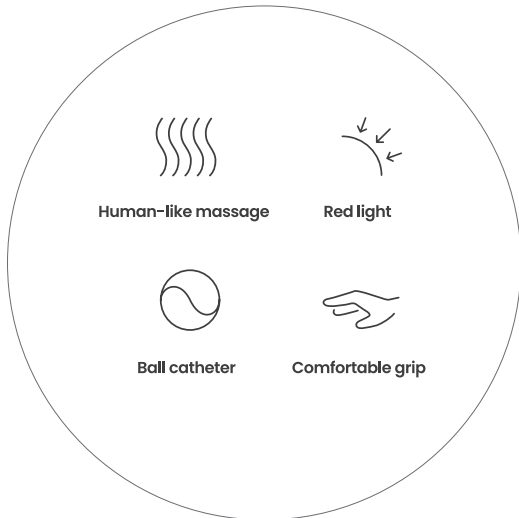

**Human-like massage**

**Red light**


**Ball catheter**

**Comfortable grip**

7



## PRODUCT OPERATION

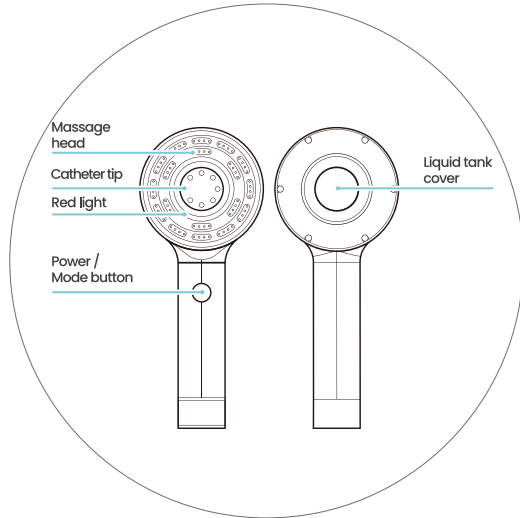

8

## HOW TO USE

This product cannot be used while charging. After charging, please unplug the power adapter before use.

### Power on/Mode Selection

Press and hold the power/mode button for about 2 seconds, and the prompt sound is long "di", which means the product is turned on; After turning on, the default mode is mode one.

After turning on the machine, short press the power button to switch the mode, and the prompt sound short "di" means switching succeeded.

Mode 1: Fast Massage + red light

Mode 2: Slow massage + red light

Mode 3: Alternate fast and slow massage + red light


### Liquid Addition

1. Open the lid of the liquid tank and add suitable products according to your own consumption. The recommended amount is NO more 3ml, so as not to waste.
2. It is recommended to clean the liquid storage tank after each use to keep the product clean.
3. Nursing solutions of different brands and concentrations may cause the uniformity of the ball discharge speed if the concentration is too low, there may be a small amount of liquid leakage, please clean up in time.
4. Special attention: Do not add unknown substances or corrosive liquids to avoid product damage.

9



## OPENING GUIDE

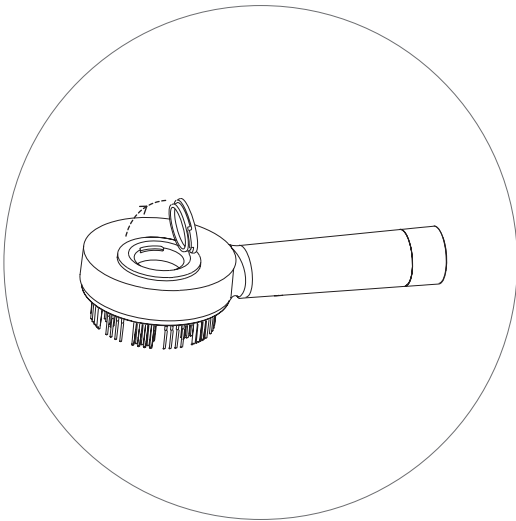

10

## INSTRUCTIONS

### Timing function

Smart timer 10 minutes.

### Shutdown

Long press the switch button, the product shuts down.

### Product Cleaning


1. If you need to clean the appearance, you can wash it with water directly. (polish surface such as plastic)
2. It is recommended to rinse the liquid storage tank and the liquid outlet hole with clean water after each use to avoid clogging. Wipe the product clean and place it in a ventilated and dry place.

### Other

1. The product may produce sounds such as "Yu" and "Crack" during normal use, which belong to sounds made by the product during normal operation.
2. When the new product is unpacked, there may be a little liquid in the liquid storage tank, which is for factory quality inspection.

11

## CHARGING

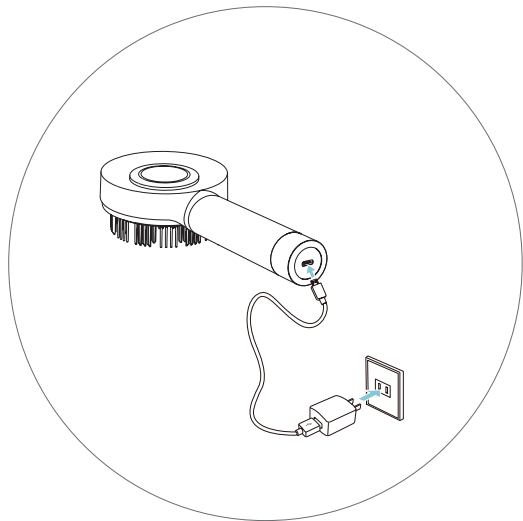

12

## CHARGING INSTRUCTIONS

1. Connect the Type C charging cable to the charging interface of the host.
2. Insert the other end of the charging cable into the rechargeable USB port or connect the 5V adapter, plug into the socket to charge the host.
3. When the product is low in power and charging, the red indicator light flashes, and the green indicator light stays on after charging is completed.
4. The charging time is about 2.5 hours.

\* Please keep the charging port dry before charging.

It is normal for the massager and the power adapter to heat up slightly during charging.

After Charging

Unplug the power adapter and unplug the charging cable from the massager.

13



## TROUBLESHOOTING

If there is a failure during operation, please check the following. If you cannot troubleshoot or the product still cannot operate satisfactorily, please cut off the power supply, and it should be repaired by professional technicians. Please do not disassemble this product by yourself!

Basic failure	Possible cause	Basic Solution
Can not power on	Is the battery low	If yes, please charge in time
Weak massage	Is the battery low	If yes, please charge in time

If you still cannot solve the problem after checking, please contact the customer service:

Email: [support@breocare.com](mailto:support@breocare.com)

Web: [us.breo.com](http://us.breo.com)

Tel: 1-866-546-3587

Working time :

PST 6:00am--11:00pm (Monday to Friday)

PST 7:00am--3:00pm ( Saturday to Sunday)

Made in China

## Notice:

- Before disposing of this product, the battery must be removed as required by regulations;
- When removing the battery, the product must be powered off; Batteries should be disposed of safely in accordance with laws and regulations.

**-Use the appliance at a temperature between 32°F and 104°F(0°C~40°C).**











**-Batteries should be charged in a location where the temperature is between 50°F and 104°F(10~40°C).**

**Non-professionals do not disassemble by themselves!**

Steps to remove the battery:

1. Cut off the external power supply and remove the product shell;
2. Unplug the power on the battery;
3. Disassemble the battery;
4. Please place waste batteries at a fixed recycling point;
5. Please abide by the local environmental protection laws and disposal guidelines, and care for the environment.

## SYMBOL DESCRIPTION

	Keep Dry		Applied part of type BF
	Refer to instruction manual		Serial number
	Manufacturer information: The manufacturer name and address		Disposal in accordance with Directive 2002/96/EC (WEEE)
	Medical device		The device complies with part 15 of the FCC rules
	This device has not been tested for use in an MR environment and should not be used exposed to MR environments while patients are wearing the device. Keep it outside the MRI scanner room.		
	The device complies with compliance with RoHS Directive 2011/65/EU and its amendment directives 2015/863/EU		

## FCC STATEMENT

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

### Disposal

Notice:

At the end of the product life cycle, do not throw this product into the normal household garbage, but bring it to a collection point for the recycling of electronic equipment and battery. Waste Electrical and Electronic Equipment may have potentially harmful effects on the environment. Improper disposal may lead to the accumulation of harmful toxins in air, water and soil, which is harmful to human health. you have obligation to dispose of device correctly. Consult your municipal authority for information about disposal.





## PRODUCT PARAMETERS

Product Name	Breo Scalp Trilogy Massager
Model	Scalp 3
Charging Input	5V  1A
Rated Power	5W
Power Supply	3.8V  1600mAh rechargeable lithium battery
Net Weight	250g
Dimension	L207*W85*H45mm
<b>Made in China</b>	



Shenzhen Breo Technology Co.,Ltd  
19/F, Hisense Southern Building, No.1777 Chuangye Road,  
Nanshan District, Shenzhen, Guangdong, CN, 518056

Version: 2.0 @202307



**breo** Scalp Trilogy Massager  
Scalp 3

SIMPLY RELAX