





 1300 054 554

 [masseusemassage.com.au](http://masseusemassage.com.au)

 [info@masseusemassage.com.au](mailto:info@masseusemassage.com.au)



## MASSEUSE HEALTH PRODUCT MANUAL

Thank you for purchasing your Masseuse Health Massage Chair.

You have purchased a chair from Australia's leading massage chair specialist!

This chair was designed to provide excellent relief from back pain, muscle fatigue and all manners of stress.

Before operating your amazing new chair, please read this instruction manual carefully so you can get the most out of your massage experience.

Please pay attention to the important safety warnings and keep this manual for future use.

## **CONTENTS**

---

<b>Safety and Maintenance</b>	02
<b>Labelled Massage Chair Diagram</b>	04
<b>Overview of Functions</b>	05
<b>Method of Usage</b>	06
<b>Product Specifications</b>	12

## Safety and Maintenance

### 1. IMPORTANT SAFETY WARNINGS

- Do not allow children to touch the moveable parts of this product.
- Please only use the ground power supply included with this product.
- Please unplug the product when not in use, or before cleaning, to avoid injury or causing damage to the product.
- Please operate this product according to this instruction manual.
- Do not use accessories not recommended by the manufacturer.
- This product is not suitable for outdoor use.
- Any other usage that is not listed in this instruction manual is forbidden.
- The recommended usage time of this product is 20 minutes.
- Please do not use the product if the seat cover or leather is broken or damaged.
- Please do not use it if any ventilation holes are covered or blocked.
- Please do not touch internal mechanisms
- Always start with the massage function on a lower strength setting before increasing it to prevent injury.
- Users must stop and turn the power off immediately if they feel ill.
- If you feel that the chair is functioning abnormally, turn off the power switch and unplug the power cord. Call our tech support department immediately. DO NOT attempt to disassemble or repair the chair.

### 2. Usage Environment

- Keep this product out of high temperature environments. Keep away from fire and long-term exposure to the sun.
- Please do not use it in a dusty or caustic environment.
- Please do not use it where there is limited space or poor ventilation.

### 3. People who are not suitable to use this product

- Persons with any other medical conditions in addition to those specified below need to consult their physician prior to use:
- People who are suffering from osteoporosis are not recommended to use this product.
- People who are suffering from heart disease or carrying electronic medical equipment such as a pacemaker are not recommended to use this product.
- Pregnant women or menstruating women are not recommended to use this product.
- People who have injuries or skin disease are not recommended to use this product.
- Children under the age of 14 or people with mental health issues or disabilities are not recommended to use this machine without supervision.
- People diagnosed with medical conditions such as spinal abnormalities are not recommended to use this product.
- People who have come into contact with water or those who are wet should not use this product.

## Safety and Maintenance

### 4. Safety

- Do not get this product wet or let it come into contact with water. This will avoid electric shock and damage to the leather or internal mechanisms.
- Do not use wet cloth or any wet material to clean the electric parts of this product such as the switch and plug.
- Stay away from this product during power interruptions to avoid injury if the power is suddenly restored.
- If this product stops working correctly, please stop using it immediately. Contact the retailer.
- If you feel unwell, stop using the product and consult a doctor immediately.
- This product is not intended for use by individuals, including children with reduced physical, sensory and mental capabilities. This product is also not intended for individuals who lack the knowledge or experience to use it, unless there is proper supervision, or they have been thoroughly instructed on how to use it.
- Children should be supervised around and while using the product.
- If the supply cord is damaged, it should be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid hazards.


### 5. Maintenance

- Product maintenance should only be done by approved service agents. Refrain from disassembling the product yourself.
- Always unplug the power cord after product use.
- Do not use this product if the socket is loose.
- If this product is not often used, fold the wires properly and store in a cool, dry and dustless place.
- Keep this product out of high temperature environments. Keep away from fire and long-term exposure to the sun.
- If the detachable cord/supply cord is damaged, only use replacements from the manufacturer, its service agent or similarly qualified individuals to avoid hazards.
- Use dry cloth to clean this product. Do not use flammable liquids such as thinner, benzene, or alcohol to clean this product.
- The mechanical components of this product are specially designed and made, it is not necessary to have it regularly serviced.
- To move this product on an uneven surface, it should be lifted. Do not pull or roll it.
- Please use this product intermittently; do not use it continuously for extended periods of time.

### 6. General Troubleshooting

- It is normal to hear the motor when the massage chair is in use.
- If the controller does not work normally, please check if the plug and socket are connected properly and make sure that the switch is turned on.
- If the set massage time is over, this product will automatically shut down. If this product is continuously used for a long period of time, the temperature protector will automatically shut down the power. When this happens, this product should be used only after 30 minutes of rest.

### Labelled Massage Chair Diagram

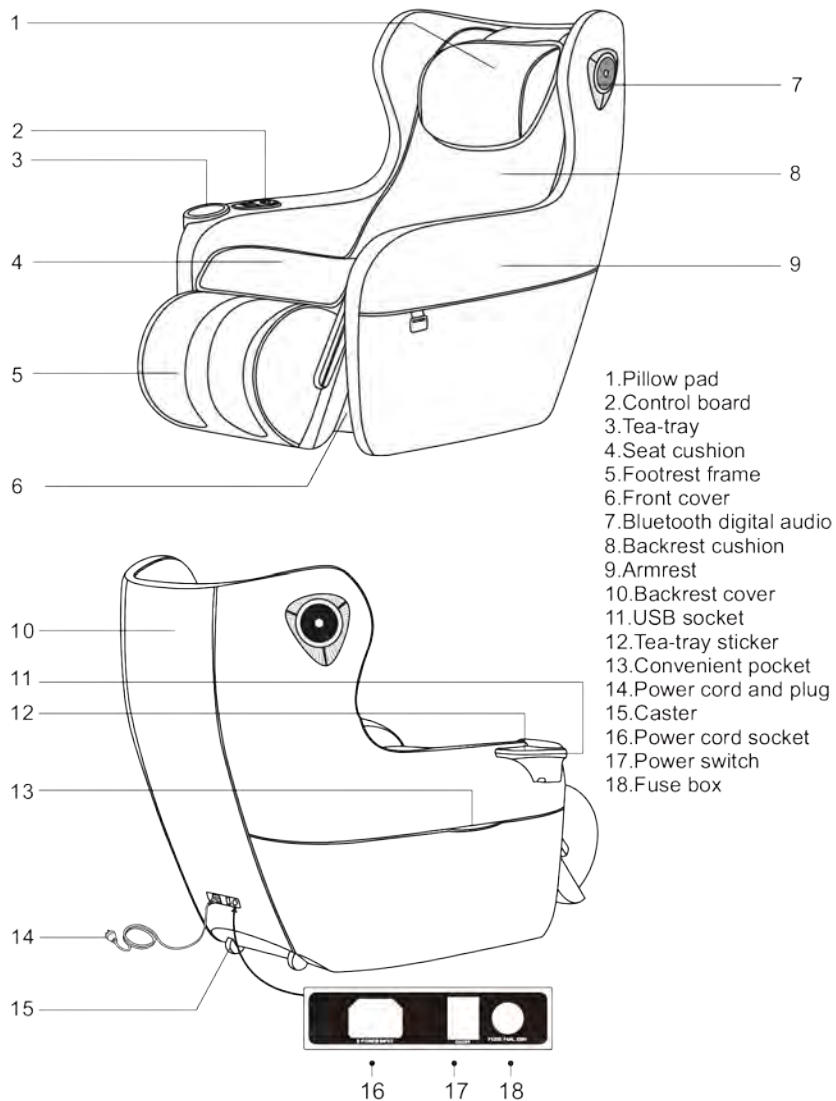

### Overview of Functions

- This product is designed with 2 sets massage mechanism hands. One set is specially designed for neck which can be adjusted up and down. The other set can move up and down, due to its four-wheel drive intelligent design.
- SL type super-long curved rail offers close massage from shoulder to thigh
- Emulates 3 professional massage techniques; kneading, tapping, and simultaneous kneading and tapping. Customise the intensity of the massage hands in kneading and tapping with 3 adjustable speed grades.
- Contains 4 automatic programs: Auto I (full body), Auto II (neck and shoulders), Auto III (lower body) and Demo.
- Fully customisable massage as the massage hands can be controlled inch-by-inch to target any acupuncture point or specific area of the body.
- Includes a specialised neck and shoulder massage program with both one-way and multi-directional massage available. Choose the perfect massage spot with adjustable massage hand positioning.
- Lower body AirSpa massage function with in-built calf airbags, containing 3 adjustable intensity levels.
- Integrated infrared carbon fibre heater will maintain temperatures of a toasty 40 - 50°C.
- Integrated USB port located at the front of the tea-tray.
- A fully reversible footrest means that it can be rotated 180° and folded-away out of sight, enabling the massage chair to double as a comfortable recliner.
- Bluetooth audio speakers are located on the upper armrest and can be accessed via a smart device. Play your favourite music while getting your favourite massage.

Method of Usage

Remote Control Display



Method of Usage

1. Turn on your Massage Chair



2. Start the Massage

- i. Press the red power button to turn on the chair.
- ii. The automatic massage program 'Auto I' will commence automatically. From here you can now customise your massage experience.

A. Automatic Massage Programs

Button	Description	Massage Features	Display
Automatic Massage Program	Auto	Auto I: Full body relaxation massage.	First Auto light on
		Auto II: Focused upper body (neck, back) remedial massage.	Second Auto light on
		Auto III: Focused lower body (waist, glutes, thighs)	Third Auto light on
		Demo: Quick 5 min massage to experience all different functions and features.	Demo light on

B. Manual Massage Techniques

Button	Description	Massage Features	Display
Manual Massage Technique	KNEAD	Kneading-style circular motion massage technique.	Kneading light on
	TAP	Tapping-style massage technique. Three adjustable intensity levels.	Tapping light on
Massage position adjustment		Move the massage hands upward.	
		Move the massage hands downward.	

## Methods of Usage

### C. Neck and Shoulder Massage

Button	Description	Massage Features	Display
Neck and Shoulder Massage		Automatic program designed especially for the neck and shoulders. One-directional and bidirectional massage available.	See relevant lights on remote display
Heater		Switch the infrared heater on or off, which will warm up to a comfortable 40 – 50°C after approx. 3 mins.	Heater light on
Massage position adjustment		Move the neck and shoulder massage hands upward.	
		Move the neck and shoulder massage hands downward.	

### D. Seat AirSpa Massage

Button	Description	Massage Features	Display
AirSpa Massage		Turns on / off the calf massage airbags. Three adjustable intensity levels.	Air light on

### E. Calf Rest Position

Button	Description	Massage Features
Calf Rest Elevation		Lifting the calf rest: Hold this button to lift the calf rest slowly. Once in a comfortable position, release to stop.
Calf Rest Decline		Reclining the calf rest: Hold this button to recline the calf rest slowly. Once in a comfortable position, release to stop.

### 3. Bluetooth Speaker Operation

- Turn on the massage chair using the power button.
- Open your smart device “Settings” app and go into the Bluetooth settings.
- Turn on your Bluetooth to locate nearby devices.
- When “IMCM-0066” is located, click to establish a connection as this is your massage chair (password: 0000 or 1234).
- Open your music player of choice on the smart device and play!

**Please Note:** Every chair has a unique model code for the Bluetooth connection (i.e. “IMCM-0066... IMCM-0061... IMCM-0088...etc.”), so don’t worry if it doesn’t perfectly match the above example.

## Methods of Usage

### 4. Turn Off / Finish your Massage

- Pressing the power button at any time during the massage will stop all massage functions immediately.



### 5. Tea Tray Installation

- Take the chair out of its packaging. *Refer to Figure 1.*
- Take out the cup holder mounting plate from the accessory bag. Align with the mounting screw hole in the right armrest, then screw in 2 x M4\*10 screws into the lock with a Phillips head screwdriver. Then pull out the male head of USB line from the right armrest of the massage chair. *Refer to Figure 2.*



Figure 1

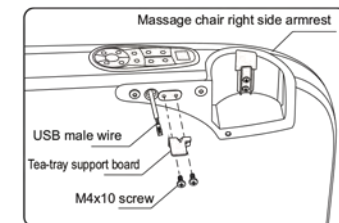


Figure 2

- Take out the upper cover of the tea-tray from the accessory bag, align the slot of the upper cover with the limit board. *Refer to Figure 3.*
- Push the upper assembly down to insert it into the limit board. Make sure it is secured in place. *Refer to Figure 4.*

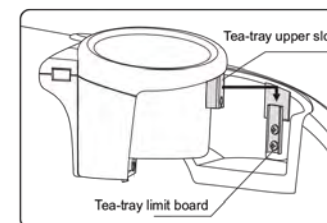


Figure 3

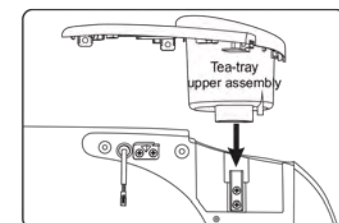


Figure 4

## Methods of Usage

- After confirming the tea-tray upper assembly is installed correctly, use a Phillips head screwdriver to 2 x M4\*10 screws in to secure the tea-tray upper assembly on the right side of the armrest. Refer to Figure 5.
- Connect the male USB wire to the female USB wire, which is on the tea-tray upper assembly. Refer to Figure 6.

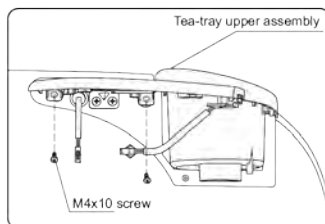


Figure 5

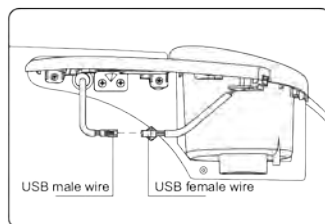


Figure 6

- Take out the tea-tray lower assembly from the accessory bag, align it with the screw hole in the upper assembly, and screw in the 3 x ST2.9\*10 self-tapping screws tightly. Refer to Figure 7.
- Once the tea-tray is installed, locate the threaded screw hold under the tray and screw in the last M4\*10 screw tightly using a Phillips head screwdriver. The installation is now complete. Refer to Figure 8.

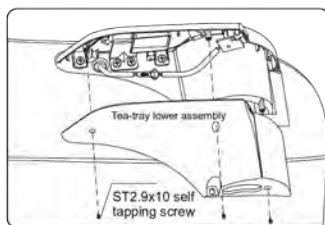


Figure 7

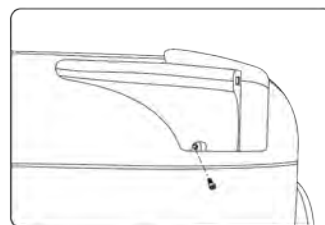


Figure 8

## Methods of Usage

### 6. Protecting your Floor

- Placing the massage chair on timber floors may cause damage after an extended period of time due to its weight, so please place carpet or other material beneath the chair to avoid damaging your floor.




CAUTION: Do not wheel the chair on wooden floors, rough surfaces or in a narrow space.

It is recommended that two people lift this chair.

### 7. Moving your Massage Chair

- Make sure all the cords are off the floor. When moving the chair, tilt the backrest down to an angle where the centre of gravity falls on the back castors (around 25°).

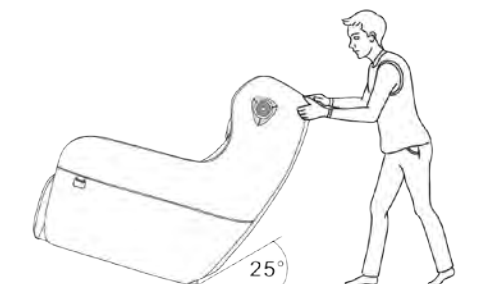

CAUTION: Please be sure to turn off the power at the chair and power point. Then disconnect the power supply and remote control before moving the chair!

## Product Specifications

Model: LW-MH

Name: Masseuse Health Massage Chair

Rated voltage: 110-120V ~ 60Hz

220-240V ~ 50Hz/60Hz

Rated power input: 150W

Rated power input: 150W

Rated time: 20min

Safety structure: Class I

Material introduction: PVC, PA, steel parts and electric & electrons parts

MATERIAL	PU	LEATHER	CLOTH	WOOD
	X			