

AMIRO

Eye Rejuvenation RF Device



User Manual

* We appreciate you buying an AMIRO product. The user manual's images of the item, its accessories, and its user interface are merely for reference.

- Carefully read these instructions before using for your safety.
- Please save this manual for future reference.
- When transferring ownership of the product to others, please include this manual.

Contents

- I. Product Information
- II. Instructions
- III. Care Advice
- IV. Safety Precautions
- V. Table of Common Malfunctions

I. Product Information

Product description

This product is a device for eye skin care. Using it with the AMIRO Eye Mask for Eye Rejuvenation RF Device can effectively improve the eye skin.

Product specifications

- Product name: Eye Care Machine
- Product model: ABF206
- Product size: 140x25x25 mm
- Net weight: 90.3g
- Size of electrode contact surface: 18x18 mm
- Rated input: 5V $\overline{=}$ 2A
- Battery capacity: 700mAh 3.7 V
- Country of origin: China


Product list

- AMIRO Eye Rejuvenation RF Device x 1
- AMIRO Eye Mask for Eye Rejuvenation RF Device x 4

Charging cable x1
User manual x1

II. Instructions


Device illustration

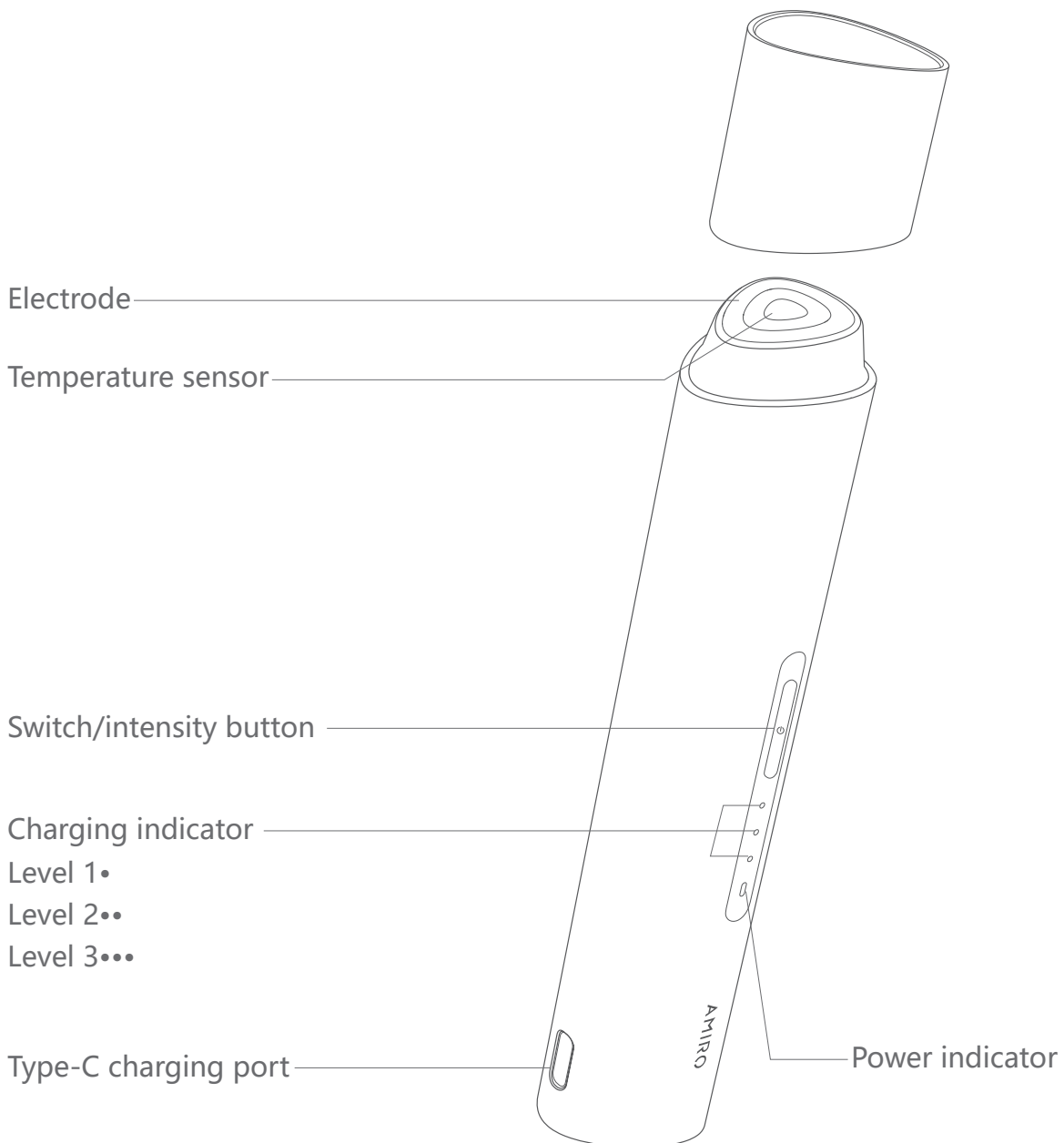


Operating procedure

1. Charging

- Connect the charging cable to a power adapter (rated inputs $\leq 15\text{W}$ and has a compliance certification such as CQC).
- Insert the Type-C connector of the charging cable completely into the charging port of the device.
- When charging starts, the mode/power indicator starts flashing
- The device is fully charged when the power indicator is lit up.
- The device cannot be used when charging.
- Lay the device down when charging to prevent it tipping over.

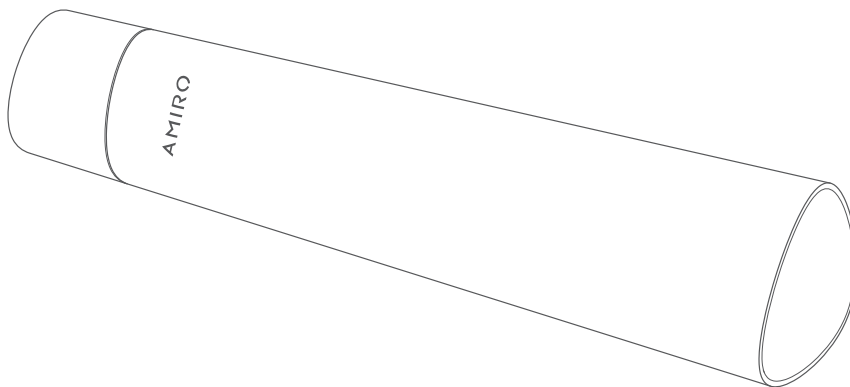

Lay down on the front side

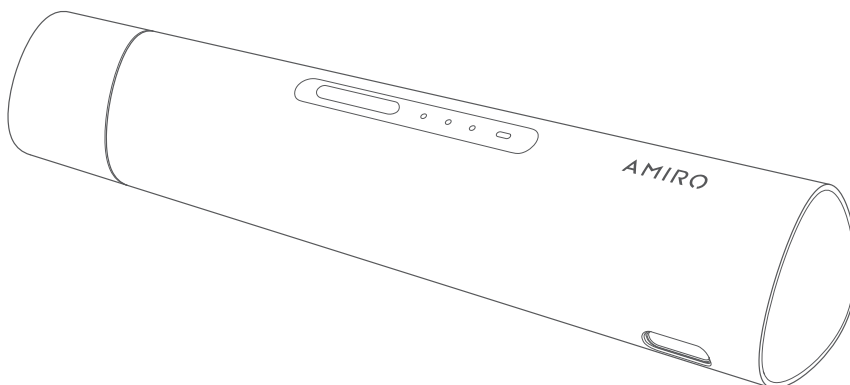

Lay down on the charging side

2. Powering on and off

- When the device is switched off, long press the switch/intensity button for 2s to turn on the device.
- The device will automatically turn off if no operation and no contact with skin is made within 60 seconds after being powered on.
- Long press the switch/intensity button for 2s to power off the device.

3. Intensity selection

- The device' s default intensity is Level 1 (1 indicator on) when first powered on; it remembers the last intensity used when it is powered on subsequently.
- After switching on, press the switch/intensity button briefly to switch intensity. The intensities are switched in the cycle of 1-2-3-2-1, and the intensity light changes accordingly.
- The highest intensity is 3 (3 lights) and the lowest is 1 (1 light).

4. Power indicator


- The white light is on when the power is higher than 21% and the red light is on when the power is lower than 20%.
- When the device is switched off, press the switch/intensity button to see the power status.
- The power indicator is on when the device is switched on.

5. Low-power protection

- The low-power protection mode is turned on when the power is less than 5%.
- Press the switch/intensity button and the red light flashes 10 times.
- Long press the switch/intensity button and the device does not power on.
- When the device enters low battery protection, charge it before using it again.

III. Care Advice

- Set intensity to Level 1 when using the device for the first time
- Use the device 2 to 3 times per week

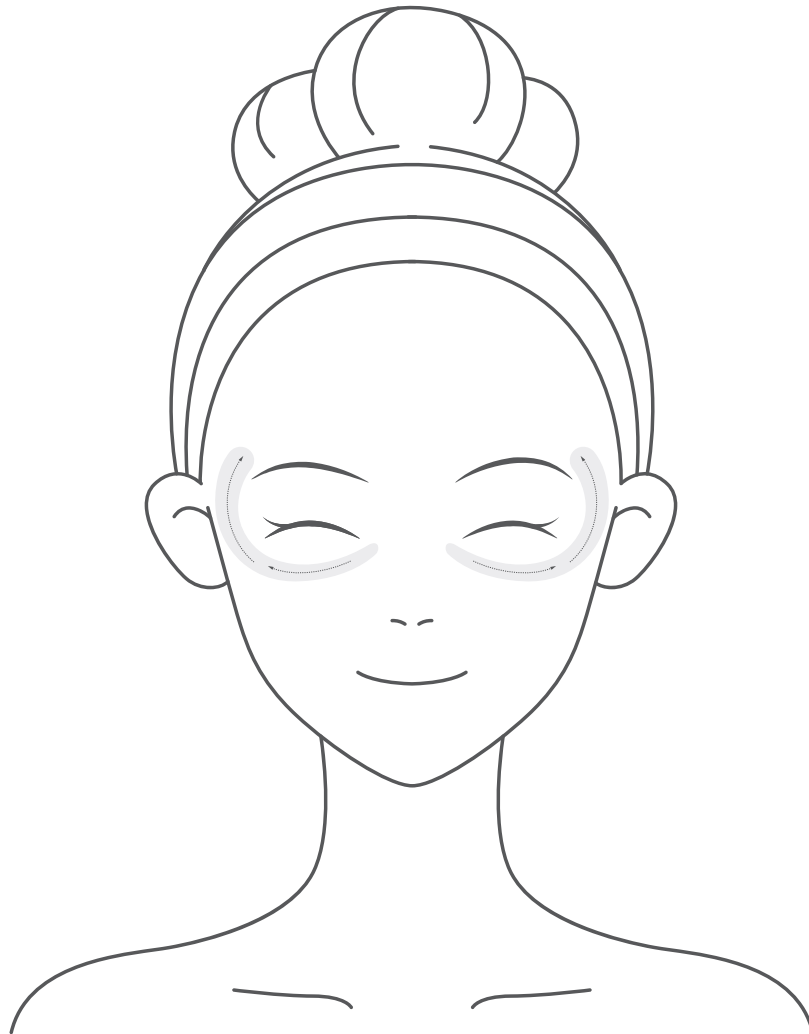

* Use the device with AMIRO Eye Mask for Eye Rejuvenation RF Device for better effect

Care steps

Step 1: Makeup removal and cleansing

- Remove makeup, cleanse your face, and apply toner to keep your face moisturized.




Step 2: Put on the eye mask

Using matching eye masks is recommended

* How to use AMIRO Eye Mask for Eye Rejuvenation RF Device

- Dry your hands, take out the freeze-dried eye mask and gently apply to your under-eye skin
- Take out the active moisture eye mask and cover it on the freeze-dried eye mask

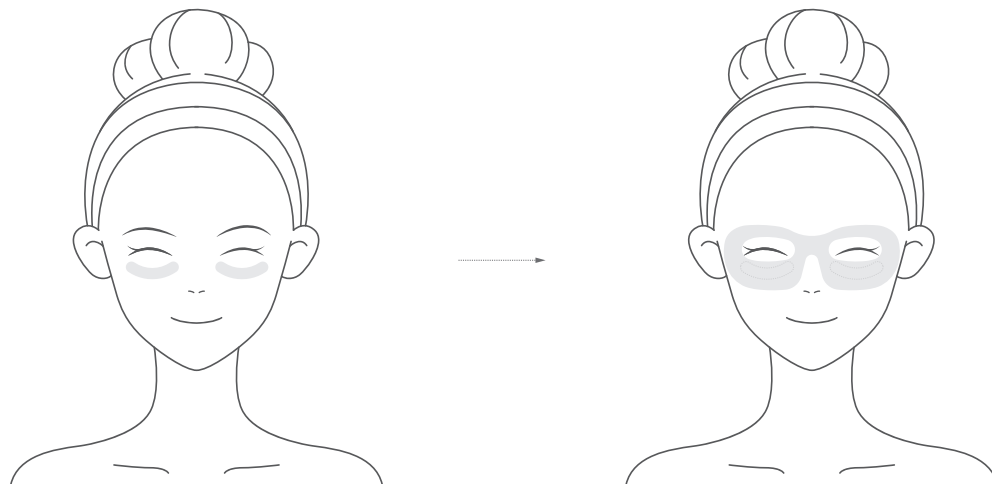

Step 3: Eye care

* The whole care procedure takes about 4 minutes and the device automatically powers off after that.

- Switch on the device and select an intensity level that suits you.
- Press the head of the device against the under-eye skin and the device will start working after a brief vibration (refer to technique illustration).
- During the care procedure, press on the left and right eyes alternatively.
- After a 4-minute care procedure, the device vibrates for 0.5s before powering off automatically.

Treatment techniques (introduction and illustrations)


- Use on the left and right eyes alternatively, for under-eye lines, corners of eyes and arch of eyebrows.
- Press each area for 1-2s, (change to the next area when you feel obvious heat and electric current), and keep the head of the device tight against the skin.

[Note: Do not use the device on the eyeballs and upper eyelids, or on the same area for too long].



Step 4: Care finished

- When using the device with an eye mask, keep the mask on for another 5-10 minutes, massage the skin around the eyes until the essence is absorbed.

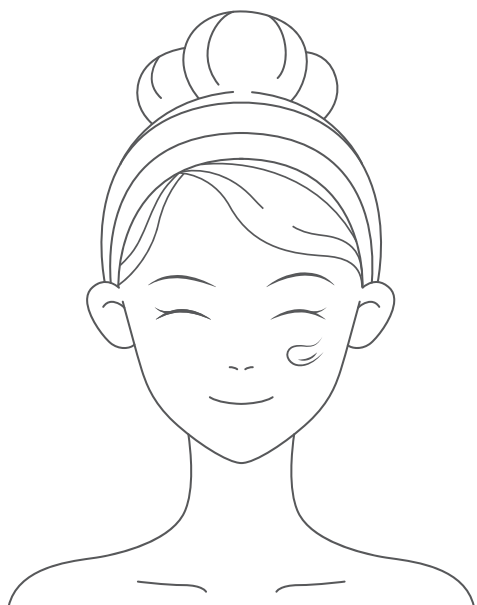

Step 5: Device maintenance

- After the care routine, wipe off the residue on the device with a slightly wet soft cloth or wet tissue before putting it back in the storage box or pouch.

IV. Safety Precautions

Tips:

Carefully read all [Safety Precautions] before using the product and adhere to all guidelines in this handbook.

* [Warning]: Failure to follow the instructions may result in adverse effects.

Operation safety measures:

- The recommended ambient temperature for this product is 10°C-30°C.
- Ensure temperature in the charging environment is between 10-30°C.
- Charge the device first and check that all of the buttons and indicators function normally before using it for the first time or after a long time off.
- Set intensity to Level 1 when using the device for the first time or before your skin gets used to this device.
- Use only disposable consumables recommended by the manufacturer.

Precautions before using the product:

- Press the Power/Mode button to check that the product is fully charged for normal use.
- Before using the product for daily care, remove any makeup and rinse off any remaining makeup remover.
- Before using, always remove any metal items, glasses, and contact lenses from the designated locations.
- Use it exclusively for cosmetic purposes.
- Do not submerge the device (and base) in water.
- Do not use or store the charging base in areas of high humidity.

Product usage precautions:

- Hold the stamp against to the skin during use.
- If you notice any unusual skin issues, stop using the device right away and get in touch with our after-sales service or a medical expert.
- When using the device, keep the RF electrodes pressed closely to the skin; if no eye mask is used, an appropriate amount of gel should be applied; make sure to glide the device slowly and evenly over your skin; if you feel overheated, slide faster or lower the intensity.
- When using the device without an eye mask, do not remain in the same position for a long time to prevent burns. Slide the device over your skin.
- To avoid eye injury, avoid staring directly at an LED for extended periods.
- Do not apply the RF electrodes to the eyelids or around the eyes to prevent damage to eyes.
- Differences in individual skin sensitivity may lead to differences in temperature perception, please adjust the device to an intensity ideal for you.
- Do not use the product on the same area for a prolonged period to avoid skin damage or infection.

After-use Precautions:

- When charging the device, please keep the charging base or charging port dry and clean.
- When cleaning the device, keep the power off and disconnect the charging cable and adapter.
- When cleaning the charging base, please detach the charging cable from the charging base and keep the charging base dry and clean after cleaning for future use.
- It is normal to have slightly red and warm skin after use, please use the device with moisturizing/repairing skincare products.
- Stop using the device right away and consult a doctor if you experience a great deal of discomfort or abnormal feelings while using it.

Precautions

- This product is for people above 18 years old. Do not use this product on minors.
- This product contains small parts. To prevent accidental intake, keep it away from children.
- Use it exclusively for cosmetic purposes.
- Do not use the device on pets, as it could harm them.
- Do not use the product for purposes other than those listed in the manual to prevent personal injury or device failure.
- It is not recommended for anyone with mental or physical disabilities to use this device as it could cause an accident or physical harm.
- Do not use the device if it is malfunctioning, as it may result in bodily harm or danger.
- Do not use the device near flammable substances (such as alcohol, thinner, gasoline, nail polish remover, or hair spray) since it could result in fire, electric shock, or bodily harm.
- Do not charge, use, or store the product close to a fire, heat source, or high-temperature area, since it may cause the device to overheat, ignite, or burst.
- Do not use the device near water sources, in restrooms, or while bathing, as it may result in electric shock, a short circuit, or a fire.
- Do not charge the device in water, as it may cause an electric shock, a

short circuit, or a fire.

- Please avoid dropping the device and the base or any collisions to prevent malfunctions or injuries due to the fall. Do not use the equipment if it has been damaged by a fall or collision.
- Avoid contact of the electrodes with metals as it may cause device malfunction and personal injury.
- Do not wipe the device with corrosive solvents like alcohol, thinner, gasoline, nail polish remover, etc., as it may damage its surface.
- When the equipment is not in use or being cleaned, turn the power off in advance to prevent an electric shock, a short circuit, or a fire.
- Please use the power adapter specified by the company and make sure that the adapter is plugged directly into the socket when using the device; do not use an extension cord to avoid electric shock, short circuit, or fire.
- Please do not disassemble or modify the device since this could result in fire, electric shock, or physical harm. Please get in touch with our after-sales service if the device needs to be fixed.
- Please unplug the power adapter immediately if any abnormalities occur during charging. Please cease operating the equipment right away if there is smoke, noise, or an odor.
- Do not use the device together with flammable cosmetics to prevent fire, electric shock, or physical harm.

Unsuitable groups

People with any of the following conditions may not use this product as it could result in accidents or skin or bodily injury:

- Heart disease, abnormal blood pressure, malignant tumors
- Infectious diseases or fever
- Hemorrhagic disease or cicatricial diathesis
- Dysalgesia
- Consult a doctor if you are receiving dental or oral treatment
- Menstruation, pregnancy or breastfeeding
- Puberty or adolescence
- Metal implants around the application area
- Diseases such as cancer (especially skin cancer and pre-cancerous moles), acute diseases, tuberculosis, febrile diseases, infectious diseases, blood diseases, visceral diseases, skin diseases, allergic dermatitis, etc.

- Impaired immune system due to immunosuppressive diseases such as AIDS, or taking immunosuppressive drugs
- Consult a doctor if you have had permanent plastic fillers within the application area, or you have temporary dermal fillers or Botox implanted within the application area
- Any disease within the application area due to thermal irritation, such as recurrent herpes simplex virus
- Possible blemishes within the application area cause by physical irritation
- Allergic, sensitive or susceptible skin
- Fatigue, illness, under the influence of alcohol or sleeping pills
- Dilated capillaries caused by prolonged use of steroids, hormones, and other drugs, or liver dysfunction
- People using these medical devices: Implanted assistive equipment, such as pacemakers, internal defibrillators, artificial lungs, or portable gadgets, like electrocardiographs

Where device should not be applied:

Do not use this product on the following body parts or if the following conditions exist to avoid malfunction of the device:

- Unspecified parts
- Throat nodules, eyeballs, eyelids, oral cavity, private parts, mucous membranes, etc.
- Lymphatic areas of the face and neck
- Severe acne in the application area
- Purulent and inflamed skin, previous scarring or abnormal wound healing
- If you have undergone surgeries such as laser peeling or chemical peeling in the last 3 months or have healing wounds within the application area
- Skin inflammation or redness due to sunburn
- Pathological pigmentation on the skin
- Protruding moles or viral warts on the skin
- Dermatological conditions such as rashes, inflammation, itchiness, etc.

Note: The list of safety measures above is not all-inclusive. Please contact us or consult a doctor if you have any questions about whether you can use this product

Instructions for removing and discarding spent batteries

Instructions for removing batteries:

1. Remove the battery first before discarding the device.
2. To safely remove the battery, the device must be power-off.

Used battery disposal:

1. Place it in special wastebin, recycling station, waste disposal facility, or any location suited for recycling spent batteries.
2. Give it to the waste battery recycling institutions for professional recycling treatment.

V. Common Malfunctions

Problem	Possible cause	Solution
Start failure	Dead battery	Charge fully before use
	End of battery life	Contact the manufacturer for repair or replacement
Automatic power off	No operation for more than 60 seconds	Restart the device
	Dead battery	Charge fully before use
	Care procedure ends	Automatic power off after a 4-minute care procedure
Charging failure	Mismatched charging device	Use the official charging device provided
Overheated eyes	High intensity	Lower intensity
Low heat	Low intensity	Increase intensity
	Fast sliding of the device	Reduce device sliding
Tingling caused by current	Device is not pressed against the skin	Press the device against the skin
	High intensity	Lower intensity
Skin redness after use	RF heat	It normally recovers soon afterward; please seek medical attention if you don't recover or have other discomforts.