


AMIRO

LumoJelly LED Light Therapy Face Mask User Manual



4-light therapy, gentle silicone for skin safety, compact and portable

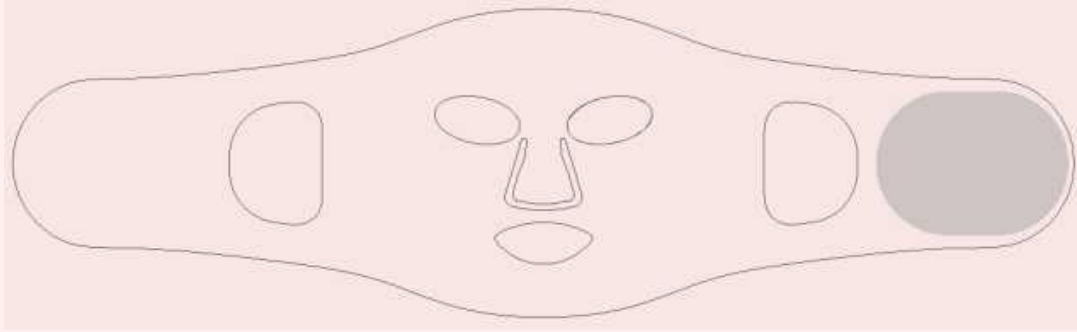

*Thank you for choosing AMIRO products. The illustrations in the manual are for reference purposes only. Please rely on the actual product for accurate information on its appearance, accessories, and user interface. For your safety, please read this manual carefully before use. Please keep this manual for future reference.

Safety Instructions: Please strictly follow the instructions in the manual when using the product. Please consult your doctor if you have any questions regarding this product or the user manual. Alternatively, you may contact us for further assistance. (Contact information can be found on the last page)
V.1.0

Product Parts

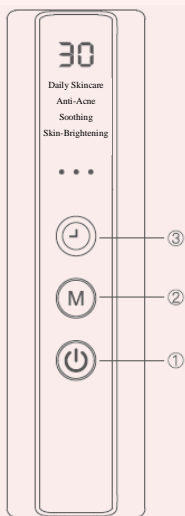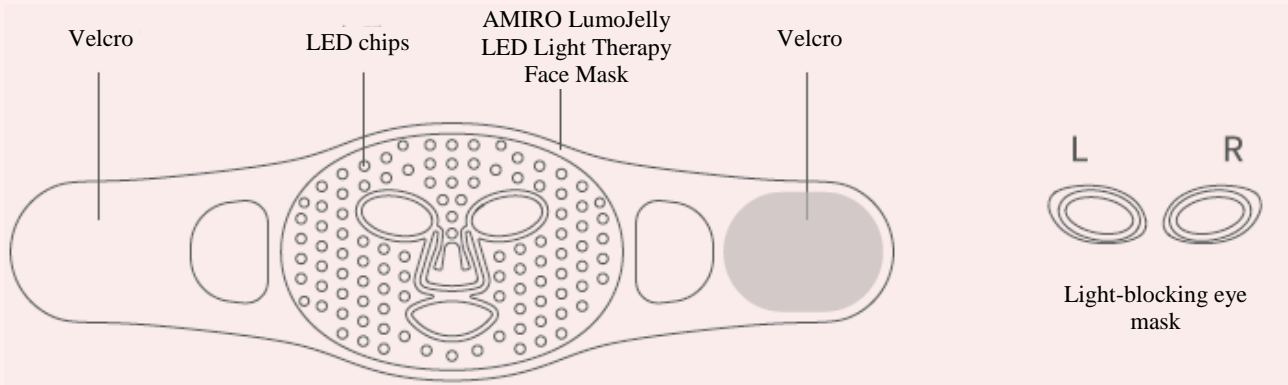

How to Use the Remote Control

01 Power Button

Press and hold the power button for 2 seconds to turn the device on/off. Powering it up will automatically restore the previous timer settings and mode. On the first use of the device, the Daily Skin Care mode will automatically be activated with a default timer set to 10 minutes.

02 Mode Button

The mode that was previously used is remembered by default. Press the button briefly to cycle through the 4 available modes: Daily Skin Care, Anti-Acne, Soothing, and Skin-Brightening.

03 Timer Button

The last timer setting that was previously used is remembered by default. Press the button briefly to adjust the timer duration. There are 3 options available: 10 minutes, 20 minutes, and 30 minutes.

Prohibited Use

Caution: This product is not suitable for individuals belonging to any of the following groups:

- Pregnant, planning to conceive, or breastfeeding
- Individuals with photosensitive skin or are taking medications that may cause photosensitivity
- People with hereditary skin conditions such as lupus, photosensitive eczema, or vitiligo
- Individuals diagnosed as ineligible for this product

Before using the product for the first time, please conduct a photosensitivity test:

- Place the device on your arm and switch it on for 10 minutes of exposure. The presence of symptoms such as redness or itching indicates that your skin may be photosensitive. To avoid side effects, discontinue use and consult a dermatologist before using the product again.
- Before use, cleanse your face and ensure it is dry (makeup-free) to prevent any photosensitive ingredients in cosmetics from affecting the skincare effect.
- It is normal for your skin to appear slightly red after using the product. Typically, the redness subsides within an hour or two and your skin will return to normal. If you feel a burning sensation or any discomfort during the treatment, please stop using the product immediately.

Instructions:

How to Use:

Step 1: Start by thoroughly cleansing your face and gently patting your skin dry.

Step 2: Place the protective eye mask on the face, keeping the mask securely in position by following the "L" and "R" markings. Press and hold the power button on the remote control to turn on the facial mask.

Step 3: The mask is available in four different modes: Daily Skin Care, Anti-Acne, Soothing, and Skin-Brightening. Choose the mode that suits your needs and preferences.

Step 4: You can set the device to run for 10 minutes, 20 minutes, and 30 minutes. The timer button can be used to adjust the time. The device will automatically shut off once the timer is up. Remove the mask to end the treatment.



Cleanse your face and gently pat it dry



Attach the protective eye masks onto the facial mask



Wear the facial mask properly and switch it on



Choose your desired skincare mode and begin the treatment




After the treatment, remove the mask and proceed with your regular skincare routine



Clean the silicone part of the mask with a soft towel or wet wipes

100 sets of LED chips with different light wavelengths:

-  Blue light (wavelength 455-465nm)  Orange light (wavelength 600-610nm)  Red light (wavelength 625-635nm)  Infrared light (wavelength 830-850nm)

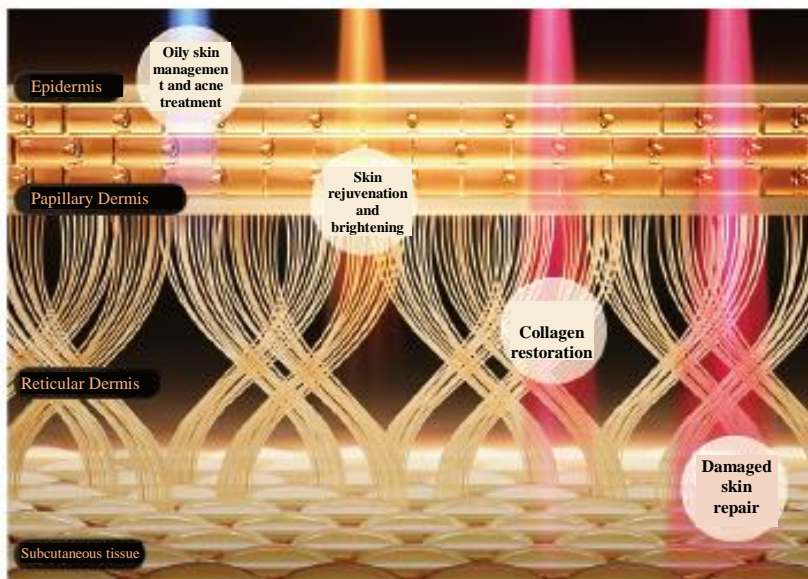

Warranty Card

Purchased from		Date of purchase	
Account/ ID		Contact information	
After-sales service options	Return <input type="checkbox"/> Exchange <input type="checkbox"/> Repair <input type="checkbox"/>	Reason for service	

Warranty Terms and Conditions

Chinese Mainland

• Our company provides a 1-year warranty for customers who purchase this product within the Chinese Mainland. The warranty period starts from the date of receipt (for online purchases, it begins from the logistics receipt date, unless proven otherwise; whereas in physical transactions, the date of physical delivery is taken as the starting date.)

• Our company provides products and services in compliance with relevant laws and regulations such as the “E-Commerce Law of the People’s Republic of China”, “Product Quality Law of the People’s Republic of China”, and “Law of the People’s Republic of China on the Protection of Consumer Rights and Interests”. For offline purchases, if quality issues arise within 7 days, you can contact the seller with a valid proof of purchase or receipt to request a return or exchange; if quality issues arise within 15 days, you can contact the seller with a valid proof of purchase to request an exchange. For online purchases, you have 7 days from the delivery date to return the product without providing any reason, but it must be in sellable condition. Customers are responsible for paying the shipping charges for returning the product within 7 days, except for returns or exchanges due to quality issues.

Countries or regions outside Chinese Mainland

• If you purchased this product in a country or region outside the Chinese Mainland, the warranty period, return and exchange policies, and repair procedures will adhere to the laws and regulations of that specific country or region. For specific matters, please contact the authorized dealer in your location for resolution.

• Upon purchasing the product, please make sure you keep the purchase documents, such as invoices. In the event that a valid invoice or purchase document is not available during the warranty period, the production date on the product body's label will be used to determine the warranty period.

• “Shenzhen Accompany technology CO., Ltd.” and “深圳市宗匠科技有限公司” are the Chinese and English names of our company, respectively.

Please be aware that the following situations are not covered by the free warranty service during the warranty period:

• Damage caused by improper use, such as incorrect power supply connection, use of inappropriate accessories, failure to comply with the instructions, incorrect use, or damage caused by negligence, etc. Repairs or modifications not authorized or approved by our company; other damages caused by force majeure (such as natural disasters, abnormal voltage, etc.).

• Product aging, wear and tear caused by normal use, not affecting normal functionality of the product.

• Free services do not include consumable parts and accessories, such as charging cables. Transportation costs are not included as well.

Safety Warnings

Please read these safety warnings carefully before using the device according to the instructions. These warnings are provided to ensure the correct and safe use of the product and to prevent harm to the user as well as others. Please strictly adhere to them.



Prohibited

Caution

- Only use this product for skincare purposes as instructed. Do not attempt to dismantle, modify, or repair it on your own.
- Do not use or store the device and its accessories in high-humidity environments, such as bathrooms.
- If you experience any pain or discomfort while using the product, please stop using it immediately.
- Do not share this product with other people.
- Do not throw the device to avoid potential collisions and damage.
- Do not touch the product with sharp or high-temperature objects.
- When using the product, please ensure that the eye mask is properly attached to effectively block out light.
- Do not look directly at the LED lights for an extended period of time, as doing so may potentially harm your eyes.
- Please strictly adhere to the recommended treatment time.
- Stop using the device immediately in the event of any malfunction

Warning


- Do not use the device in ways other than those described in the instructions, as this may cause malfunctions or result in accidents.
- If the user suffers from any physical, sensory, or mental illnesses, or is undergoing treatment for dermatological conditions, it is important to ensure that this product is used under the supervision of others or with guidance provided by a professional, such as a dermatologist.
- Do not wrap the power cord around the power plug.
- Do not damage the power cord by bending, twisting, stretching, or placing heavy objects on it.
- Do not place the device in high-temperature environments to avoid the risk of fire and electric shock.
- Do not put the device into fire.
- Do not use the product on children and ensure it is kept out of their reach.
- Do not immerse the device in water or wash it with water.
- When charging the mask, make sure to use an adapter with an output power of less than 15W and that meets the requirements of the GB 4706.1-2005 standard.



Prohibited

Parameter

Product name	AMIRO LumoJelly LED Light Therapy Face Mask	Model	ABM004A	Number of LED Chips	100 sets
Charging Voltage	5 V = 2A	Operating Current	1.6A	Rated Power	6W
Light Wavelengths	Blue light: 455-465nm, Orange light: 600-610nm, Red light: 625-635nm, Infrared: 830-850nm				
List of contents in packaging	AMIRO LumoJelly LED Light Therapy Face Mask, user manual and guide, eye masks, charging cable, remote control				
Battery capacity	2000mAh	Product Dimensions	665×7×193mm		

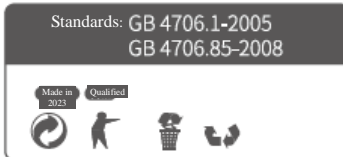

Customer Service Hotline: 4000480888