

# AMIRO

**\* We appreciate your purchase of an AMIRO product. Kindly be aware that the images in the user manual are for reference purposes only and may differ slightly from the actual item, its accessories, or its user interface.**

- For your safety, please read this manual carefully before use.
- Please keep this manual for future reference.
- Please include this manual when transferring ownership of the product to others.

V1.0

# Contents

---

- I. Product Information
- II. Instructions
- III. Usage
- IV. Safety Precautions
- V. Table of Common Malfunctions
- VI. Product Warranty

## I. Product Information

---

### Product Description

This product is an LED face mask designed for carrying out light therapy at home for effective treatment of a range of skin conditions using different wavelengths.

### Product Specifications

- Name: AMIRO Spectra 5-in-1 LED Light Therapy Mask
- Model: ABM502
- Size: 225×192×76 mm
- Input: 5 V=2 A
- Battery Capacity: 3200 mAh
- Country of Origin: China

### Product List

---

AMIRO Spectra 5-in-1 LED Light Therapy Mask ×1

Cable x1

Charging Cable ×1

Controller x1

User Manual x1

Instructions Card x1


Cleaning Cloth x1

## II. Instructions


---

### Device illustration

Spectra 5-in-1 LED  
Light Therapy Mask

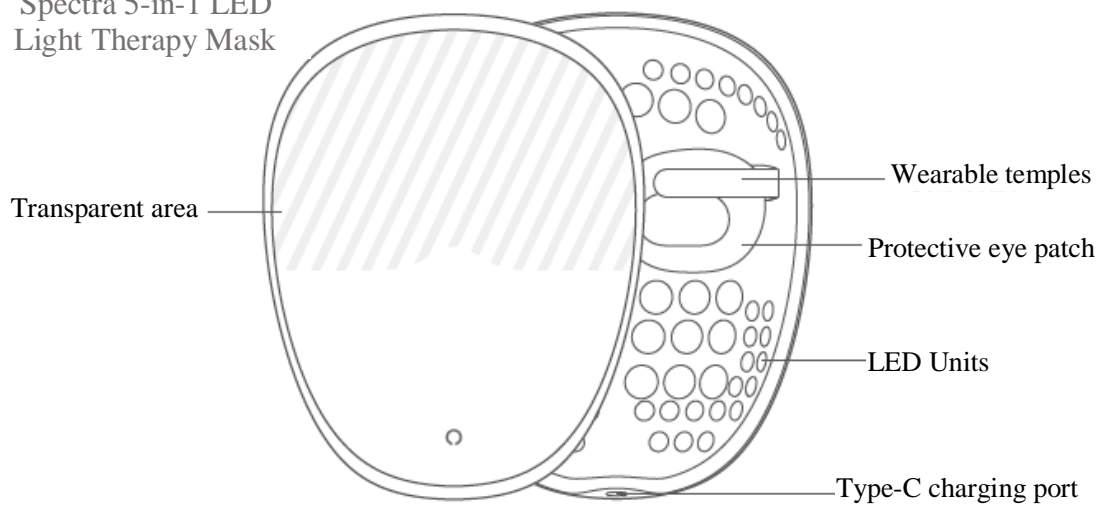

Controller (Power)

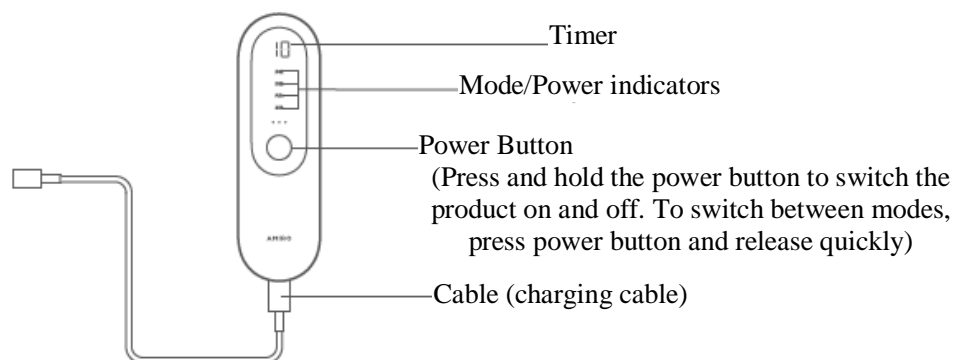

## Operation Procedure

---

### 1. Powering on/off and mode switching

- Connect the controller and the mask using the cable that is included in the packaging
- Press and hold the controller's power button for 2s to turn the device on. Then, press the power button to switch between the four modes: Anti-Acne, Skin-Soothing, Skin-Brightening, and Skin-Tightening
- In each mode, the device will operate for 10 minutes before turning off automatically. Hold down the power button for two seconds to turn off the device before the time is up.



Connect the controller and  
the mask


### 2. Charging the controller

- You only need to charge the controller because the mask is powered by the controller.
- Connect the Type-C charging cable to the adapter. Using an adapter with CQC certification and an output power of 5-10 W is highly recommended. Plug the adapter into the power source and make sure the charging voltage is below 20 V.

- When charging starts, the mode/power indicators starts flashing
- When the mask is fully charged, all three mode/power indicators will be lit up.

### 3. Checking the battery level

- The battery level can be checked by pressing or holding the controller if the mask is not connected.
- The battery level can be checked by pressing the power/mode button of the controller.
- The highest battery level is indicated when all three indicators are illuminated. One illuminated indicator indicates low battery power.
- A low battery level will be indicated by the orange indicators flashing. Please recharge the mask before using it.

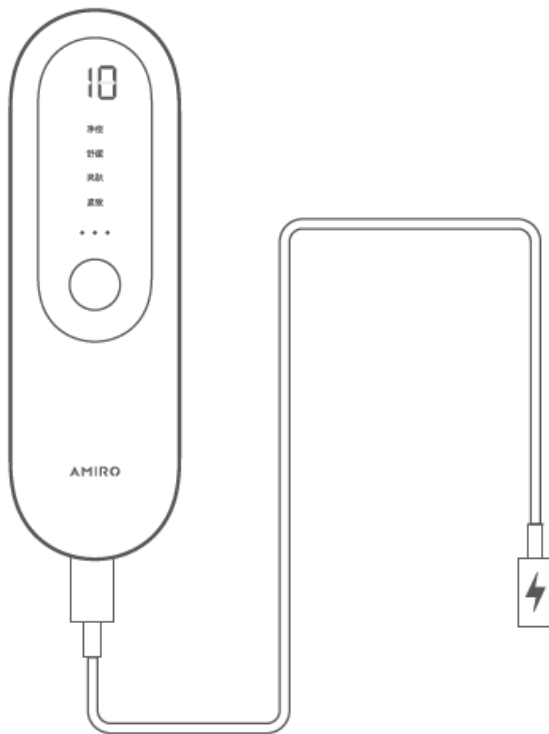


Controller (power)

### III. Method

---

#### Light therapy duration

- The device should only be used 3-5 times per week and for no longer than 30 minutes at a time
- The device will operate in each mode for 10 minutes before automatically turning off. Several modes can be used simultaneously.

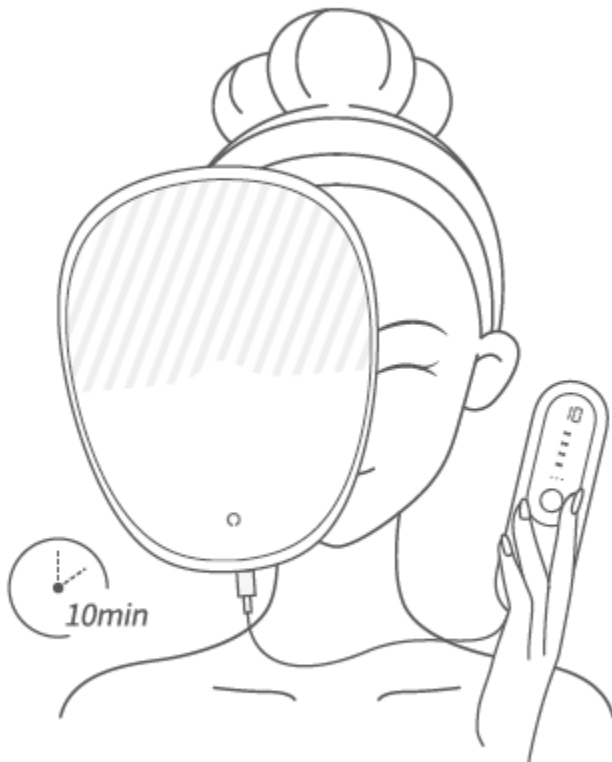

## Procedures

---

### Step 1: Makeup removal and cleansing

- Remove all makeup and cleanse your face before using the device. You may use a toner to moisten the skin, but avoid products that contain photosensitive ingredients or whiteners



### Step 2: Connecting the mask and the controller

- Connect the cable included in the packaging to both the ports on the controller and the mask to ensure proper connection.




### Step 3: Put on the mask and choose your desired mode

- Please select the mode that best suits your skin condition according to the treatment advice. Different modes are suitable for different skin types.
- Applying a facial sheet mask after phototherapy will increase post-treatment effectiveness.

#### [Anti-acne mode]

**Treatment Advice:** After cleansing your face, switch to anti-acne mode and use a moisturizing toner or lotion with simple ingredients. Apply more moisturizing products if your skin still feels dry after using the device. Maintaining a healthy lifestyle is also helpful to lower the likelihood of acne recurrence (including exercise, a healthy diet, and enough sleep)

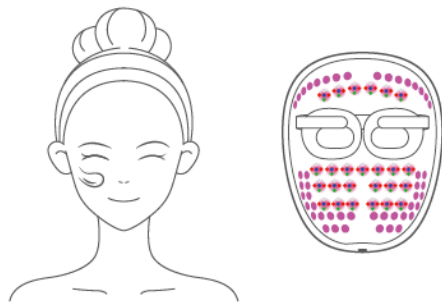

#### [Skin-Soothing mode]


**Treatment Advice:** Set the device to skin-soothing mode and use in conjunction with skin care products made with gentle, non-irritating ingredients. The skin barrier is fragile during the sensitive period, hence daily sun protection is advised






### [Skin-Brightening mode]

**Treatment Advice:** Set the device to skin-brightening mode after removing all makeup and cleansing your face. To double the effect, use skin-brightening products after using the device. Please avoid using any products that contain vitamin C, retinol, or other photosensitive ingredients prior to using the device.

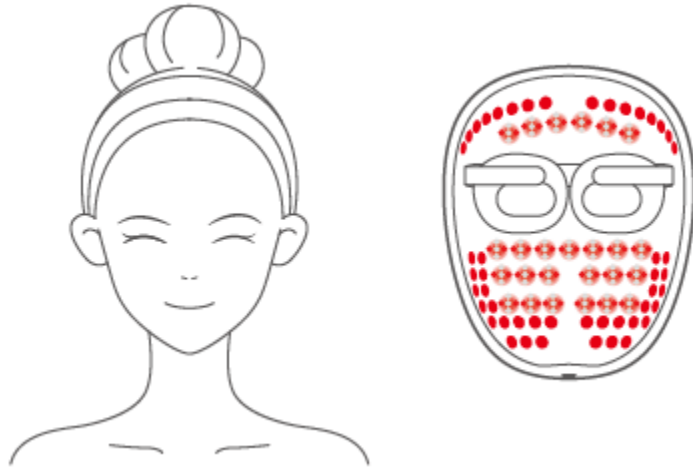

### [Skin-Tightening mode]

**Treatment Advice:** For cleansing, use gentle, non-irritating products, and for moisturizing, use a hydrating product that's free of photosensitive ingredients. After that, switch to skin-tightening mode to begin the treatment.

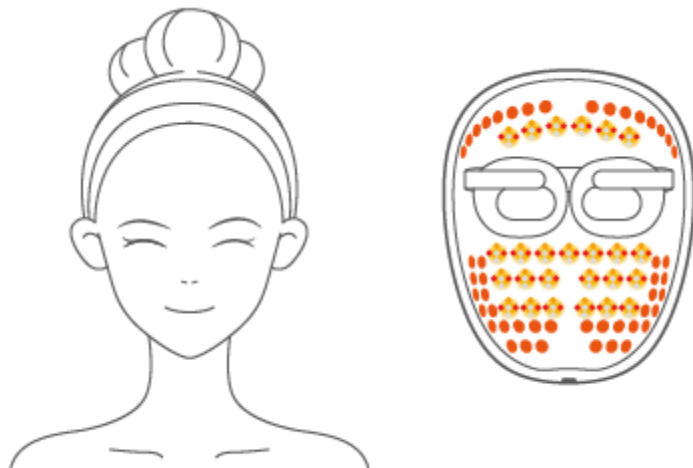

### Step 4: Device maintenance

- Place the device back into its packaging box or storage bag after use. Gently wipe it down with a wet cloth or tissue if it gets dirty during treatment.

## IV. Safety Precautions

---

Precautions:

- This product is intended for people over the age of 18 and should not be used by minors.
- This product contains small parts. To prevent accidental intake, keep it away from children.
- Do not use the device for non-cosmetic purposes.
- Do not use the device on pets as it may cause injury to them.
- Do not use the product for purposes other than those listed in the manual to prevent personal injury or device failure.
- People with physical or mental disabilities should not use this device as it may result in unwanted accidents or bodily injuries.
- Do not use the device if it is malfunctioning as doing so may result in bodily injury or danger.
- To prevent fire, electrocution, or bodily harm, please refrain from using the device near flammable substances such as alcohol, thinners, gasoline, nail polish remover, or hair spray.
- Do not charge, use, or leave the product near a fire, heat source, or high temperature area, as doing so may cause the device to overheat, catch fire, or explode.
- Do not use this device near water sources, in bathrooms or while bathing, as this may result in electrocution, short circuit, or fire.
- Do not charge the device in water, as it may cause electrocution, a short circuit, or a fire.
- The mask or controller should not be dropped, or any impact on them should be avoided in order to prevent malfunctions. Do not use the device if it has been damaged after a collision or being dropped.
- Do not wipe the device with any corrosive solvents (alcohol, thinner, gasoline, nail polish remover, etc.) as doing so may damage its surface.
- Turn the power off when the device is not in use or being cleaned to prevent electrocution, short circuit, or fire.
- Make sure to use the power adapter specified by the company and plug the adapter directly into the socket when using the device. Do not use an extension cord to prevent electrocution, short circuit, or fire.
- Please do not disassemble or modify the device as doing so could result in fire, electrocution, or physical harm. Please contact our after-sales service if the device needs to be repaired.

- Please unplug the power adapter immediately if any abnormalities occur during charging. Please stop operating the equipment immediately if there is any smoke, noise or odor.
- Do not use the device near combustible cosmetics to prevent fire, electrocution, or physical harm.

### **Tips:**

Read all [Safety Precautions] carefully before using the product and adhere to all guidelines in this handbook.

- \* [Warning]: Failure to follow the instructions for use may result in adverse effects.
- \* [Warning]: The device could cause damage to the skin if it is used repeatedly or for a long period of time on the same area.
- \* [Warning]: Avoid bypassing or removing the built-in safety mechanism.
- \* [Device]: Charge the device using an adapter with a maximum measured output power of less than 10 W and ensure that it complies with the requirements of the GB4706.1-2005 standard.

### **Attention:**

- The charging environment should be maintained at a recommended temperature between 10°C and 30°C.

Ensure the temperature in the charging environment is between 10-30°C.

- Before using the device for the first time or after a long period of time, please charge it and make sure everything is working as it should.

As soon as you receive the product, it is important to conduct a photosensitivity test. The use of this product should be discontinued immediately if you experience any discomfort due to photosensitive skin.

### **Caution:**

- Use of this product may cause discomfort or adverse effects if you suffer from

allergies, skin conditions, allergic dermatitis, or sensitive skin. Stop using the product immediately if you experience any adverse effects.

- Avoid applying it to areas with compromised skin integrity (e.g. with burns, blisters, cuts, wounds, or active or infectious skin diseases)
- Do not use it on a child
- Check that the device is fully charged before using it
- Before using the product for daily care, remove any makeup and rinse off any remaining makeup remover.
- Make sure to always remove any metal items, glasses, and contact lenses from the designated locations before using the device.
- Do not use the device for non-cosmetic purposes.
- Avoid immersing the device (and controller) in water
- Avoid using or keeping the controller and mask in damp places

### **Precautions when using the product:**

- When using the device, make sure to wear the mask as instructed to ensure the eye patch fits properly and there is no light entering the eyes.
- Contact the after-sales service or seek medical attention immediately if abnormal skin conditions develop.
- Avoid looking directly at the LED light for a prolonged period of time to prevent damage to the eyes.
- Select an appropriate mode based on your skin conditions.
- To prevent skin injury, do not use the device for longer than 30 minutes in a single day.

### **Precautions after using the product:**

- Please keep the controller or charging port dry and clean when the device is charging.

- When cleaning the controller, disconnect the charging cable and adapter while ensuring the power is off.
- Please remove the charging cable from the mask before cleaning it. Keep the mask dry and clean after cleaning it for future use.
- It is normal for the skin to become slightly warm and dry after using the device. The device can be used with moisturizing skin care products
- Stop using the device right away and consult a doctor if you experience excessive discomfort or any abnormal symptoms while using it.

## **The device is unsuitable for the following groups of people:**

---

People with any of the following conditions may not use this device as it could result in accidents or skin and bodily injury:

- Heart disease, abnormal blood pressure, malignant tumors
- Infectious diseases or fever
- Hemorrhagic disease or cicatricial diathesis
- Sensation disorders
- If you are undergoing dental or oral treatment, consult a doctor before using the device
- Menstruation, pregnancy, or breastfeeding
- Photosensitivity, photosensitive skin conditions, or under photosensitive medications
- Skin conditions, such as lupus erythematosus, dermatomyositis, porphyria, and albinism
- Presence of metal implants in the targeted area
- Suffering from, or having previously suffered from, cancer (especially skin cancer and precancerous moles), acute illnesses, tuberculosis, febrile diseases, infectious diseases, blood diseases, visceral diseases, skin diseases, or allergic dermatitis
- Impaired immune system due to immunosuppressive diseases such as AIDS or the consumption of immunosuppressive drugs
- Please consult a doctor before using the device if you have had permanent dermal fillers injected or have received temporary dermal fillers/Botox injections in the facial areas of application
- Diseases caused by thermal irritation within the application area, such as recurrent herpes simplex
- Physical irritation within the application area may result in blemishes
- Experiencing fatigue, illness, under the influence of alcohol or on sleeping pills
- Dilated capillaries caused by prolonged use of steroids, hormones, and other drugs,

or liver dysfunction

- Individuals using these medical devices: Implanted assistive devices like pacemakers, internal defibrillators and artificial lungs, as well as portable devices like electrocardiographs

## **Areas on the body where the device should not be applied:**

---

To avoid malfunctions, do not use this product on the following parts of the body or if any of the following conditions exist:

- Unspecified parts
- Throat nodules, eyeballs, eyelids, inside the mouth, private parts, mucous membranes, etc.
- Face and neck lymphatic areas
- Purulent and inflamed skin, previous scarring, or abnormal wound healing
- You should not use this product if you have had laser peeling or chemical peeling in the last three months or have open, healing wounds near the application site
- Pathological pigmentation on the skin
- Protruding moles or viral warts on the skin
- Dermatological conditions such as rashes, inflammation, itchiness, etc.

Safety Instructions: The list of safety measures above is not exhaustive. Should you have any questions regarding the use of this product, please do not hesitate to contact us or consult a physician.



## Instructions for removing and discarding used batteries

---

Instructions for removing batteries:

1. Do not replace the battery during routine use.
2. Battery leakage hazards: When the electrolyte of a battery comes into contact with your skin or clothing, immediately wipe it off with a dry cloth as it is extremely corrosive.
3. Remove the battery before discarding the device. as per the steps listed below:
  - ① Remove the product casing
  - ② Remove the internal wires, then take out the battery
4. The device must be switched off and unplugged from the power source before battery removal.
5. Safely dispose of the battery.

Used battery disposal:

1. Place it in special waste bin, recycling station, waste disposal facility, or any location suited for recycling spent batteries.
2. Hand them over to waste battery recycling institutions for professional recycling treatment.

## V. Common Malfunctions

Problem	Possible causes	Solution
Unable to switch on the device	Battery depleted	Charge to full before use
	End of battery life	Contact the manufacturer for repair or replacement
Automatic power off	Battery depleted	Charge to full before use

Problem	Possible causes	Solution
Automatic power off	Light therapy ends	Automatic power off after a 10-minute light therapy
Charging failure	Mismatched charging device	Use the official charging device provided
Warm skin	The effects of light therapy on the skin	Normal
Dry skin	Skin dehydration after light therapy	Use moisturizing skin care products

## VI. Warranty

---

### Chinese Mainland

- Customers who purchase this product within Chinese Mainland are covered by a 1-year warranty. Warranty periods begin with the date of receipt for online purchases (unless proven otherwise, logistics receipt dates are taken for online purchases; for physical transactions, physical delivery dates are taken for physical transactions).
- Our company delivers products and services in compliance with relevant laws and regulations such as the “E-Commerce Law of the People’s Republic of China”, “Product Quality Law of the People’s Republic of China”, and “Law of the People’s Republic of China on the Protection of Consumer Rights and Interests”. When quality issues arise with in-store purchases within 7 days, contact the seller with a valid receipt or proof of purchase to request a return or exchange. An exchange can be requested from the seller if quality issues arise within 15 days of purchase. For online purchases, you have 7 days from the delivery date to return the item without providing a reason, but it must be in sellable condition. Customers are responsible for the shipping charges for returning items within 7 days, except for returns or exchanges due to quality issues.

## **Countries or regions outside Chinese Mainland**

- If you purchased this product in a country or region outside the Chinese Mainland, the warranty period, return and exchange policies, and repair procedures will be in accordance with the laws and regulations of that specific country or region. For specific matters, please contact the authorized dealer in your location for resolution.
- Be sure to keep all purchase documentation, including invoices, after purchasing the product. The warranty period will be determined by the production date on the label of the product body if a valid invoice or purchase document is not available during the warranty period.
- Accompany Tech and Shenzhen Accompany Technology Co., Ltd. refer to the same company.

## **Please be aware that the following situations are not covered by the free warranty service during the warranty period:**

- Damage caused by improper use, such as incorrect power supply connection, use of inappropriate accessories, failure to follow the instructions, incorrect use, or damage caused by negligence, etc. Repairs or modifications not authorized or approved by our company; other damages caused by force majeure (such as natural disasters, abnormal voltage, etc.).
- Product aging and wear and tear caused by normal use of the product, which does not affect its normal functionality.
- Free servicing does not cover consumable components and add-ons, such as adapters. Delivery charges are not included as well.

# Repair information

**Purchase platform:**

**Purchase date:**

**Account/ID:**

**Contact number:**

**After-sales service type:** **Return** **Exchange**  
**Repair**

**Reason:**

**AMIRO**

**Spectra 5-in-1 LED Light Therapy Facial Mask**



**User Manual**