


AMIRO L1 LED Light Facial Mask

- User Manual

AMIRO
LED Light Facial Mask

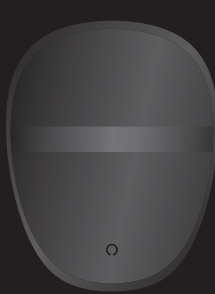

User Manual

* We appreciate you buying an AMIRO product. The user manual's images of the item, its accessories, and its user interface are merely for reference.

- Carefully read these instructions before usage for your safety.
 - Please save this manual for future reference
 - When transferring ownership of the product to others, please include this manual.

AccompanyTech
Customer Service Tel: 400-048-0888
Manufacturer: Shenzhen Accompany Technology Co., Ltd.

1

Contents

I. Product Information
II. Instructions
III. Usage
IV. Safety Precautions
V. Table of Common Malfunctions
VI. Warranty

I. Product Information

Product description
This home LED light facial mask uses lights with various wavelengths to repair and shield facial skin

Product specifications

- Product Name: AMIRO L1 LED Light Facial Mask
- Product model: ABM502
- Product size: 225x192x76mm
- Rated input: 5V 2A
- Battery capacity: 2500mAh
- Country of origin: China

Product list


AMIRO L1 LED Light Facial Mask x1
Cable x1

2


II. Instructions

Device illustration

LED Light Facial Mask

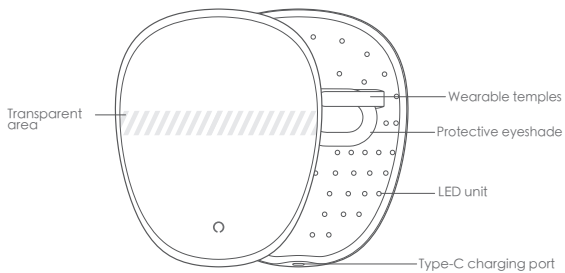

Controller (power)

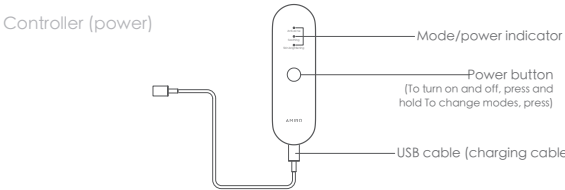


3

Operating procedure

1. Powering on and off and mode switching

- Connect both the controller and the mask port using the official cable provided
- To switch between the three modes, press and hold the controller's power button for 2s to power on, then press the power button: Anti-acne, soothing, and skin-brightening modes
- For each mode, the device will operate for 10 minutes before automatically turning off. Press and hold the power button for 2s to turn the device off in advance.

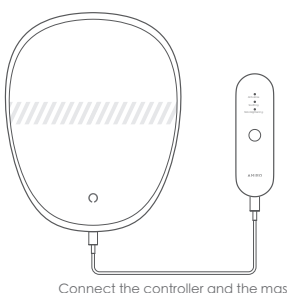

Connect the controller and the mask

2. Charging the controller


- Only the controller needs to be charged individually because it supplies the electrical energy for the mask's lighting modules
- Plug in a power adapter (one that supports a rated output $\leq 10W$ and has a compliance certification like CQC) after connecting the Type-C charging cable to it

4

- When charging starts, the mode/power indicator starts flashing
- All 3 mode/power indicators are lit up when fully charged.

3. Checking the power status

- You can press or press and hold the controller to check the power status even if the mask is not connected
- You can use the controller switch/mode button to verify the power status while the mask is connected and powered off
- 1 indicator represents low power, while 3 indicators represent high power.
- When only one indicator is left flashing rapidly 5 times before going off, the device has insufficient power and enters low battery protection, and shall be charged before using it again.

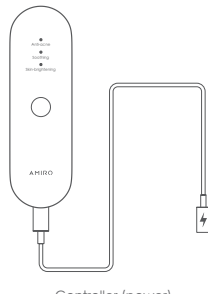


Controller (power)

5

III. Usage

Care duration

- The device should only be used 3-5 times per week and for no longer than 30 minutes at a time
- For each mode, the device will operate for 10 minutes before automatically turning off. More than one mode can be used at a time

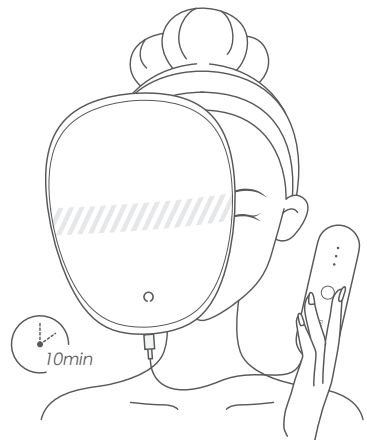


6

Care procedure


Step 1: Makeup removal and cleansing

- Remove any makeup and cleanse your face before using the device. You can then apply toner for basic moisturizing, but avoid using any products that include photosensitive or whitening components



Step 2: Connecting the mask and the controller

- Plug the officially provided cable into the controller and the mask port, and verify the successful connection




7


Step 3: Put on the mask and start the care mode

- Please select the mode that best suits your skin condition using the care recommendations. Different modes are suitable for different types of skin status
- After the phototherapy care, applying a sheet facial mask will bring greater caring efficacy


[Anti-acne mode]
 Suitable for excessively-oily skin, and acne or pimple breakouts
 Care advice: Turn on the anti-acne mode after face cleansing and use a simple-ingredient moisturizing toner or lotion. If your skin still feels dry after using the device, apply more moisturizing products. Maintaining a healthy lifestyle is also helpful to lower the likelihood of acne recurrence (including exercise, a healthy diet, and enough sleep)

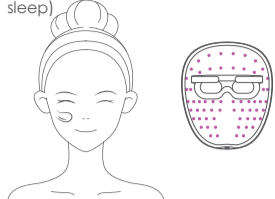

[Soothing mode]
 Suitable for sensitive skin during seasonal changes, sunburn aftercare, healing from allergic reactions, or cosmetic surgery scars recovery
 Care advice: Set the device to soothing mode and use skin care products with gentle, non-irritating ingredients. The skin barrier is fragile during the sensitive period, hence daily sun protection is advised

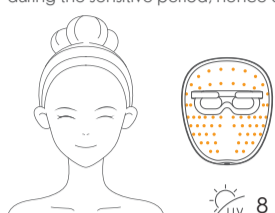

[Skin-brightening mode]
 Suitable for everyday skin-brightening or skin rejuvenation after staying up late
 Care advice: Set the device to skin-brightening mode after face cleansing and makeup removal. To double the effect, use skin-brightening products after using the device. Before using the device, please avoid using any products containing vitamin C, retinol, or other photosensitive or irritating components

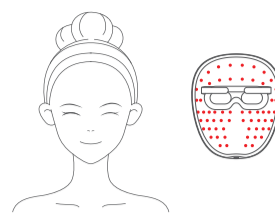

Step 4: Device maintenance

- After completing the care routine, you can place the device back in its storage box or pouch. If there are any smudges, wipe them off with a wet tissue or a soft, slightly damp cloth

9

IV. Safety Precautions

Tips

Carefully read all [Safety Precautions] before using the product and adhere to all guidelines in this handbook

* [Warning]: Failure to follow the instructions for use may result in adverse effects

* [Warning]: Skin damage could develop by using the device on one area repeatedly or for an extended period of time

* [Warning]: Avoid bypassing or removing the built-in safety mechanism

Operation safety measures:

- The recommended ambient temperature for this product is 10°C-30°C
- Be sure to charge the device in an ambient temperature between 10-30°C
- Charge the device first and check that all of the buttons and indicators function normally before using it for the first time or after a long time off
- When you use it for the first time, perform a photosensitive test on your arm

Precautions before using the product:

- For those with allergies, skin conditions, allergic dermatitis, sensitive skin, etc., using this product may cause discomfort or adverse effects. In case of adverse effects, immediately cease its use
- Avoid applying it to areas with compromised skin integrity (e.g. with burns, blisters, cuts, wounds, or active or infectious skin diseases)
- Do not use it on a child
- Verify that the device is fully charged before using it

10

- Before using the product for daily care, remove any makeup and rinse off any remaining makeup remover
- Before using, always remove any metal items, glasses, and contact lenses from the designated locations
- Use it exclusively for cosmetic purposes
- Avoid immersing the device (and controller) in water
- Avoid using or keeping the controller in damp places

Product usage precautions:

- Wear the mask as directed when using the device to ensure it fits the eyes without light leakage
- If you notice any unusual skin issues, stop using the device right away and get in touch with our customer support team or a medical expert
- To avoid eye injury, avoid staring directly at an LED for extended periods
- Choose an appropriate setting based on your skin conditions
- To prevent skin injury, do not use the device for longer than 30 minutes in a single day

After-use Precautions:

Please keep the controller or charging port dry and clean when charging the device

- Disconnect the charging cable and adapter and keep the power off while cleaning the controller
- Please remove the charging cable from the mask before cleaning it, and keep the mask dry and clean after cleaning it for future use
- After using the device, it is normal for the skin to become slightly warm and dry. Use moisturizing skin care products with the device
- Stop using the device right away and see a doctor if you experience a great deal of discomfort or abnormal feelings while using it

11

Precautions

- Do not use this product on anyone under the age of 18; it is intended for adults only
- This product contains small parts. To prevent accidental intake, keep it away from minors
- Use it exclusively for cosmetic purposes
- Do not use the device on pets, as it could harm them
- Do not use the product for purposes other than those listed in the manual to prevent personal injury or device failure
- It is not recommended for anyone with mental or physical limitations to use this device as it could cause an accident or physical harm
- Do not use the device if it is broken, as it may result in bodily harm or danger
- Do not use the device near flammable substances (such as alcohol, thinner, gasoline, nail polish remover, or hair spray) since it could result in fire, electric shock, or bodily harm
- Do not charge, use, or store the product close to a fire, heat source, or high-temperature area, since it may cause the device to overheat, ignite, or burst
- Do not use the device near water sources, in restrooms, or while bathing, as it may result in electric shock, a short circuit, or a fire
- Do not charge the device in water, as it may cause an electric shock, a short circuit, or a fire
- To prevent malfunctions, avoid dropping the mask or the controller, or any impacts on them. Do not use the equipment if it has been damaged by a fall or collision
- Do not wipe the device with corrosive solvents like alcohol, thinner, gasoline, nail polish remover, etc., because doing so could harm the device's surface
- When the equipment is not in use or being cleaned, turn the power off in advance to prevent an electric shock, a short circuit, or a fire
- Please use the power adapter specified by the company and make sure that the adapter is plugged directly into the socket when using the device; do not use an extension cord to avoid electric shock, short circuit, or fire
- The device should not be disassembled or modified since this could result in fire, electric shock, or physical harm. Please get in touch with our after-sales support if the device needs to be fixed

12

- Please unplug the power adapter immediately if any irregularities are observed during charging. Please cease operating the equipment right away if there is smoke, noise, or an odor
- Do not use the device near combustible cosmetics to prevent fire, electric shock, or physical harm

Unsuitable groups

This product should not be used by anyone who has any of the following conditions as it could result in accidents or harm to the skin or other parts of the body:

- Heart disease, abnormal blood pressure, and malignant tumors
- Infectious diseases or fever
- Hemorrhagic disease or cicatricial diathesis
- Dysalgasia
- If you are undergoing dental or oral therapy, consult a doctor before using
- Menstruation, pregnancy, or breastfeeding
- Photosensitivity, photosensitive skin conditions, or taking photosensitive medications
- Skin conditions, such as lupus erythematosus, dermatomyositis, porphyria, and albinism
- Metal implants around the application area
- Current or past illnesses such as cancer (particularly skin cancer and pre-cancerous moles), acute illnesses, TB, febrile diseases, infectious diseases, blood diseases, visceral diseases, skin diseases, allergic dermatitis, etc.
- Immune system impairment brought on by immunosuppressive conditions like AIDS or using immunosuppressive medications
- If you've had temporary dermal fillers, Botox, or permanent plastic fillers placed within the application area, consult a doctor before using
- Any condition that causes thermal sensitivity in the application area, such as recurrent herpes simplex virus

13

- Possible physical irritation-related blemishes in the application area
- Fatigue, illness, drunk, or using sleeping medicines
- Hepatic dysfunction or dilated capillaries brought on by extended use of steroids, hormones, and other medications
- Individuals that use the following medical devices: Implanted assistive equipment, such as pacemakers, internal defibrillators, artificial lungs, or portable gadgets, like electrocardiograms

Locations where the device shouldn't be used:

To prevent the device from malfunctioning, avoid using this product on the following body parts or if the following circumstances exist:

- Non-designated body parts
- Adam's apple, eyeballs, eyelids, oral cavity, private parts, mucous membranes, etc.
- Face and neck lymphatic areas
- Puerulent and inflamed skin, previous scarring, or abnormal wound healing
- If you have been undergoing a surgical procedure such as a laser or chemical peel within the last 3 months, or if you have healing wounds close to the application area
- Pathological pigmentation on the skin
- Protruding moles or viral warts on the skin
- Dermatological disorders, such as rashes, inflammation, itching, etc.

Note: The list of safety measures above is not all-inclusive. Please contact us or consult a doctor if you have any questions about whether you can use this product

14

Instructions for removing and discarding spent batteries

Instructions for removing batteries:

1. During routine use, avoid replacing the battery
2. Battery leakage hazards: Battery electrolyte is extremely corrosive; if it comes in contact with your skin or clothing, wipe it off immediately with a dry cloth
3. Remove the battery first before discarding the device, as per the steps listed below:
 - ① Remove the product casing
 - ② Remove the internal wires, and take out the battery
4. The device must be turned off and unplugged from the power source before removing the battery
5. Safely dispose of the battery

Used battery disposal:

1. Place it in a trash can, recycling container, waste disposal facility, or any location suited for recycling spent batteries
2. Give it to a reputable company that recycles spent batteries

V. Common Malfunctions

Problem	Possible cause	Solution
Start failure	Dead battery	Charge fully before use
	End of battery life	Contact the manufacturer for repair or replacement
Automatic power off	Dead battery	Charge fully before use

15

Problem	Possible cause	Solution
Automatic power off	Care procedure ends	Automatic power off after a 10-minute care
Charging failure	Mismatched charging device	Use the official charging device provided
Warm skin	Effect after phototherapy	Normal
Dry skin	Dehydrated skin after phototherapy	Use moisturizing skin care products

VI. Warranty

Chinese Mainland

- Customers who buy the product in the Chinese mainland are eligible for a warranty that lasts for 1 year starting from the day of receipt (in the absence of evidence to the contrary, the date of receipt is the date of signing for the package for online transactions and the date of delivery for offline transactions)
- The Company complies with all applicable laws and regulations, including the E-commerce Law of the People's Republic of China, the Product Quality Law of the People's Republic of China, the Law of the People's Republic of China on the Protection of Consumers' Rights and Interests, and other laws and regulations. For offline transactions, if quality issues arise within 7 days, you can get in touch with the seller and present a valid proof of purchase or receipt to have the item returned or exchanged; if quality issues arise within 15 days, you can get in touch with the seller and present a valid proof of purchase to have the product returned or exchanged. For online purchases, you have 7 days from the delivery date to return the item

16

without giving a reason, but it must be in sellable condition. Customers are responsible for paying the shipping charges for returning items within 7 days, except for returns or exchanges due to quality issues

Other countries and regions

- If you acquire the goods outside of the Chinese Mainland, the warranty period, return policy, maintenance, and other relevant issues will be governed by the laws and regulations of that nation or region. If you have any special questions, please get in touch with your local authorized dealer

- When buying the item, kindly request and keep the receipt or other proof of purchase. The production date found on the product body will be used as the starting time of the warranty period if the receipt or other evidence of purchase is not presented to determine the warranty period
- Accompany Tech and Shenzhen Accompany Technology Co., Ltd. refer to the same company

Please be aware that the following situations will not qualify for free warranty services throughout the warranty period:

- Damage arising from inappropriate usages, such as improper access to the power supply, improper use of accessories, failure to follow instructions, improper use or damage caused by negligence, etc.; unapproved repairs or modifications; or force majeure (such as natural disasters, and voltage abnormalities)
- The device's typical use-induced deterioration and aging that has no impact on how it performs
- Free servicing does not include consumable components and add-ons, such as the adapter, npr does it include delivery charges

17