


AMIRO

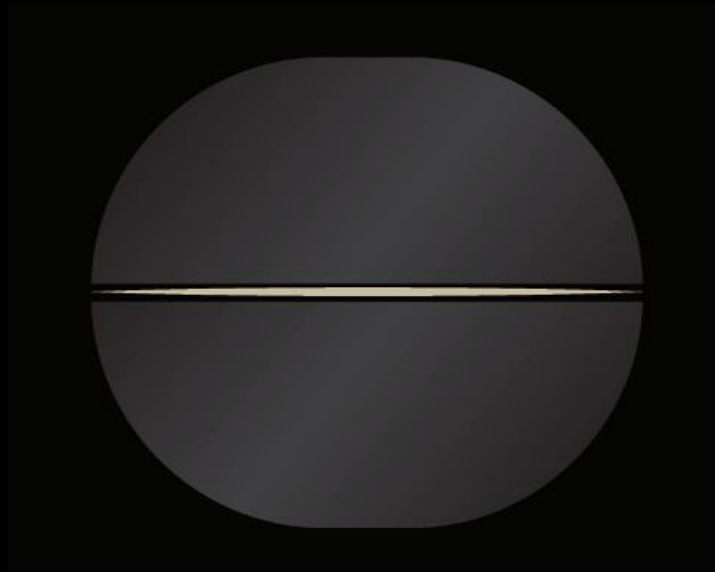

Amiro Essence Boosting Beauty Device User Manual

We appreciate your purchase of an AMIRO product.

Please read this section carefully and closely follow the instructions in order to ensure correct and safe use of the product.

Please save this manual for future reference.

Please include this manual when transferring ownership of the product to others.



Instructions for Use:

Dear customer, thank you for choosing Amiro Essence Boosting Beauty Device by Shenzhen Accompany Technology Co., Ltd. (“our company”). Before using the device, please ensure that you have carefully read and fully understood the user manual. To use the device properly, make sure the instructions are followed carefully. If you have any questions while using the device, feel free to reach out to our company for assistance. We promise to offer you the best of our service at all times.

Safety Warnings:

Our company guarantees the safety, reliability, and functionality of the device only when all specified operations are strictly adhered to. Our company is not liable for any other circumstances.

- Device adjustment and maintenance should be carried out by designated and trained personnel of our company.
- The electrical equipment and work environment must meet national standards, industry norms, and the manual requirements.
- Strictly adhere to the instructions for the use, charging, and storage of the device.
- Use only the accessories provided by our company

Contents

| | |
|---|----|
| Device and Manufacturer Information..... | 1 |
| Safety Instructions..... | 1 |
| People who should not use the device..... | 6 |
| Applicable scenarios..... | 8 |
| Operating principle..... | 8 |
| Packaging List..... | 9 |
| Product Illustration..... | 10 |
| Main Technical Parameters..... | 11 |
| Product parameters..... | 12 |
| Charging..... | 13 |
| Preparation before first-time use..... | 14 |
| Usage..... | 16 |
| To switch off the device:..... | 21 |
| Device maintenance tips:..... | 21 |
| Post-treatment skincare tips:..... | 21 |
| Product pairing..... | 22 |
| Product maintenance..... | 22 |
| Product troubleshooting..... | 23 |
| Other Specifications..... | 25 |
| Explanation of graph and symbol..... | 25 |
| Waste disposal declaration..... | 26 |

Device and Manufacturer Information

Product Information

Name: Amiro Essence Boosting Beauty Device

Model / Specification ABS001

Production Date: Check the nameplate for details

Manufacturing Information

Manufacturer: Shenzhen Accompany Technology Co., Ltd.

Address: F26, High-Tech Zone Union Tower, 63 South 10th Avenue, High-Tech Zone, Binhai Community, Yuehai sub-District, Nanshan District, Shenzhen, China

Postal Code: 518063



Contact number: 400 048 0888

After-sales Service Provider: Shenzhen Accompany Technology Co., Ltd.

Address of the After-sales Service Provider: F26, High-Tech Zone Union Tower, 63 South 10th Avenue, High-Tech Zone, Binhai Community, Yuehai sub-District, Nanshan District, Shenzhen, China

Safety Instructions

To use this device safely and to prevent bodily harm to the user or others, please read the safety precautions carefully and follow the precautions outlined.

| | |
|---|---|
|  | Warning: Operations that may result in bodily injury to the user or damage to the device. |
|  | Safety Instructions: Additional instructions for better usage of the device |

⚠Warning: Use, Charging and Storage Safety

To prevent injury or discomfort, and to ensure the device's life span, please ensure that the environment where the product is used and stored is dry, free from flammable or explosive materials, not exposed to strong direct light, and kept at an appropriate indoor temperature.

- Keep away from fire sources or flammable materials during use, charging, or storage.
- Avoid proximity to strong electromagnetic fields, shortwaves, or microwaves to prevent interference and unstable energy output.
- Refrain from using, charging, or storing device in damp or immersion-prone areas to prevent electric leakage or malfunction.
- Avoid using, charging, or storing device in environments with high levels of corrosive substances (such as hydrogen sulfide, sulfur dioxide, ammonia, nitrogen, hydrochloric acid, sulfuric acid, etc.) to prevent corrosion of the main unit and accessories
- Avoid using, charging, or storing the device in areas with significant temperature fluctuations (near heaters, hot pots, air conditioning, etc.) to prevent potential risks due to varying temperature.
- To prevent dropping or damaging the device, avoid using, charging, or storing it in tilted, unstable, or vibrating conditions.
- To protect the environment and prevent skin or body discomfort, please prevent the device from coming into contact with pollutants or waste during use.

⚠Warning: Fire Risks

- Keep away from fire sources or flammable materials during use, charging, or storage.
- Do not use the device near flammable and explosive materials (such as alcohol, hair gel, gasoline, nail polish, nail polish remover, high-pressure sprays).
- Do not use this device together with flammable beauty or hair devices, such as essential oils, hair gel, hair wax, foam mousse, and similar items.
- Avoid using, charging, or storing the device in high-temperature areas like near open flames, heaters, or direct sunlight to prevent overheating, explosions, or fires.
- Avoid direct contact of the electrode head with flammable materials to prevent fire risks caused by short circuits.
- Avoid using, charging, or storing the device in oxygen-rich environments to prevent fire hazards and other dangers.
- Remember to turn off the power and unplug the charger when the device is idle or during cleaning to prevent short circuits, electrocution, or fires.
- To avoid short circuits, electrocution, or fires, use the power adapter and its accompanying charging cable provided by our company for charging. When using the device, plug the adapter directly into the power outlet instead of using an extension cord to prevent short circuits, electrocution, or fires.
- For your safety and to prevent electrocution, fires, or injuries, do not attempt to disassemble, repair, or modify the device on your own. Please contact customer service for any necessary repairs.

⚠Warning: Biological Hazards and Contact Safety

·The components of this device that come into contact with the human body are crafted using 316 stainless steel and have been approved through nickel release (GB/T 19719-2005 & GB/T 28485-2012) and biocompatibility tests (GB/T 16886.5-2017/GB/T 16886.10-2017). This guarantees that they are safe for use. If you have a history of allergies to these materials, please avoid direct contact with the device or stop using it if you experience any allergic reactions. Consult a doctor if necessary.

·This device is for personal use on the face only. Please refrain from sharing it with others. If the device is to be used on multiple people, please clean the electrode head before each use to prevent cross-contamination.

·Due to the presence of bacteria on facial skin, it is important to clean the electrode head thoroughly before and after each use to prevent bacterial growth.

·For hygiene purposes and to avoid cross-infection of bacteria or viruses, do not use this device on non-human animals, such as pets.

⚠Warning: Hazards of Water Damage


This device has passed IPX6 waterproof testing.

·Prevent liquid from entering the main unit or charging port. After cleaning, ensure the device is completely dry before storing it in a dry environment to avoid potential short circuits, electrocution, or fires caused by water damage.


·To prevent electrocution, short circuits, or fires, refrain from using, charging, or storing the device near water sources, in bathrooms, or while bathing.

·Do not immerse any part of the device in liquid to prevent electrocution, short circuits, or fires.

·The device's bottom charging port features a waterproof plug design. When cleaning, make sure the waterproof plug is securely in place to prevent liquid from entering the charging port and potentially causing electrocution, short circuits, or fires.

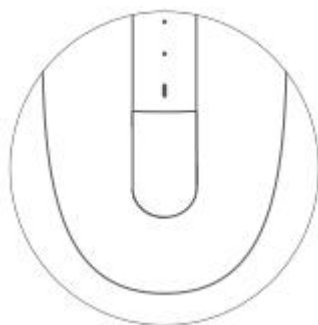


Waterproof Plug Properly Installed

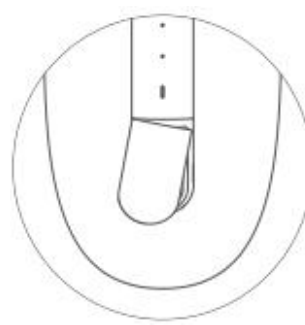

Waterproof Plug Not Properly Installed

⚠ Warning: Hazards Associated with Incorrect Usage and Wrong Settings

- Do not use this device on the eye area to avoid potential harm from excessive energy.
- Do not keep the device in the same position for an extended period to avoid discomfort and injury from accumulated heat.

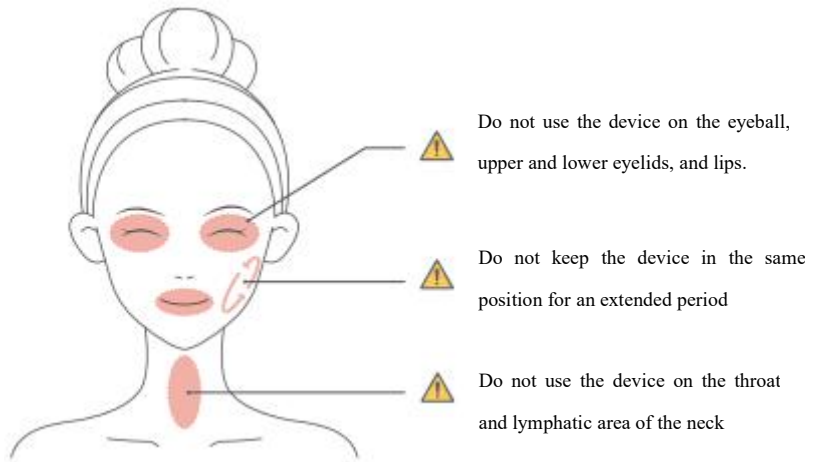

⚠ Warning: Operational Safety

Before using the device, ensure all functions and components are intact for effectiveness and safety.

- Do not use the device for non-cosmetic purposes.
- Do not use it for pet care to prevent injuries to your pet.
- To prevent physical discomfort and skin damage, avoid overusing it or similar devices.
- Remember to switch off the power after using the device or when not in use.
- Use the device after you have removed your makeup.
- Avoid looking directly at the LED light for a prolonged period of time to prevent damage to the eyes.
- Avoid using the device on the eyes to prevent eye injury.
- Feeling different temperatures while using the device is normal, as individual skin sensitivity differs. Adjust the temperature to your preference as necessary.
- To prevent skin damage or discomfort, it is recommended to avoid prolonged use on the same area.
- To prevent personal injury or accidents, it is important to follow the instructions in the user manual and refrain from unauthorized operations.
- Avoid using the device near strong magnetic or high-frequency devices such as microwave ovens.
- To prevent device malfunctions and personal injury, ensure that the electrode pads do not touch metal while in use or during charging.
- To prevent component failures or injuries from dropping the device, make sure to keep the main unit and base secure and stable to avoid any impact or falls. Do not use the device if it has been damaged by a fall or collision.

- Do not use the device if it is malfunctioning as doing so may result in bodily injury or danger.
- To avoid overheating hazards, do not cover the main unit with materials that can hinder heat dissipation, like fabric or blankets.
- Keep the device away from children as it contains small parts that could be swallowed or cause choking.
- Properly dispose of plastic packaging bags after unpacking to prevent risks of children suffocating or pets playing with them.
- When using, charging, or storing the device, make sure to keep it out of reach of infants, toddlers, and individuals with mental disorders to prevent any potential injuries.

Warning: Electrical Safety

- Do not disassemble, repair, modify, or damage the main unit or charging cable. Do not bend, pull, or twist the power cord.
- Do not expose the main unit or charging cable to water or high humidity to prevent electrocution or leakage. Avoid charging, using, or storing the device in areas prone to moisture, such as bathrooms or sinks.
- Use only the original charging cables provided with the device.
- This device does not come with an adapter. When charging, make sure to use an adapter that complies with national certification standards and provides a 5V=3A output.
- Do not connect the device to unstable power sources to prevent inconsistent charging efficiency or damage to the device.
- Do not connect this device in parallel or series with other electrical appliances. When the rated capacity is exceeded due to improper usage such as connecting multiple devices in parallel, it may lead to a fire hazard.
- If there is a power outage during use, stop charging the device, disconnect it from the charging cable, and store the main unit properly.
- Only use the charging cable provided by our company. Using other cables may cause damage to the main unit.
- Avoid using improper adapters to prevent safety risks like overheating, electrocution, and leakage.

Warning: Risk of Electrocution

- Do not use the device if you are wearing contact lenses or metal accessories such as glasses, earrings, necklaces, watches, or rings.
- To avoid risking electrocution, do not attempt to open the main unit's casing by yourself.
- Do not handle the charger with wet hands. Ensure that the device, charging cable, and adapter are dry at all times.
- Do not use the charging cable if it is loose or damaged when inserting it into the adapter.
- After use, unplug the charging cable and hold it properly when pulling it out. To

avoid the risk of electrocution, do not plug or unplug the charging cable with wet hands After removing the charging cable, ensure that the waterproof plug is secured in place.

- Avoid using improper adapters to prevent safety risks like overheating, electrocution, and leakage.
- If the charging cable is damaged, contact the manufacturer or dealer for a replacement.

▲Caution: Improper storage may reduce the lifespan of the device.

Ensure to store this device in a safe location and protect the electrode heads properly. If the device hasn't been used for an extended period, ensure to verify that it starts up correctly.

- Keep the main unit and charging cable out of reach of children, infants, or pets.
- Avoid putting excessive weight on the main unit or pulling on the charging cable.
- Do not place the device in a microwave or high-pressure container.
- Disconnect the charging cable after using the device or while cleaning the main unit, during power outages, lightning, or similar situations to avoid any potential damage.
- Check the electrode heads for cracks, deformities, or damage before using the device. Do not use the device if any problems are found.
- Do not submerge the device in liquid-filled containers.
- Do not clean the device with solvents like alcohol, diluents, gasoline, nail polish remover, diesel, makeup remover, etc.

▲Caution: Transportation Safety

- When transporting the device, avoid shaking or placing heavy objects on top of it. Otherwise, it may lead to damage or malfunction of the device.
- When transporting the device, use proper packaging and handle it carefully to prevent the main unit from falling or experiencing heavy impacts that may cause damage.
- Avoid exposing the device to rain or submerging it into water during transportation to prevent the risk of short circuits or electrocution.

People who should not use the device

People with any of the following conditions may not use this device as it could result in accidents or skin and bodily injury:

- Individuals with abnormal blood pressure, heart disease, cardiac functional disorders, or tendencies towards functional disorders.
- Infectious diseases or fever
- Hemorrhagic disease or cicatricial diathesis
- Individuals with abnormal pain perception, inability to express consciousness, or lack of temperature and stimuli perception.
- If you are undergoing dental or oral treatment, consult a doctor before using the

device

- Individuals experiencing menstruation (which can lead to hormonal instability and significant skin changes), pregnancy, or breastfeeding.
- Infants, young children, or individuals under 18 years old (in adolescence).
- Presence of metal implants in the targeted area.
- Individuals using implantable electronic devices such as pacemakers.
- Individuals using life-sustaining electronic devices such as heart-lung machines.
- Individuals using wearable medical devices such as electrocardiographs, infusion pumps.
- Suffering from, or having previously suffered from, cancer (especially skin cancer and precancerous moles), acute illnesses, tuberculosis, febrile diseases, infectious diseases, blood diseases, visceral diseases, skin diseases, or allergic dermatitis.
- Impaired immune system due to immunosuppressive diseases such as AIDS or the consumption of immunosuppressive drugs.
- Please consult a doctor before using the device if you have had permanent dermal fillers injected or have received temporary dermal fillers/Botox injections in the facial areas of application.
- Diseases caused by thermal irritation within the application area, such as recurrent herpes simplex.
- Individuals with allergies, sensitive skin, or fragile skin.
- Experiencing fatigue, illness, under the influence of alcohol or on sleeping pills.
- Dilated capillaries caused by prolonged use of steroids, hormones, and other drugs, or liver dysfunction.
- Individuals with skin diseases, genetic eczema, allergic dermatitis, scar-prone skin, or sensitive and fragile skin.
- Avoid using the device if you are using Isotretinoin, alpha hydroxy acids (AHAs), beta hydroxy acids (BHAs), or salicylic acid. Individuals using sedatives, local anesthetics, or anesthetic painkillers. Avoid using the device if you are taking aspirin medication, using immunosuppressants, or using retinoic acid and similar medications.

Consult a doctor before use if you have high myopia (above 600 degrees).

·Use of this product may cause discomfort or adverse effects if you suffer from allergies, skin conditions, allergic dermatitis, or sensitive skin. Stop using the product immediately if you experience any adverse effects.

Do not use the device on areas of skin that are injured or covered with burns, blisters, cuts, wounds, active skin diseases, and/or infections.

⚠ Warning: Failure to follow the instructions for use may result in adverse effects;

⚠ Warning: The device could cause damage to the skin if it is used repeatedly or for a long period of time on the same area.

·Avoid bypassing or removing the built-in safety mechanism.

Applicable scenarios

·Target Consumers Recommended for individuals aged 18 and above seeking to diminish facial wrinkles, enhance skin tone, and address skin repair needs (please check contraindications before use).

·Applicable Area: Apply to the face, avoiding the eye area.

Operating principle

Micro-current

Microcurrent therapy can boost cellular ATP production, stimulate the synthesis of collagen types I, IV, and XVII, and elastic fibers. It also utilizes medium and low-frequency pulse currents to induce muscle contractions, elevating the skin and tightening the tissues and muscles.

Light Therapy

Photodynamics is the use of a light-emitting device to convert electrical energy into controlled light energy of specific wavelengths for cosmetic and beautifying purposes. Utilizes visible light to stimulate tissues through photothermal, photochemical, and biological effects, enhancing skin cell metabolism, reducing cell aging, increasing skin elasticity, and revitalizing the skin.

·Red light: Enhances cell activity, stimulates collagen regeneration, diminishes fine lines, shrinks pores, increases skin elasticity, and revitalizes dull skin.

·Yellow light: Regulates skin pigmentation, fixes hyperpigmentation, evens out skin tone, and enhances skin brightness and smoothness.

·Infrared light: Accelerates tissue repair and regeneration, relieves skin inflammation, promotes tissue recovery and healing, improves skin elasticity and tightness, reducing the appearance of fine lines and wrinkles.

Thermotherapy:


Utilizes infrared radiation and high-energy heating plates to circulate heat within the skin, boosting blood flow, activating skin cells, and enhancing pore function.

Iontophoresis


Generates an electric field on the stratum corneum, the outermost layer of the epidermis, to facilitate electroporation, improving skin permeability for increased absorption of active ingredients from skincare products into deeper skin layers, maximizing skincare benefits.

Packaging List

- ① Amiro Essence
Boosting Beauty
Device × 1

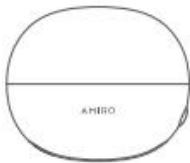

- ② Charging
Cable × 1

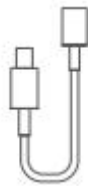

- ③ Quick Start
Guide Card × 1



- ④ User Manual × 1




- ⑤ After-sales
Service Card × 1



Product Illustration

Device illustration

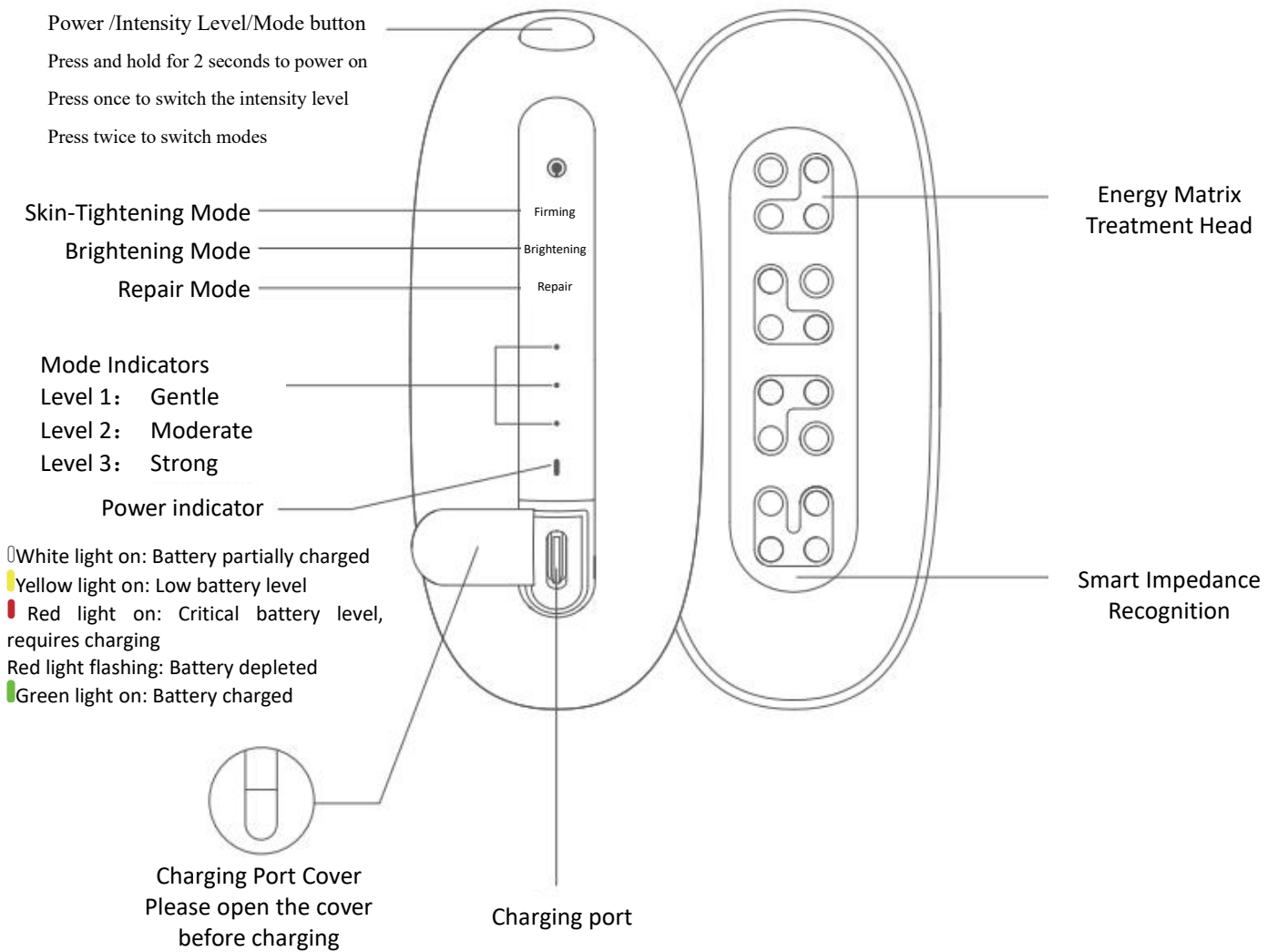

Main Technical Parameters

Normal working conditions

Input: Type-C 5 V = 3 A

*It is advisable to use a charger that complies with the national regulations for charging

Product parameters


| | |
|----------------------------|--|
| Product name | Amiro Essence Boosting Beauty Device |
| Product Model | ABS001 |
| Product mode | Skin-Tightening Mode / Brightening Mode / Repair Mode |
| Device weight | ≈ 170 g |
| Device dimensions | 90.6 mm*38.5 mm*75.6 mm |
| Red light wavelength | 633 nm |
| Yellow light wavelength | 590 nm |
| Infrared light wavelength | 850 nm |
| Battery capacity | 1300 mAh |
| Electrode Size | 135.96 mm ² |
| Maximum Treatment Duration | 5 min |
| Compliance Standards | GB4706.1-2005, GB4706.15-2008 GB4343.1-2018, GB/T 4343.2-2020 |

Charging

Verify that the device is fully charged during the first use

Charging time: Approximately 60 minutes (charging time may vary due to different charging conditions)

Battery life: 50 minutes (may vary depending on usage environment and intensity level)

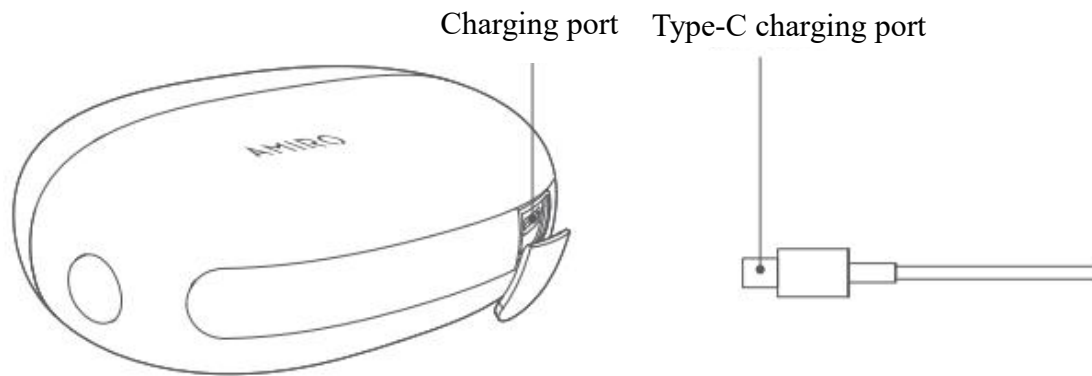

Preparation before first-time use

Checking

Open-box checking: Unpack the main unit and accessories from the box, verify the contents against the packaging list, and retain the packaging materials for future use.

Check the device and accessories for any damage. If there are signs of mechanical damage or equipment malfunction, do not use the device and contact the distributor or after-sales service immediately.

Charging


Make sure the device is fully charged before using it for the first time.

Skin adaptation test

This device treats the skin with heat, microcurrents, and light. To prevent skin discomfort, conduct a skin sensitivity test before using it for the first time. Follow the illustrated steps below:


·Cleanse the inner arm with soap, body wash, or facial cleanser. Dry the area using a towel or tissue. Then, an appropriate amount of gel in the center of the electrode head. Turn on the device and select Skin-Tightening mode 1. Gently glide the device over the cleaned skin area for 2 minutes.

·After 2 minutes, gently massage any leftover gel until it is fully absorbed. If there are no signs of redness, itching, or heat on the skin within 24 hours, you can proceed with using the device.

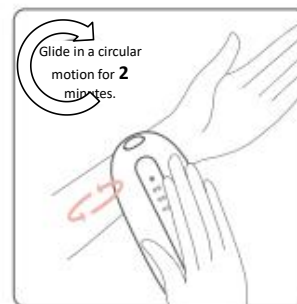

⚠ Warnings

- If you notice any abnormal conditions such as rashes, redness, or itching on the test area, please discontinue the test immediately.
- If you experience any skin conditions before, during, or after the test, please discontinue the use of this device.
- Avoid using the device on the same spot for excessive periods of time and remember to glide it over your skin during use.
- Make sure to glide the electrode head smoothly if the essence, gel, or other commonly used device dries out during treatment.
- If you experience a strong tingling sensation on your skin, reduce the intensity level or switch to a different area.
- Do not apply the device to the sensitive areas indicated in the illustration below.

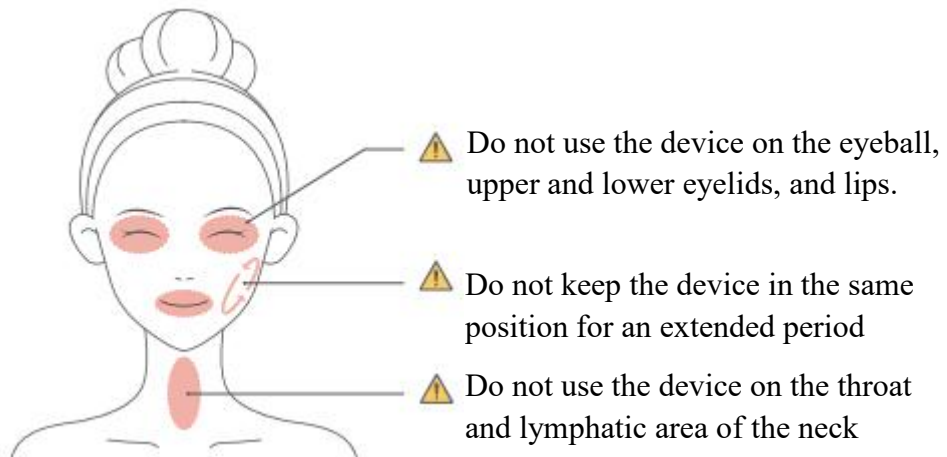

Usage

Device Instructions


1. Remove makeup and cleanse your face.

Before use, remove all makeup and cleanse the facial, eye, jawline areas. Make sure to also remove contact lenses and any metal jewelry and from the face, ears, neck, hands, etc.




2. Apply a sufficient amount of essence to your face or insert the essence cartridge into the device.

1. Apply an appropriate amount of essence or gel, then switch on the device to start your skincare routine. Reapply the essence or gel immediately if your skin feels dry or resistant during the therapy.

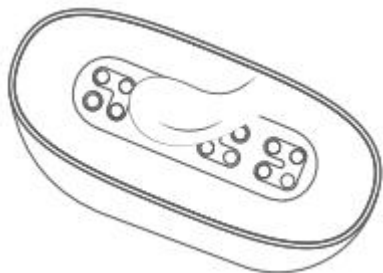


2. Make sure the essence cartridge is attached securely before turning on the device.

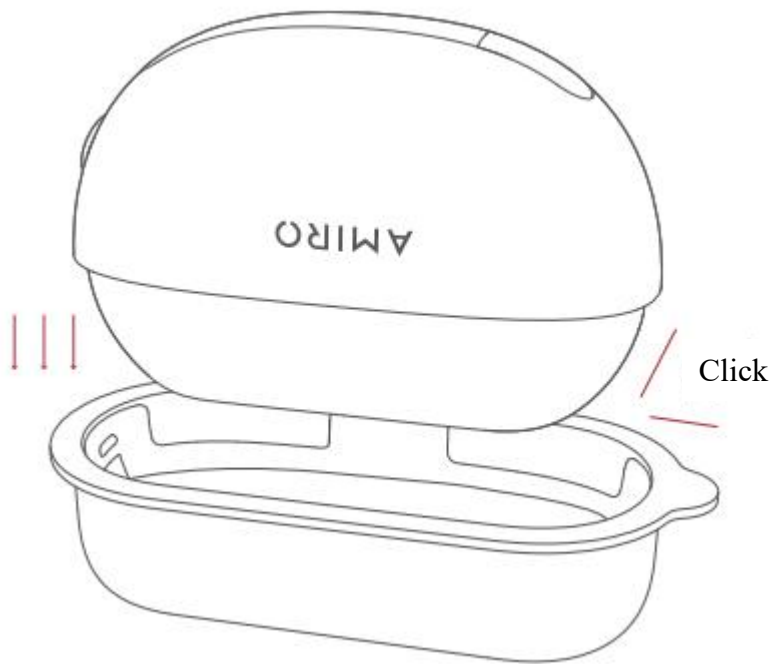

3. To switch on the device, press and hold the power button until you feel a prolonged vibration.

After switching on the device, you may monitor the battery level by checking the battery indicator.

Regardless of whether the device is on or off, simply plug the charging cable into the device's charging port to connect it to the power source and begin charging.

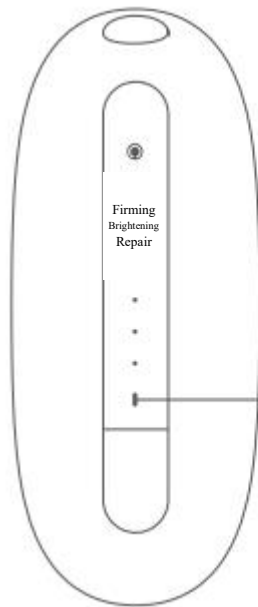

The indicator light will turn on when the device is connected to the power source.

■ Yellow light flashing - Charging

■ Green light on – Battery charged

Status of Battery Indicator

| Not charging | Charging |
|--|---|
| Red light on: Critical battery level, please charge before use | Red light flashing: Charging, critical battery level |
| Yellow light on: Low battery level, please charge after use | Yellow light flashing: Charging, low battery level |
| White light on: Battery partly charged | White light flashing: Charging, battery partly charged |
| Green light on: Battery charged | Green light on: Battery fully charged, unplug the charging cable and secure the waterproof charging plug cover. |

4. Modes and Levels

Press twice to select the desired mode


Select the mode that matches the main effect of the essence you are using. For instance, if you are using an essence for skin tightening, select the “Skin-Tightening mode”.

When you select a mode, the indicator light will switch to white, and the electrode panel will light up, indicating that the treatment is ready to begin. The default mode is “Skin-Tightening mode”.

When the device is on, press the power button twice to cycle through the modes in the following order:

Skin-Tightening Mode → Brightening Mode → Repair Mode

After each mode, the device will switch off after long vibration. To switch to another mode, restart the device and make a new selection.

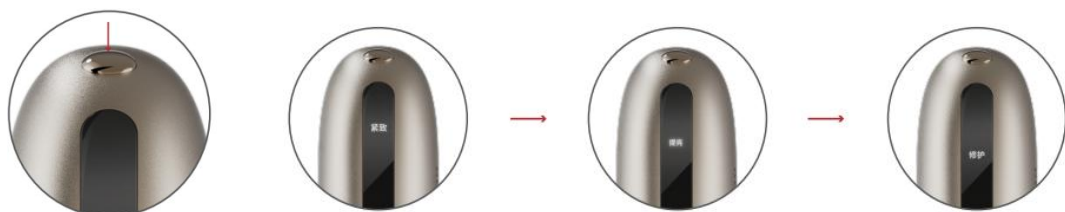

Press the power button twice to change modes while the device is on,

Skin-Tightening Mode

Brightening Mode

Repair Mode

Select intensity level

To set an intensity level, press the gold button to switch intensity level. The indicator lights will light up to display the intensity level. While the device is on, press the button once to cycle through the in tensity levels in this order:

Level 1→Level 2→Level 3→Level 2→Level 1



Press the button once to switch modes while the device is on.



Level 1



Level 2

Level 3

Level 2

Level 1

5. End of Treatment

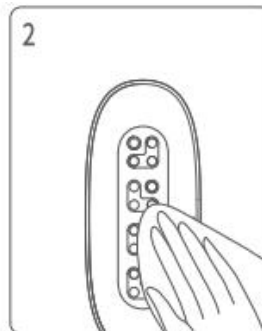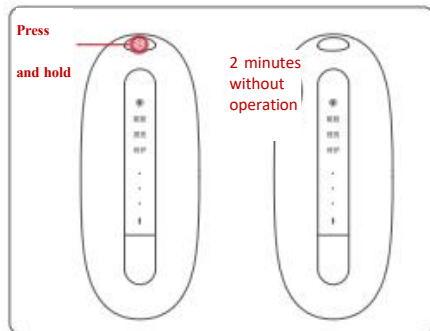

To switch off the device:

Automatic power off

- The Skin-Tightening mode will be switched off automatically after 5 minutes of treatment.
- The Brightening Mode will be switched off automatically after 5 minutes of treatment.
- The Repair Mode will be switched off automatically after 5 minutes of treatment.

The device will be switched off automatically if it is not in use for a certain amount of time.

- The device will be switched off automatically after 120 seconds of inactivity.

To manually switch the device off,

- press and hold the power button for 2 seconds.

Device maintenance tips:

After powering off, wipe the treatment head with a dry cloth, damp tissue, etc.

*Note: Do not clean the device with any corrosive liquids such as alcohol, alkaline detergents, or nail polish remover.

Post-treatment skincare tips:

After the treatment, you can apply skincare devices or sunscreen without needing to wash your face.

Safety Instructions: Increase device usage in progressive stages to build tolerance. Avoid harsh skincare products, such as acids, retinol creams, or whitening products, at the beginning to prevent skin irritation. Protect your skin from sun damage and premature aging by consistently applying sunscreen and moisturizer during treatment.

Product pairing

Recommended Essence:

AMIRO's essence cartridge, gel, serums, etc.

Essence not recommended to use:

The following types of skincare products are not suitable for use with the device due to their low electrical conductivity and viscous texture.

| Types | Examples |
|--|---|
| Facial oils | Oil-based serums, Essential oils |
| Acid-based facial essence | Containing salicylic acid, mandelic acid, alpha hydroxy acids, etc. |
| Facial essence with thick, viscous texture | Thick, viscous facial essence |

Product maintenance

- To prevent electrocution, keep the device dry and avoid immersing it or any other parts in water.
- Store the device and its accessories properly after each treatment to prevent damage or loss. Remember to clean the device before putting it away.
- When the device is not in use for an extended period of time, store it in a dry, shaded, and cool place to prevent damage. Avoid contact with hard objects.
- Do not store the device and its accessories in a damp environment for extended periods of time as it can affect the device's performance.
- Before using the device after a long break, make sure to charge it to full, check its functionality, and start at the lowest setting to gradually build tolerance.

Product troubleshooting







| FAQ | Possible cause | Recommended Solutions |
|------------------------------|--|--|
| Unable to start the device | Battery depleted | Use the device after it has been charged |
| | End of battery life | Contact the manufacturer for repair or replacement |
| | Incorrect operation | Press and hold the power button until the device vibrates briefly and the treatment head panel light is on |
| Automatic shutdown | Not in use for a specific duration | Restart the device |
| | Dead battery | Charge to full before use |
| Device suddenly malfunctions | Not being used for an extended period of time | Device will be switched off automatically after 120 seconds of inactivity |
| | Battery depleted | Charge the device fully or use it while it is connected to the power source |
| | Device malfunction | Please contact the store where you purchased the device or call the customer service hotline at 400-048-0888 |
| Device overheating | Using the device for too long | Turn off the device and allow it to rest for at least half an hour before using it again |
| Heated sensation on the face | Selecting the wrong intensity level | Switch to a gentler level. Start with Level 1 to build tolerance |
| | Using the device on the same spot repeatedly | Avoid leaving the device on the same spot for an extended period of time |
| | Using the device in combination with compatible products | Do not use devices with similar functions simultaneously |

| FAQ | Possible cause | Recommended Solutions |
|--|--|---|
| Feeling a tingling sensation on the skin | Device not touching the skin correctly | Press the device against the skin and glide it evenly across |
| | Choosing the wrong level and mode | Switch to a gentler level. Start with Level 1 to build tolerance |
| Sensation is too mild | The intensity selected is too low | Switch to higher intensity levels |
| | Insufficient amount of gel/essence | Use more gel/essence to allow the device to glide across skin more easily |
| | The essence cartridge has dried up | Stop using device |
| | The conductivity of the skin care products used with the device is too low | Contact customer service to inquire about the conductivities of typical skincare devices |
| | Gliding too fast | Decrease the gliding speed and make sure the device glides smoothly and evenly. |
| | The skin has good tolerance | It is normal for different parts of the face to have varying sensitivity to the current |
| Skin discomfort | Redness | If your skin has a temporary reaction to the warmth of the device, the condition should improve soon. Seek medical attention immediately if it doesn't improve or if you experience other discomfort. |
| | Skin dryness | Apply sunscreen and use hydrating skincare products |

Other Specifications

- 1). The device has a waterproof rating of IPX6.
- 2) Voltage range: 100-240 V~, 50 Hz / 60 Hz, input power: Type-C 5 V = 3 A, 15 W

Explanation of graph and symbol

| Symbol | Explanation | Symbol | Explanation |
|---|---------------------------------|---|----------------|
|  | Warning |  | Keep dry |
|  | Power/Standby |  | Direct current |
|  | Fragile items, handle with care |  | Up |

Waste disposal declaration

Disposal of used battery:

- Do not attempt to replace the battery on your own while the device is operating normally.
- Risks of battery leakage Do not come into direct contact with battery liquid since it is highly corrosive. In the event that you come into contact with battery liquid, please wipe it immediately with a dry cloth, and seek medical attention if you experience any discomfort.
- Make sure to remove the battery before disposing of the device.
 - ① Remove the product casing
 - ② Remove the internal wires and take out the battery
- Remember to turn off and unplug the device before removing the battery.
- Dispose of the battery properly and safely.