
AMIRO

AMIRO Electrical Stimulation Device For
Skin Lifting In Skincare

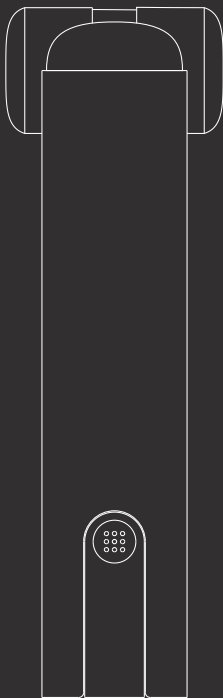USER MANUAL

AMIRO

AMIRO Electrical Stimulation Device For Skin Lifting In Skincare USER MANUAL

*Thank you for choosing AMIRO.

Illustrations in the user manual of the product,
accessories, and user interface are for reference only.

For your safety, please read this manual carefully before use.

Please keep this manual for future reference.

Please include this manual when transferring the product to others.

Product and manufacturer information

PRODUCT INFORMATION

Product Name: AMIRO Electrical Stimulation Device For Skin Lifting In Skincare

Model/Specification: ABS004

product manufacturing information

Manufacturer: Shenzhen Zongjiang Technology Co., Ltd

Address: 25th Floor, High tech Zone Joint Headquarters Building, No.63 Gaoxin
South 10th Road, Binhai Community, Yuehai Street, Nanshan District, Shenzhen
Postal Code: 518063

Contact information: 400 048 0888

After sales service unit: Shenzhen Zongjiang Technology Co., Ltd

After sales service unit address: 25th floor, High tech Zone Joint Headquarters
Building, No. 63 Gaoxin South 10th Road, Binhai Community, Yuehai Street,
Nanshan District, Shenzhen

Safe usage instructions

Warning: Instructions that may cause harm to the user or damage to the instrument

Attention: Please provide additional instructions for better use of the product

Warning: Risks arising from the environment

To prevent injury or discomfort and ensure the longevity of the instrument, please ensure that the environment is dry, free of flammable and explosive materials, without direct sunlight, and at a suitable temperature when using and storing the product.

- Do not use, charge, or store in areas near sources of fire or flammable materials.
- Do not use, charge, or store in areas close to strong radio waves, strong magnetic fields, short waves, or microwaves to avoid equipment output interference and unstable energy causing danger.
- Do not use, charge, or store in damp or hazardous environments to avoid electrical leakage or malfunction.
- Do not use, charge or store in corrosive environments (such as environments with high levels of hydrogen sulfide, sulfur dioxide, ammonia, nitrogen, hydrochloric acid, sulfuric acid, etc.) to avoid corrosion of the host and accessories.
- Do not use, charge or store in areas with large temperature changes (such as heating, hot pot, air conditioning, etc.) to avoid risks caused by individual differences in temperature perception.
- Do not use, charge, or store in a tilted, unstable, or vibrating environment to prevent the product from falling or being impacted.
- When using this product, please avoid pollutants and waste to prevent environmental pollution from causing discomfort to the skin or body.


Warning: Risks associated with contact

- This product comes into direct contact with the skin and is intended for personal use. Please do not share it with others. If multiple people need to use or experience the product, please clean the electrode head area before use to avoid cross infection.
- This product is intended for facial use, and the facial skin has a rich microbiota. Please clean the electrode tip promptly before or after use to avoid bacterial growth.
- Do not use it on non-human bodies (such as pets), as it may cause cross infection of bacteria/viruses.


Warning: Risk of liquid ingress protection

- Do not allow any liquid to enter the host and charging port. After cleaning, please dry the equipment and store it in a dry environment to avoid short circuits, electric shock, and fire caused by liquid ingress.
- Do not use, charge or store near water sources, in bathrooms or while bathing, as it may cause electric shock, short circuit, or fire.
- Do not immerse any part of the instrument in liquid, as it may cause electric shock, short circuit, or fire.

The bottom charging port of this product is designed with a waterproof plug. When cleaning, please ensure that the waterproof plug is installed in place to prevent liquid from entering the charging port, otherwise it may cause electric shock, short circuit, and fire.

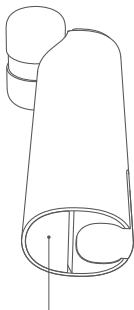✓ Waterproof plug
installed correctly

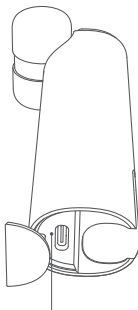✗ Waterproof plug not installed

⚠ Warning: Risks caused by incorrect use of regions and parameter settings

· Do not use this product on the eyeball, upper and lower eyelids, and lip areas to avoid human injury caused by excessive energy.

Please do not use this product in the areas of the eyeball, upper and lower eyelids, and lips

Please do not use this product in the areas of the eyeball, upper and lower eyelids, and lips

Warning: Risks caused by incorrect operation

To ensure effectiveness and safety, please ensure that all functions and components of this product are intact before use.

- Do not use for non cosmetic purposes.
 - Do not use for pet care to avoid injury.
 - Do not overuse this product or similar products to avoid physical discomfort and skin damage. Please turn off the power promptly after using the product or when not in use after turning it on.
 - Please use this product after removing makeup.
 - Do not stare directly at LED lights for a long time, as it may cause eye damage.
 - Do not use the guide head around the eyelids and eyeballs, as it may cause eye injury.
 - When using a beauty device, it is normal to feel different temperatures due to individual differences in skin perception. Please adjust the appropriate level according to your personal situation.
 - Do not use it continuously on the same area for a long time to avoid skin damage or discomfort.
 - Do not bypass or cancel the inherent safety devices of the equipment.
- Failure to follow the instructions may result in adverse effects.

Do not use it for purposes other than this user manual or perform the following operations, otherwise it may cause personal injury or accidents.

- Do not use this product with strong magnetic or high-frequency devices (such as microwave ovens, etc.).
- Do not contact the electrode pads of this product with metal during use or charging, as it may cause product malfunction and personal injury.
- Please avoid dropping or colliding the host and base, otherwise it may cause mechanical failure or injury due to the host falling. Do not use if the instrument has been damaged due to falling or collision.
- Do not use the machine in a malfunctioning state, as it may have adverse effects on the body and pose a danger.
- Do not wrap the host with items that affect heat dissipation, such as fabrics, blankets, etc., to avoid danger caused by overheating.
- This product contains small components. Please keep it out of reach of children to avoid the risk of ingestion or suffocation.
- After disassembling the packaging, please handle the plastic packaging bag properly to avoid the risk of suffocation caused by children and pets playing.
- When using, charging, or storing this product, please be sure to prevent infants, young children, and mentally ill individuals from touching it to avoid injury.
- Do not bypass or cancel the inherent safety devices of the equipment.
- Failure to follow the instructions may result in adverse effects.

Warning: Risk of power supply

- Do not disassemble, repair, modify, or damage the host or charging cable by yourself. Do not bend, pull, or twist the power cord arbitrarily.
- Do not charge, use or store in high humidity places such as water, bathrooms, washbasins, etc; Prevent the host and charging cable from getting wet and causing electric shock and leakage.
- Do not use non original charging cables.
- This product does not come with an adapter. When charging, please choose an adapter that meets national certification standards and has a measured output power of no more than 15W.
- If the power supply is unstable, it may cause unstable charging efficiency or damage to the instrument.
- Do not connect this product in parallel or series with other electrical appliances. If multiple parallel connections are used beyond the rated range, it may cause a fire.
- If there is a power outage during use, please stop charging, unplug the charging cable, and properly place the host.
- Please use the charging cable provided with this product. Using a charging cable that does not meet the requirements may cause damage to the host or reduce charging efficiency.
- Do not use unauthorized adapters as they may pose safety risks such as overheating, electric shock, leakage, etc.

Warning: Risk of electric shock

- Do not use metal jewelry such as glasses and contact lenses, earrings, necklaces, watches, rings, etc.
- To avoid electric shock, do not open the casing of the host by yourself.
- Do not operate charging with wet hands. Please keep this product, charging cable, and adapter dry.
- Do not use if the charging cable plug is loose or damaged when inserted into the adapter.
- Do not plug or unplug the charging cable with wet hands, as it may cause electric shock. After unplugging the charging cable, please install the waterproof plug in place.
- Do not use unauthorized adapters as they may pose safety risks such as overheating, electric shock, leakage, etc.
- Please contact the manufacturer or distributor when the charging wire is damaged.

Attention: Risk of unreasonable storage and reduced service life

Reasonable storage is important for the normal use of this product, and it should be placed on a safe surface while paying attention to the protection of the electrode tip. When not in use for a long time, please check if the instrument can start normally before using it again.

- Do not store the host and charging cable within the reach of children, infants, or pets.
- Do not apply heavy pressure or pull the charging cable on the host.
- Do not place in microwave ovens or high-pressure containers.
- Do not place the product and its supporting components in the exposed area.
- Do not leave the charging cable for a long time. When cleaning the main unit after using the instrument, or in case of power outage, lightning, etc., be sure to disconnect the power supply.
- Before use, please confirm if the electrode tip is cracked, deformed, or damaged. If so, do not use.
- Do not immerse this product in containers containing liquids.
- Do not use solvents such as alcohol, diluents, gasoline, nail polish remover, diesel, makeup remover, etc. to wipe or soak the instrument.

Attention: Risks of transportation

- During transportation, avoid shaking or being covered by heavy objects. Otherwise, it may cause damage or malfunction to the host.
- During transportation, it should be properly packaged and handled with care to avoid damage to the main unit caused by falling or heavy impact.
- Do not expose the product to rain or soak it during transportation, as it may cause a short circuit or electric shock, posing a safety risk.

Contraindications and Scope of Application

Do not use this product if any of the following situations occur, as it may cause accidents or other issues:

- Do not use if you have thyroid disease or thyroid abnormalities, and avoid the thyroid area when using.
- People with abnormal blood pressure, heart disease, heart dysfunction, or tendencies towards dysfunction.
- People with infectious diseases or fever
- People with bleeding disorders or scar constitution.
- People with abnormal perception of pain, unable to express their consciousness, or unable to perceive temperature and stimulation.
- People who are undergoing dental or oral treatment should consult a doctor before use.
- During menstruation (hormonal instability, significant changes in skin condition), pregnancy or lactation.
- Infants and young children or individuals under the age of 18 who are currently in the developmental or adolescent stage.
- There is a metal implant in the area of use.
- People who are currently using implantable medical electronic devices such as pacemakers.
- People who are using electronic devices such as artificial heart and lung machines to sustain their lives.
- People who are using wearable medical electronic devices such as electrocardiographic recorders and infusion sets.
- Currently or previously suffering from the following diseases, such as cancer (especially skin cancer, precancerous nevi), acute diseases, tuberculosis, febrile diseases, infectious diseases, blood diseases, visceral diseases, skin diseases, allergic dermatitis.
- Immune system damage caused by immunosuppressive diseases, such as taking immunosuppressive drugs.
- If permanent plastic fillers have been implanted within the scope of use, or temporary subcutaneous fillers or botulinum toxin have been implanted within the scope of use, please consult a doctor's opinion first.
- Any disease caused by thermal stimulation within the scope of use, such as recurrent herpes simplex virus.
- Allergic constitution and skin during the allergic period.
- The body is in a state of excessive fatigue, feeling unwell, under the influence of alcohol, or taking sleeping pills.
- Capillary dilation caused by long-term use of steroids, hormones, and other medications, as well as liver dysfunction.
- People with skin diseases, hereditary dermatitis, allergic dermatitis, scar constitution, or sensitive and fragile skin.
- External use of isotretinoin, alpha hydroxy acid (AHA), salicylic acid (BHA), or rhododendron acid is prohibited. Do not use sedatives, local blockers, or anesthetic painkillers during the period of taking them. Prohibited during the use of aspirin, immunosuppressants, isotretinoin, and other medications.
- People with high myopia (over 600 degrees), please consult a doctor.
- Do not use this device in areas where skin integrity is damaged (such as burns, blisters, incisions, wounds, active skin diseases, and/or infections).

Attention: The above safety precautions are not listed in detail. If you are unsure whether you can use this product, please consult a doctor or contact our company before making a decision.

scope of application

- Target audience: Suitable for individuals aged 18 and above (please refer to the contraindications section)
- Applicable areas: human face (avoid eyeballs, upper and lower eyelids, lips, throat, and neck lymphatic areas)

Expected usage environment: home environment

Working principle

Photodynamic

Convert electrical energy into controlled light energy of specific wavelengths for cosmetic purposes through luminescent devices. Mainly utilizing visible light to undergo photothermal, photochemical, and/or biological stimulation with human tissues, increasing cellular metabolism, slowing down cellular aging, increasing skin elasticity, and achieving skin tenderness.

Red light: effectively enhances cell activity, promotes collagen regeneration, reduces fine wrinkles, shrinks pores, increases skin elasticity, and improves skin dullness.

Infrared light: Deep heating activates cell regeneration at the root.

EMS microcurrent import technology

simulates microbial currents to mimic human bioelectricity, activating the vitality of skin cells at the bottom layer. Positive and negative ions are introduced in both directions to promote penetration, and the combined ingredients reach the dermis layer directly.

The multi frequency respiratory electric pulse hybrid technology simulates the brain nerve electrical signal, mobilizes the muscle to move actively, realizes the facial promotion of skin care products and essence penetration, and promotes the absorption of essence; Relaxing muscles, physical massage, and soothing facial muscles.


Product composition



① AMIRO Skin Lifting Device x1

② User Manual x1

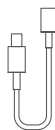③ Charging Cable x1

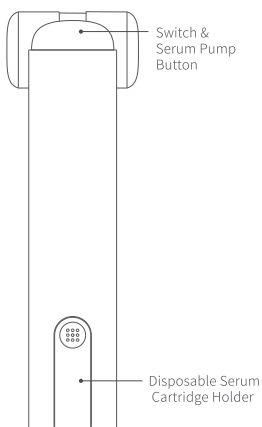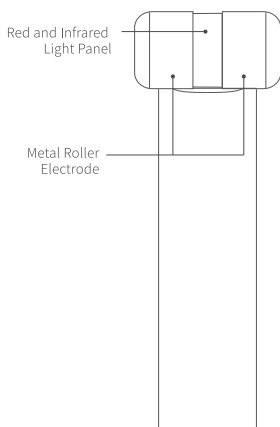Normal Operating

Ambient Temperature: +5°C~+35°C

Relative Humidity: 15%~80% RH



Atmospheric Pressure: 70 kPa - 106 kPa

Product Parameters

Material Composition	Roller Electrode	Stainless Steel
	Device Housing	Aluminium Alloy

Product Name	AMIRO Electrical Stimulation Device For Skin Lifting In Skincare
Product Model	ABS004
Net Weight	≈103g
Dimensions	≈126mm x 50mm x 40mm
Battery Capacity	350mAh
Electrical Impulse Frequency	100Hz~1000Hz
Red Light Parameters	625nm±10nm
Infrared Light Parameters	855nm±10nm


Product Buttons and Functions

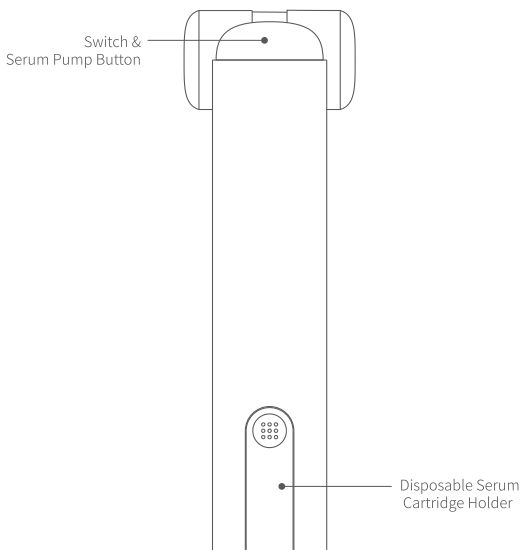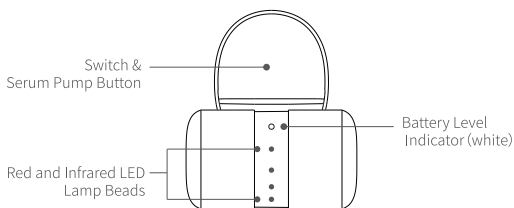Power On/Off

Press the Switch & Serum Pump Button– the red lights turns on, and the device is ready for treatment.

**Short-press to power on
press and hold to power off**

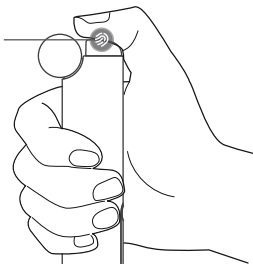When powered on, press and hold the Power/Dispense button for 2 seconds – the indicator light turns off and the device shuts down.


Powered on
(light illuminated)

Powered off
(light off)


Short-press the Power/Dispense button to dispense the serum.

**Press once to
dispense serum**


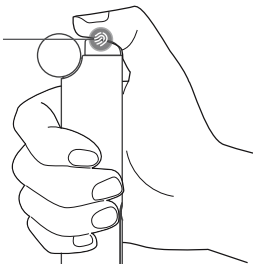
The use and replacement of secondary warehouses

Press the slot with your thumb and slide down to remove the serum cartridge holder..

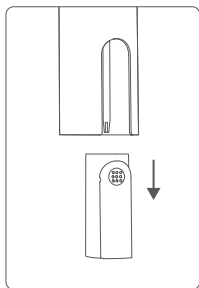Insert the twisted-open single-dose serum into the cartridge holder, align it with the catheter port, and insert it into the device.

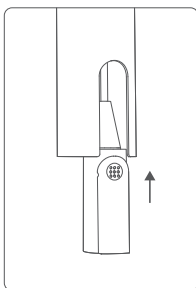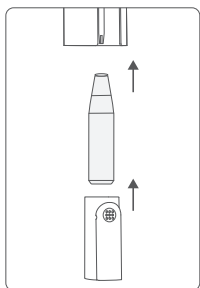
Usage

After loading the single-dose serum and powering on, briefly press the [Power/Serum Pump Button] to dispense the serum, then roll the device over your face.

Battery Level & Charging Methods

Checking battery level during operation:

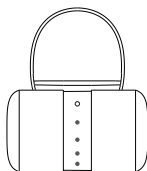Sufficient battery	Low battery (<20%)
○ White light stays on	○ White light flashes


*When battery drops below 10%, the device will auto power off.

Charging instructions:

Before initial use, fully charge the device.

Charging time: 40 mins (varies by environment).

Battery life: Approx. 300 mins (varies by usage)

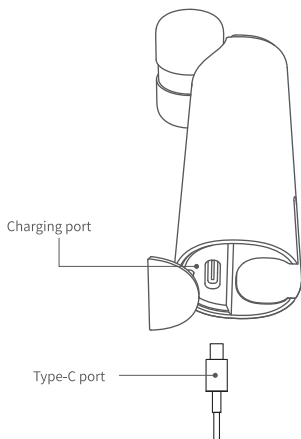Pre-Use Preparatio

Unboxing Inspection

Remove the main unit and accessories from the packaging. Verify all items against the packing list. Retain packaging materials for future transport/storage.



Inspect the main unit and accessories for mechanical damage. Do not use if functional defects are observed; contact the agent or after-sales service immediately.

Charging

Fully charge the device before initial use

Skin Compatibility Test

To avoid adverse reactions (redness/itching/burning) from multi-frequency electrical pulses and light energy:


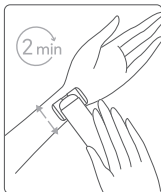
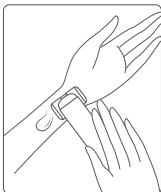
Clean inner forearm with soap, then dry thoroughly.

Power on the device and dispense single-dose serum onto the forearm.

If no skin abnormality occurs within 24 hours, proceed to normal use.

Precautions After Using Instruments

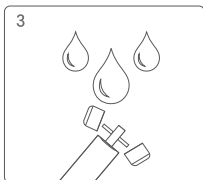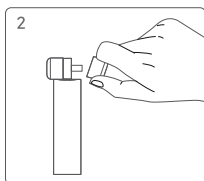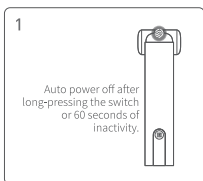Device Maintenance:

After power off, detach the metal guide head. Clean with a damp cloth or rinse under water, then wipe dry.

WARNING:

Clean the device only with water. Do not use alcohol, alkaline detergents, nail polish remover, or other corrosive liquids.

Ensure the power rubber seal is fully closed before rinsing to prevent water ingress into the power port.

Product maintenance and upkeep

- Keep the instrument dry; Do not immerse the device or any of its components in water to avoid the risk of electric shock.
To avoid damage or loss of the instrument and its accessories, please store them properly after each use. Before storage, please make sure to clean the instrument first;
 - Please keep the instrument properly when not in use for a long time, store it in a dry, dark, and cool place, avoid contact with hard objects to prevent damage to the instrument;
 - Do not store the instrument and its accessories in a humid environment for a long time, as it may have adverse effects on the performance of the instrument;
 - If the instrument has not been used for a long time and needs to be used, please fully charge the instrument and check its status.
- Product troubleshooting

Malfunction and maintenance

Problem	Potential cause	Solution
Start failure	No power left	Charge fully before use
	End of battery life	Contact the manufacturer for repair or replacement
Automatic power off	No operation for more than 60 seconds	Restart the device
	No power left	Charge fully before use
Unable to charge	Incompatible adapter used	Use the officially recommended adapter
Tingling caused by current	Device not pressed against the skin	Press the device against the skin
Redness of skin after use	Caused by RF heat	Normally recovers soon afterwards, please seek medical attention if you don't recover or have other discomfort.

This appliance contains batteries that are non-replaceable.
The battery must be removed from the appliance before it is scrapped.
The appliance must be disconnected from the supply mains when removing the battery.
The battery is to be disposed of safely.

Warranty

1 year limited warranty

For more information about the consumer warranty benefits, please refer to AMIRO' s official website www.amirobeauty.com. Except as prohibited by laws or otherwise promised by AMIRO, the after-sales services shall be limited to the country or region of the original purchase. Under the consumer warranty, to the fullest extent permitted by law, AMIRO, will, at its discretion, repair, replace or refund your product. Normal wear and tear, force majeure, abuse or damage caused by the user's negligence or fault are not warranted. If in doubt please contact the relevant person as AMIRO may identify.

Products which were not duly imported and/or were not duly manufactured by AMIRO and/or were not duly acquired from AMIRO or a AMIRO' s official seller are not covered by the present warranties. As per applicable law you may benefit from warranties from the non-official retailer who sold the product. Therefore, AMIRO invites you to contact the retailer from whom you purchased the product.

Manufacturer: Shenzhen Accompany Technology Co., Ltd.

Address: 25th Floor, Gaoxin Community United Headquarters Building, No.63, Gaoxin South 10th Road, Binhai Community, Yuehai Street, Nanshan District, Shenzhen, China.

For detailed e-manual, please go to www.amirobeauty.com.

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

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Fabricante: Shenzhen Accompany Technology Co., Ltd.