

Contents

- I. Product Information
- II. Instructions
- III. Usage
- IV. Safety Precautions
- V. Table of Common Malfunctions
- VI. Warranty

I. Product Information

Product description

The home-use Amiro LumoMax High-energy Light Therapy Device (LumoMax) is designed to address different skin issues through the use of four types of narrow-spectrum light.

Product specifications

- Product name: Amiro LumoMax High-energy Light Therapy Device
- Product model: ABM003
- Product size: 275×197×13 MM
- Rated input: 15 V===2A
- Country of origin: China


Product list

Amiro LumoMax High-energy Light Therapy Device x1
Cable x1
Adapter x1
Protective goggles x1
User manual x1
Instructions card x1

II. Instructions

Device illustration


High-energy Light Therapy Device

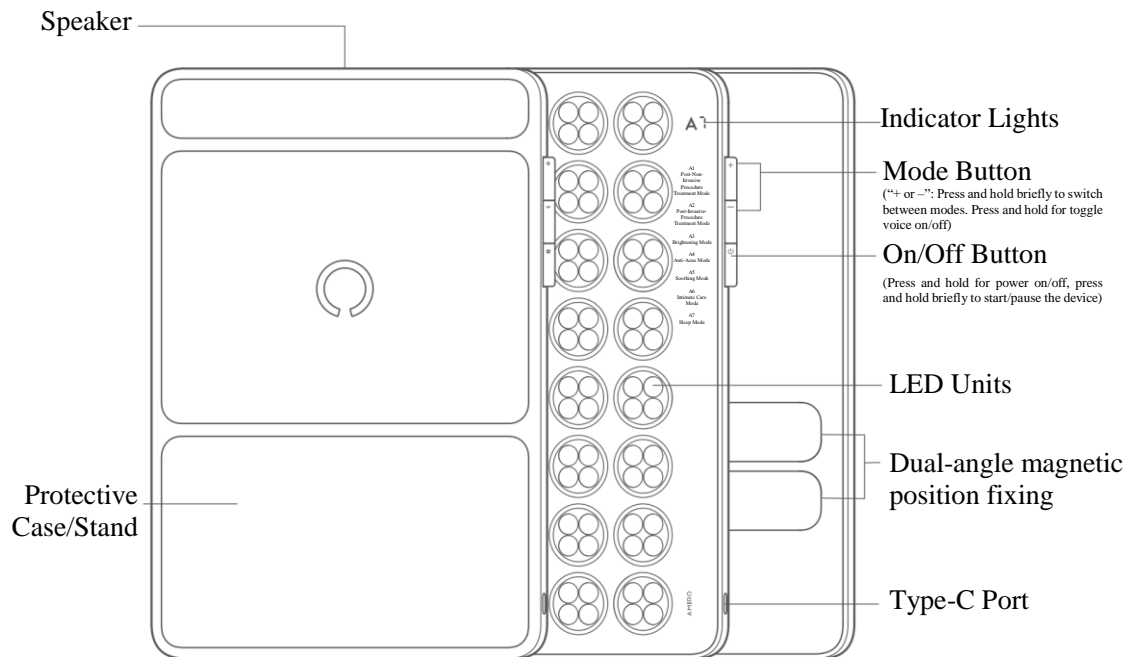

Operating procedure

1. Powering on and off and mode switching

- Connect the device to the power adapter that came with the unit. If the mode light in the upper right side of the panel lights up, it means the device has been successfully powered on.
- Press and hold the "+" or "-" button for 0.5 seconds to switch between different modes.
- Press and hold the power button for 0.5 seconds to begin treatment. Remember to wear protective goggles.

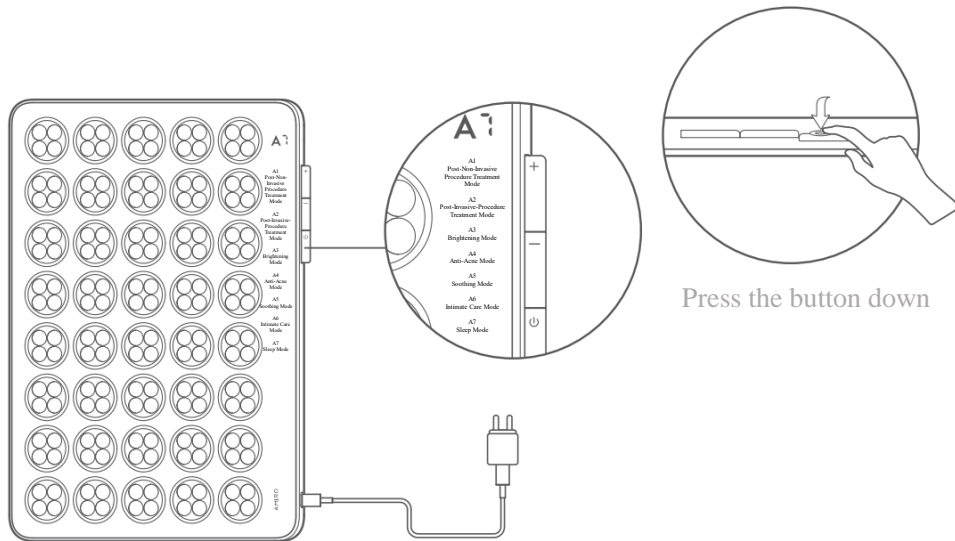


Connect the device to power supply

2. Pause and Resume Mode


- During treatment, press and hold the power button for 0.5 seconds to pause.
- To resume the light therapy after putting it on pause, simply press and hold the power button for 0.5 seconds.

*Note: When the device has been paused for more than 10 minutes, it will be automatically shut down and you will need to switch it back on before the next use.

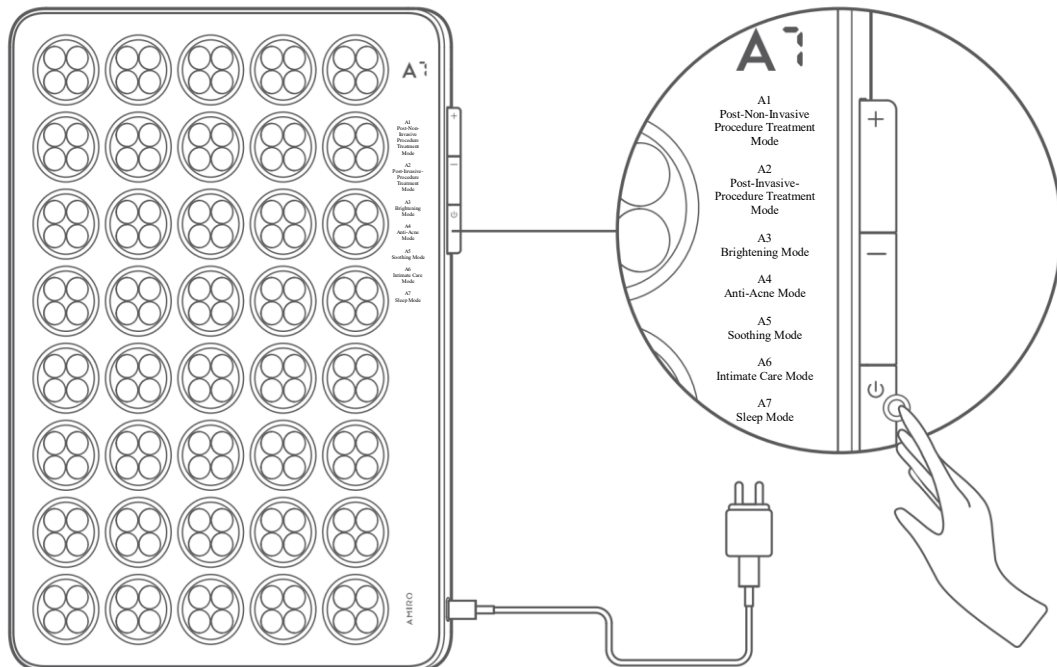

3. To Turn On/Off Voice Reminder


- When the device is powered on, press the "+" and "-" buttons simultaneously and hold for 2 seconds to turn off the voice reminder.
- When the device is powered on, press the "+" and "-" buttons simultaneously and hold for 2 seconds to turn on the voice reminder again.

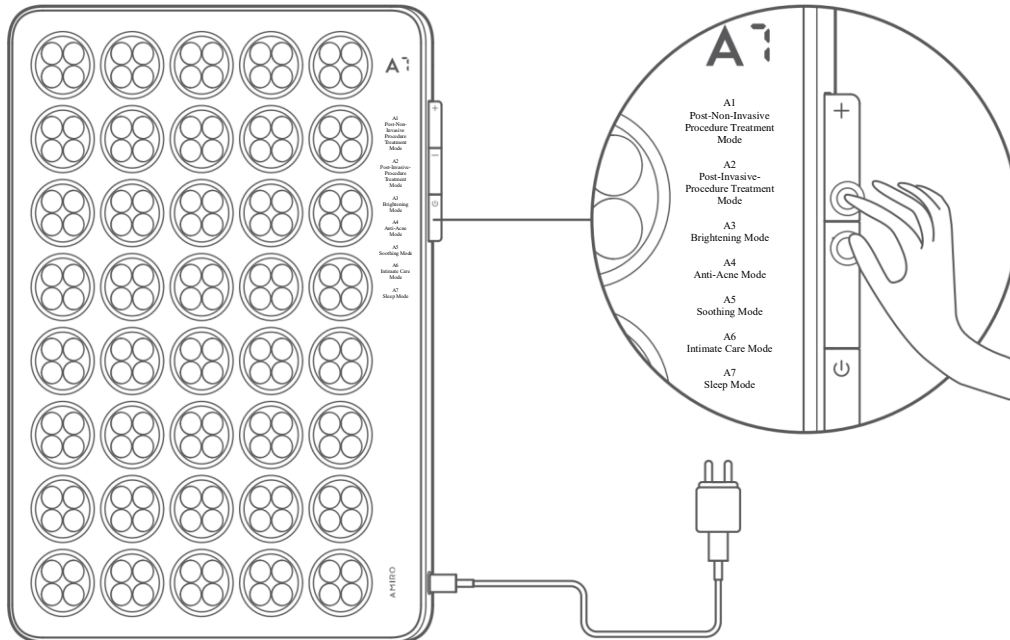

III. Usage

Care duration

- The device should only be used 3-5 times per week and for no longer than 30 minutes at a time
- For modes A1-A6, the device will automatically stop after 10 minutes of use. For mode A7, it automatically stops after 30 minutes of use.

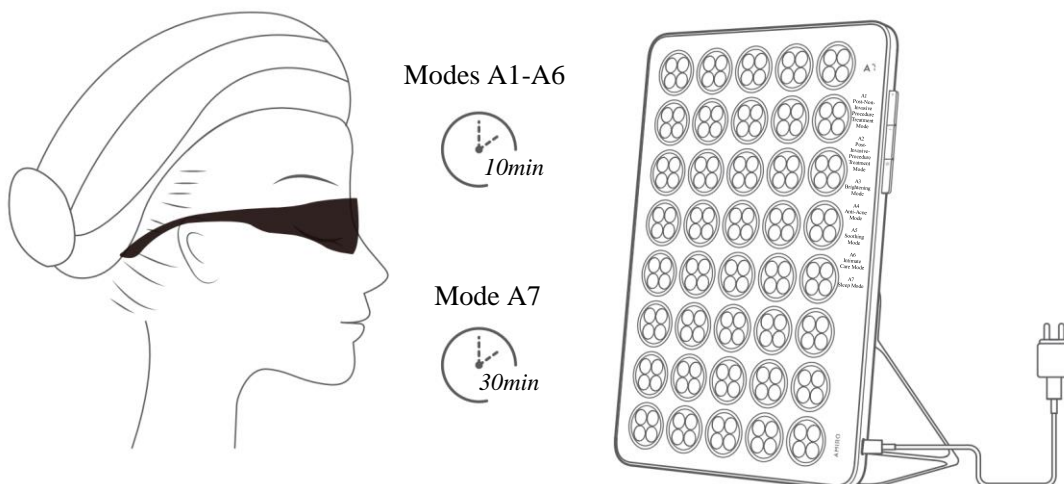

Care steps


Step 1: Makeup removal and cleansing

- Remove all makeup and cleanse your face thoroughly before using the product. Avoid using skincare products with photosensitive/whitening ingredients.



Step 2: Connecting the device and adapter

- Use only the official power adapter that comes with the product.

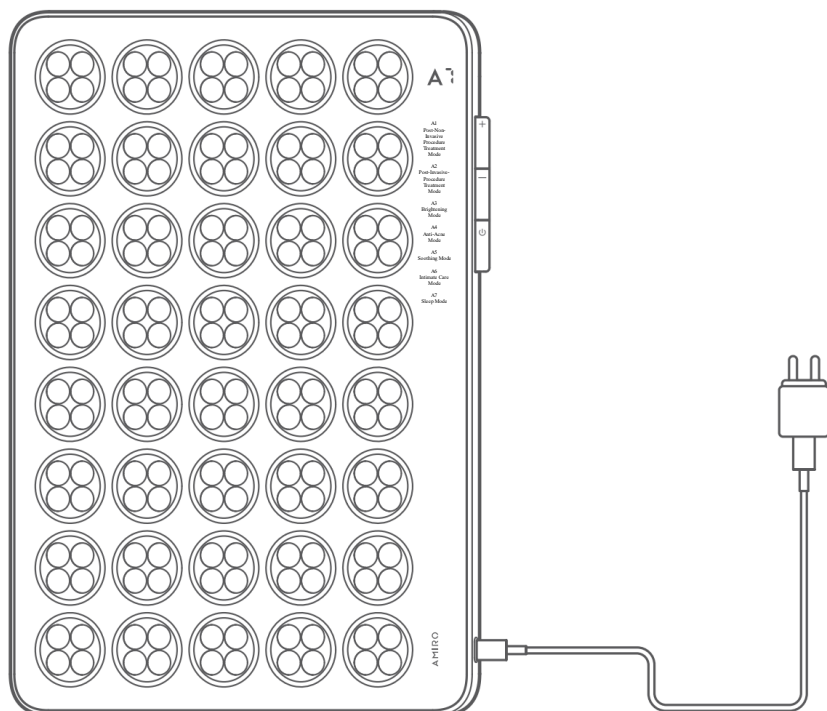


Adapter (Power Supply)

Step 3: To activate a mode

- Press and hold the "+" or "-" button for 0.5 seconds to switch between modes.
- Press and hold the power button for 0.5 seconds to begin treatment. Remember to wear protective goggles.

* Note: Please select the mode that best suits your skin condition using the treatment advice. Different modes are suitable for different types of skin status.

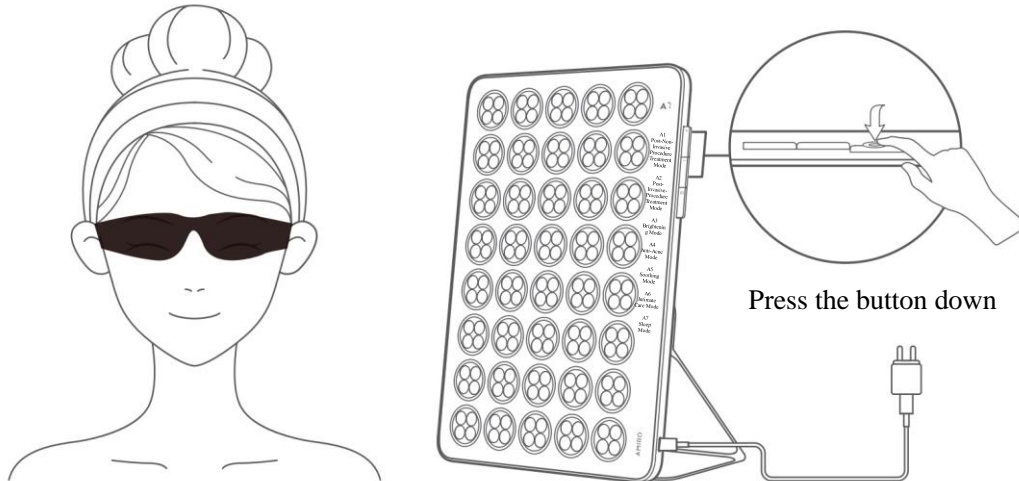


A1 [Post-Non-Invasive Procedure Treatment Mode]

This mode is designed for light therapy after non-invasive aesthetic procedures such as photo-rejuvenation, Thermage, chemical peels, etc.

Treatment Advice: Follow the post-procedure cleansing instructions and use skincare products without photosensitive ingredients for moisturizing and hydration. Switch to Post-Non-Invasive Procedure Treatment Mode to begin the therapy. Make sure to moisturize and protect your skin from the sun on a daily basis.

Wavelength: ● Yellow Light: 590nm ● Red Light: 633nm ● Infrared Light: 850nm

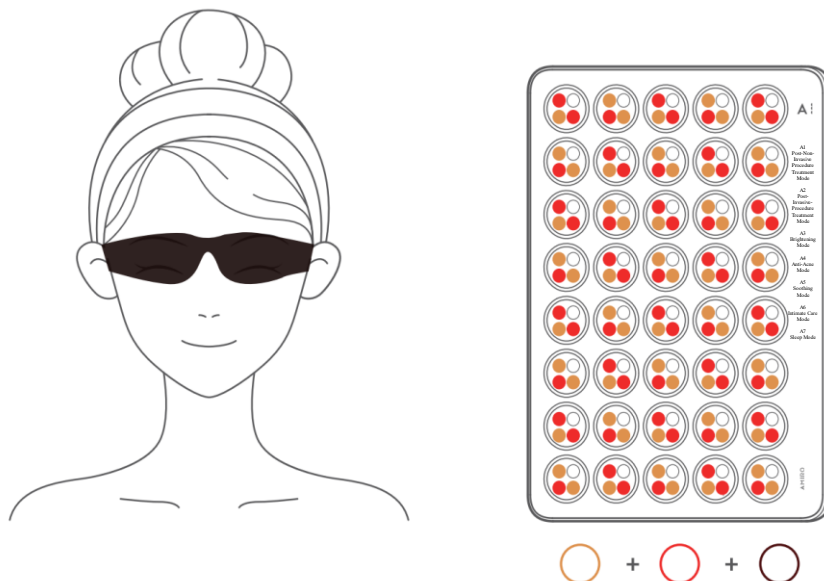


A2 [Post-Invasive-Procedure Treatment Mode]

This mode is designed for post-procedure treatment and care after invasive aesthetic procedures like mesotherapy, RF micro-needling, and acne extraction.

Treatment Advice: Follow the post-procedure cleansing instructions, then moisturize and hydrate the skin with skincare products that do not contain photosensitive ingredients. Next, switch to the Post-Invasive-Procedure Treatment Mode to start the therapy. It is recommended to moisturize and protect the skin from the sun daily after the therapy.

Wavelength: ● Yellow Light: 590nm ● Red Light: 633nm ● Infrared Light: 850nm

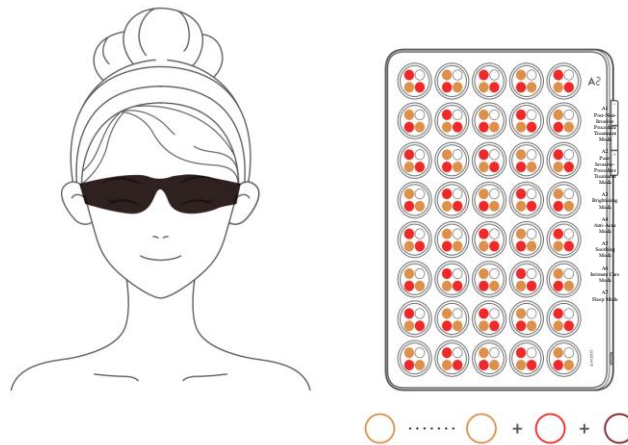


A3 [Brightening Mode]

This mode is designed for daily skin brightening and rejuvenation in order to help you achieve a brighter complexion and tender skin texture.

Treatment Advice: Cleanse your face with gentle and non-irritating cleansing products for cleansing, and moisturize and hydrate with a photosensitive ingredient-free skincare product. Then, switch to Brightening Mode to begin the therapy. After light therapy, you may double the brightening effects by using skincare products containing whitening ingredients.



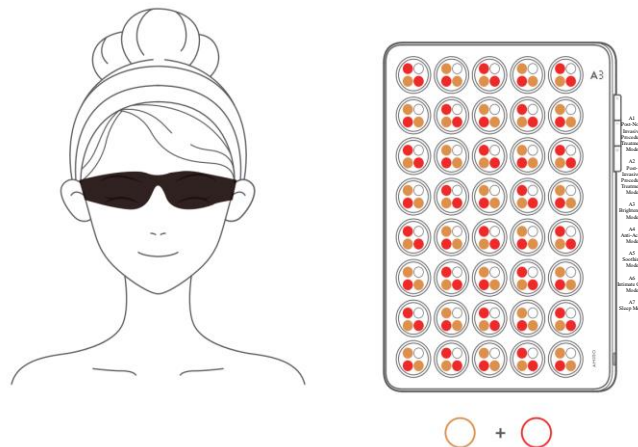

Wavelength: ● Yellow Light: 590nm ● Red Light: 633nm




A4 [Anti-Acne Mode]

This mode is suitable for treating acne and pimples.

Treatment Advice: Cleanse your face with gentle and non-irritating cleansing products for cleansing and moisturize and hydrate with a photosensitive ingredient-free skincare product. Then, switch to Anti-Acne Mode to begin the therapy. Maintaining healthy lifestyle habits such as regular exercise, adopting a healthy diet, and getting sufficient sleep can help reduce the likelihood of acne recurrence.

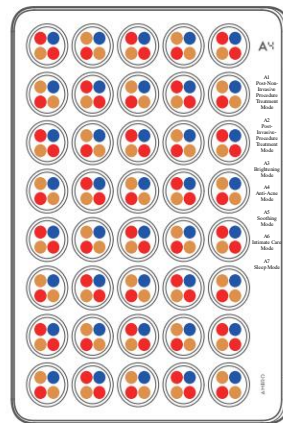

Wavelength:  Blue Light: 415nm  Red Light: 633nm  Yellow Light: 590nm  Infrared Light: 850nm




A5 [Soothing Mode]

This mode is designed for soothing sensitive skin and after-sun care.

Treatment Advice: Cleanse your face with gentle and non-irritating cleansing products for cleansing and moisturize and hydrate with a photosensitive ingredient-free skincare product. Then, switch to Soothing Mode to begin the therapy. The skin barrier is vulnerable when suffering from sensitivity, so it is important to moisturize and protect your skin from the sun on a daily basis.

Wavelength:  Yellow Light: 590nm  Red Light: 633nm

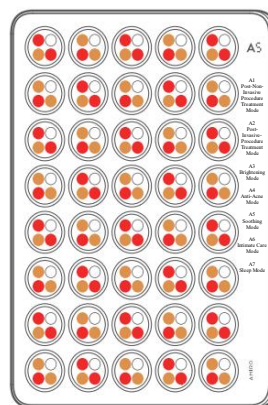


A6 [Intimate Care Mode]

This mode is designed for light therapy on intimate areas.

Treatment Advice: After proper cleansing, adjust to Intimate Care Mode to begin the therapy. If any abnormal symptoms such as redness, swelling, spots, or ulcers appear on intimate areas after therapy, it is necessary to consult a doctor for examination and treatment.

Wavelength: ● Blue Light: 415nm ● Red Light: 633nm ● Infrared Light: 850nm

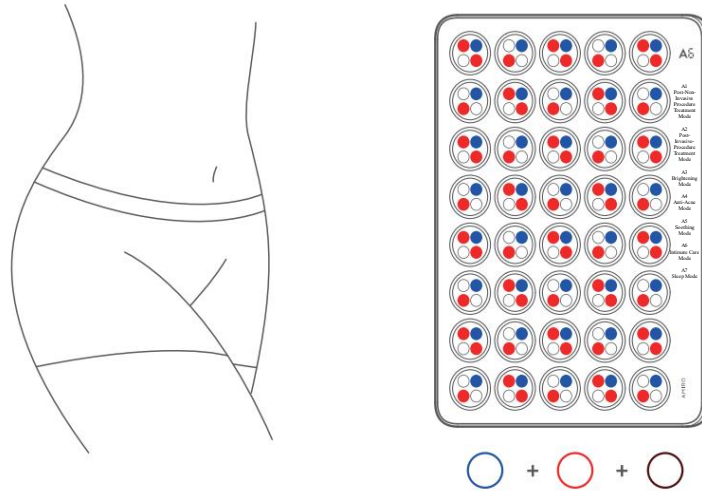


A7 [Sleep Mode]

This mode is designed for overall skincare before sleep.

Treatment Advice: Cleanse your face with gentle and non-irritating cleansing products for cleansing and moisturize and hydrate with a photosensitive ingredient-free skincare product. Then, switch to Sleep Mode to begin the therapy.

Wavelength: ● Yellow Light: 590nm ● Red Light: 633nm ● Infrared Light: 850nm ● Blue Light: 415nm

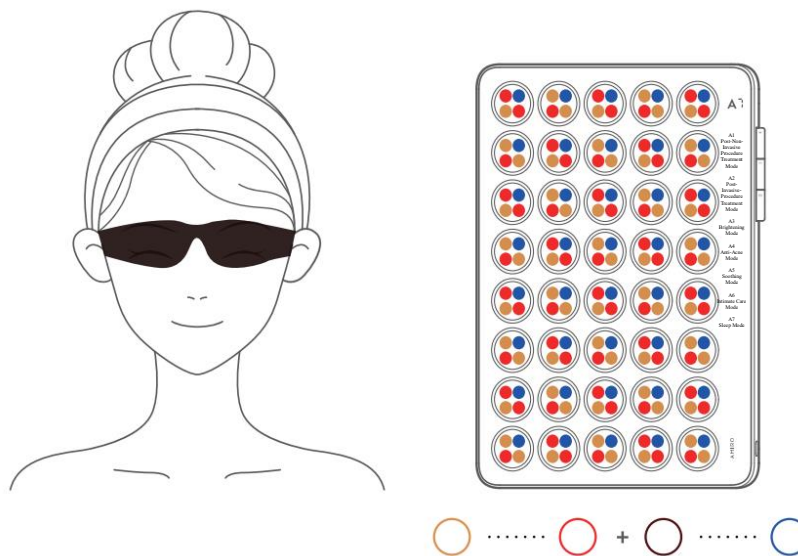

Step 4: Device maintenance

- After using the device, place it back into the packaging box or storage bag. If it gets dirty during therapy, gently wipe it with a slightly wet cloth or wet tissue.

IV. Safety Precautions

Tips

Before using the device, please read all the [Safety Precautions] carefully and strictly abide by the instructions in this manual.

- * [Warning]: Failure to follow the instructions for use may result in adverse effects.
- * [Warning]: Skin damage could develop by using the device on one area repeatedly or for an extended period of time.
- * [Warning]: Do not bypass or remove any built-in safety features of the product.

Operation safety measures:

- The recommended ambient temperature for this device is 10°C-30°C.
- When using the device for the first time or after a long period of storage, please ensure that the buttons and indicator lights are functioning properly.
- Perform a photosensitivity test on your skin upon receiving the product. Stop using immediately if you have photosensitive skin.

Precautions before everyday use

- Individuals who are prone to allergies, skin diseases, allergic dermatitis, sensitive skin, etc., may experience discomfort or adverse effects when using the product. If this occurs, stop using the device immediately.
- Do not use the device on damaged skin, such as areas with burns, blisters, incisions, wounds, active skin diseases, or infections.
- Use only the official power adapter that comes with the product.
- Before using the device, make sure to remove all makeup and cleanse the area you want to treat.
- Remove any metal products, glasses, and contact lenses before using the device.
- Do not use it for non-beauty purposes.
- Avoid immersing the device, including the adapter, in water.
- Do not use or store the adapter in an environment with high humidity.
- This device is not intended for use on children.

Product usage precautions:

- Do not look directly at the LED for a prolonged period to prevent damage to eyes.
- If you experience abnormal skin conditions, stop using the device immediately and contact our after-sales service or consult a professional doctor.
- Please select the appropriate mode based on your skin conditions when using this device.
- Avoid using the device on the same skin area for more than 30 minutes in one day to prevent skin damage.

After-Use Precautions

- After use, ensure that the device and adapter are kept in a dry and clean condition.
- When cleaning the device, make sure it is disconnected from the power cable. Keep the device dry and clean after cleaning for future use.
- It is normal for the skin to feel slightly warm and dry after use. Use the device with moisturizing skincare products to hydrate the skin.
- If you experience strong discomfort or abnormal feelings during use, stop using the device immediately and consult a doctor.

Precautions

- This product contains small parts. Keep it out of the reach of children to avoid accidental ingestion.
- Do not use it for non-beauty purposes.
- Do not use the device on pets as it may cause injury to them.
- Do not use the product for purposes other than those specified in the manual so as to avoid personal injury or device failure
- Do not allow individuals with physical or mental disabilities to use the device without supervision, as it may result in accidents or bodily injuries.
- Do not use the device if it is malfunctioning as it may result in bodily injury or danger.
- Do not use the device near flammable substances (such as alcohol, thinner, gasoline, nail polish remover, or hair spray) as it may result in fire, electric shock, or bodily injury
- Do not charge, use or store the product near a fire, heat source, or high temperature area, as it may cause the device to overheat, catch fire, or burst.
- Do not use near water sources, in bathrooms or while bathing, as it may result in electric shock, short circuit, or fire.
- Do not charge in water, as it may result in electric shock, short circuit, or fire.
- Avoid dropping or colliding the device and controller, as it may result in malfunctions. Do not use the device if it has been damaged by a fall or collision.
- Do not use corrosive solvents (alcohol, thinner, gasoline, nail polish remover, etc.) to wipe the device, as it may damage its surface
- Turn the power off when the device is not in use or being cleaned to prevent electric shock, short circuit, or fire.
- Please use power adapter specified by the company and make sure that the adapter is plugged directly into the socket when using the device; do not use an extension cord to avoid electric shock, short circuit or fire.
- Please do not disassemble or modify the device, as it may cause fire, electric shock, or physical injury. If the device needs repair, please contact our after-sales service.
- Please unplug the power adapter immediately if any abnormalities occur during use. Please stop using the device immediately if there is any smoke, noise or odor.
- Do not use the device after applying flammable cosmetics, as it may lead to fire, electrocution, or bodily harm.
- Do not use the device while using pacemakers, electrocardiogram scanners, and other electronic medical instruments, as it may cause instability in the medical instruments.

Situations where the device should not be applied:

People with any of the following conditions may not use this device as it could result in accidents or skin and bodily injury:

- The device is not recommended for individuals with extensive skin diseases or photosensitivity.
- Individuals with severe heart diseases, liver problems, blood disorders, nervous system diseases, malignant tumors, facial neuralgia, hemophilia, febrile diseases, and epilepsy triggered by light should strictly refrain from using this device.
- Pregnant women are advised not to use this device. If necessary, avoid using on the abdominal area.
- If you are taking medication that may cause photosensitivity, use this device with caution as advised by a doctor.
- If itching, swelling, erythema, rash, blisters, or any other abnormal symptoms occur during or after therapy, it could be a result of photosensitivity. It is important to stop using the device immediately and consult a doctor.
- If you feel uncomfortable when using the device, immediately stop the operation and consult a doctor.

Note: The safety warnings and instructions provided are not exhaustive. Please consult a doctor or contact us before deciding to use it if you are unsure about using the product.

V. Common Malfunctions

| Problem | Possible causes | Solution |
|---------------------|------------------------------------|---|
| Start failure | Using an incompatible adapter | Use only the official adapter that comes with the unit. |
| Automatic power off | Care procedure ends | After 10 minutes of treatment |
| | Pause time exceeds 10 minutes | When the device has been paused for more than 10 minutes, it will be automatically shut down. |
| Warm skin | Effect after phototherapy | Normal |
| Dry skin | Dehydrated skin after phototherapy | Use moisturizing skin care products |