




PORTABLE WATER FLOSSER

User Manual

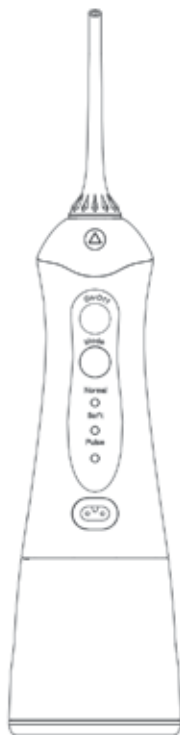

CONTENTS

Contents	01
ULStandard	02-05
Product Description	06
Product Accessories	07
Operating Procedures	08
Usage Instruction	09
Recommended Techniques	10
Disassembly Example	11
Battery Maintenance and Handling	12
Common Problems	13
After-sales Services	14

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

READ ALL INSTRUCTIONS BEFORE USING.

DANGERS:

- Do not handle charger with wet hands.
- Do not immerse in water or other liquid.
- Do not use in the bath or shower.
- Do not place or store product where it can fall or be pulled into a bath, shower or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the charger cord for damage before the first use and during the life of the product.


WARNING:

WARNING – To reduce the risk of burns, electrocution, fire, or injury to persons.


- A product should never be left unattended when plugged in.
- Close supervision is necessary when the product is being used by children.
- Use this product only for its intended use as described in this manual.
- Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug. If it is not working properly or is damaged, or has been dropped into water, return the product to the 360PRO Service Centre for examination and repair.
- Keep the cord away from heated surfaces.
- Never block the air openings of the product or place it on a soft surface such as a bed or couch. Keep the air openings free of lint and hair.
- Do not use outdoors or operate where aerosol products are being used or where oxygen is being administered.
- Connect this product to a properly grounded outlet only.
- Fill reservoir with water or mouth-rinse only. Never use salt or baking soda.

SAVE THESE INSTRUCTIONS.

PRODUCT DESCRIPTION

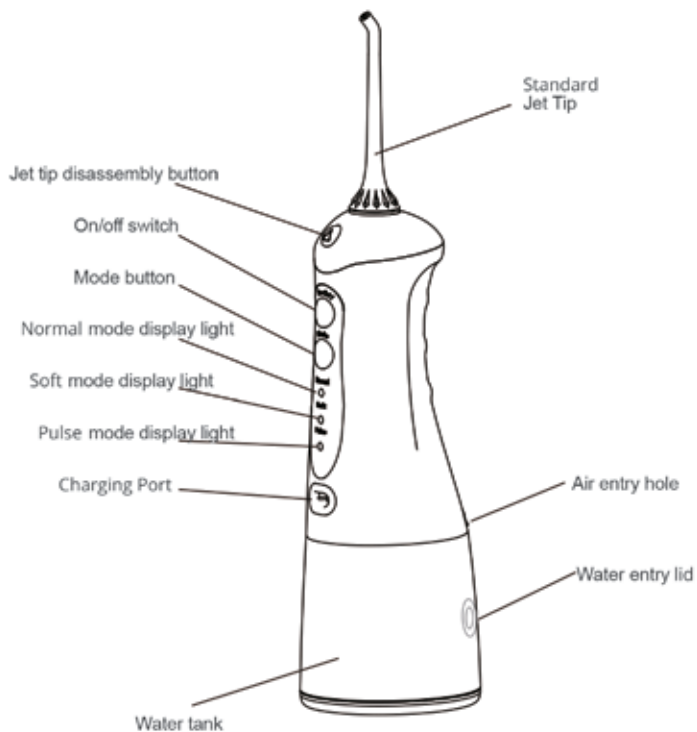

PRODUCT ACCESSORIES

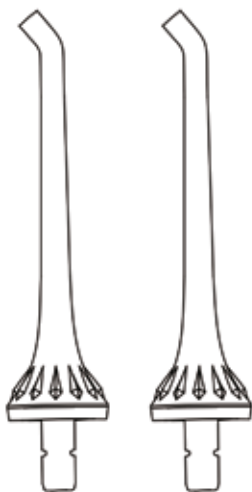

Standard jet tips x 2

The jet tips and other changable components can be purchased at www.my360pro.co.nz

It is recommended that you replace your standard jet tips every 3 months or less if they show signs of wear.

OPERATING PROCEDURES



SET-UP


Extremely important!! Prior to using for the first time, please fully charge the device, and insert the plug into the power socket, as well as connect to a power supply. During charging, the product cannot be used.

INSTALLATION AND DISASSEMBLY OF TIPS

Insert the standard jet tip into the spray head hole until you hear a click.



Open the water tank lid and pour water into the tank. Note: please ensure that when pouring water into water tank, the device is switched off.

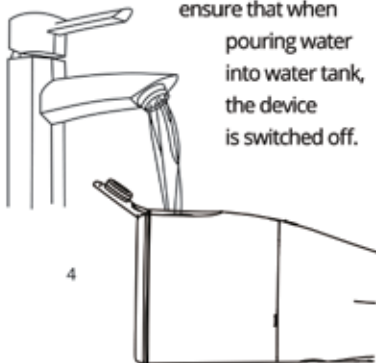

If changing or removing the tip, please press on the disassembly button on the jet tip handle before removing the jet tip.




INSTRUCTIONS FOR USE

- The unit has 3 modes: Normal, soft and pulse.
- The modes may be changed when the device is switched off, and you can select the mode you like before switching the device back on.
- Irrespective of the mode, after operating for 120 seconds the product will automatically switch off.
- During recharge, the three display lights will flash, when all three lights are constantly on, this indicates that the unit is recharged.
- During usage, if all three display lights are flashing, this indicates that the products power is insufficient, and recharging is required.
- During recharge, this product cannot be used.
- Rotate the jet tip dial to adjust the direction of the water flow.




RECOMMENDED TECHNIQUES

Holding the device at a downward slant, position the jet tip in the mouth, and direct the jet tip at teeth before switching on the device.

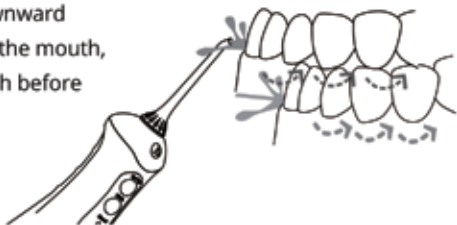

When the jet stream is at 90 degree angle to the gum line, slightly close your mouth to avoid the water spraying out. Slide the tip along the gum line as well as pausing briefly between teeth until the area around the gums and the spaces between the teeth are fully cleaned.

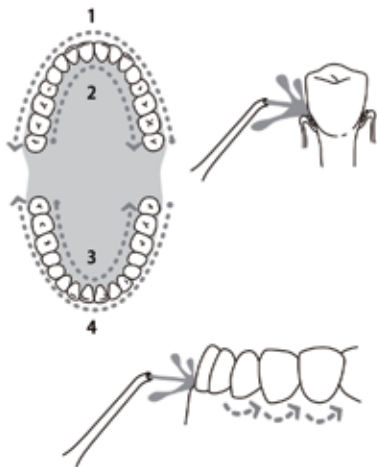


DISASSEMBLY EXAMPLE

AFTER USE

Please switch off the power, as well as emptying the water tank.



Using both hands as indicated in Diagram 2 – use your left hand to hold the product body, and your right hand to hold the water tank. At the same time rotate in an anti-clockwise direction, in order to remove the water tank. Empty.

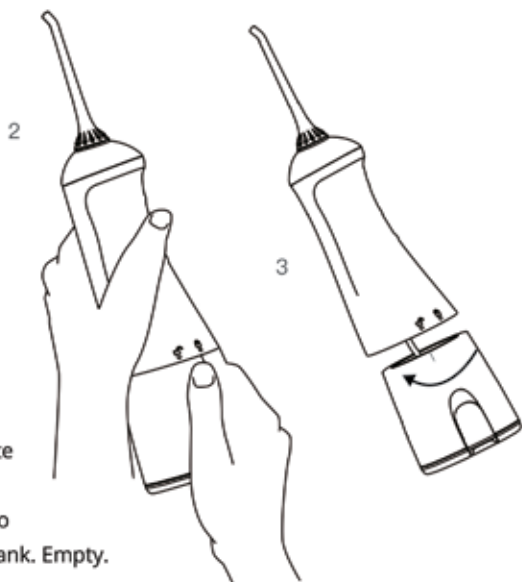

BATTERY MAINTENANCE AND HANDLING

BATTERY MAINTENANCE

- Prior to initial usage the device must be fully charged, for approximately four hours.
- If to be stored for a long period of time, please fully charge the device in advance.
- In order to maximize the usage life of the battery, please promptly recharge when power is insufficient.

BATTERY HANDLING

- The device includes a lithium battery which cannot be changed. Towards the end of the battery life, please dispose of responsibly.

COMMON PROBLEMS AND SOLUTIONS

Problems	Reason	Method of handling
Leakage of water between jet tip and screw part	The seal ring on the fixed socket for the jet tip is damaged	Contact the 360PRO Service Centre for instruction
Water tank leaks water	Seal ring is loose	Re-install correctly
Inadequate pressure	Low voltage	Recharge
Device cannot be switched on	Low voltage	Recharge

If the above methods cannot resolve these problems, please consult the 360PRO Service Centre at www.my360pro.co.nz

Gentle Reminder:

- At the time of initial usage, please select "gentle mode", in order to avoid gums bleeding due to lack of adaption of usage. After consistent usage for two weeks there should be improvements. If bleeding is excessive or bleeding continues after two weeks, chronic gum inflammation could be an issue, and it is recommended that professional dental treatment be sought.
- It is not recommended that the product be used by people with oral ulcers or open oral sores.

AFTER-SALES SERVICES

THANK YOU FOR PURCHASING THE 360PRO PORTABLE WATERFLOSSER- PLEASE TAKE NOTE OF THE FOLLOWING ISSUES:

- Please retain the sales invoice issued by the retailer as a record.
- One-year product warranty for non-human malfunctions is provided with the purchase of this product.
- Where the following circumstances arise during the warranty period, costs and services fees shall be charged:
 - Damage due to inappropriate usage or maintenance including lack of cleaning.
 - Damage caused by a third party repairer.
 - Damage caused by weather events (such as storms, earthquake or fire).
- The unit must only be repaired by the 360PRO Service Team or the warranty will be voided.
- Replacement tips can be purchased at www.my360pro.co.nz

Thank you for purchasing the 360PRO Portable water flosser.

Water flossers pulse water to clean the teeth, and remove harmful bacteria below the gumline, places which cannot be cleaned by brushing the teeth alone. Manual toothbrushes are unable to effectively clean all areas of the mouth. Dental floss can only make contact with a limited area beneath the gum. The 360PRO portable Water flosser can clean out harmful bacteria hidden between teeth and under the gum line.

Extra replacement tips can be purchased at www.my360pro.co.nz



360PRO

PO Box 101 613, NSMC, Auckland, NZ 0632

P: +64 9 477 0285 | www.my360pro.co.nz

Designed in New Zealand, made in China.

