


waterpik®

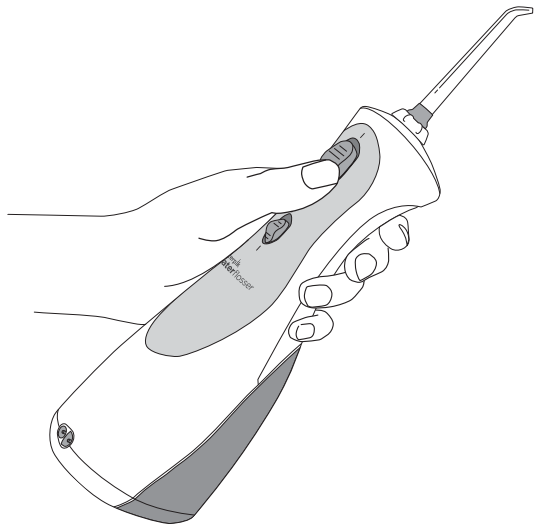

Waterpik® Water Flosser Model WP-450

.....
www.waterpik.com

TABLE OF CONTENTS

Important Safeguards	3
Product Description	5
Getting Started	7
Cleaning and Troubleshooting Your Waterpik® Water Flosser	8
Limited Two-Year Warranty	10

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Do not handle charger with wet hands.
- Do not place in or drop into water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the charger cord for damage before the first use and during the life of the product.

WARNING:


To reduce the risk of burns, electrocution, fire, or injury to persons:





- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Contact Waterpik International, Inc. or your local distributor which can be found at www.waterpik.com/intl.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use tips and accessories that are recommended by Water Pik, Inc.


- Do not drop or insert any foreign object into any opening on the product.
- Keep charging unit and cord away from heated surfaces.
- Do not use this product outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use iodine, bleach, or water insoluble concentrated essential oils in this product. Use of these can reduce product performance and will shorten the life of the product.
- Fill reservoir with water or other dental professional recommended solution only.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- Closely instruct and supervise children 6 years and above and individuals with special needs in the proper use of this product.
- Children should be supervised to ensure that they do not play with the appliance.

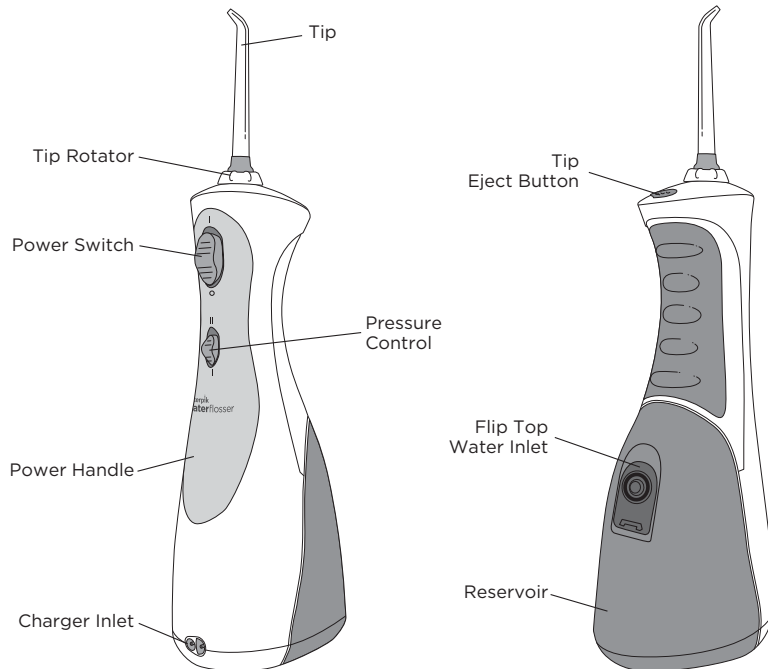
SAVE THESE INSTRUCTIONS.

PRODUCT DESCRIPTION



Replace Every 3 Months	
Orthodontic Tip	BEST FOR
	<input checked="" type="checkbox"/> Braces <input checked="" type="checkbox"/> General Use
Pik Pocket™ Tip	BEST FOR
	<input checked="" type="checkbox"/> Periodontal Pockets <input checked="" type="checkbox"/> Furcations
Plaque Seeker® Tip	BEST FOR
	<input checked="" type="checkbox"/> Implants <input checked="" type="checkbox"/> Bridges <input checked="" type="checkbox"/> Crowns <input checked="" type="checkbox"/> General Use
Toothbrush Tip	BEST FOR
	<input checked="" type="checkbox"/> General Use

Replace Every 6 Months	
Classic Jet Tip	BEST FOR
	<input checked="" type="checkbox"/> General Use
Tongue Cleaner	BEST FOR
	<input checked="" type="checkbox"/> Fresher Breath

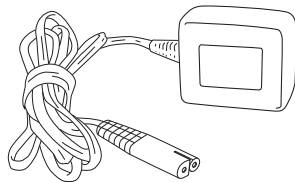

Charger
Waterpik® p/n 20016188

NOTE: All tips not included in all models.


Preparation and Use

Refer to the Quick Start Guide for proper use instructions.


When Finished

Turn the unit OFF (I). Empty any liquid left in the reservoir.

If desired, the reservoir may be removed from the power handle by pulling the reservoir release latch while sliding the reservoir down towards the base of the unit.

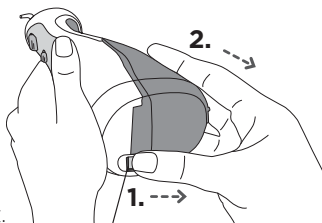

Reservoir
Release Latch

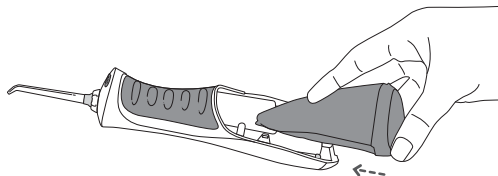

To replace the reservoir on the power handle, simply slide it up towards the top of the unit. The reservoir will snap into place when fully engaged.

TIP USAGE

Pik Pocket™ Tip

The Pik Pocket™ Tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets.

To use the Pik Pocket™ Tip, **set the unit to the lowest pressure setting**. Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket. Turn the unit on and continue tracing along the gumline.




Tongue Cleaner Tip

To use the tongue cleaner tip, set to lowest pressure setting, place tip in the center/middle of your tongue about half way back and start the Water Flosser. Pull forward with light pressure. Increase pressure as you prefer.



Orthodontic Tip

To use the Orthodontic Tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

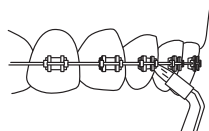

CLEANING AND TROUBLESHOOTING YOUR WATERPIK® WATER FLOSSER

Toothbrush Tip

Place the toothbrush tip in mouth with the brush head on the tooth at the gumline. The toothbrush tip can be used with or without toothpaste. Turn the Water Flosser on so the water flows through the tip. Using a **light pressure** (bristles should not bend) massage the brush back and forth with very short strokes-much like you would with a manual toothbrush.



Plaque Seeker® Tip

To use the Plaque Seeker® Tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.



Using Mouthwash and Other Solutions

Your Waterpik® Water Flosser can be used to deliver mouthwash and antibacterial solutions. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

Cleaning

Never immerse unit in water. Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser to wipe exterior.

The water reservoir is removable for easy cleaning and is top rack dishwasher safe. It is recommended that the reservoir is removed and cleaned weekly.

Removing Hard Water Deposits

Add 1 tablespoon of white vinegar to a full reservoir of warm water. Point the handle and tip into sink. Turn unit ON (I) and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water.

Battery Care and Disposal

The unit will not charge if the switch is in the ON (I) position. If you use your Water Flosser once a day or less, it is recommended you charge it overnight, once a week. If you use it twice a day or more, charge the unit more frequently. Do not run the batteries down completely as this may shorten the lifespan of the batteries and will require charging the unit for 24 hours to obtain a full charge.

If the Water Flosser is to be stored for an extended period of time (six months or more), be sure to fully charge it prior to storage. The unit

contains a Nickel Metal Hydride battery. At the end of the battery life, please recycle the unit at your nearest recycle center per local or state requirements. The battery is not replaceable.

Service Maintenance

Waterpik® Water Flossers have no consumer-serviceable electrical items and do not require routine service maintenance. For all your service needs contact your local distributor which can be found at www.waterpik.com/intl.

Some replacement parts and accessories for the Water Flosser are also available for sale at the customer service numbers shown below. Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

Distributed By:

Brandfarm Pty Ltd
Lvl 5 Nexus Building
4 Columbia Crt
Norwest Business Park
Norwest NSW 2153
+61 2 8860 6463

LIMITED TWO-YEAR WARRANTY

Water Pik, Inc. warrants to the original purchaser/owner of this product that it is free from defects in materials and workmanship for two years from date of purchase. Save your receipt as your proof of purchase date. Water Pik, Inc. reserves the right to request the product to be returned for analysis to confirm the claim falls within our warranty policy. We will replace any part of the product, which in our opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase and was used according to these instructions. This limited warranty excludes accessories or consumable parts such as tips, etc. Water Pik, Inc. will not be responsible for any expenses incurred as a result of any warranty claim accepted by Water Pik, Inc. This warranty gives you specific additional legal rights to those which you may have depending on the law in your place of residence. Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and for compensation for any other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure.

To make a warranty claim and before sending the product, please contact:

Brandfarm Pty Ltd
Lvl 5 Nexus Building
4 Columbia Crt
Norwest Business Park
Norwest NSW 2153
+61 2 8860 6463

This warranty is given by:

Water Pik, Inc.
1730 East Prospect Road
Fort Collins, CO 80553-0001
USA



Water Pik, Inc.

1730 East Prospect Road
Fort Collins, CO 80553-0001 USA

www.waterpik.com

Distributed By:

Brandfarm Pty Ltd
Lvl 5 Nexus Building
4 Columbia Crt
Norwest Business Park
Norwest NSW 2153
+61 2 8860 6463
www.waterpik.com.au

Waterpik® is a trademark of Water Pik, Inc. registered in Argentina, Australia, Benelux, Canada, Chile, China, Hong Kong, India, Japan, Korea, Mexico, New Zealand, Norway, Pakistan, Russian Federation, Singapore, South Africa, Switzerland, Taiwan, Ukraine, and the United States.

Waterpik® (stylized) is a trademark of Water Pik, Inc. registered in Australia, Brazil, Canada, China, EU, Hong Kong, Japan, Mexico, Korea, Russian Federation, Switzerland, Trinidad & Tobago, and the United States.

Pik Pocket™ is a trademark of Water Pik, Inc.

Plaque Seeker® is a trademark of Water Pik, Inc. registered in Canada, China, EU, and the United States.

Designed in USA. Made in CHINA.

Form No. 20016977-F AA

©2012 Water Pik, Inc.