



## QUICK START GUIDE

### 1 Charge Unit for 24 Hours Prior to First Use.


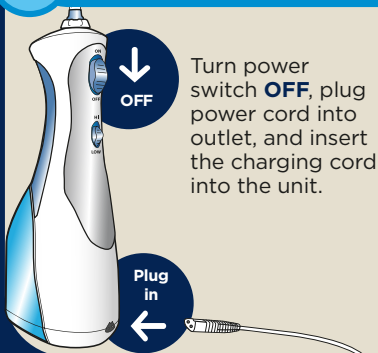


### 2 Insert Tip:

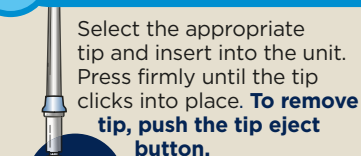


**Caution:** Do not press tip eject button while unit is running.

### 3 Fill the Reservoir:

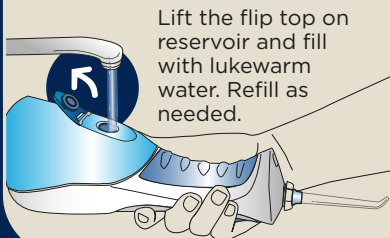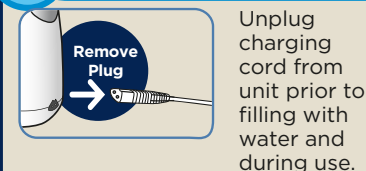

### 4 Usage:

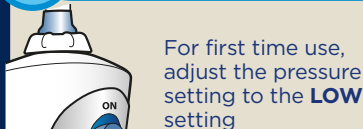

#### Pressure Setting




To increase water pressure move to **HI** setting.

*continued on back*

## 5 On/Off Procedure:

**Always place tip in mouth before turning unit on.**


Lean over sink and close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink. Turn unit **ON**.

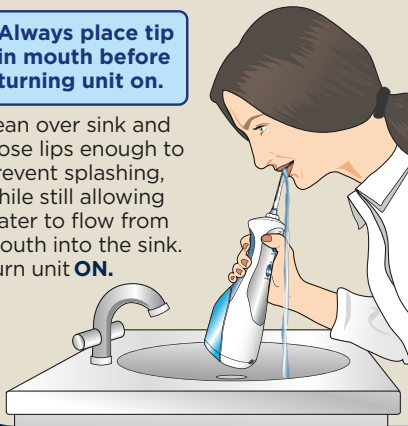


Always turn unit **OFF** before removing tip from mouth.

## 6 Technique Tips for Best Results:

Aim the water stream just above the gumline at a 90 degree angle.

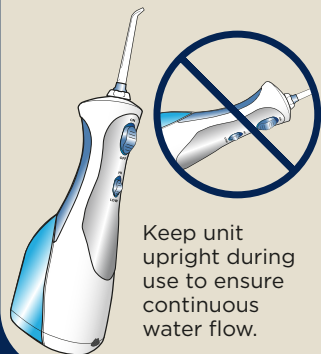

Keep unit upright during use to ensure continuous water flow.

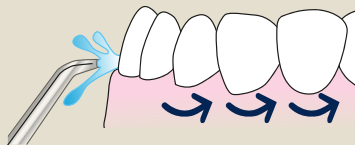

Starting with the back teeth, follow the gumline and pause briefly between teeth. Continue until all areas around and between teeth have been cleaned.