


# MIO

Pore Extraction  
& Diamond  
Microderm  
Skin  
Resurfacing  
System



## A Closer Look MIO

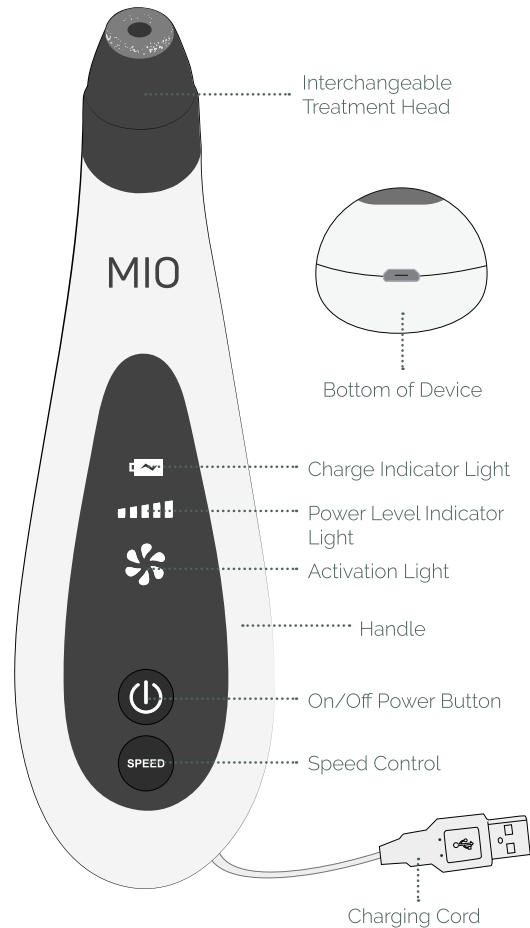

Test a small area first on forearm in a way that you feel suction. Move in a straight line and wait one hour to evaluate results.

## Usage Guide MIO



**PRE-TREAT.** Thoroughly cleanse skin with warm water, steam or a warm wash cloth to open pores. Pat dry with a clean towel.


1. Attach the **Diamond Microdermabrasion Tip** and follow below "**Usage Directions**" to expertly remove dead skin and rough dry patches with professional grade exfoliation. Remove Diamond Microdermabrasion Tip after use.



**Usage Directions:** Press the "ON/OFF POWER" button to activate or deactivate the device. Select the desired power level by pressing the "Speed Control" button. **We suggest beginning on Low power to start.** Begin gliding MIO over 1"-2" treatment sections while holding your skin taut with the opposite hand. \*See Illustrations under Expert Advice, other side. **Do not stop and hold MIO over an area of skin. Keep MIO moving at all times.** When treatment is complete, rinse your face clean and pat dry. **Discontinue use if you experience any discomfort or pain.**




## Usage Guide MIO

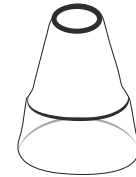


2. Extracting Suction Tip - Medium Circular Opening

2. Attach the **Extracting Suction Tip** with the medium or large circular hole at tip to remove congestion and detoxify pores. These tips are designed to be used in easy to reach places like the chin, cheeks and forehead. Begin gliding MIO over 2" sections of skin while holding your skin taut with the opposite hand, avoiding the eye and lip area.

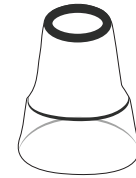

2. Extracting Suction Tip - Large Opening

3. Attach the **Extracting Suction Tip** with the small oval opening to remove congestion and detoxify pores in hard to reach places. The angled opening makes this tip ideal to purify pores around the nose area and in the T-Zone. Begin gliding MIO over 1"-2" sections of skin while holding your skin taut with your opposite hand. Do not stop and hold MIO over an area of skin. **Keep MIO moving at all times**, to avoid too much suction in one area.

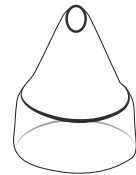

3. Extracting Suction Tip - Oval Opening

**Note:** MIO will shut off automatically after 15 minutes.

## Simple Charging / Filters



Filter




Plastic Cap




Plastic Rings

- MIO is cordless and rechargeable. Before using for the first time, charge fully for 24 hours. MIO can run up to 40 minutes on a single charge.
- To recharge, insert the small end of the Charging Cord into the port on the bottom of the MIO. Then connect the USB end of the Charging Cord to (i) a USB port in a computer or other USB power supply, or (ii) a wall adapter and then plug into an electrical outlet.
- When the MIO's battery is in need of charging, the Power Indicator Light will turn to red. If the device is fully discharged the device will power off.
- As the battery charges, the flashing red light will turn solid red, and then will turn green when the device is fully charged.
- Recharge the MIO only as needed. DO NOT overcharge which can shorten battery life. Disconnect the Charging Cord when not charging the device.

## Expert Advice MIO

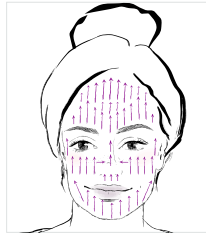

Vertical Pattern for Diamond Microdermabrasion Tip

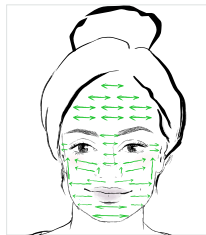

Horizontal Pattern for Diamond Microdermabrasion Tip

**Keep Moving**



**Test a small area first on forearm in a way that you feel suction. Move in a straight line and wait one hour to evaluate results.**

- Diamond Microdermabrasion Tip.** Aestheticians recommend using it 1 time per week. It is best to use it in gentle, short upward movements of about 1"-2" long. **Avoid thin skinned areas around the eyes and lips. Hold skin taut at all times.**
- After the first Vertical pass with the **Diamond Microdermabrasion Tip**, you may decide to make a second pass in the opposite direction. Aestheticians refer to this as cross hatching. It is best to make only one pass as your skin is getting used to microdermabrasion.

- Keep MIO Moving.** It's important to keep MIO moving at all times. Never keep MIO stationary over one area of skin. Moving in short 1"-2" sections is recommended and always hold skin taut with opposite hand.

## Warnings & Warranty

### BASIC SAFETY PRECAUTIONS.

When using electrical products, especially when children are present, basic safety precautions should always be taken. IF YOU THINK THE POWER IS OFF WHEN THE SWITCH IS OFF, YOU ARE WRONG. If the plug is in, the power is on. Unplug electrical appliances when you are not using them.

### DANGERS.

To reduce the risk of electrocution: DO NOT get the Product or Charging Cable wet or immerse in water. Keep this product away from open flame or flammable aerosol products. The Device is not suitable for use in the presence of a flammable anesthetic mixture with air or oxygen or nitrous oxide.

### WARNINGS.

To reduce the risk of burns, electrocution, fire or injury to persons: Close supervision is necessary when this Product is used by, on, or near children or physically challenged persons. Use this product only for its intended use as described in this User Guide. DO NOT use this product where oxygen is being administered. The product has no serviceable parts. Never attempt to open the Device Handle for any reason. DO NOT disassemble, heat above 600°C ( 1400°F) or incinerate.

### IMPORTANT CAUTIONS.

This is not a medical Device and is not intended to treat medical conditions. If you are under a doctor's care for any skin condition, if you think you may have a skin condition or if you have very sensitive skin, consult your physician before using this product. This product uses magnets which may affect internal medical devices such as pacemakers. Keep away from internal medical Devices. This Device is electrical and contains small parts. Keep out of reach of children. THIS PRODUCT IS NOT DESIGNED OR INTENDED PRIMARILY FOR CHILDREN 12 YEARS OF AGE OR YOUNGER.

**Warranty:** Your warranty is not effective unless you register the Product within 90 days of purchase. To register go online to [www.warrantyspasciences.com](http://www.warrantyspasciences.com) and complete the information required for registration.

For a period of one (1) year from the date of purchase, Spa Sciences (the "Company") warrants the MIO Microdermabrasion & Pore Extraction Skin Resurfacing System (the "Product") against damage, defects or failure resulting from ordinary consumer use under normal home conditions according to the instructions included with the Product (a "Covered Defect"). This warranty DOES NOT include damage,

## Warranty

defects or failure caused by misuse or improper care of the Product; impact or insertion of foreign objects into or on the Product; application of harmful chemicals or abrasives to the Product; failure to follow manufacturer's instructions; normal wear and tear; alteration or overuse of the Product; or any professional or commercial use of the Product. If there is a Covered Defect during the one (1) year warranty period the Company agrees to exchange the Product with a new equivalent Product. This is the sole remedy available for breach of this limited warranty. Any claim for breach of this limited warranty must be commenced within one (1) year following the original purchase. This limited warranty applies only to a purchase of the Product from or an authorized reseller, is valid for the original end user of the Product only and cannot be transferred under any circumstances.

THE PROVISIONS OF THE FOREGOING WARRANTIES ARE IN LIEU OF ANY OTHER WARRANTY, WHETHER EXPRESS OR IMPLIED, WRITTEN OR ORAL (INCLUDING ANY WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE). IN NO EVENT SHALL THE COMPANY BE LIABLE FOR ANY SPECIAL, INCIDENTAL, CONSEQUENTIAL, PUNITIVE, OR INDIRECT DAMAGES BASED UPON ANY LEGAL THEORY. SOME STATES DO NOT ALLOW LIMITS ON WARRANTIES, OR ON REMEDIES FOR BREACH IN CERTAIN TRANSACTIONS. IF APPLICABLE LAW DOES NOT PERMIT THE DISCLAIMER OF IMPLIED WARRANTIES, THEN THE DURATION OF SUCH IMPLIED WARRANTIES, ARE EXPRESSLY LIMITED TO THE DURATION OF THE EXPRESS LIMITED WARRANTY SET FORTH HEREIN.

## Additional Info

The entire operating manual including safety, care and warranty information is available at [www.miomanual.com](http://www.miomanual.com)

Email customer service at [customerservice@spasciences.com](mailto:customerservice@spasciences.com)