





A Closer Look

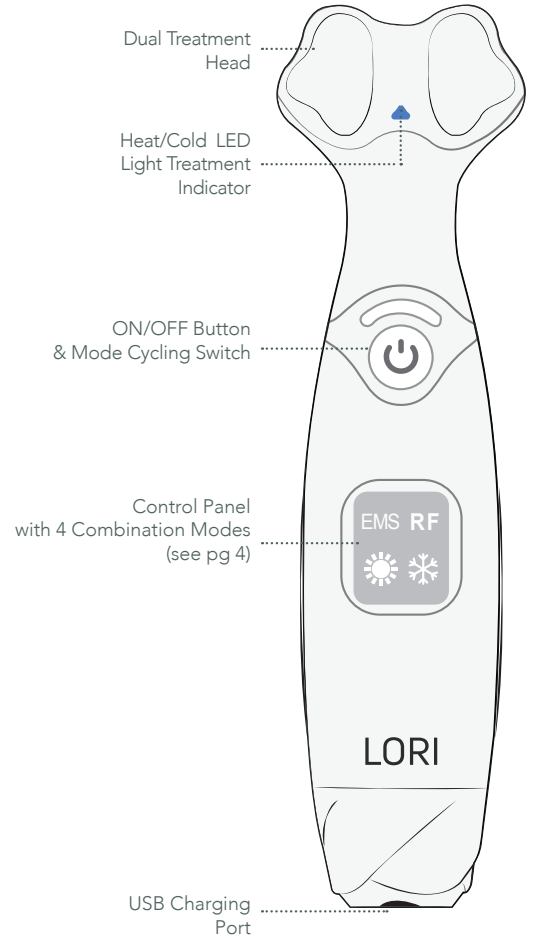

Charging


BEFORE FIRST USE

LORI is cordless and rechargeable. Before use, fully charge for 4 hours, even if the indicator lights turn green sooner than 4 hours.

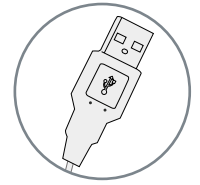


CHARGING

To access the USB charging port, open the waterproof seal cover located at the base of the handle. Connect the USB end of the provided charging cord to a USB port of a computer, or a USB wall adapter, and plug into an electrical outlet. As the LORI charges all the lights on its control panel will be red (figure a.) The indicator light will turn green when the device is fully charged.

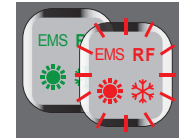

Recharge the LORI only as needed and up to 4 hours only. DO NOT overcharge which can shorten the battery life. Disconnect the charging cord when not charging the device. LORI should be powered OFF and cannot be used while charging.

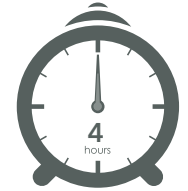

POWERING OFF AND ON

Long Press the power button to turn the device ON/OFF. LORI will shut off automatically after about 10 minutes of continuous use.



CYCLING THROUGH MODES

Short Press the same main button to cycle between the different modes of LORI.

Usage Guide

OPERATIONAL MODES:

Once you have assured of the full charge of the LORI you may proceed with the treatment of your choice:



MODE 1. = EMS + heat + red indicator light.

* Hot temperature is about 45°C (±3-5°C) = 113°F (±37-41°F)

MODE 2. = RF + heat + red indicator light.

MODE 3. = EMS + cold + blue indicator light.

* Cold temperature is about 9-15°C (48-59°F)

MODE 4. = RF + cold + blue indicator light.

Usage Guide

Expert Advice


Warnings & Warranty

Warnings & Warranty

HOW TO USE THE LORI:

1. Apply an appropriate amount of skincare product to the area of your face you wish to treat.
2. Select the appropriate mode with the heat/cold and EMS/RF combination of your choice (refer to the next page for tips and expected benefits.)
3. After use, wipe the dual treatment massage head with a damp cloth or paper towel.

WHERE TO USE LORI:

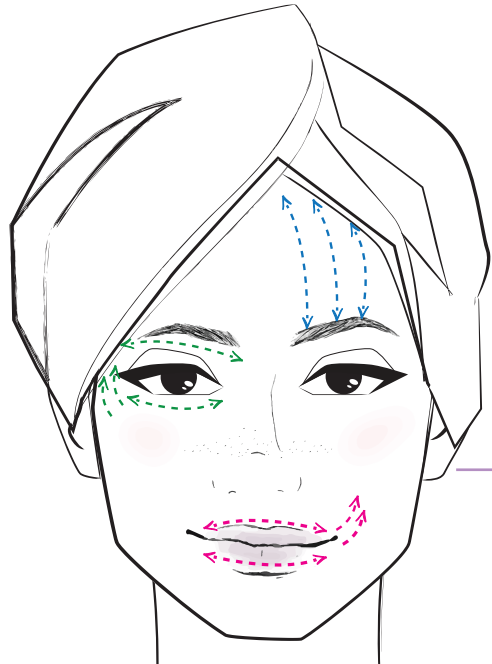

ADDITIONAL TIPS ON HOW TO USE THE LORI:

- LORI is designed to rejuvenate and promote blood circulation to the smaller contours of the face; revitalizing the delicate skin around the eyes and lips, and on your forehead.
- EMS (Electric Muscle Stimulation) vibrations stimulate collagen regeneration and skincare ingredient absorption.
- Thermal Heat promotes blood circulation to reduce eye bags and dark circles.
- For hygiene reasons, we do not recommend sharing LORI with other people.
- Exercise caution when using LORI around the eyes, ensuring that it does not touch the eye balls.
- Remember to turn off the device after use. Store LORI in a clean, dry place.

EYE AREA: Slowly move LORI massage head around your eye and in an upward motion at the outside end of the eye and crows feet.

LIPS AREA: Slowly move LORI massage head around your mouth and in an upward motion at the outside crease of the lips.

FOREHEAD: Slowly move LORI massage head in an upward/downward direction from your eyebrows to the top of your forehead.

BASIC SAFETY PRECAUTIONS.

When using electrical products, especially when children are present, basic safety precautions should always be taken. IF YOU THINK THE POWER IS OFF WHEN THE SWITCH IS OFF. YOU ARE WRONG. If the plug is in, the power is on. Unplug electrical appliances when you are not using them.

DANGERS.

To reduce the risk of electrocution: DO NOT get the Product or Charging Cable wet or immerse in water. Keep this product away from open flame or flammable aerosol products. The Device is not suitable for use in the presence of a flammable anesthetic mixture with air or oxygen or nitrous oxide.

WARNINGS.

To reduce the risk of burns, electrocution, fire or injury to persons: Close supervision is necessary when this Product is used by, on, or near children or physically challenged persons. Use this product only for its intended use as described in this User Guide. DO NOT use this product where oxygen is being administered. The product has no serviceable parts. Never attempt to open the Device Handle for any reason. DO NOT disassemble, heat above 600°C (140°F) or incinerate.

IMPORTANT CAUTIONS.

This is not a medical Device and is not intended to treat medical conditions. If you are under a doctor's care for any skin condition, if you think you may have a skin condition or if you have very sensitive skin, consult your physician before using this product. Keep away from internal medical devices. This Device is electrical and contains small parts. Keep out of reach of children. THIS PRODUCT IS NOT DESIGNED OR INTENDED PRIMARILY FOR CHILDREN 12 YEARS OF AGE OR YOUNGER.

If you have any medical concerns about using the LORI, please consult your dermatologist before use.

If any of the components are damaged (LORI Device Handle, or USB Charging Cord) DO NOT use. DO NOT attempt to open or services components.

Email customer service at
customerservice@spasciences.com

© 2022 Spa Sciences. All rights reserved.

Warranty:

Your warranty is not effective unless you register the Product within 90 days of purchase. To register go online to www.warrantyspasciences.com and complete the information required for registration.

For a period of one (1) year from the date of purchase, Spa Sciences LP (the "Company") warrants the LORI - LED Sonic Thermal Eye & Lip Massager (the "Product") against damage, defects or failure resulting from ordinary consumer use under normal home conditions according to the instructions included with the Product (a "Covered Defect"). This warranty DOES NOT include damage, defects or failure caused by misuse or improper care of the Product; impact or insertion of foreign objects into or on the Product; application of harmful chemicals or abrasives to the Product; failure to follow manufacturer's instructions; normal wear and tear; alteration or overuse of the Product; or any professional or commercial use of the Product. If there is a Covered Defect during the one (1) year warranty period the Company agrees to exchange the Product with a new equivalent Product. This is the sole remedy available for breach of this limited warranty. Any claim for breach of this limited warranty must be commenced within one (1) year following the original purchase. This limited warranty applies only to a purchase of the Product from or an authorized reseller, is valid for the original end user of the Product only and cannot be transferred under any circumstances.

THE PROVISIONS OF THE FOREGOING WARRANTIES ARE IN LIEU OF ANY OTHER WARRANTY, WHETHER EXPRESS OR IMPLIED, WRITTEN OR ORAL (INCLUDING ANY WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE). IN NO EVENT SHALL THE COMPANY BE LIABLE FOR ANY SPECIAL, INCIDENTAL, CONSEQUENTIAL, PUNITIVE, OR INDIRECT DAMAGES BASED UPON ANY LEGAL THEORY. SOME STATES DO NOT ALLOW LIMITS ON WARRANTIES, OR ON REMEDIES FOR BREACH IN CERTAIN TRANSACTIONS. IF APPLICABLE LAW DOES NOT PERMIT THE DISCLAIMER OF IMPLIED WARRANTIES, THEN THE DURATION OF SUCH IMPLIED WARRANTIES, ARE EXPRESSLY LIMITED TO THE DURATION OF THE EXPRESS LIMITED WARRANTY SET FORTH HEREIN.

Additional Info

Entire operating instructions including safety, care and warranty information is available at www.spasciences.com/products/lori or hover camera over the QR CODE. NO APP NEEDED

