
A Closer Look





Function Guide

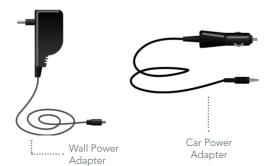

ON/OFF

Quick and easy to turn on and off. Automatic timer turns off after 15 minutes.



Easily adjust speed of kneading according to your requirements.



DIRECTION

Easily adjust the direction of the kneading movement to target key areas of tension



HEAT MODE

Easily toggle heat function on and off.

TESA is therapeutic all-body massage wrap, that features an accessible control panel and optional heat to relax and relieve tense muscles.

POWERING FROM AC OUTLET

For use inside a home or office, use the power adapter. Insert the DC adapter into the round power jack on the side of the device, and plug the AC adapter into a 120V outlet on the wall.

POWERING FROM VEHICLE

For use inside a vehicle, use the car power adapter. Insert the DC adapter into the round power jack on the side of the device, and plug the car power adapter into the car charger socket or lighter outlet.

CONTROL BUTTONS

The Control Buttons can be found on the right handle of the device. Press the ON/OFF button once to power the device on or off.


AUTO SHUT OFF

The device will automatically shut off after 15 minutes of continuous use.

Usage Guide

- 1. Press the Power Button once to turn the device ON, the massaging modes will be in continuous mode.
- 2. Place the massager on the targeted massage area and use the handle straps to hold the massager in place.
- 3. Press the Speed Button to adjust between the three (3) speeds of the massage modes.
- 4. Change the direction of movement of the massage modes by pressing the Direction Button once, change the direction back by pressing the Direction Button again.
- 5. Optionally, press the Heat Button to turn on the gentle heat. The modes will shine red light when the heat is on. Press the heat button again to turn the heat off.
- 6. Use the handle straps to vary pressure according to your needs. If the kneading is too rough, pull less at the straps.

MASSAGE AREAS

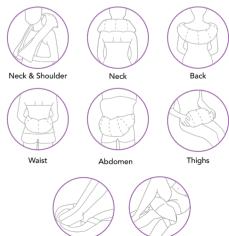





Usage Guide

Warnings

Warnings & Warranty

Problem	Solution
Unable to start the massager	The massager will automatically finish after 15 minutes. Check if the power plug has pulled out. Check for any bad connections and cease use immediately if any are found. If one of the above work, turn off the main power supply and wait for one hour so that the product can cool down. The massager has an automatic overheating protection feature.
Suddenly stops working during use	

IMPORTANT SAFETY INSTRUCTIONS

The following basic safety guidelines must always be followed when using electrical appliances to reduce the risk of fire, electric shock and/or injury to any person:

WARNING

- 1. Read all instructions & keep all instructions for reference.
- 2. Not suitable for children
- 3. Do not operate the device while driving.
- 4. Do not use the device if you are pregnant.
- Do not use if you have any of the following health conditions: Heat sensitivity, a circulation disorder, varicose vein, bruises, broken skin, skin irritation, muscle or tissue inflammation, vein inflammation, injury, or any open wound.
- 6. If you experience any pain or discomfort during the massage, please stop use immediately and consult your doctor.
- 7. Do not use if you have any of the following: Pacemaker, artificial joints, or electronic implants.
- 8. Do not use this product if you have applied any oils or creams to your skin.
- 9. Power OFF device and disconnect from the power source after use.

- 10. Take precautions to make sure you do not fall asleep during use.
- 11. Do not carry this product by the power cord.
- 12. Never clean this product unless it is unplugged. Use a dampened cloth to wipe the device and dry immediately.
- 13. Do not plug the device into an outlet until it is completely dry.
- 14. Do not operate under a blanket or pillow. Excessive heating can occur or cause fire, electric shock, or injury.
- 15. Never use this product anywhere near water or high humidity.
- 16. Never operate with wet hands. Always make sure to be completely dry. Do not immerse or expose the device to liquid of any type.
- 17. This appliance is suitable for personal household use only. Commercial or professional use will void the warranty. The supplier cannot be held responsible for injury or damage if the device has been used for anything other than its intended use.
- 18. Keep device out of reach of children.
- 19. Do not leave the device unattended while plugged in.
- 20. Do not allow the cord to rest in any position that may cause trips or pull hazards.
- 21. Do not allow the power cord near hot surfaces.
- 22. Do not use the device if it has been dropped or damaged in any way.
- Check the cord and the plug regularly for any damage. If damage is present, it must be repaired by the manufacturer or qualified professional.
- 24. Use only original parts to repair or maintain the device.
- 25. Unplug and stop using the device immediately if it begins to malfunction. Never try to inspect or open the device.
- 26. Avoid using before bed. The device may have a stimulating effect and may delay sleep.
- 27. Keep all packaging provided with this device away from the reach of children.

BASIC SAFETY PRECAUTIONS.

When using electrical products, especially when children are present, basic safety precautions should always be taken. IF YOU THINK THE POWER IS OFF WHEN THE SWITCH IS OFF. YOU ARE WRONG. If the plug is in, the power is on. Unplug electrical appliances when you are not using them.

DANGERS

To reduce the risk of electrocution: DO NOT get the Product or Charging Cable wet or immerse in water. Keep this product away from open flame or flammable aerosol products. The Device is not suitable for use in the presence of a flammable anesthetic mixture with air or oxygen or nitrous oxide.

WARNINGS.

To reduce the risk of burns, electrocution, fire or injury to persons: Close supervision is necessary when this product is used by, on, or near children or physically challenged persons. Use this product only for its intended use as described in this User Guide. DO NOT use this product where oxygen is being administered. The product has no serviceable parts. Never attempt to open the Device Handle for any reason. DO NOT disassemble, heat above 600°C (1400°F) or incinerate.

IMPORTANT CAUTIONS.

This is not a medical device and is not intended to treat medical conditions. If you are under a doctor's care for any skin condition, if you think you may have a skin condition or if you have very sensitive skin, consult your physician before using this product. This product uses magnets which may affect internal medical devices such as pacemakers. Keep away from internal medical devices. This Device is electrical and contains small parts. Keep out of reach of children. THIS PRODUCT IS NOT DESIGNED OR INTENDED PRIMARILY FOR CHILDREN 12 YEARS OF AGE OR YOUNGER.

If you have any medical concerns about using TESA please consult your dermatologist before use.

If any of the components are damaged (TESA Device, handles or Adapter Cord) DO NOT use. DO NOT attempt to open or service components.

Email Customer Service at customerservice@spasciences.com

© 2021 Spa Sciences. All rights reserved.

Warranty:

Your warranty is not effective unless you register the Product within 90 days of purchase. To register go online to www.warrantyspasciences. com and complete the information required for registration.

For a period of TWO (2) years from the date of purchase, Spa Sciences LP (the "Company") warrants the TESA 3D Shiatsu Kneading Massager (the "Product") against damage, defects or failure resulting from ordinary consumer use under normal home conditions according to the instructions included with the Product (a "Covered Defect"). This warranty DOES NOT include damage, defects or failure caused by misuse or improper care of the Product; impact or insertion of foreign objects into or on the Product; application of harmful chemicals or abrasives to the Product; failure to follow manufacturer's instructions; normal wear and tear; alteration or overuse of the Product; or any professional or commercial use of the Product.

If there is a Covered Defect during the two (2) years warranty period, the Company agrees to exchange the Product with a new equivalent Product. This is the sole remedy available for breach of this limited warranty. Any claim for breach of this limited warranty must be commenced within two (2) years following the original purchase. This limited warranty applies only to a purchase of the Product from an authorized reseller, is valid for the original end user of the Product only and cannot be transferred under any circumstances.

THE PROVISIONS OF THE FOREGOING WARRANTIES ARE IN LIEU OF ANY OTHER WARRANTY, WHETHER EXPRESS OR IMPLIED, WRITTEN OR ORAL (INCLUDING ANY WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE). IN NO EVENT SHALL THE COMPANY BE LIABLE FOR ANY SPECIAL, INCIDENTAL, CONSEQUENTIAL, PUNITIVE, OR INDIRECT DAMAGES BASED UPON ANY LEGAL THEORY. SOME STATES DO NOT ALLOW LIMITS ON WARRANTIES, OR ON REMEDIES FOR BREACH IN CERTAIN TRANSACTIONS. IF APPLICABLE LAW DOES NOT PERMIT THE DISCLAIMER OF IMPLIED WARRANTIES, THEN THE DURATION OF SUCH IMPLIED WARRANTIES, ARE EXPRESSLY LIMITED TO THE DURATION OF THE EXPRESS LIMITED WARRANTY SET FORTH HEFEIN

Additional Info

Entire operating instructions including safety, care and warranty information is available at www.spasciences.com/products/tesa or hover camera over the QR CODE.

NO APP NEEDED

