

The NuLifeBeauty Light Therapy Mask provides you with a salon-quality light treatment at the comfort of your own home.

STEP 1

Before using the mask, prepare your skin by washing or cleansing your face.



STEP 2

Plug in the mask (with neck attachment and/or microcurrent plugs if desired) to the controller using the USB cables. Then plug in the power cord and secure the mask on using the strap.



STEP 3

Select your desired color and intensity level using the remote control. Adjust the timer to your liking and begin your treatment session.

NOTE: Selecting "Color 8" will alternate between all the seven different colors. We recommend starting off with a 10 minutes session and increasing as you feel comfortable.



STEP 4

After your session, turn off and remove the mask. Then apply your moisturizer and/or serum. Once the timer runs out, the mask controller will beep to notify you that your session is complete. We recommend using it for 4-6 weeks and then continued maintenance thereafter for all natural, radiant-looking skin!



STEP 5

Gently clean treatment surface area as needed with a damp cloth or a nonirritating water-based wet wipe. For deeper cleaning you can wipe gently with 70% alcohol. Store your mask in a safe place.




Any questions?

Reach out to our support team and they will get back to you usually within 1 business day or visit <https://nulifebeauty.com/pages/faq>

help@nulifebeauty.com

Device Parts & Description

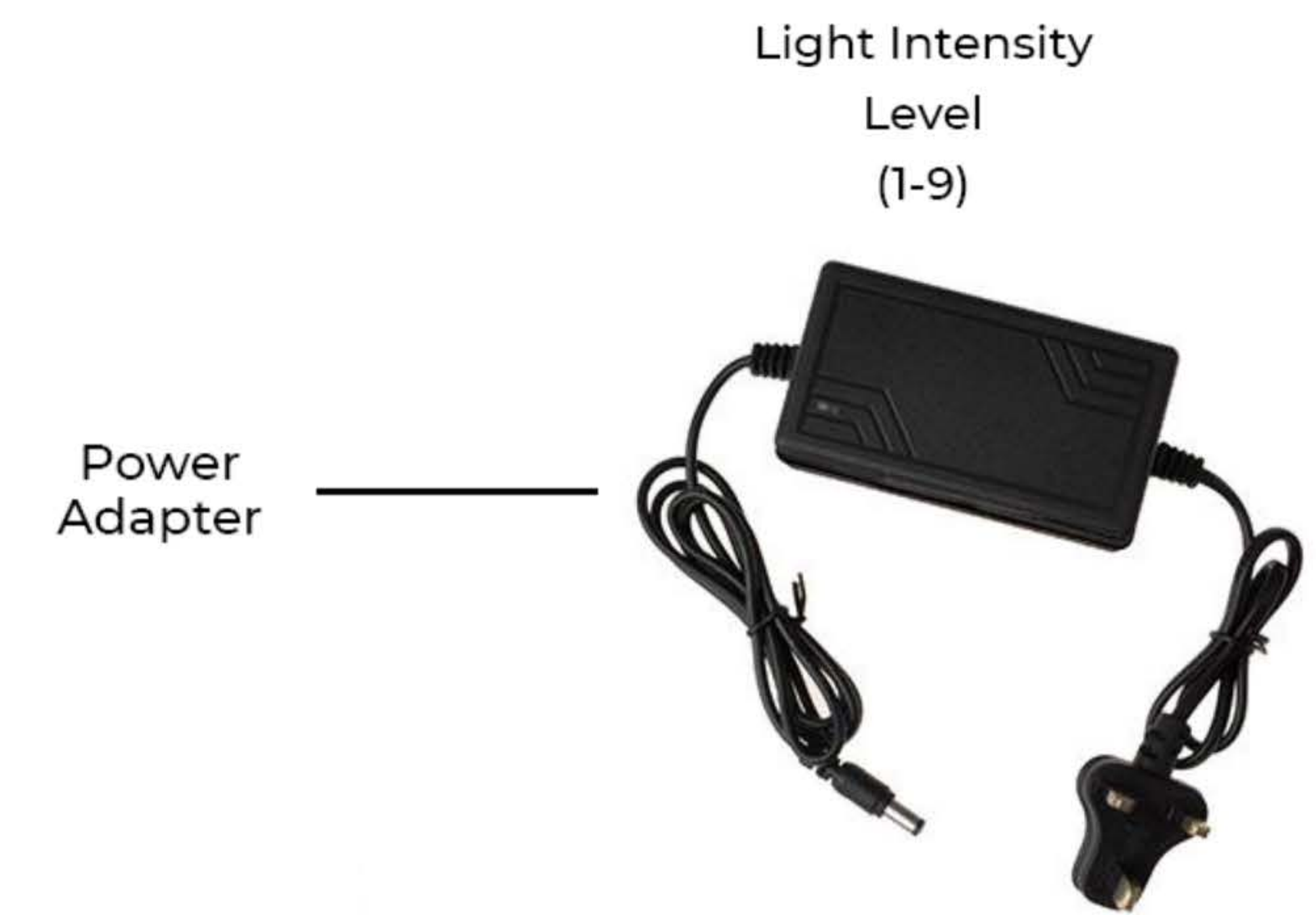

Light Intensity Level (1-9)