


waterpik® waterflosser®


QUICK START GUIDE

WP-100

1 Unit Prep

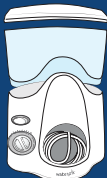

A. Remove the reservoir from base and fill with lukewarm water. Replace the reservoir firmly on the base.

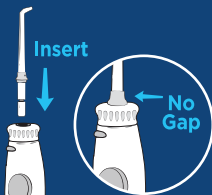


B. Press tip firmly into handle. There should be no gap between handle and colored tip ring.



C. Turn pressure to high, point tip into sink, turn unit **ON** until water flows. Turn unit **OFF**.

2 Water Floss

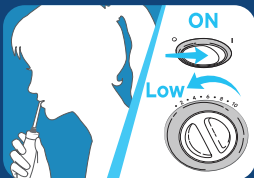

A. Start on low pressure and increase to your liking. With unit **OFF**, place tip in mouth. Turn power **ON**.



B. Close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink.



C. Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.

This device is only compatible with a voltage system of 220-230VAC, 50Hz.
Do not plug this device into a voltage system that is different.