

waterpik® waterflosser

QUICK START GUIDE

WP-660


1 Unit Prep



A. Fill reservoir with warm water and replace it firmly on the base.




B. Insert tip until it clicks into place.

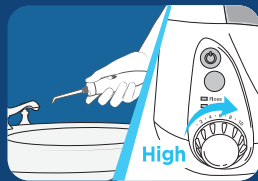

C. Turn pressure to high, point tip into sink, turn power and water switches ON until water flows. Turn switches OFF.


2 Water Floss

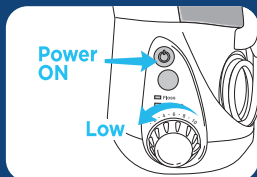

A. Start on low pressure and increase to your liking. Turn power ON.



B. Turn water ON. Close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink.

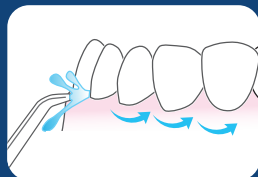

C. Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.

This device is only compatible with a voltage system of 230VAC, 50Hz. Do not plug this device into a voltage system that is different.