

QUICK START GUIDE

1 Prepare the Reservoir:




Remove the reservoir from base and fill with lukewarm water.

Replace the reservoir onto the base, **PRESSING DOWN FIRMLY**

Replace



2 Select and Insert Tip:

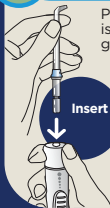

Press tip firmly into handle. When tip is properly inserted there should be no gap between control knob and tip.



No Gap




3 Prime Unit (First Time Use Only):

Prior to Very First Use of Unit: Set pressure control slide switch on handle to **high**, point the tip into sink. Turn unit **ON** until water flows from the tip, turn unit **OFF**.



4 Adjust Pressure and Pause Control:

For first time users, adjust pressure control to **lowest position**. Gradually increase pressure as you prefer.

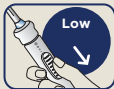

You can temporarily stop water flow by pressing the pause button on the handle.



continued on back

This device is only compatible with a voltage system of 220-230VAC, 50Hz.
Do not plug this device into a voltage system that is different.

5

On/Off Procedure:

Always place tip in mouth before turning unit on.

Lean over sink and close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink. Turn unit **ON**.

Always **PAUSE** or turn unit **OFF** before removing tip from mouth.



6

Technique Tips for Best Results:

Aim the tip just above the gumline at a 90 degree angle.

Pointer l'embout juste au-dessus de la ligne gingivale à un angle de 90 degrés.

Starting with the back teeth, follow the gumline and pause briefly between teeth.

Continue until all areas around and between teeth have been cleaned.



7

Tip Removal:

To remove tip, simultaneously hold control knob down and lift tip up.

