


waterpik®

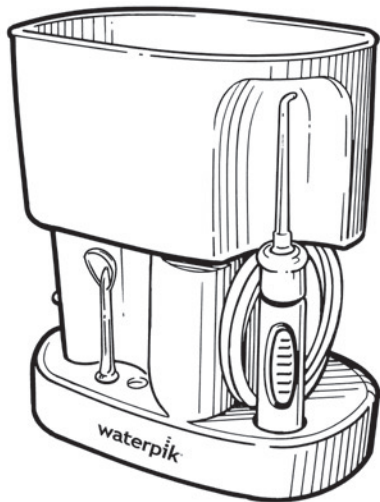

**Waterpik® Water Flosser
Model WP-70**

.....
www.waterpik.com


WELCOME

Congratulations! By selecting the Waterpik® Water Flosser you are now on your way to better gum health, and a healthier smile!


With a clinically proven combination of water pressure and pulsations, the Waterpik® Water Flosser removes harmful bacteria deep between teeth and below the gumline – where brushing alone won't reach.

Thank you for choosing the Waterpik® Water Flosser.

Traditional Oral Care

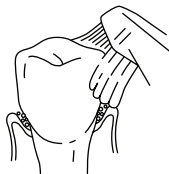


Traditional brushing can't easily reach hidden food debris and harmful bacteria.

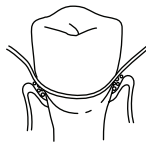

Flossing has limited access below the gumline.

Oral Care With The Waterpik® Water Flosser

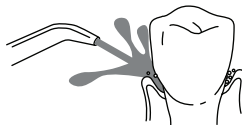

The Waterpik® Water Flosser removes trapped food debris and harmful bacteria.

TABLE OF CONTENTS

Important Safeguards	4
Product Description	6
Getting Started	8
Care of Your Waterpik® Water Flosser	14
Troubleshooting Guide	15
Limited One-Year Warranty	16

IMPORTANT SAFEGUARDS

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Always unplug product after using.
- Do not handle plug with wet hands.
- Do not immerse in water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the power supply cord for damage before the first use and during the life of the product.

WARNING:


To reduce the risk of burns, electric shock, fire, serious injury or damage to tissue:


- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact Water Pik, Inc. or your local distributor which can be found at www.waterpik.com.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas. See operating instructions for correct usage.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use tips and accessories recommended by Water Pik, Inc.

- Fill reservoir with water or other dental professional recommended solution only.
- Do not drop or insert any foreign object into any opening or hose.
- Keep electrical cord away from heated surfaces.
- Do not operate where oxygen or aerosol sprays are being used.
- Do not use mouthwash/rinses that contain iodine, bleach, or tea tree oil.
- Do not use while wearing any oral jewelry. Remove prior to use.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- For household use only. Do not use outdoors.
- Closely instruct and supervise children and individuals with special needs in the proper use of this product.
- Do not remove the Pik Pocket™ flexible tip from the clear shaft. If it separates for any reason, discard both the tip and shaft and replace with a new Pik Pocket™ Tip. Pik Pocket™ Tip not included in all models.
- Do not use this product for more than 5 minutes in each two-hour period. Longer usage can overheat the motor and cause product failure as well as voiding the warranty.
- Children should be supervised to ensure that they do not play with the appliance.
- This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug is intended to fit into a polarized outlet only one way. If the plug does not fit fully into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to modify the plug in any way.

SAVE THESE INSTRUCTIONS.

PRODUCT DESCRIPTION

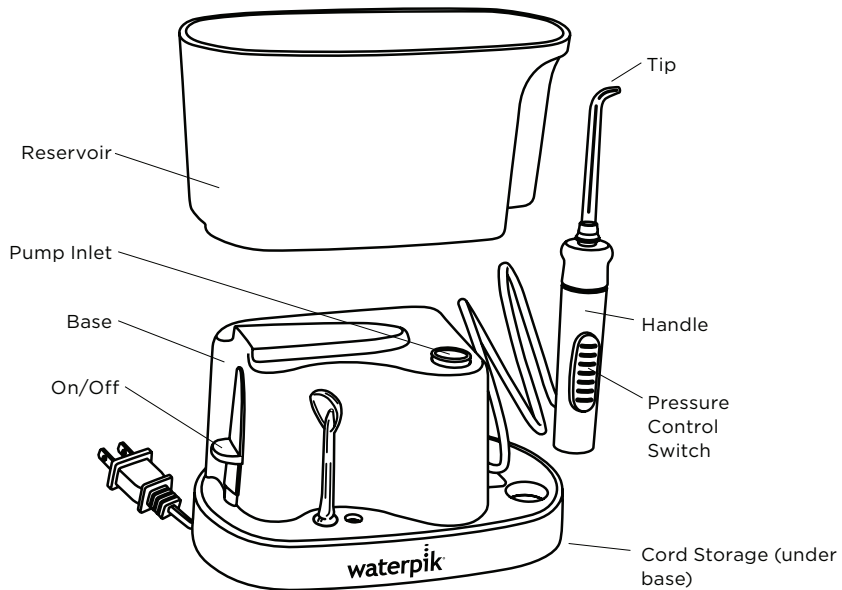


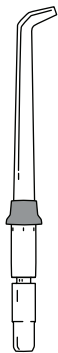

Classic
Jet Tip



Tongue
Cleaner Tip

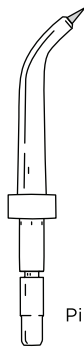

Pik Pocket™ Tip


NOTE: all tips not included in all models

GETTING STARTED

Plug the cord into an electrical outlet. If the outlet is controlled by a wall switch, be sure it is turned on.


Preparing the Reservoir

Remove reservoir from base. Fill the reservoir with lukewarm water. Replace the reservoir over the base and press down firmly.

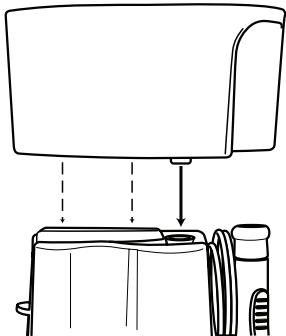


Inserting and Removing Tips

Pull the knob on the top of the handle toward the hose and insert tip into the center of the knob. The colored ring will be flush with the end of the knob if the tip is correctly locked into place. To remove tip from the handle, pull the knob on the top of the handle toward the hose and pull tip from handle.

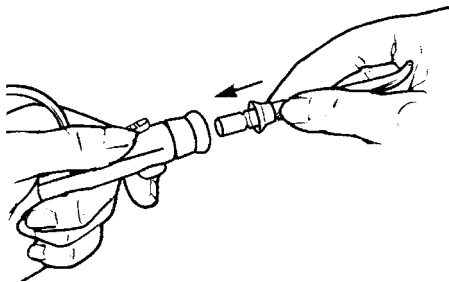

Adjusting the Pressure Setting


Slide the pressure control switch located on the handle to the lowest setting (away from the tip as far as possible) for first-time use. Gradually increase pressure over time to the setting you prefer – or as instructed by your dental professional.

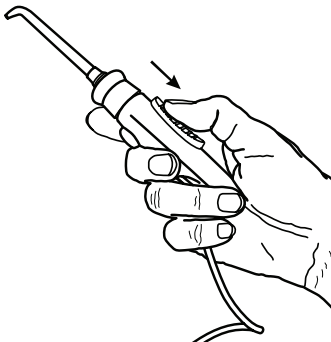


Placement in Mouth

Lean low over sink and place the jet tip in mouth. Aim the tip toward teeth.

With your free hand, slide the on/off switch on the base of the unit to the “up” position. Set pressure at desired setting.


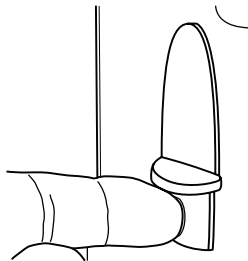


Direct the stream by turning the knob at the top of the handle.

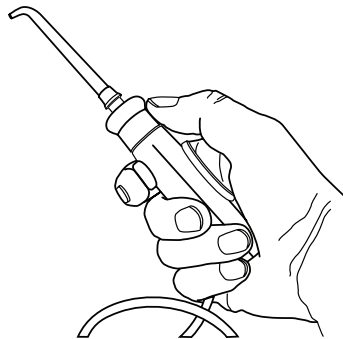


Recommended Technique

Direct the jet stream at a 90-degree angle to your gumline. **Close lips to avoid splashing** but allow water to flow freely from mouth into the sink. For best results, start in the molar area (back teeth) working toward the front teeth. Glide jet tip along gumline and pause briefly between teeth. Continue until all areas around and between the teeth have been cleaned.

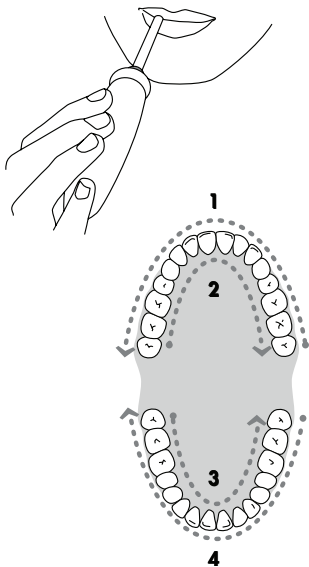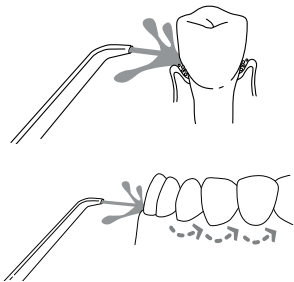

Pause Control

You can temporarily stop the flow any time by pressing the pause button on the handle.

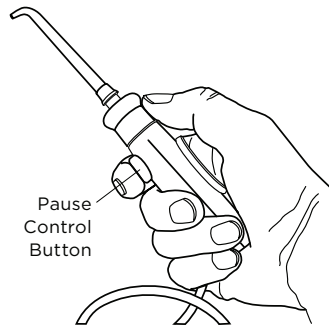

When Finished

Turn the unit off. Empty any liquid left in the reservoir. Unplug the unit.

TIP USAGE

Pik Pocket™ Tip

The Pik Pocket™ Tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets. It is not intended for general use.


To use the Pik Pocket™ Tip, lean low over your sink and place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket. Turn the unit on to the lowest pressure setting and clean by gently following gumline, inserting tip in between teeth. Continue tracing along the gumline until all areas of the mouth are clean.



NOTE: This tip must be used only with the Pressure Control set to the lowest setting. Use on a higher setting will shorten life of tip and handle.

Tongue Cleaner Tip

Lean over the sink. Place the tongue cleaner on your tongue and start the Water Flosser – use the lowest pressure setting, increase as comfortable. Gently place the tongue cleaner in the center/midline of your tongue about half way back and with light pressure pull forward.

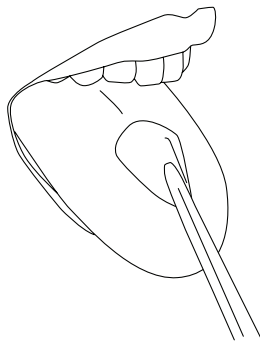

Replacement tips/attachments may be purchased directly from the local distributor.


GETTING STARTED

Using Mouthwash and Other Solutions

Your Waterpik® Water Flosser can be used to deliver mouthwash and antibacterial solutions. Ask your dental professional to recommend the appropriate solution for you. Note: use of some antibacterial solutions could shorten the life of your Waterpik® Water Flosser.

After using any special solution in the reservoir, take a moment after every use to rinse the unit. Partially fill the reservoir with warm tap water. Aim the tip into your sink. Turn on the unit and let it run until empty. This will help keep the tip from clogging and help your unit to last longer.

RECOMMENDED TIP REPLACEMENT INTERVALS

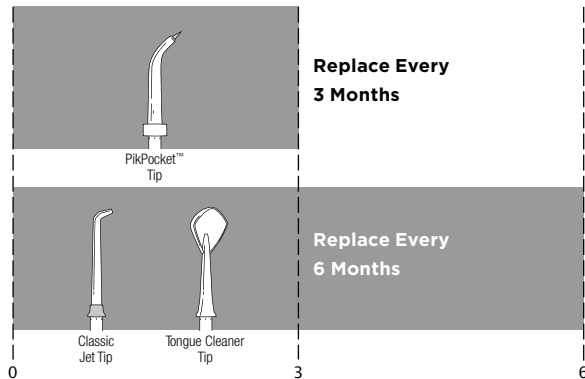

CARE OF YOUR WATERPIK® WATER FLOSSER

Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser.

Before cleaning, unplug from the electrical outlet.

Before exposing the Water Flosser to freezing temperatures, remove the reservoir and run the unit until completely empty.

Service Maintenance

Waterpik® Water Flossers have no consumer-serviceable electrical items and do not require routine service maintenance.

Some replacement parts and accessories for the Water Flosser are also available from Water Pik, Inc. or your local distributor which can be found at www.waterpik.com.

Refer to serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

PROBLEM	CAUSE	CORRECTIVE ACTION
Tips and Tongue Cleaner		
Leakage between the tip and handle knob.	1) Tip not fully engaged.	1) Remove tip and reinstall.
Pik Pocket™ Tip		
Tear in soft tip.	1) Tip is worn.	1) Replace with new tip.
Water Flosser Base		
Hose to the handle leaks.	1) Damage to the hose.	1) Replace the Handle & Hose kit. Contact your local distributor for more information.
Reservoir leaks.	1) Missing reservoir valve.	1) Replace reservoir valve. Contact your local distributor for more information.
	2) Reservoir valve is upside down.	2) Turn valve right-side up.
Inadequate pressure.	1) Reservoir not positioned properly.	1) Re-position reservoir. Press down firmly on the reservoir.
	2) Tip is clogged with debris.	2) Flush tip under faucet or replace tip.
Unit doesn't start.	1) Motor not running.	1) Be sure outlet is functioning; test with another appliance (check electrical reset button and/or wall switch).

Still have questions? Go to www.waterpik.com. If none of these corrective actions are successful, call the Water Pik, Inc. Customer Service Center. Toll-Free Number: 1-800-525-2774 (USA) or 1-888-226-3042 (Canada).

LIMITED ONE-YEAR WARRANTY

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 1 year from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will repair or replace at its discretion any part of the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as tips, etc.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.