




HOOGA WAND

USER MANUAL

What's Included:

- 1 Hooga Wand
- 1 USB Type C Charging Cable
- 1 User Manual

Product Introduction:

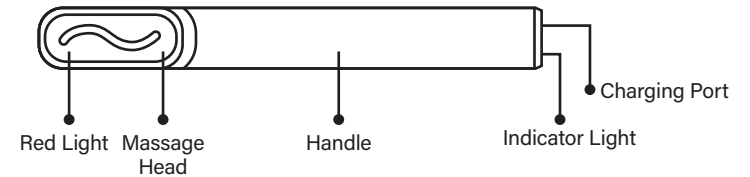

Product Features

- 1. Microcurrent:** Microcurrent stimulates the muscles and deeper layers of the skin, helping reduce signs of aging and promoting smoother, healthier skin.
- 2. Light Therapy:** Red light deeply rejuvenates the skin and helps combat signs of aging. Red light also stimulates collagen and elastin production for smoother, younger looking skin.
- 3. Facial Massage:** Provides a gentle massage, helping decrease the appearance of puffiness and dark circles.
- 4. Therapeutic Warmth:** Helps increase the absorption of skin care products, reduce redness, and support healthy, glowing skin.



Smart Touch Function

Our Hooga Wand is equipped with a smart sensor chip so that when you hold the handle of the wand and the massage head touches your skin, the wand automatically turns on and begins to work and stops automatically when the massage head leaves your skin. This results in easy to use functionality as well as optimal power savings, reducing the frequency with which you need to charge the wand.

Note: If the wand isn't working when you are holding the handle and touching the massage head to your face, it is likely that your skin is too dry. The wand is meant to be used on well hydrated skin and works best in conjunction with a hyaluronic acid or aloe vera infused serum.

Rotating Massage Head

Our Hooga Wand is designed to be used in multiple different positions, with the massage head in either the vertical position or horizontal T-Bar position. The massage head rotates 90 degrees, allowing for customized treatments that can reach any area of your face.

In the vertical position, the massage head can easily treat small corners or otherwise hard to reach areas. The vertical position is best for treating under the eyes.

In the horizontal T-Bar position, the massage head can treat larger areas such as the forehead, cheeks, or neck.

How To Use?

1. Wash and dry your face so you start with clean, hydrated skin.
2. Apply a priming serum of your choice infused with hyaluronic acid and/or aloe vera to prep your skin for treatment.
3. Glide the Wand across your face in an upward and outward moment, slowly covering your forehead, cheeks, under eyes, jaw, neck and upper lip. The Wand will automatically turn on when making contact with hydrated, clean skin. If the wand does not turn on, either your skin needs to be more hydrated or the Wand is low on charge.
4. Use 3 to 5 times per week for at least 5 minutes per day.
5. We recommend using the bottom-up method to achieve a lifting effect.



Eye Area



Face



Neck

Note: When the wand is automatically activated by holding the handle and touching the massage head to your skin, the microcurrent, light therapy, facial massage and therapeutic warmth will all be activated at the same time. It may take 30 seconds or so of continued use until the warmth is felt.

Benefits



Fade blemishes & smooth fine lines



Amplify serums & moisturizers



Boost skin's glow



Reveal a youthful radiance

Charging & Battery Life

To charge: Insert the USB Type C end of the charging cable gently into the charge port at the bottom of the wand. Connect the other end of the charging cable into a 5V/1A power adapter or a computer USB port. When successfully charging, the indicator light will be flashing. For battery health, we don't recommend leaving the wand charging for more than 24 hours. The Hooga Wand cannot be used while charging.

- When the battery level is low, the indicator light will flash quickly.
- When charging, the indicator light will flash.
- When fully charged, the indicator light will remain illuminated.

Charge time from zero to full battery is approximately 100 minutes. The battery life is approximately 60 minutes of continuous use. This equates to 12 five minute treatments.

Precautions & FAQ

- Do not disassemble or modify this product.
- Store in a cool, dry place out of reach of children.
- We recommend testing the product on your hand before use on the face to ensure there is no negative reaction. Remember the skin on your hand or other test area needs to be hydrated for the massage head to turn on properly.
- If you find any abnormality or variability such as unusual vibration during use, please check whether the battery is low (indicator light flashing quickly).

How do I clean the wand?

Clean the Hooga Wand by spraying a dry cloth with water or disinfecting spray and wiping down the massage head. Do not spray the wand directly or submerge your wand in any way.

Who should not use the wand?

Out of an abundance of caution, we don't recommend the wand for anyone who is pregnant, under the age of 18, has an active or recent cancer diagnosis, or a pacemaker.

Is the wand waterproof?

No. Do not submerge the wand in water and do not use it in the shower.

What is a recommended serum to use with the wand?

We recommend incorporating products that contain conductive ingredients (e.g. aloe vera or hyaluronic acid) to fully benefit from the microcurrent technology. If you use products without conductive ingredients, you will still benefit from the wand's other features. Please note that for the massage head to be activated when you touch it to your skin, your skin needs to be well hydrated.

Product Specifications

Model Name: Hooga Wand
Product Material: Aluminum Alloy/PC
Power: 1.5W
Voltage: 3.7V
Battery Capacity: 350mAh
Charging Input: 5V/1A
Dimensions: 5.9" x 0.59" x 0.47"