


Red Light Therapy Panel Wheeled Rack Instruction Manual

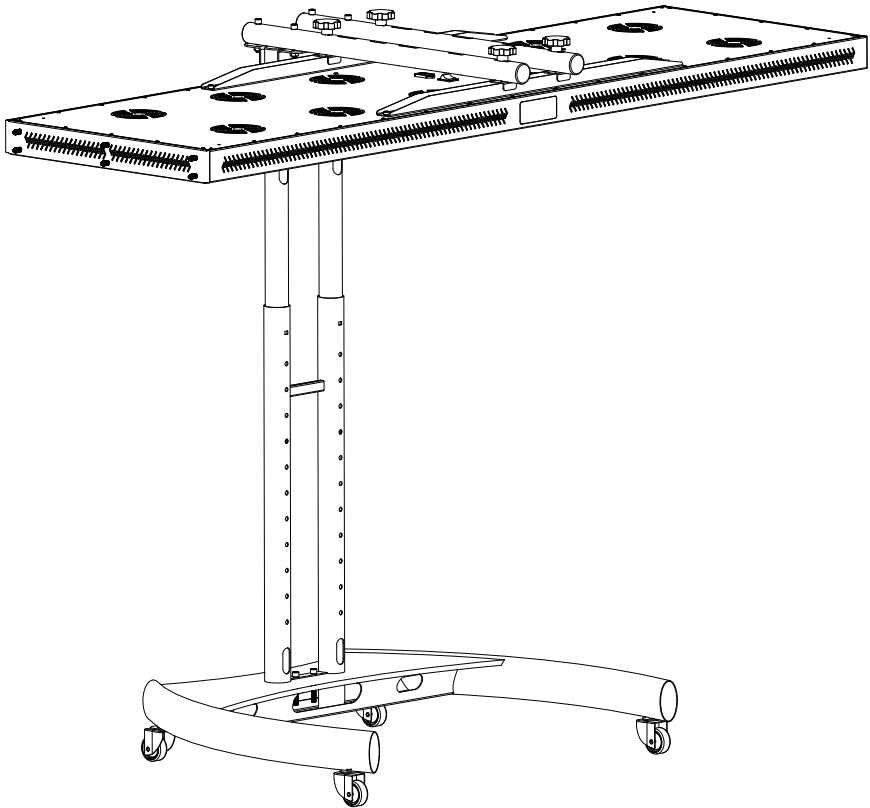

Stand-180

Multi-functional and professional mobile stand

For LED red light therapy panels: 55" - 80"

Maximum load capacity: 200 lbs / 90.9 kg

Shelf max load: 25 lbs | 11.4 kg

Tray max load: 10 lbs | 4.5 kg

Product Catalog





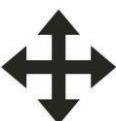


Table of Contents and Warning	2/12
Parts list	3/12
Assembly and Installation	
1. Castor Wheels Assembly To Base Installation	4/12
2. Telescopic Mounting Poles Installation.	5/12
3. Interface Installation	6/12
4. Selection of Panel Mounting Patterns	7/12
5. Product Specifications	8/12
6. Horizontal Stand Assembly Instruction	9/12
7. Screw Horizontal Metal Structure To Panel	10/12
8. Contact us	11/12

WARNING!


Serve personal injury and property damage can result from improper installation or assembly. Read this user manual carefully before beginning.

- If you do not understand the instructions or have any concerns or questions, please contact a qualified local installer (page 11).
 - Do not install or assemble if the product or hardware is damaged or missing, if you require replacement parts, please contact your local distributor for assistance.
 - This product fits most 21"x71"; Maximum weight of red light therapy panel: 60kg; Shelf max load: 10 lbs/4.5kg; Tray max load: 10 lbs/4.5kg.
 - Do not use this product for anything other than the original design.
 - This product contains moving parts, please use it with caution.
 - The manufacturer disclaims any liability for the modifications, improper installation, or installation over the specified weight range. The manufacturer will not be liable for any damage arising out of the use of, or inability to use the product.
- All images are for reference, to prevail in product.

Manual

						
Tightening Direction	Loosening Direction	Phillips Screwdriver	Metre Rule	Adjustment	By Hand	Hex Wrench

Parts List

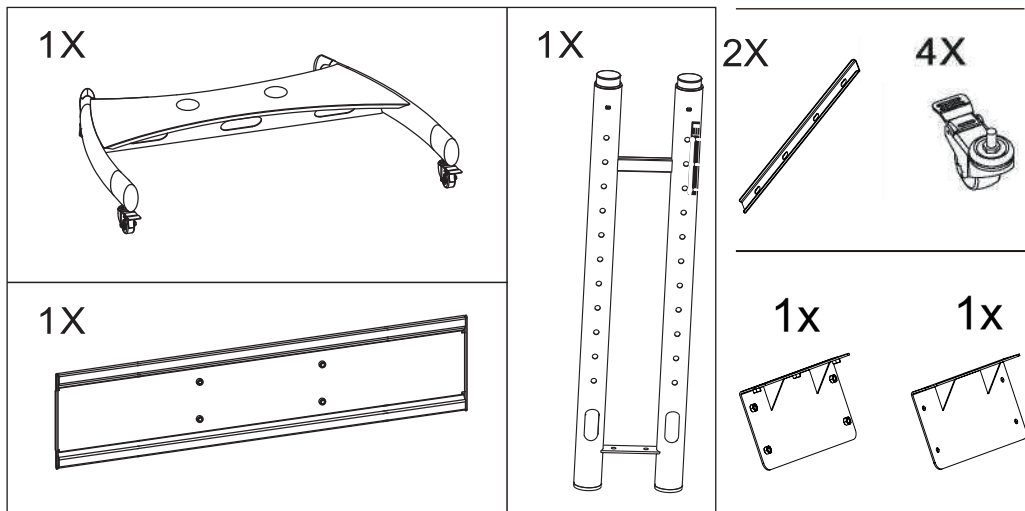

2X




4X



2X



1X

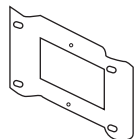

A

6X



M8x20mm

B

8X



M8x60mm

C

4X



M8x50X35mm

D

26X



M6x8mm

E

1X



M6

F

1X



M8

1. Castor Wheels Assembly To Base Installation


Step 1

Step 2

Step 3

Note:

Make sure cables management windows are facing the rear of base.

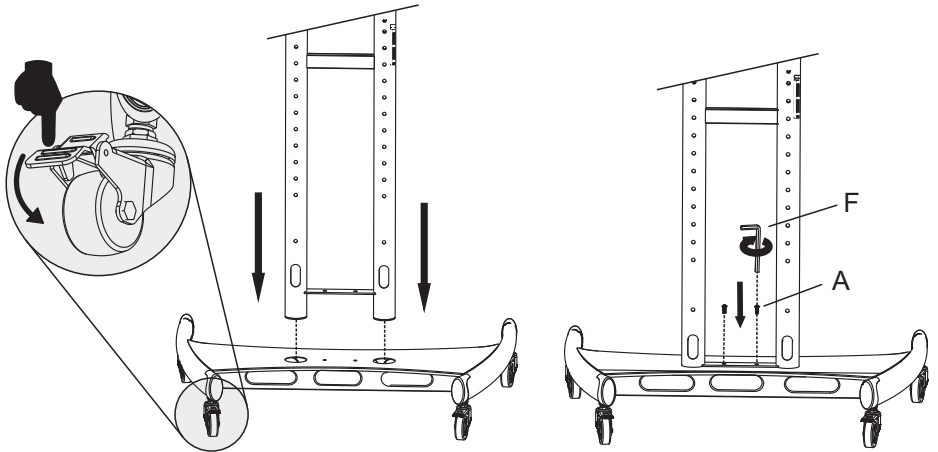


Caution:

Ensure screws are secured firmly.



Step 4

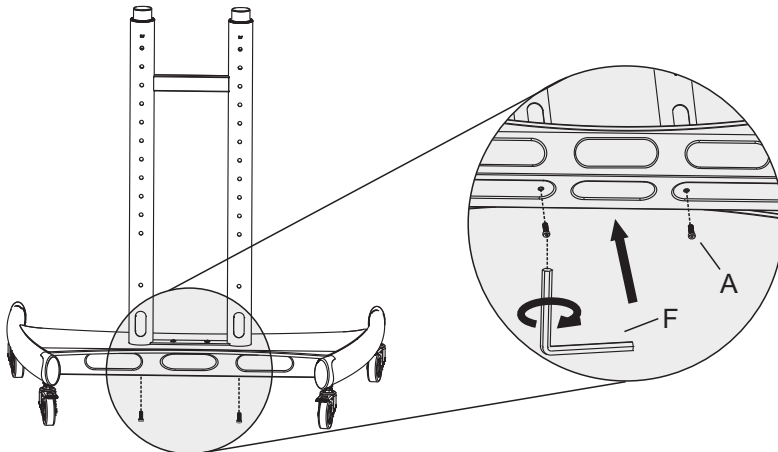

2. Telescopic Mounting Poles Installation.

- Step 1 : Unscrew cap screw on right lower poles assembly by using J (6mm Allen key) from hardware pouch. Put screw aside. As shown in Fig-1.
- Step 2 : Use both hands to grip upper pole right, twist and pull upper pole upward until your desired height. Then twist upper pole again toward plunger holes direction and engaged spring ball into the plunger hole and set. Secure screw firmly as shown in fig-2.
- Step 3 : Repeat step 2 on left poles assembly.

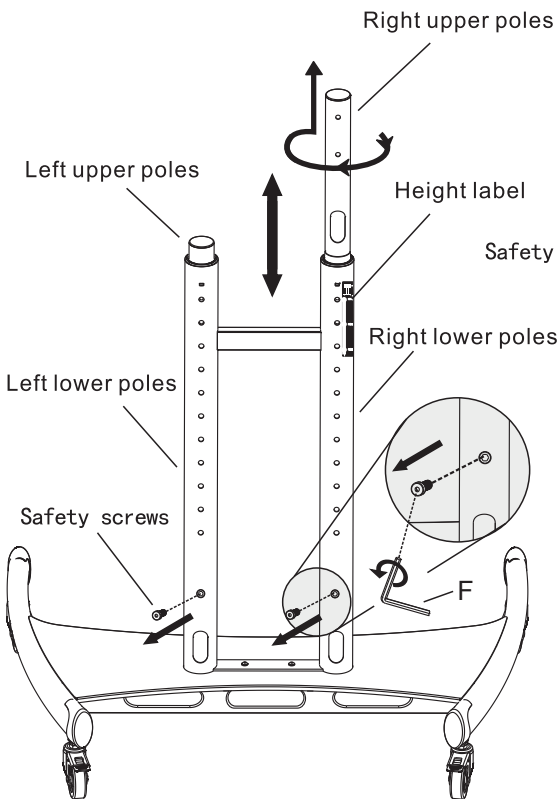


Figure 1


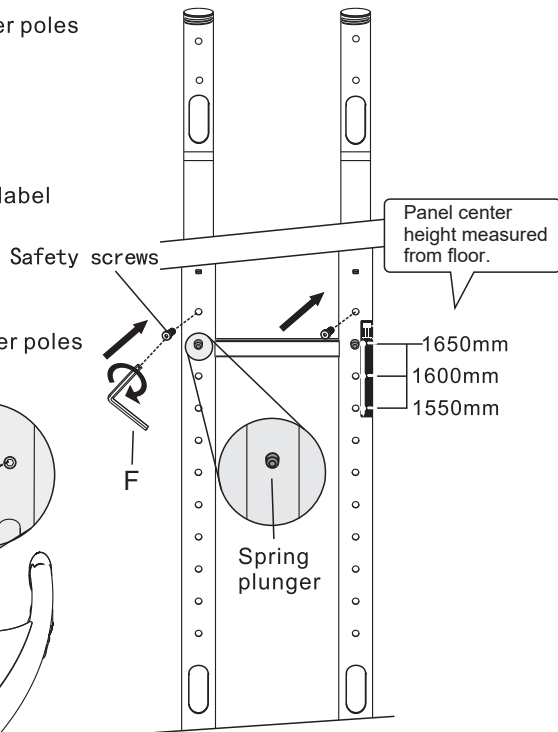


Figure 2

3. Interface Installation

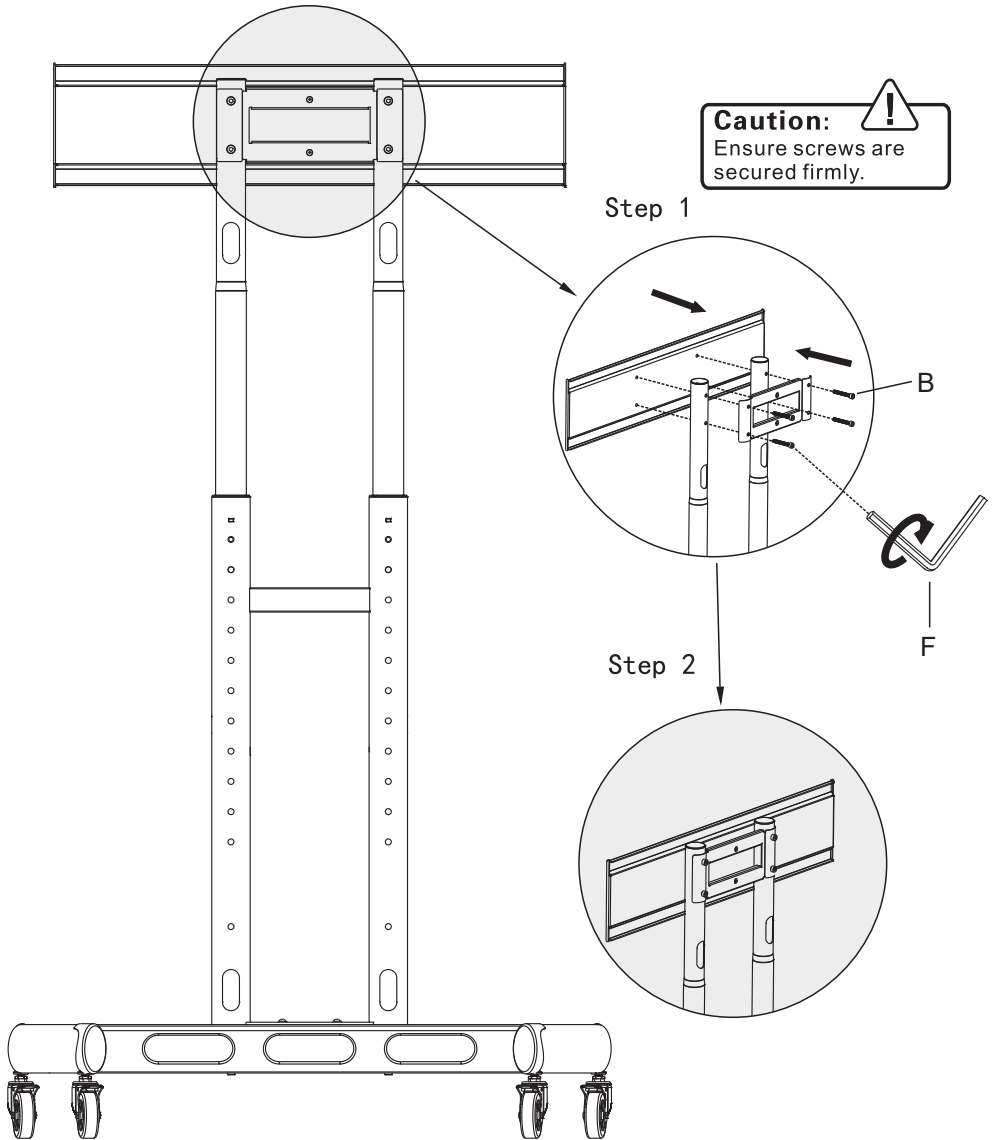

4. Selection of Panel Mounting Patterns

Step 1. Screw the vertical stand mounting kits to panel as following

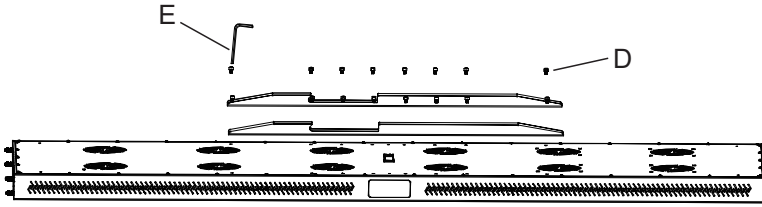


fig-1

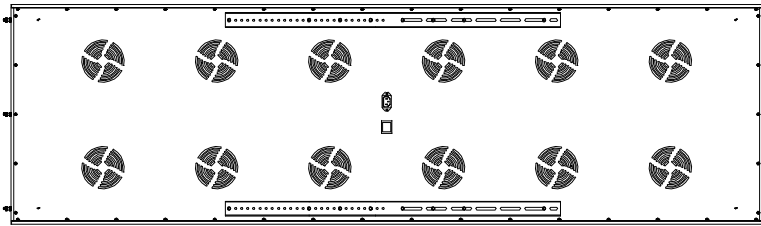


fig-2

Step2. Hang panel (with vertical mounting kit) on the telescopic mounting parts, as following

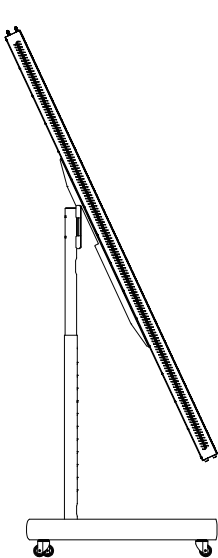


fig-3






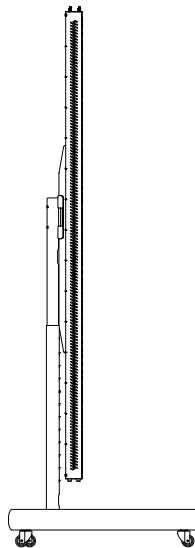



fig-4

5. Product Specifications

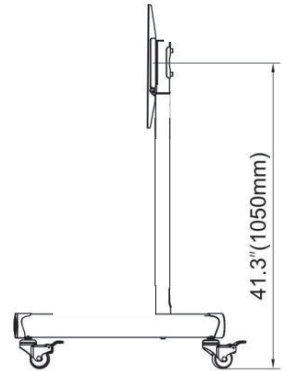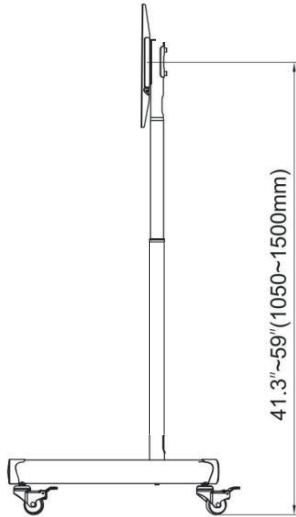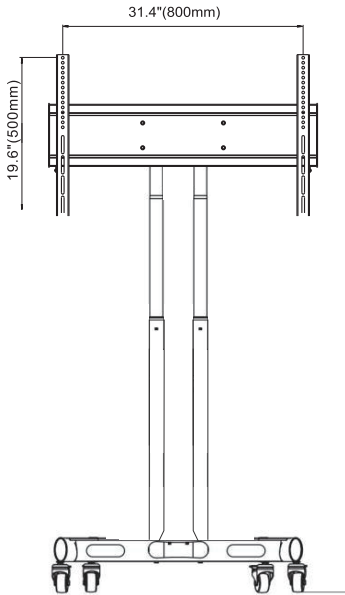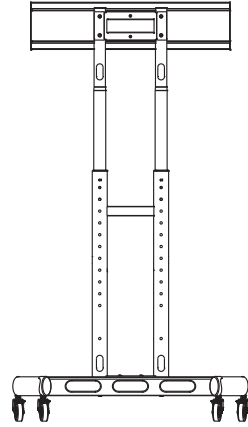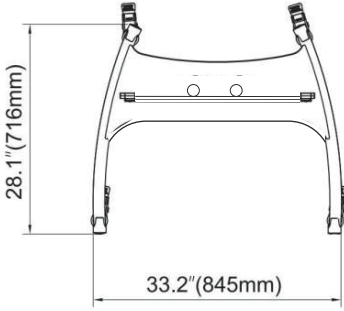


6. Horizontal Stand Assembly Instruction

- Step1 : 1.The base and telescopic mounting installation is same as previous example
2.Connect and fixed Horizontal L-form mounting kit 'A' with the telescopic mounting, as

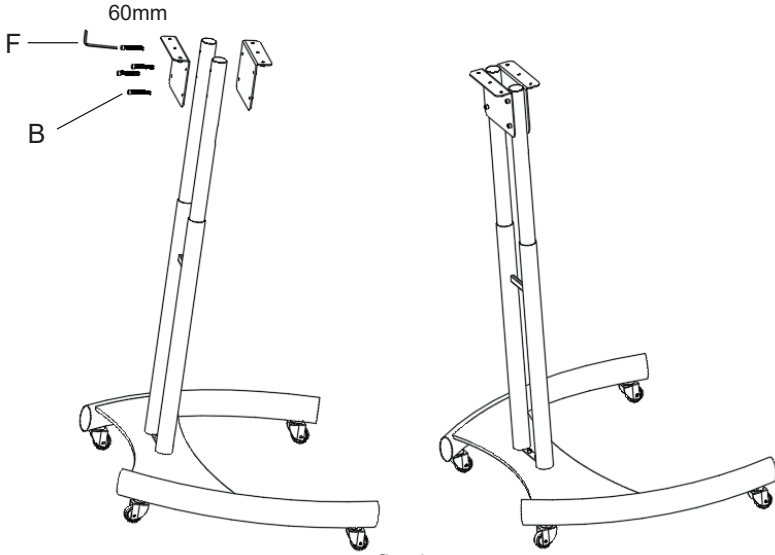

Step2

- 1.Fix panel(with horizontal metal structure screwed to) with bolt as following

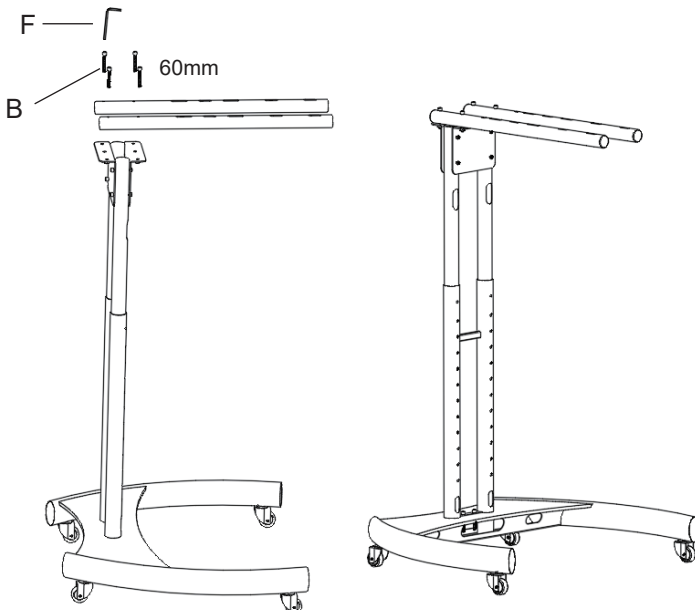

7. Screw Horizontal Metal Structure To Panel

1. Screw horizontal metal structure to panel as following

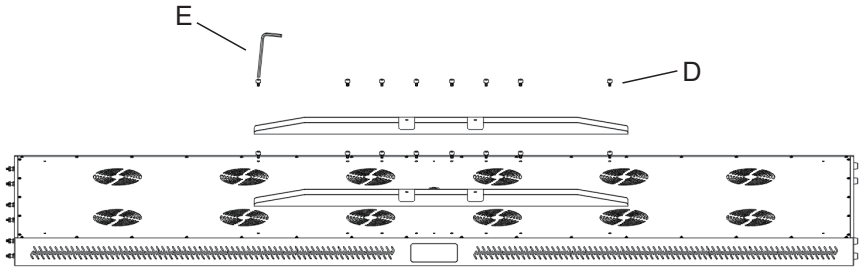


fig-1

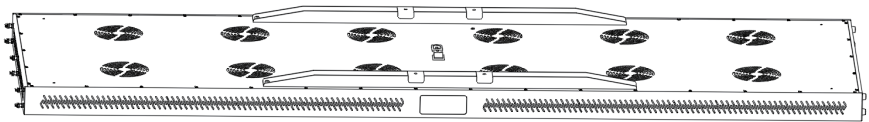


fig-2

2. Fix panel(with horizontal metal structure screwed to) with bolt as following

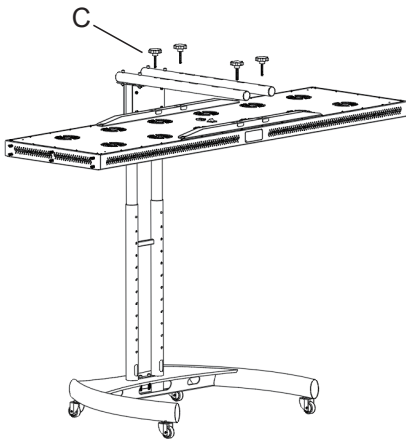


fig-3

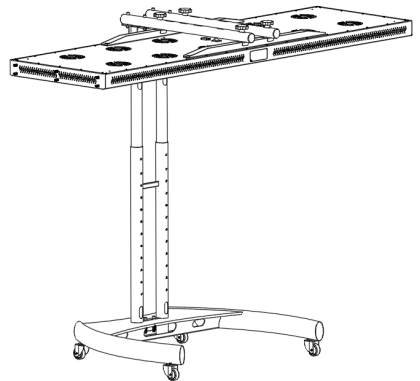


fig-4

8. Contact us

**We reserve the right to modify or alter instructions.
No modification or alteration without formal notice.**

- This mobile stand is used for our Red Light Therapy panels only.
- Please ensure the stand is placed on level floor.
- Please ensure the wheels are locked after installation.
- Children should not use unless approved by physician.

If you have any questions about assembly for this mobile stand, please feel free to contact us for online installationservice.

