# MAGNET-PERFECT™ SMART CERVICAL MASSAGER

USER MANUAL

INTELLIGENT IMPULSE MASSAGE

magnetinnov  
MEDICAL INNOVATION

# TABLE DES MATIÈRES

FR

EN

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FR

EN

Precautions.....	13
Technical features .....	16
Product Diagram (Rechargeable model).....	17
Product Diagram (Battery operated model).....	18
Operational Mode .....	19
Tips .....	20
Daily maintenance .....	20
Troubleshooting .....	20
Warning .....	22

# Precautions

FR  
EN

## Danger

- Never use the product with the following electronic medical devices:
  1. Pacemakers and other implantable electronic medical devices;
  2. Artificial hearts, lungs, and other electronic medical devices that maintain vital functions
  3. Electrocardiographs and other electronic medical devices.

\* These electronic medical devices can cause malfunctions and turn out to be fatal.

## Warning

- If a high-frequency surgical device is applied at the same time as the massager, the massager's electrode could cause burns and damages to the device.
- If the product is used near a shortwave or microwave therapy (1M), its output may be unstable.
- The electrode should not be in direct contact with wounds or used on scars during treatment.
- To avoid the increasing the risk of cardiac fibrillation not place the electrode near your heart or chest while you use the product. Do not use it over your neck, head, mouth, genitals, or on rashes, and do not try to place it on the soles of your feet.

## Warning

- Product use is recommended for only one user. If the product is being used by multiple people, clean it before use.
- Attach the electrode evenly and firmly on your skin while using the massager.
- If you don't feel well during or after use, stop using the product immediately and consult your doctor.
- To avoid accidents or physical illness, consult your doctor before using the product if any of the following cases apply to you:
  1. If you are receiving medical treatment and are dizzy.
  2. If you have heart or nerve problems.
  3. If your body temperature is over 100.4° F
  4. Infectious diseases.
  5. Rashes or skin sensitivity.
  6. DO NOT use while driving.
  7. Elderly people should use the product under the instruction of their doctor.
  8. DO NOT use the product after drinking or taking sleeping pills.
  9. NEVER try to inspect or open the product.
  10. DO NOT put heavy objects the product.
  11. DO NOT continue to use the product if it's damaged.
  14. DO NOT drop the product from a high place.
  12. DO NOT submerge or drop the device in water or any other liquid.
  13. NEVER use pins or other metal clips with this device.

# Precautions

---

FR

## Notes

EN

- If your skin or body feel strange, consult your doctor and follow their instructions.
- To avoid electric shocks, cut off the power supply before moving the device if you want to move the massager or put it on someone else during the massage.
- To avoid accidents, keep out of reach of children.
- Do not use while sleeping. Use while sleeping could cause bodily dysfunction, or the massager could move to other parts of the body and cause a loss of consciousness.
- Do not put the electrode in the same place for a long time. Doing so can cause skin inflammation.
- To avoid electric shocks, do not allow metal objects to come into contact with the electrode during the massage (necklaces, belts, etc.)
- Do not use the massager in the same place for an extended period of time (more than 30 minutes). Prolonged use can cause muscle fatigue and discomfort in the massaged area.
- Wear the product correctly by putting the device on the cervical spine in a symmetrical manner. Improper use of this product could cause accidents or physical illness.


# Technical Features

FR  
EN


Model	ALE-20
Rated voltage	5V ===
Rated power	5W
Automatic timer	15min
Administrative regulation	GB4706.1-2005

# Product diagram (rechargeable product)



FR  
EN

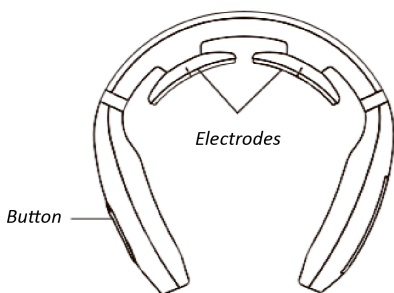Patch connection Cable




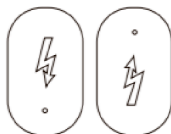
The external pulsing patch can be used to massage other parts of the body

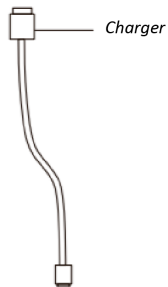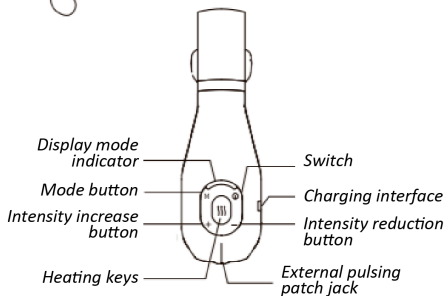
# Product Diagram (battery-operated)

FR

EN
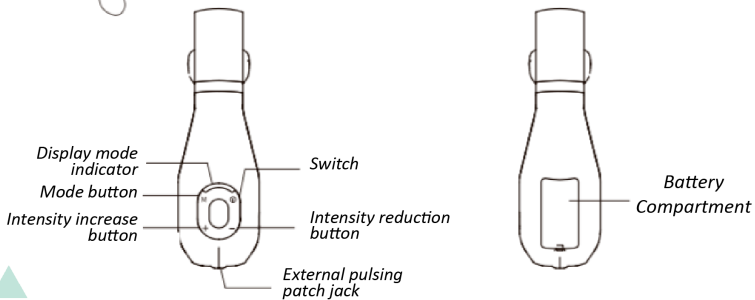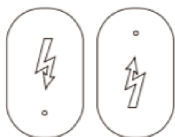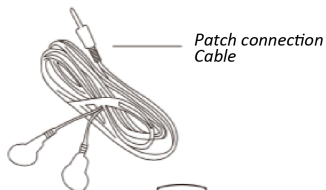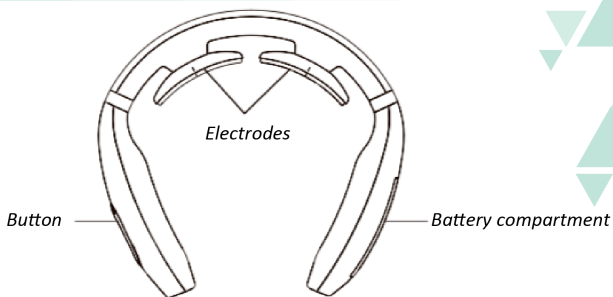
## Massager application

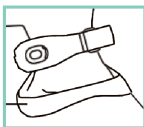Magnetinnov  
MEDICAL INNOVATION

# Operational Mode

FR  
EN

1. Dampen the skin on your cervical spine with a moist towel or a water-based gel before each use. Then put the product around your neck, bring the electrode to the skin on your neck, adjust the product position, and apply power to the device, then adjust the device to the appropriate intensity and massage mode.

2. Press the "0" button to turn the product on (the light indicator will turn on). You'll be in default mode with the device in default position (on speed level 0.)

3. "M" represents the mode button. Press the button to change modes. Each time you press the M Mode button, the button will beep. There are a total of six modes to choose from that will cyclically rotate.

4. Press the "+" button to increase the intensity. Every time the button is pressed, the intensity increases by one level. There are a total of 15 levels. (0, 1, 2, ..., 15)

5. Press the "-" button decrease the intensity. Every time the button is pressed, the intensity decreases by one level. There are a total of 15 levels. (0, 1, 2, ..., 15) (note : 15, 14, 13..., 1).

6. Press the "///" button to activate the thermal function. Select 1 for low heat, 2 for medium heat, and 3 for high heat, then press the thermal function button to turn the function off

7. You can use the patches on painful areas of your body.

8. The machine is programmed to provide massages for 15 minutes. The system will turn off 15 minutes after it is turned on.

9. The perfect magnet neck massager modes are:

**Beating mode (Red LED) :** This mode uses magnetic frequency to reach the joints and soothe recurring pain.

**Acupuncture mode (blue LED) :** The neck massager will massage acupuncture points.

**Kneading mode (green LED) :** allows for better recovery and better blood circulation sanguine.

**Automatic mode (multi-color LED) :**

Automatic adaptation of different massager features.

We recommend using the automatic mode on loop for best results

## Tips

---

FR  
EN

- The product will automatically turn off after 15 minutes of operation. We recommend no longer than 30 minutes of use per day.
- The product should be stored in a dry, ventilated place after use.
- Do not put other weight on the product.

## Daily maintenance

---

- To clean and maintain the product, turn it off and use a moist dish towel to clean the electrode. Use a bit of medical-grade rubbing alcohol if the product is particularly dirty (75% alcohol concentration) to wipe and disinfect it.
  - Use a neutral detergent and soft cloth for daily maintenance; do not use corrosive detergent, and do not put the device in water to soak or clean it.
  - To prevent short-circuiting and electric shocks during use, make sure that water and liquids do not penetrate the inside of the product while washing.

## Troubleshooting

---

- If you encounter problems during product use, please refer to the following suggestions and solutions.
  - If the product is still out of order or cannot be used, cut off the power supply.

Situation	Likely cause/ solution
<b>Device does not turn on</b>	<p>The internal rechargeable battery hasn't been charged enough.</p> <p>Check to see if the on button has been pressed.</p>
<b>Device shows normal status but no feeling</b>	<p>Is the electrode sheet in direct contact with your skin? If not, adjust the sheet's position so it's in direct contact with your skin.</p> <p>Is there any hair, jewelry, or clothing, etc. caught between the electrode patch and your skin? If so, remove it.</p> <p>Is the resistance at level 0? You can increase the massage intensity.</p>
<b>No battery or battery status display</b>	<p>Check the USB connection to make sure it's properly connected to the device and the adaptor.</p>
<b>Tingling</b>	<p>Is the electrode patch properly adjusted on your skin? Adjust the position to hold your skin taut. Is the massage duration too long? We recommend 15 minute massages at a time, twice a day.</p> <p>Is the massage intensity too strong? The intensity can be adjusted. Are you allergic to the electrode? If you suffer from severe allergies, you may need to desensitize yourself to it at first. Is it cold outside or is your skin dry? Use a damp towel to hydrate your skin.</p> <p>Do you have any open wounds? Do not place the electrode on open wounds.</p>
<b>Very weak stimulus</b>	<p>Has the massager electrode been properly applied to the neck? Adjust the electrode's position so it's adapted to your skin and adjust the intensity accordingly.</p> <p>* Is the electrode sheet dirty or oily? Clean the electrode surface.</p> <p>Adjust to increase the massage intensity.</p>

**Avertissement:**

This symbol on a device or its packaging means: do not throw away electronic devices in unapproved municipal landfills: use separate collection facilities. Contact your local authority to obtain information on available collection systems in your area. If electronic devices are thrown away in landfills, dangerous substances can seep into underground water and end up in the food chain which is dangerous to your health and wellbeing.

## Contact Us

If you have any questions, feel free to contact us, and we'll do the best we can to address your questions or concerns.

[contact@magnet-innov.com](mailto:contact@magnet-innov.com)

Magnet-Innov

[magnet\\_innov](https://www.instagram.com/magnet_innov)

**magnet**innov  
MEDICAL INNOVATION

*Thank you  
for purchasing the  
Magnet-perfect™ smart  
cervical massager*

Discover now our best products

**Magnet-Relax 2.0**  
**The Smart Eye**  
**Massager**

**Magnet-Zen™**  
**The Shiatsu**  
**Massage Pillow**

**Magnet Relax™**  
**The Massage**  
**Mask**

 Magnet-Innov

 magnet\_innov

**magnet**innov  
MEDICAL INNOVATION  
[www.magnet-innov.com](http://www.magnet-innov.com)

 Scan me  
for a gift