
# MAGNET-ZEN™ SHIATSU MASSAGE CUSHION (PROFESSIONAL EDITION)

## USER MANUAL

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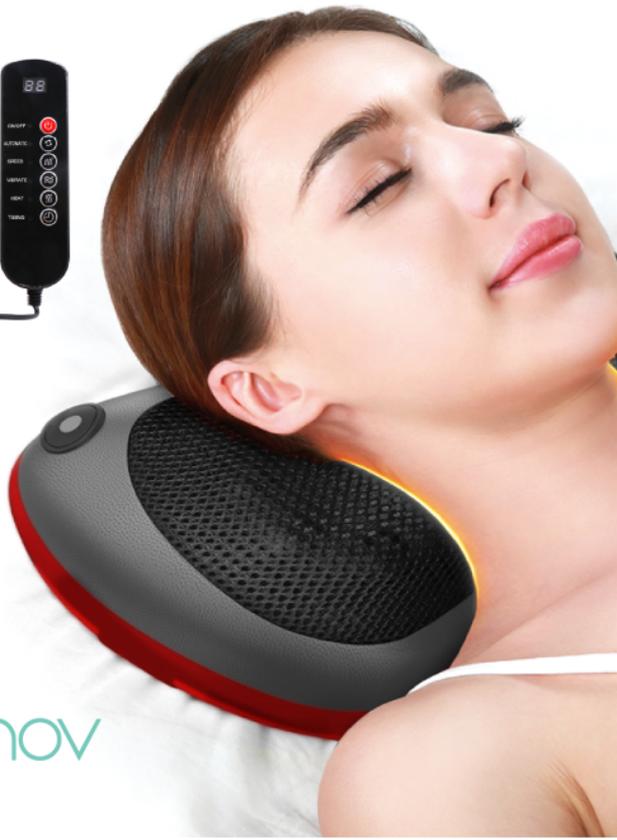**magnet**innov  
MEDICAL INNOVATION

# TABLE OF CONTENTS

Precautions and Recommendations .....	11
Product maintenance and care .....	12
Use and Warnings .....	13
Product Features .....	17
Technical Specifications .....	18
Warning .....	19

# Precautions and Recommendations

FR  
EN

**To keep the massage cushion in good condition, use the product properly and make sure to be careful of the following:**

**To avoid electric shocks, follow the instructions below :**

- 1.** Unplug the device after use or before washing.
- 2.** Do not use the massage cushion in the water or in a wet place.
- 3.** If the massage cushion falls into water, disconnect the power supply immediately and unplug the device before retrieving it.
- 4.** Do not put the massager in water or any other liquids.
- 5.** Never use the device in the bath or shower.
- 6.** Never put pins or any other metal in the device.

## Product maintenance and care:

FR

EN

1. Do not use the massager continuously for an extended period of time at the beginning. You can increase the time of use after you get used to the device.
2. In order to avoid overworking the massage cushion, do not stand or sit on it.
3. Do not exceed 30 minutes of use.
4. Never use detergent or abrasive products to clean the device. Never submerge the device in liquid cleaner.
5. Turn off the system by putting the button into the "OFF" position then unplug the device. Do not forcefully pull on the power cord.
6. Use a soft cotton cloth and neutral detergent to clean the massager.
7. Do not use or clean the massage cushion in a wet or hot place.
  8. Do not stomp on, drop, or destroy the cushion.
9. Do not forcefully pull the power cord out from the massage device.
10. Do not use the massage cushion near a stove or other heat source such as under a pillow or mat.

**We recommend that you drink at least two cups of hot water before and after using the massage cushion to boost your metabolism.**

# Use and Warnings

FR

## **Do not use the massage cushion in the following cases:**

EN

1. If you're suffering from trauma or an ulceration on a large area of your body.
2. If you're under medical care or are suffering from cancer or severe heart disease, consult your doctor and use the massage cushion according to the doctor's instructions.
3. Children or pregnant women
4. Suffering from mental illness
5. Never use the device for any other purpose other than its intended use. This massage cushion is intended for solely for indoor home use. The device is not intended to replace medical consultations.

## **Massager assembly and care of the massager**

1. Use accessories or pieces made by the manufacturer of this massage cushion .
2. Do not detach or disassemble the massage tool and do not attempt to repair it yourself without prior authorization from us. Attempts to repair it yourself will make you ineligible to benefit from our warranty and customer service.
3. For your safety and the lifespan of this massage cushion, do not open the zipper. We will not be held responsible for the consequences of this behavior.

## Plugging in instruction

1. Check if the voltage is correct and use the AC power supply.
2. Do not plug in the device with a wet hand.
3. Keep the panel, the plug, and the cord dry.
4. If the AC outlet isn't stable, do not plug in.
5. If the cord is broken, send it to our customer to receive a new one.
6. Do not wrap the cord around the massage cushion.
7. In order to avoid an electric accident, make sure that the plug is plugged in correctly. Never leave the device unattended while it's plugged in.

## Remote Control Use

1. Press the ON/OFF button once to turn the device on or off.
2. Press the automatic button once to change the massage head rotation direction.
3. Press the speed button once to increase or decrease the massage speed.
4. Press the vibration button once to modify the vibration intensity.
5. Press the heat button once to increase or decrease the heat diffusion.
6. Press the time button once to choose the duration of your massage session.

# Use and Warnings

FR  
EN

**We suggest using the massage cushion for 15 minutes at a time in order to keep it in the best conditions of use. Let the device rest for fifteen minutes between each use. In order to prevent the device from overheating, do not use continuously. The device will stop automatically after 15 minutes of use.**

## **Appropriate product use :**

1. At home, at the office, or in the car.
2. Unplug the device and keep it in a dry place if you aren't using it for a long time.
3. Do not use the massage cushion an hour before or after each meal.
4. The massage cushion does not replace medical equipment.

## **Problems and solutions :**

1. The massage cushion does not work check the outlet and plug in the power supply.
2. The massage cushion stops working: the device overheats after use that lasts too long and shuts off automatically: turn of the device and plug it in, then use it again after at least an hour.
3. The device makes noise during use: this is normal due to the cushion massage head.

4. The intensity is too strong or weak: ensure that the massage cushion is at a good angle and choose the right massage intensity.
5. The intensity isn't correctly adapted at the beginning: you should feel better after a certain time of use.
6. The fabric is damaged: cease use and change it out for new netting.
7. Cease use and turn off the massage cushion immediately if you do not feel well during use.

 **Attention :**

1. Adult supervision is required if children want to use the massage cushion.
2. The massage cushion is not a toy.
3. The device has a hot surface. Those sensitive to heat must be careful while using it. This warning is also applicable to those who suffer from hypersensitivity or sensory impairments.
4. Do not use on any body part that is burned, swollen, or suffering from inflammation, open wounds, rashes, or body parts with recent scarring without consulting your doctor before use.

FR

EN

- Home use: plug in the adapter and pillow then plug the adaptor into the outlet.
- Use in the car: plug the pillow into the plug, then the massage cushion will turn itself on.
- **How to use the control switch:**  
Press the first speed button to activate the kneading and heating function, press the second speed to turn the heating function off, press the third speed to power off.  
Unplug the device after every use in order to maintain its high performance,
- In order to the massage cushion in the best conditions for use, we recommend that you use the product for no more than 15 minutes at a time and then use it again at least 30 minutes later.

## Product Features

---

This product's original design with several massage heads can produce a sensation of warmth and the effects of a massage at the same time to deeply soothe muscle fatigue. The cushion also fosters blood circulation, pain relief, and relaxation, which allows you to enjoy complete comfort.

# Technical Specifications

FR

EN

Rated voltage : DC12V FR

Rated current : 2A.

Rated power : 24W EN.

Rated duration : 15 minutes.

Security level. 

Rated frequency : 50Hz.

AC adaptor voltage: AC110-240V

**Avertissement:**

This symbol on a device or its packaging means: do not throw away electronic devices in unapproved municipal landfills: use separate collection facilities. Contact your local authority to obtain information on available collection systems in your area. If electronic devices are thrown away in landfills, dangerous substances can seep into underground water and end up in the food chain which is dangerous to your health and wellbeing.

## Contact Us

If you have any questions, feel free to contact us, and we'll do the best we can to address your questions or concerns.

[contact@magnet-innov.com](mailto:contact@magnet-innov.com)

Magnet-Innov

[magnet\\_innov](https://www.instagram.com/magnet_innov)

Magnet<sup>+</sup>innov  
MEDICAL INNOVATION

*Thank you  
for purchasing the  
Magnet-Ten™  
shiatsu massage cushion*

Discover now our best products

**Magnet-Relax 2.0  
Smart Eye  
Massager**

**Magnet-Perfect™  
The Therapeutic  
Neck Massager**

**Magnet-Innov©  
The Eye  
Massager**

 Magnet-Innov

 magnet\_innov

**MagnetInnov**  
MEDICAL INNOVATION  
[www.magnet-innov.com](http://www.magnet-innov.com)

 Scan me  
for a gift